


Assay Class: High Sensitivity DNA Assay
Data Path: C:\...r1_High Sensitivity DNA Assay_DE34903152_2026-01-05_001.xad

Created: 1/5/2026 9:41:04 AM
Modified: 1/5/2026 10:31:04 AM

Electrophoresis File Run Summary

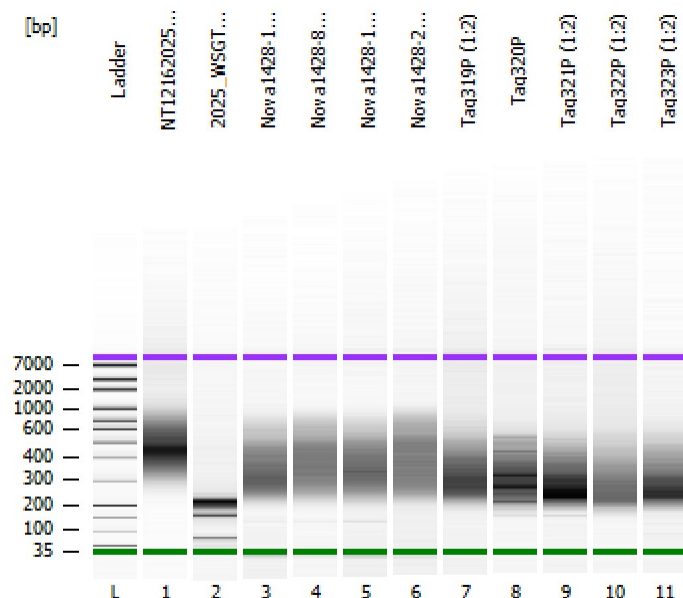

Instrument Information:












Instrument Name: DE34903152 Firmware: C.01.069
Serial#: DE34903152 Type: G2938C

Assay Information:

Assay Origin Path: C:\Program Files (x86)\Agilent\2100 bioanalyzer\2100 expert\assays\dsDNA\High Sensitivity DNA.xsy
Assay Class: High Sensitivity DNA Assay
Version: 1.03
Assay Comments: Copyright © 2003-2010 Agilent Technologies

Chip Information:

Chip Lot #:
Reagent Kit Lot #:
Chip Comments:

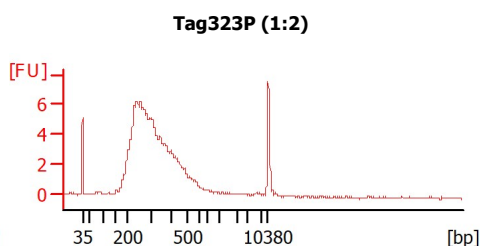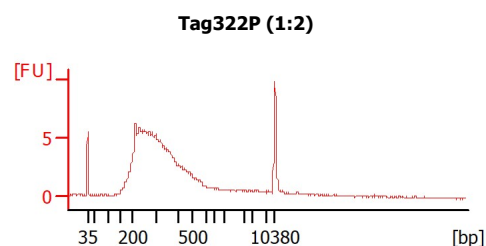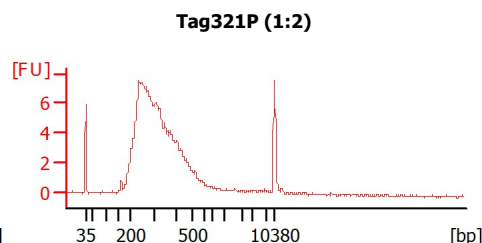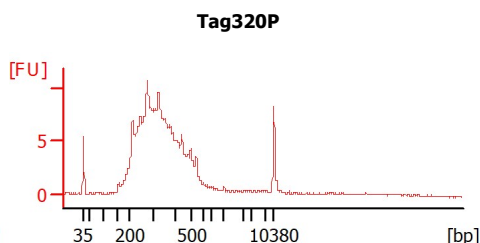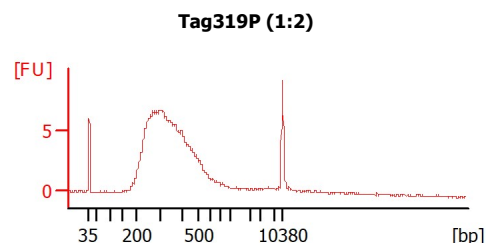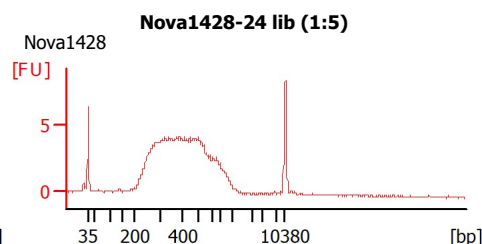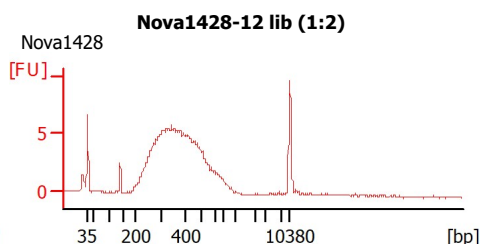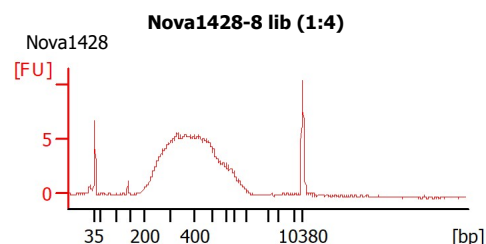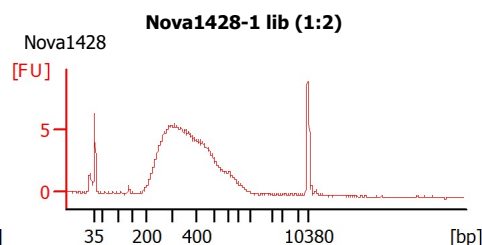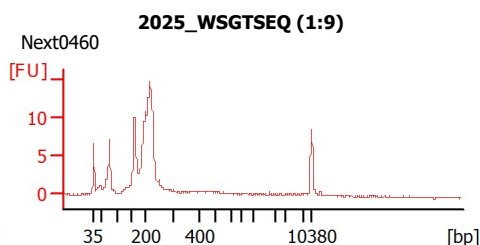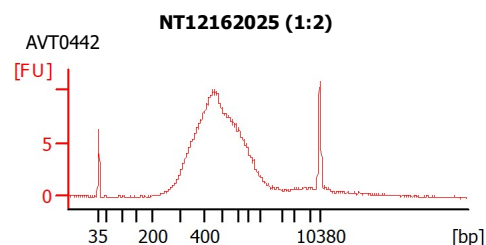

Assay Class: High Sensitivity DNA Assay
Data Path: C:\...r1_High Sensitivity DNA Assay_DE34903152_2026-01-05_001.xad

Created: 1/5/2026 9:41:04 AM
Modified: 1/5/2026 10:31:04 AM

Electrophoresis File Run Summary (Chip Summary)

Sample Name	Sample Comment	Rest. Digest	Status	Observation	Result Label	Result Color
NT12162025 (1:2)	AVT0442	☐	✓			
2025_WSGTSEQ (1:9)	Next0460	☐	✓			
Nova1428-1 lib (1:2)	Nova1428	☐	✓			
Nova1428-8 lib (1:4)	Nova1428	☐	✓			
Nova1428-12 lib (1:2)	Nova1428	☐	✓			
Nova1428-24 lib (1:5)	Nova1428	☐	✓			
Tag319P (1:2)		☐	✓			
Tag320P		☐	✓			
Tag321P (1:2)		☐	✓			
Tag322P (1:2)		☐	✓			
Tag323P (1:2)		☐	✓			
Ladder		☐	✓			

Chip Lot #**Reagent Kit Lot #****Chip Comments :**

Assay Class: High Sensitivity DNA Assay
Data Path: C:\...r1_High Sensitivity DNA Assay_DE34903152_2026-01-05_001.xad

Created: 1/5/2026 9:41:04 AM
Modified: 1/5/2026 10:31:04 AM

Electrophoresis Assay Details

General Analysis Settings

Number of Available Sample and Ladder Wells (Max.) : 12
Minimum Visible Range [s] : 32
Maximum Visible Range [s] : 138
Start Analysis Time Range [s] : 33
End Analysis Time Range [s] : 137.5
Ladder Concentration [pg/μl] : 1950
Uses Standard Area for Ladder Fragments
Lower Marker Concentration [pg/μl] : 125
Upper Marker Concentration [pg/μl] : 75
Used Upper Marker for Quantitation
Standard Curve Fit is Point to Point
Show Data Aligned to Lower and Upper Marker

Integrator Settings

Integration Start Time [s] : 33.05
Integration End Time [s] : 137
Slope Threshold : 0.8
Height Threshold [FU] : 5
Area Threshold : 0.1
Width Threshold [s] : 0.6
Baseline Plateau [s] : 0.5

Filter Settings

Filter Width [s] : 0.5
Polynomial Order : 4


Ladder

Ladder Peak	Size	Area
1	35	160
2	50	210
3	100	208
4	150	221
5	200	242
6	300	270
7	400	305
8	500	306
9	600	336
10	700	321
11	1000	366
12	2000	413
13	3000	411
14	7000	400
15	10380	214

Assay Class: High Sensitivity DNA Assay
 Data Path: C:\...r1_High Sensitivity DNA Assay_DE34903152_2026-01-05_001.xad

Created: 1/5/2026 9:41:04 AM
 Modified: 1/5/2026 10:31:04 AM

Electropherogram Summary

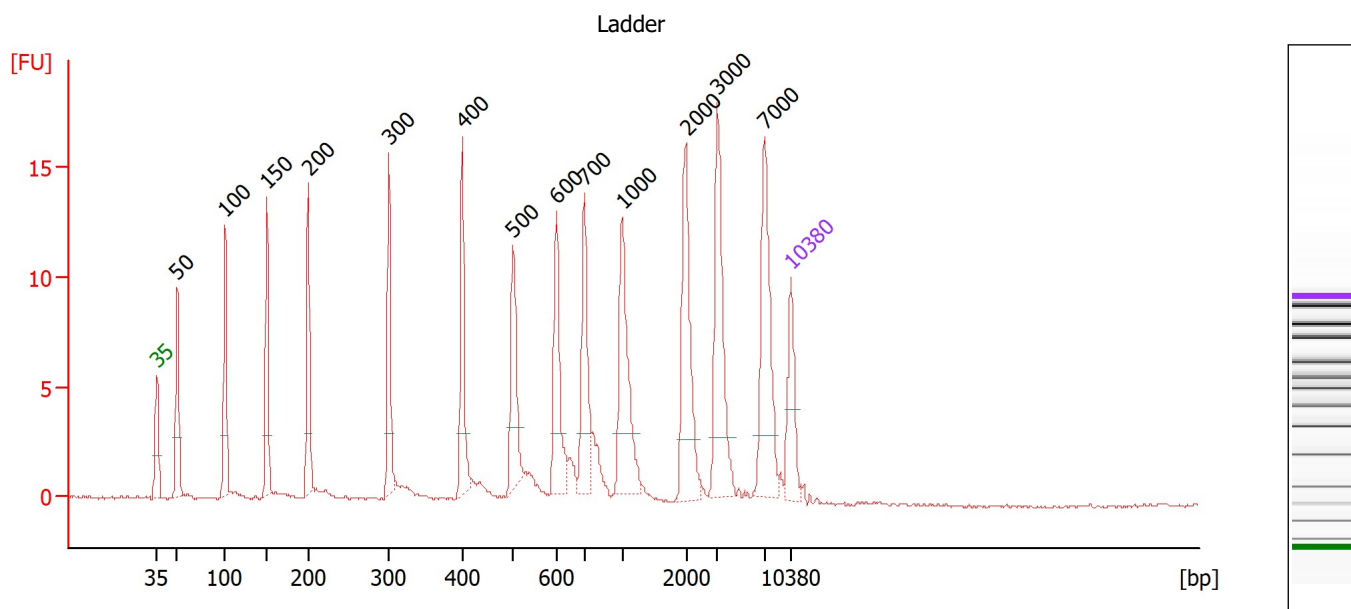

Overall Results for Ladder

Noise: 0.1


Peak table for Ladder

Peak	Size [bp]	Conc. [pg/μl]	Molarity [pmol/l]	Observations
1	35	125.00	5,411.3	Lower Marker
2	50	150.00	4,545.5	Ladder Peak
3	100	150.00	2,272.7	Ladder Peak
4	150	150.00	1,515.2	Ladder Peak
5	200	150.00	1,136.4	Ladder Peak
6	300	150.00	757.6	Ladder Peak
7	400	150.00	568.2	Ladder Peak
8	500	150.00	454.5	Ladder Peak
9	600	150.00	378.8	Ladder Peak
10	700	150.00	324.7	Ladder Peak
11	1,000	150.00	227.3	Ladder Peak
12	2,000	150.00	113.6	Ladder Peak
13	3,000	150.00	75.8	Ladder Peak
14	7,000	150.00	32.5	Ladder Peak
15	10,380	75.00	10.9	Upper Marker

Assay Class: High Sensitivity DNA Assay
Data Path: C:\...r1_High Sensitivity DNA Assay_DE34903152_2026-01-05_001.xad

Created: 1/5/2026 9:41:04 AM
Modified: 1/5/2026 10:31:04 AM

Electropherogram Summary Continued ...

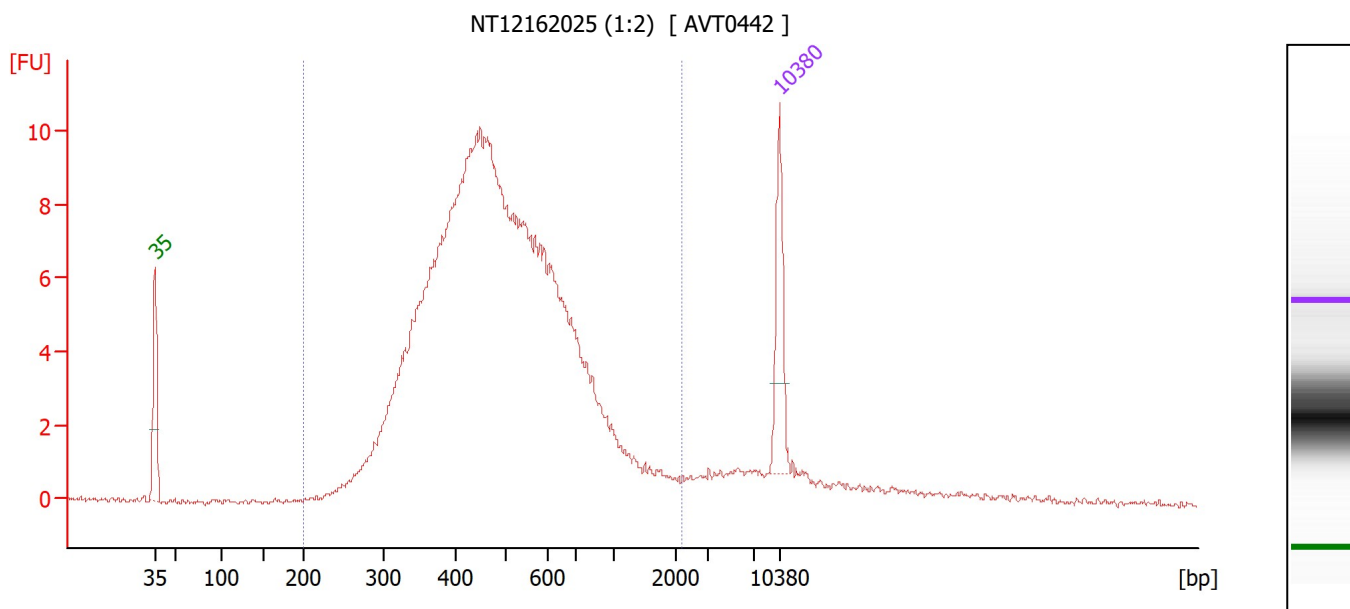

Overall Results for sample 1 : NT12162025 (1:2)

Number of peaks found: 0 Corr. Area 1: 198.3
Noise: 0.1

Peak table for sample 1 : NT12162025 (1:2)

Peak	Size [bp]	Conc. [pg/μl]	Molarity [pmol/l]	Observations
1	35	125.00	5,411.3	Lower Marker
2	10,380	75.00	10.9	Upper Marker


Region table for sample 1 : NT12162025 (1:2)

From [bp]	To [bp]	Corr. Area	% of Total	Average Size [bp]	Size distribution in CV [%]	Conc. [pg/μl]	Molarity [pmol/l]	Color
200	2,213	198.3	93	533	46.6	2,716.75	9,078.4	Blue

Assay Class: High Sensitivity DNA Assay
Data Path: C:\...r1_High Sensitivity DNA Assay_DE34903152_2026-01-05_001.xad

Created: 1/5/2026 9:41:04 AM
Modified: 1/5/2026 10:31:04 AM

Electropherogram Summary Continued ...

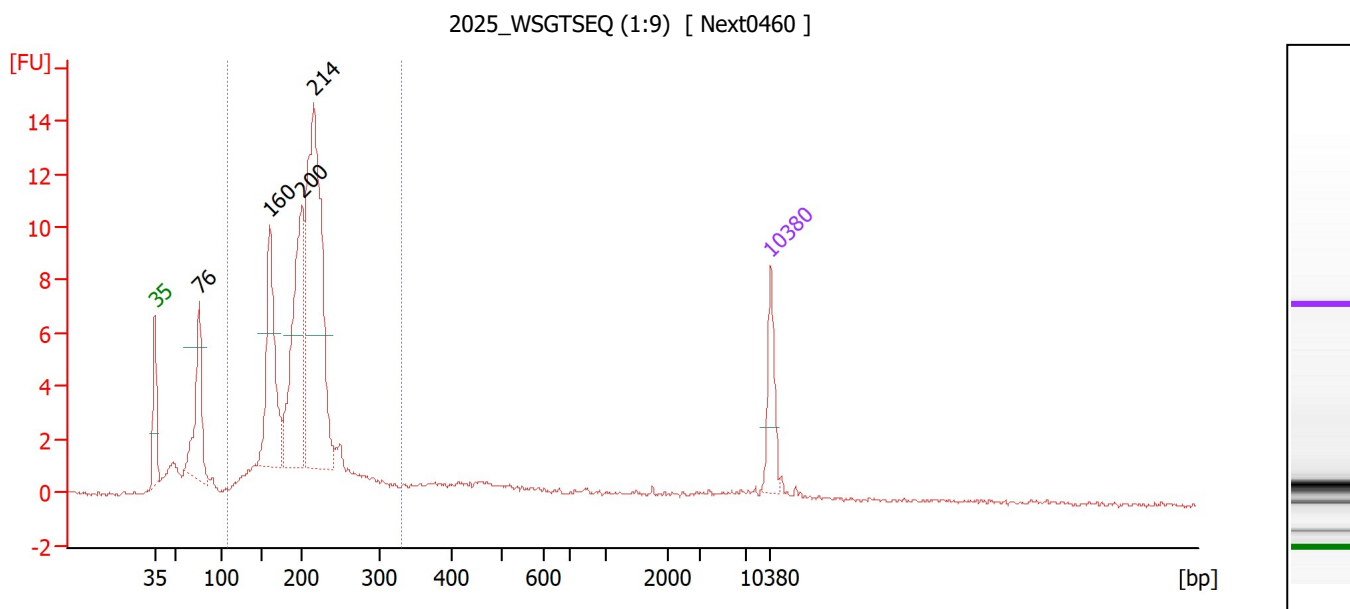

Overall Results for sample 2 : 2025_WSGTSEQ (1:9)

Number of peaks found: 4 Corr. Area 1: 107.4
Noise: 0.1

Peak table for sample 2 : 2025_WSGTSEQ (1:9)

Peak	Size [bp]	Conc. [pg/μl]	Molarity [pmol/l]	Observations
1	35	125.00	5,411.3	Lower Marker
2	76	249.33	4,947.5	
3	160	356.29	3,377.6	
4	200	456.65	3,458.0	
5	214	827.65	5,848.9	
6	10,380	75.00	10.9	Upper Marker


Region table for sample 2 : 2025_WSGTSEQ (1:9)

From [bp]	To [bp]	Corr. Area	% of Total	Average Size [bp]	Size distribution in CV [%]	Conc. [pg/μl]	Molarity [pmol/l]	Color
107	331	107.4	73	203	16.9	2,230.93	17,236.6	Blue

Assay Class: High Sensitivity DNA Assay
 Data Path: C:\...r1_High Sensitivity DNA Assay_DE34903152_2026-01-05_001.xad

Created: 1/5/2026 9:41:04 AM
 Modified: 1/5/2026 10:31:04 AM

Electropherogram Summary Continued ...

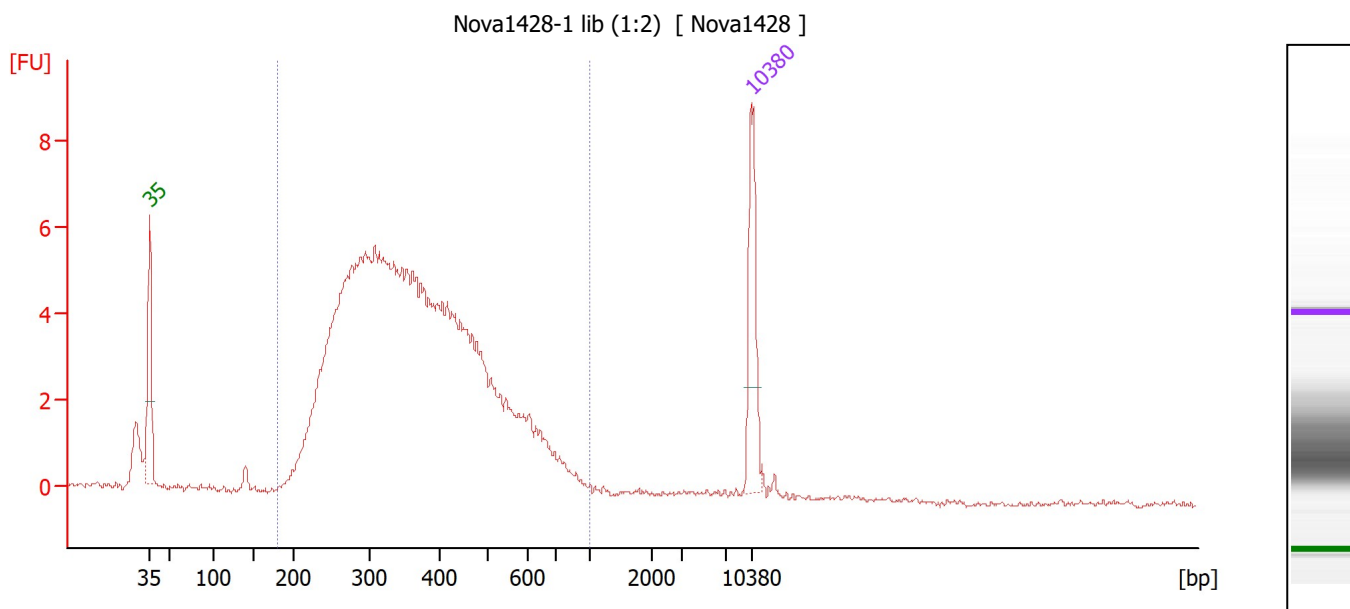

Overall Results for sample 3 : Nova1428-1 lib (1:2)

Number of peaks found: 0 Corr. Area 1: 135.7
 Noise: 0.1

Peak table for sample 3 : Nova1428-1 lib (1:2)

Peak	Size [bp]	Conc. [pg/μl]	Molarity [pmol/l]	Observations
1	35	125.00	5,411.3	Lower Marker
2	10,380	75.00	10.9	Upper Marker


Region table for sample 3 : Nova1428-1 lib (1:2)

From [bp]	To [bp]	Corr. Area	% of Total	Average Size [bp]	Size distribution in CV [%]	Conc. [pg/μl]	Molarity [pmol/l]	Color
181	1,000	135.7	96	380	31.6	2,036.03	9,131.9	Blue

Assay Class: High Sensitivity DNA Assay
 Data Path: C:\...r1_High Sensitivity DNA Assay_DE34903152_2026-01-05_001.xad

Created: 1/5/2026 9:41:04 AM
 Modified: 1/5/2026 10:31:04 AM

Electropherogram Summary Continued ...

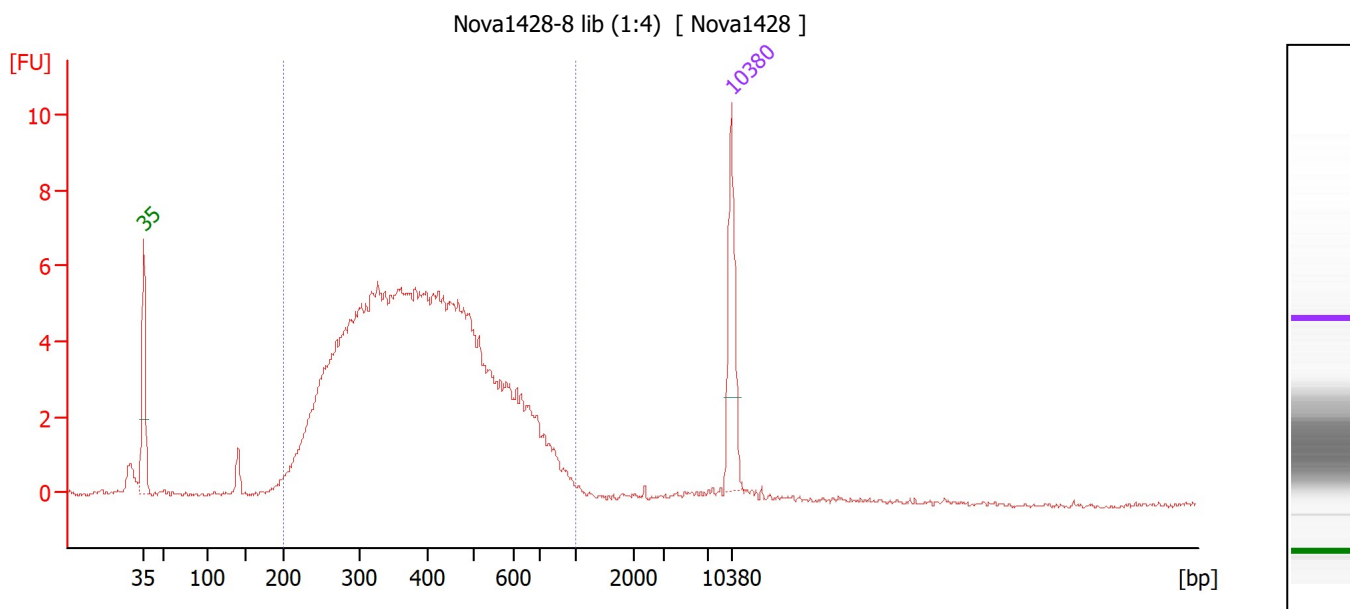

Overall Results for sample 4 : Nova1428-8 lib (1:4)

Number of peaks found: 0 Corr. Area 1: 151.8
 Noise: 0.1

Peak table for sample 4 : Nova1428-8 lib (1:4)

Peak	Size [bp]	Conc. [pg/μl]	Molarity [pmol/l]	Observations
1	35	125.00	5,411.3	Lower Marker
2	10,380	75.00	10.9	Upper Marker


Region table for sample 4 : Nova1428-8 lib (1:4)

From [bp]	To [bp]	Corr. Area	% of Total	Average Size [bp]	Size distribution in CV [%]	Conc. [pg/μl]	Molarity [pmol/l]	Color
200	1,000	151.8	95	413	32.7	2,159.40	9,045.8	Blue

Assay Class: High Sensitivity DNA Assay
Data Path: C:\...r1_High Sensitivity DNA Assay_DE34903152_2026-01-05_001.xad

Created: 1/5/2026 9:41:04 AM
Modified: 1/5/2026 10:31:04 AM

Electropherogram Summary Continued ...

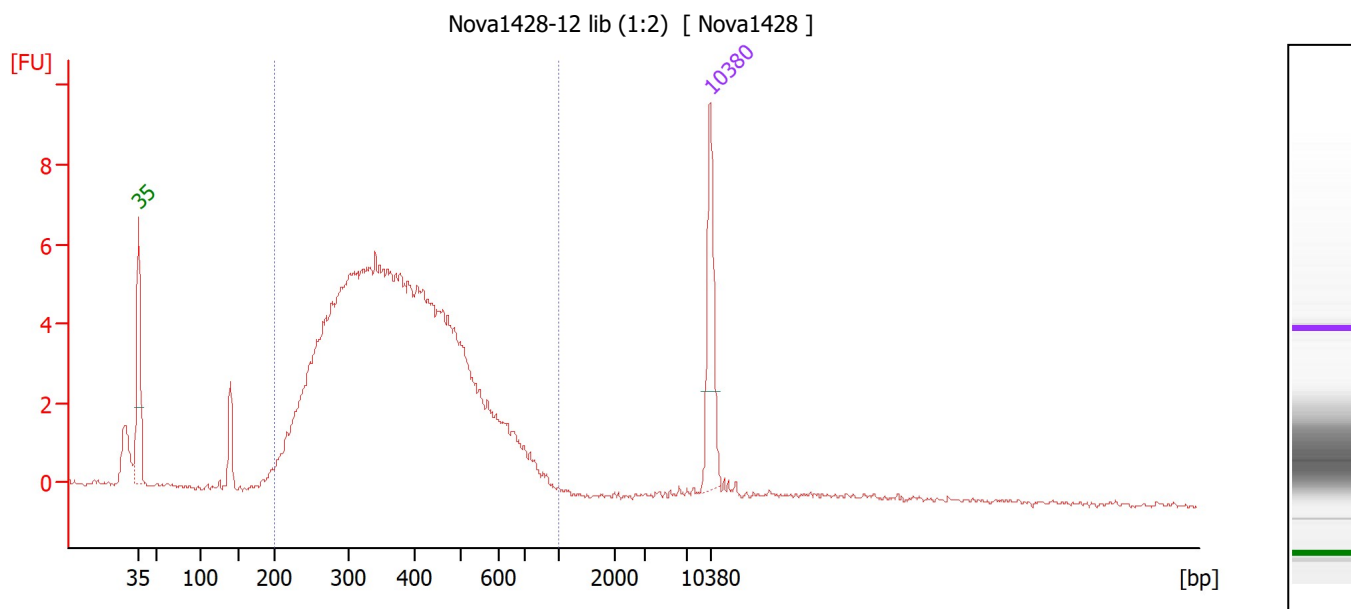

Overall Results for sample 5 : Nova1428-12 lib (1:2)

Number of peaks found: 0 Corr. Area 1: 141.6
Noise: 0.1

Peak table for sample 5 : Nova1428-12 lib (1:2)

Peak	Size [bp]	Conc. [pg/μl]	Molarity [pmol/l]	Observations
1	35	125.00	5,411.3	Lower Marker
2	10,380	75.00	10.9	Upper Marker


Region table for sample 5 : Nova1428-12 lib (1:2)

From [bp]	To [bp]	Corr. Area	% of Total	Average Size [bp]	Size distribution in CV [%]	Conc. [pg/μl]	Molarity [pmol/l]	Color
200	1,000	141.6	95	392	30.5	2,079.08	9,026.7	Blue

Assay Class: High Sensitivity DNA Assay
Data Path: C:\...r1_High Sensitivity DNA Assay_DE34903152_2026-01-05_001.xad

Created: 1/5/2026 9:41:04 AM
Modified: 1/5/2026 10:31:04 AM

Electropherogram Summary Continued ...

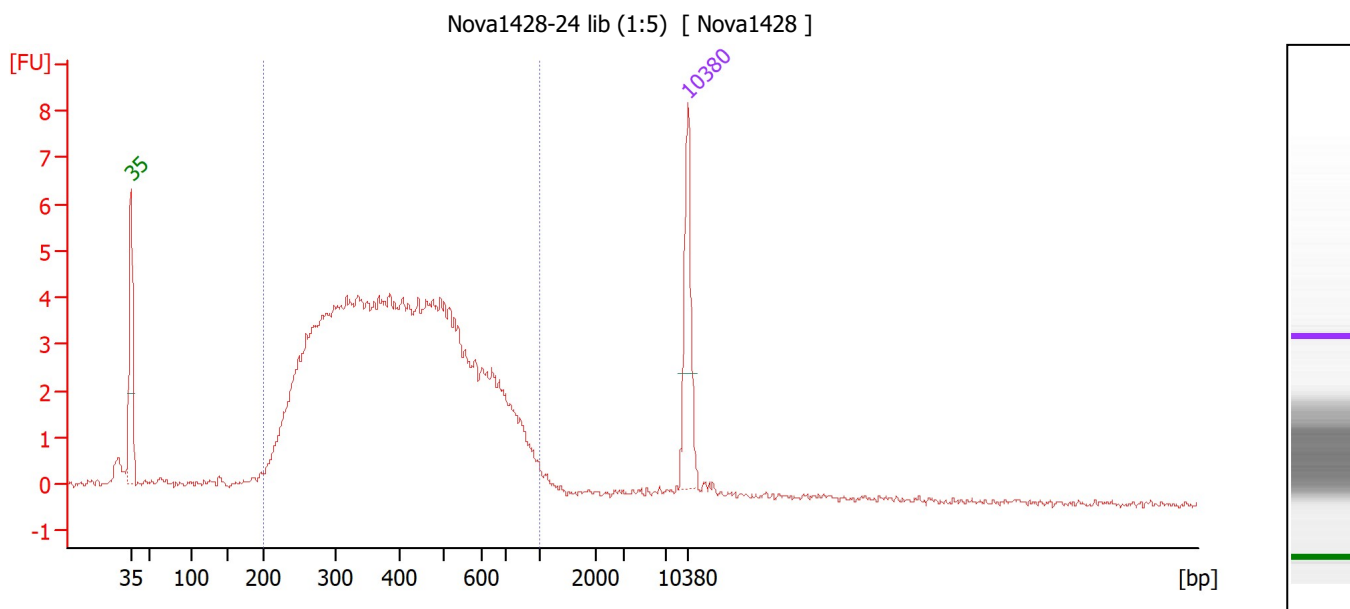

Overall Results for sample 6 : Nova1428-24 lib (1:5)

Number of peaks found: 0 Corr. Area 1: 125.2
Noise: 0.1

Peak table for sample 6 : Nova1428-24 lib (1:5)

Peak	Size [bp]	Conc. [pg/μl]	Molarity [pmol/l]	Observations
1	35	125.00	5,411.3	Lower Marker
2	10,380	75.00	10.9	Upper Marker


Region table for sample 6 : Nova1428-24 lib (1:5)

From [bp]	To [bp]	Corr. Area	% of Total	Average Size [bp]	Size distribution in CV [%]	Conc. [pg/μl]	Molarity [pmol/l]	Color
200	1,000	125.2	95	427	34.9	2,069.17	8,565.7	Blue

Assay Class: High Sensitivity DNA Assay
 Data Path: C:\...r1_High Sensitivity DNA Assay_DE34903152_2026-01-05_001.xad

Created: 1/5/2026 9:41:04 AM
 Modified: 1/5/2026 10:31:04 AM

Electropherogram Summary Continued ...

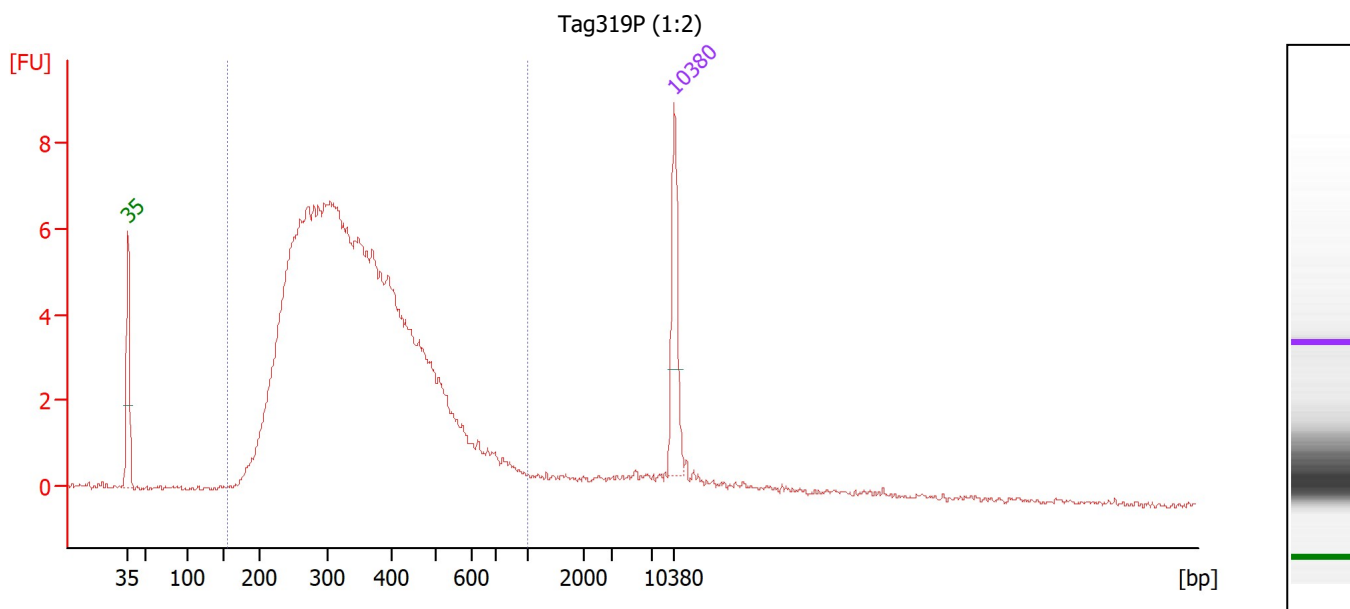

Overall Results for sample 7 : Tag319P (1:2)

Number of peaks found: 0 Corr. Area 1: 157.9
 Noise: 0.1

Peak table for sample 7 : Tag319P (1:2)

Peak	Size [bp]	Conc. [pg/μl]	Molarity [pmol/l]	Observations
1	35	125.00	5,411.3	Lower Marker
2	10,380	75.00	10.9	Upper Marker


Region table for sample 7 : Tag319P (1:2)

From [bp]	To [bp]	Corr. Area	% of Total	Average Size [bp]	Size distribution in CV [%]	Conc. [pg/μl]	Molarity [pmol/l]	Color
157	1,000	157.9	93	362	33.6	2,796.37	13,264.0	Blue

Assay Class: High Sensitivity DNA Assay
 Data Path: C:\...r1_High Sensitivity DNA Assay_DE34903152_2026-01-05_001.xad

Created: 1/5/2026 9:41:04 AM
 Modified: 1/5/2026 10:31:04 AM

Electropherogram Summary Continued ...

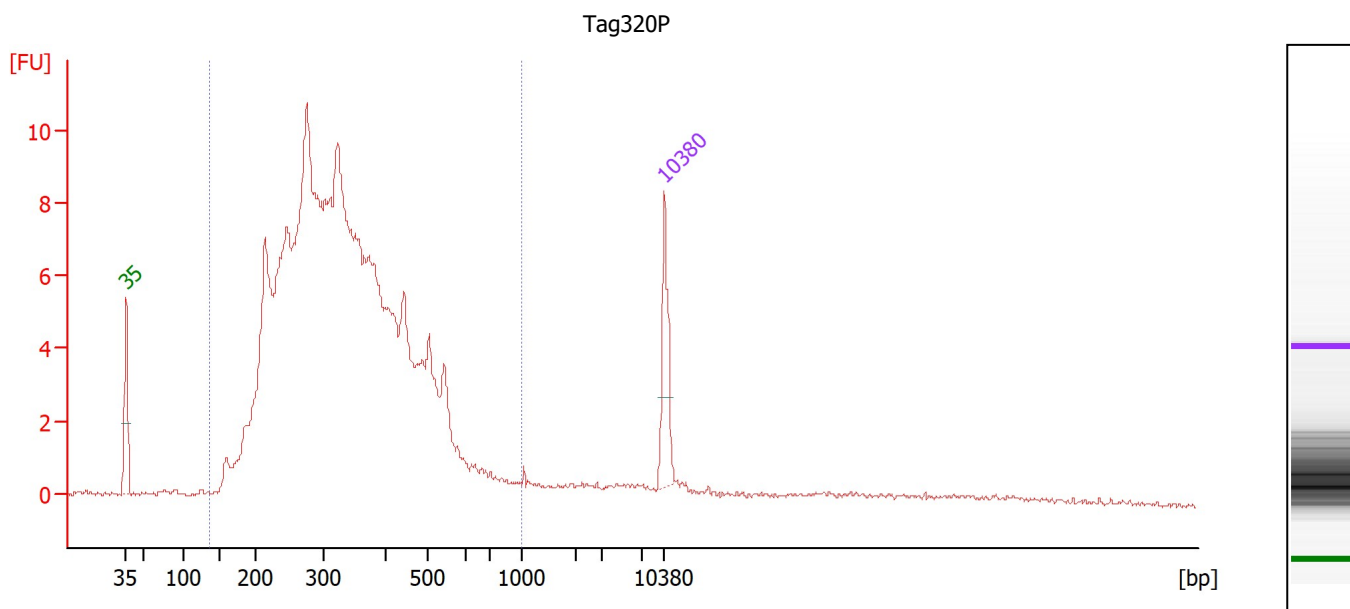

Overall Results for sample 8 : Tag320P

Number of peaks found: 0 Corr. Area 1: 209.1
 Noise: 0.1

Peak table for sample 8 : Tag320P

Peak	Size [bp]	Conc. [pg/μl]	Molarity [pmol/l]	Observations
1	35	125.00	5,411.3	Lower Marker
2	10,380	75.00	10.9	Upper Marker


Region table for sample 8 : Tag320P

From [bp]	To [bp]	Corr. Area	% of Total	Average Size [bp]	Size distribution in CV [%]	Conc. [pg/μl]	Molarity [pmol/l]	Color
137	1,000	209.1	94	345	32.9	4,049.08	20,217.7	Blue

Assay Class: High Sensitivity DNA Assay
Data Path: C:\...r1_High Sensitivity DNA Assay_DE34903152_2026-01-05_001.xad

Created: 1/5/2026 9:41:04 AM
Modified: 1/5/2026 10:31:04 AM

Electropherogram Summary Continued ...

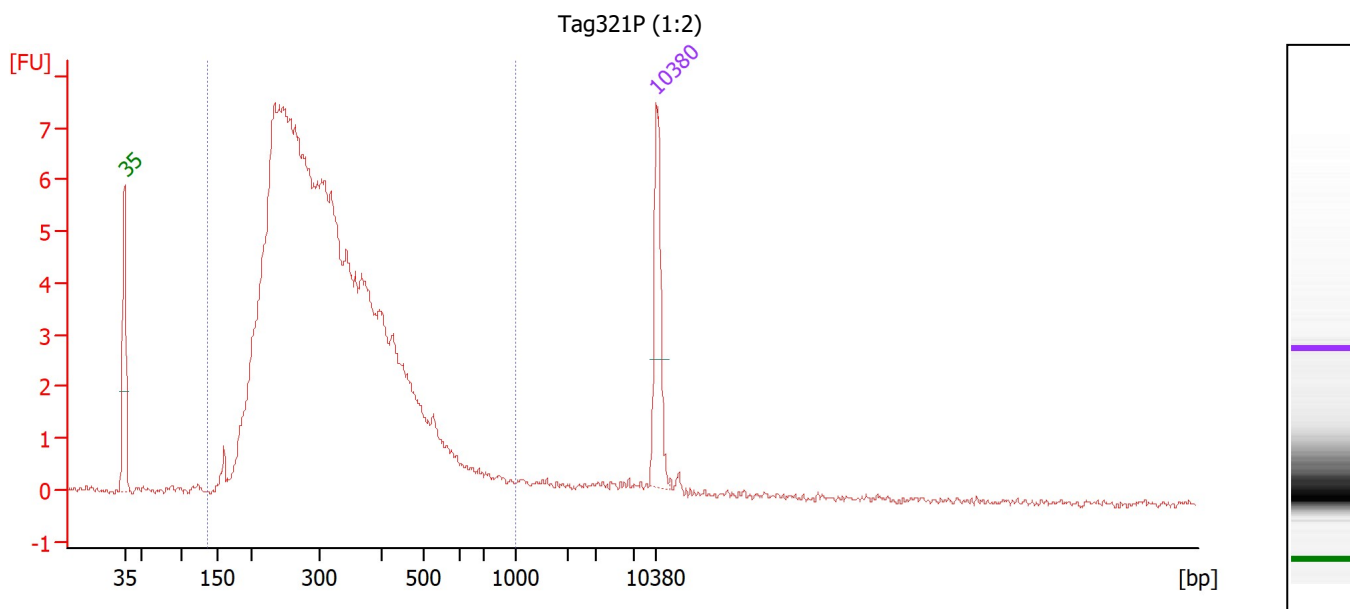

Overall Results for sample 9 : Tag321P (1:2)

Number of peaks found: 0 Corr. Area 1: 150.0
Noise: 0.1

Peak table for sample 9 : Tag321P (1:2)

Peak	Size [bp]	Conc. [pg/μl]	Molarity [pmol/l]	Observations
1	35	125.00	5,411.3	Lower Marker
2	10,380	75.00	10.9	Upper Marker


Region table for sample 9 : Tag321P (1:2)

From [bp]	To [bp]	Corr. Area	% of Total	Average Size [bp]	Size distribution in CV [%]	Conc. [pg/μl]	Molarity [pmol/l]	Color
137	1,000	150.0	96	327	33.9	2,787.87	14,571.4	Blue

Assay Class: High Sensitivity DNA Assay
 Data Path: C:\...r1_High Sensitivity DNA Assay_DE34903152_2026-01-05_001.xad

Created: 1/5/2026 9:41:04 AM
 Modified: 1/5/2026 10:31:04 AM

Electropherogram Summary Continued ...

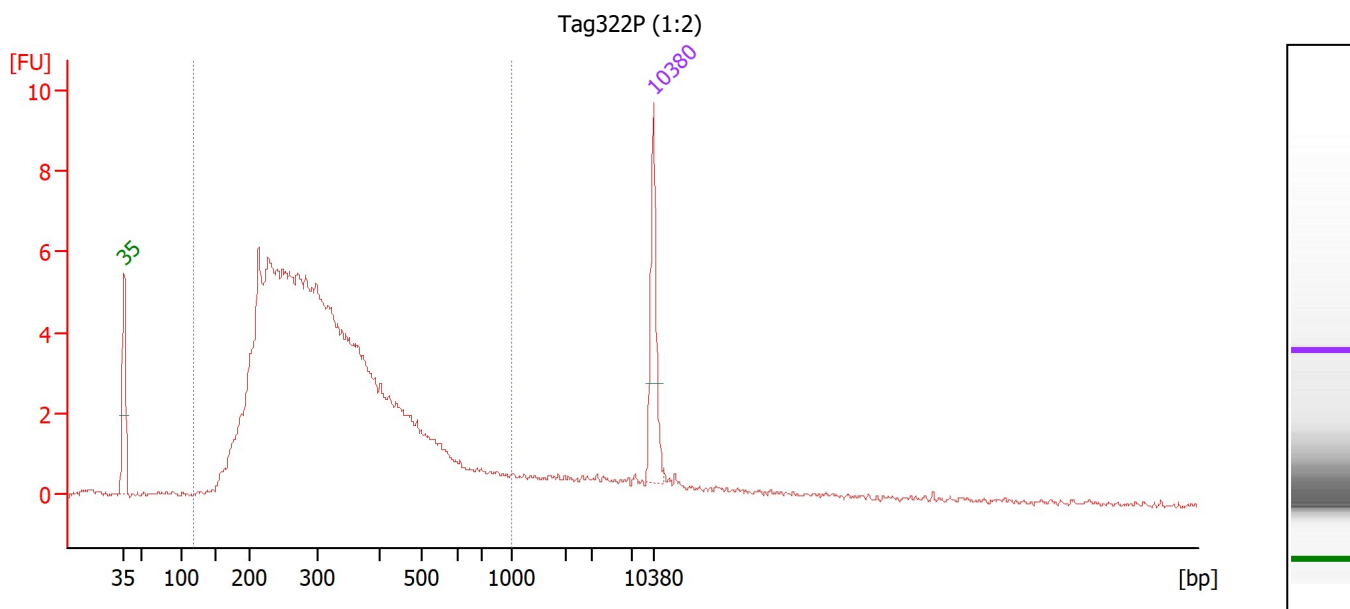

Overall Results for sample 10 : Tag322P (1:2)

Number of peaks found: 0 Corr. Area 1: 136.2
 Noise: 0.1

Peak table for sample 10 : Tag322P (1:2)

Peak	Size [bp]	Conc. [pg/μl]	Molarity [pmol/l]	Observations
1	35	125.00	5,411.3	Lower Marker
2	10,380	75.00	10.9	Upper Marker


Region table for sample 10 : Tag322P (1:2)

From [bp]	To [bp]	Corr. Area	% of Total	Average Size [bp]	Size distribution in CV [%]	Conc. [pg/μl]	Molarity [pmol/l]	Color
117	1,000	136.2	91	334	39.2	2,527.30	13,427.4	Blue

Assay Class: High Sensitivity DNA Assay
 Data Path: C:\...r1_High Sensitivity DNA Assay_DE34903152_2026-01-05_001.xad

Created: 1/5/2026 9:41:04 AM
 Modified: 1/5/2026 10:31:04 AM

Electropherogram Summary Continued ...

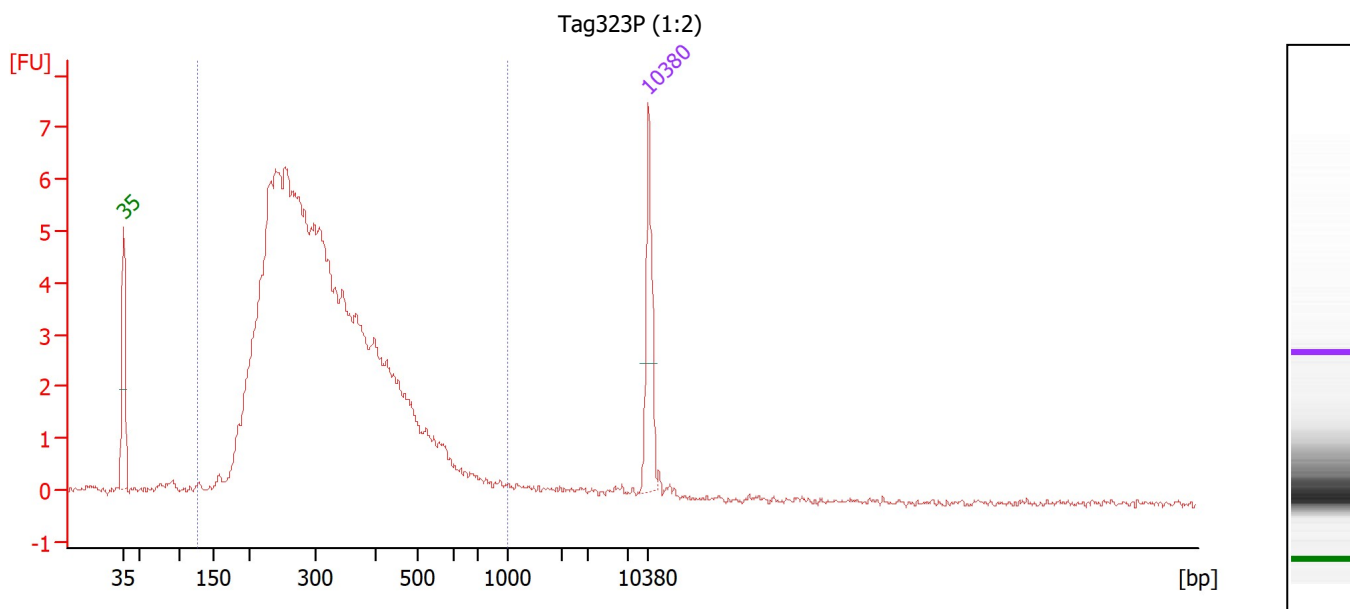

Overall Results for sample 11 : Tag323P (1:2)

Number of peaks found: 0 Corr. Area 1: 124.9
 Noise: 0.1

Peak table for sample 11 : Tag323P (1:2)

Peak	Size [bp]	Conc. [pg/μl]	Molarity [pmol/l]	Observations
1	35	125.00	5,411.3	Lower Marker
2	10,380	75.00	10.9	Upper Marker


Region table for sample 11 : Tag323P (1:2)

From [bp]	To [bp]	Corr. Area	% of Total	Average Size [bp]	Size distribution in CV [%]	Conc. [pg/μl]	Molarity [pmol/l]	Color
125	1,000	124.9	97	326	33.5	2,706.94	14,224.5	Blue

Assay Class: High Sensitivity DNA Assay
 Data Path: C:\...r1_High Sensitivity DNA Assay_DE34903152_2026-01-05_001.xad

Created: 1/5/2026 9:41:04 AM
 Modified: 1/5/2026 10:31:04 AM

Gel Image

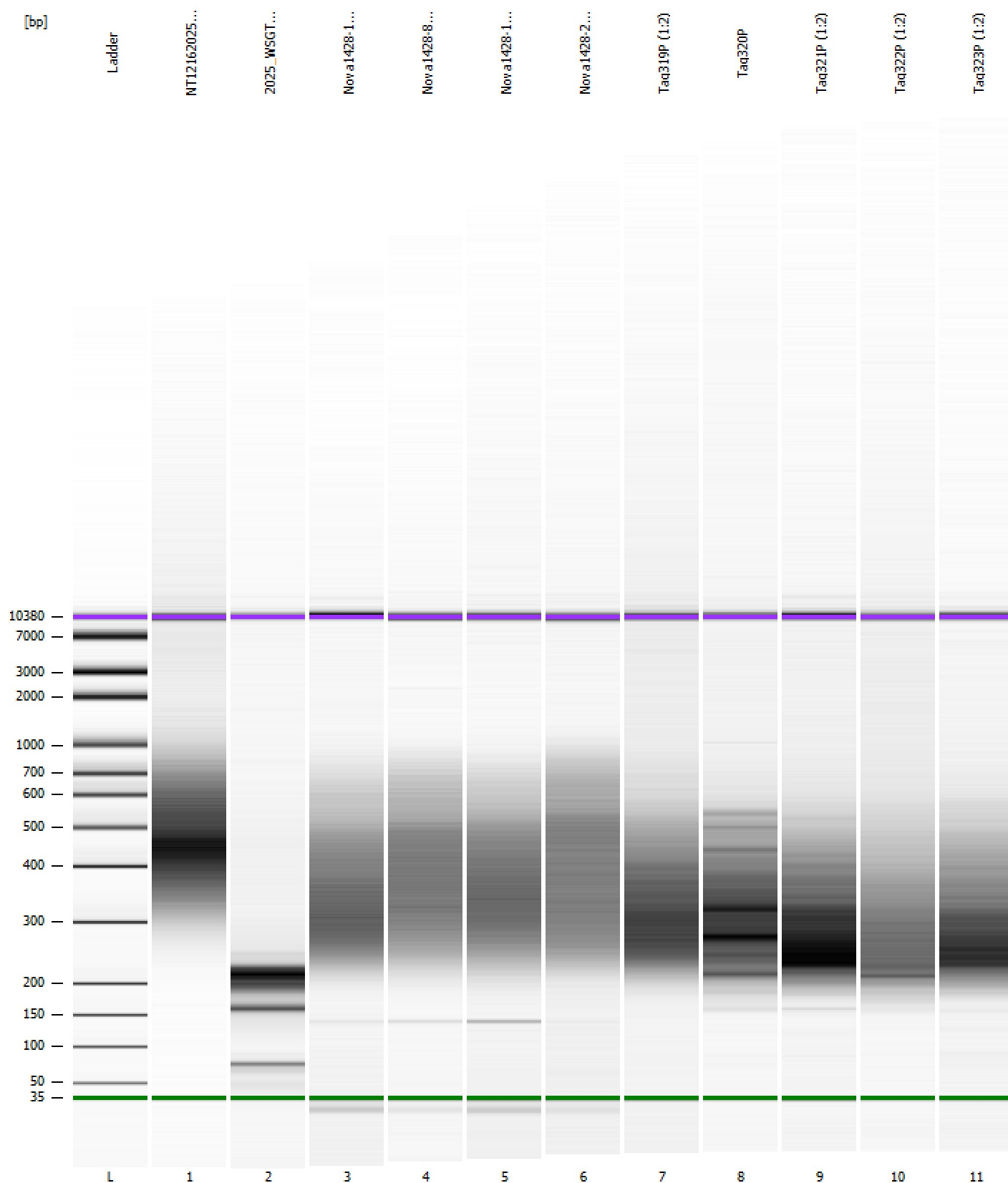

Assay Class: High Sensitivity DNA Assay
 Data Path: C:\...r1_High Sensitivity DNA Assay_DE34903152_2026-01-05_001.xad

Created: 1/5/2026 9:41:04 AM
 Modified: 1/5/2026 10:31:04 AM

Run Logbook

Description	Number	Source	Category	Sub Category	Time	Time Zone	User	Host
Run ended on port 1 (Number of wells acquired: 12)		Instrument	Run		1/5/2026 10:22:27 AM	(GMT --08:00) Pacific Standard Time	sbsuser	DESKTOP-4UNV VOB
Run started on port 1 (File: C:\Program Files (x86)\Agilent\2100 bioanalyzer\2100 expert\data\2026-01-05\Bioanalyzer1_High Sensitivity DNA Assay_DE34903152_2026-01-05_001.xad)		Instrument	Run		1/5/2026 9:41:13 AM	(GMT --08:00) Pacific Standard Time	sbsuser	DESKTOP-4UNV VOB
Product Number : G2938C		Instrument	Run		1/5/2026 9:41:13 AM	(GMT --08:00) Pacific Standard Time	sbsuser	DESKTOP-4UNV VOB
Name :		Instrument	Run		1/5/2026 9:41:13 AM	(GMT --08:00) Pacific Standard Time	sbsuser	DESKTOP-4UNV VOB
Vendor : Agilent Technologies		Instrument	Run		1/5/2026 9:41:13 AM	(GMT --08:00) Pacific Standard Time	sbsuser	DESKTOP-4UNV VOB
Serial# : DE34903152		Instrument	Run		1/5/2026 9:41:13 AM	(GMT --08:00) Pacific Standard Time	sbsuser	DESKTOP-4UNV VOB
Firmware : C.01.069		Instrument	Run		1/5/2026 9:41:13 AM	(GMT --08:00) Pacific Standard Time	sbsuser	DESKTOP-4UNV VOB
Cartridge : Electrode		Instrument	Run		1/5/2026 9:41:13 AM	(GMT --08:00) Pacific Standard Time	sbsuser	DESKTOP-4UNV VOB