

Assay Class: High Sensitivity DNA Assay
Data Path: C:\...pert\data\2025-06-02\10XSC0429_10XSC0430_ATAC_Libraries.xad

Created: 6/2/2025 2:37:26 PM
Modified: 6/2/2025 3:21:59 PM

Electrophoresis File Run Summary

Instrument Information:





Instrument Name: DE34903152 Firmware: C.01.069
Serial#: DE34903152 Type: G2938C

Assay Information:

Assay Origin Path: C:\Program Files (x86)\Agilent\2100 bioanalyzer\2100 expert\assays\dsDNA\High Sensitivity DNA.xsy
Assay Class: High Sensitivity DNA Assay
Version: 1.03
Assay Comments: Copyright © 2003-2010 Agilent Technologies

Chip Information:

Chip Lot #:
Reagent Kit Lot #:
Chip Comments:

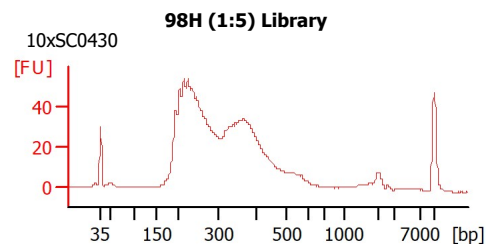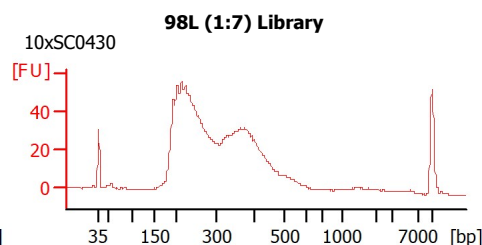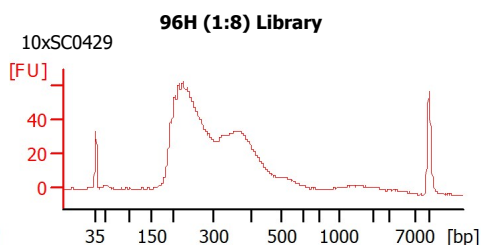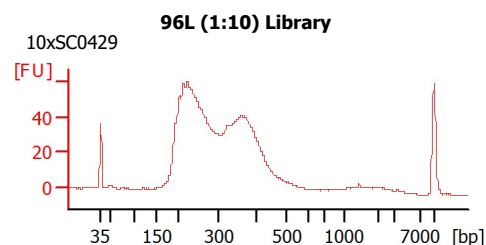

Assay Class: High Sensitivity DNA Assay
Data Path: C:\...pert\data\2025-06-02\10XSC0429_10XSC0430_ATAC_Libraries.xad

Created: 6/2/2025 2:37:26 PM
Modified: 6/2/2025 3:21:59 PM

Electrophoresis File Run Summary (Chip Summary)

Sample Name	Sample Comment	Rest. Digest	Status Observation	Result Label	Result Color
96L (1:10) Library	10xSC0429	☐	✓		
96H (1:8) Library	10xSC0429	☐	✓		
98L (1:7) Library	10xSC0430	☐	✓		
98H (1:5) Library	10xSC0430	☐	✓		
Ladder		☐	✓		

Chip Lot #**Reagent Kit Lot #****Chip Comments :**

Assay Class: High Sensitivity DNA Assay
Data Path: C:\...pert\data\2025-06-02\10XSC0429_10XSC0430_ATAC_Libraries.xad

Created: 6/2/2025 2:37:26 PM
Modified: 6/2/2025 3:21:59 PM

Electrophoresis Assay Details

General Analysis Settings

Number of Available Sample and Ladder Wells (Max.) : 12
Minimum Visible Range [s] : 32
Maximum Visible Range [s] : 138
Start Analysis Time Range [s] : 33
End Analysis Time Range [s] : 137.5
Ladder Concentration [pg/ μ l] : 1950
Uses Standard Area for Ladder Fragments
Lower Marker Concentration [pg/ μ l] : 125
Upper Marker Concentration [pg/ μ l] : 75
Used Upper Marker for Quantitation
Standard Curve Fit is Point to Point
Show Data Aligned to Lower and Upper Marker

Integrator Settings

Integration Start Time [s] : 33.05
Integration End Time [s] : 137
Slope Threshold : 0.8
Height Threshold [FU] : 5
Area Threshold : 0.1
Width Threshold [s] : 0.6
Baseline Plateau [s] : 0.5

Filter Settings

Filter Width [s] : 0.5
Polynomial Order : 4


Ladder

Ladder Peak	Size	Area
1	35	160
2	50	210
3	100	208
4	150	221
5	200	242
6	300	270
7	400	305
8	500	306
9	600	336
10	700	321
11	1000	366
12	2000	413
13	3000	411
14	7000	400
15	10380	214

Assay Class: High Sensitivity DNA Assay
 Data Path: C:\...pert\data\2025-06-02\10XSC0429_10XSC0430_ATAC_Libraries.xad

Created: 6/2/2025 2:37:26 PM
 Modified: 6/2/2025 3:21:59 PM

Electropherogram Summary

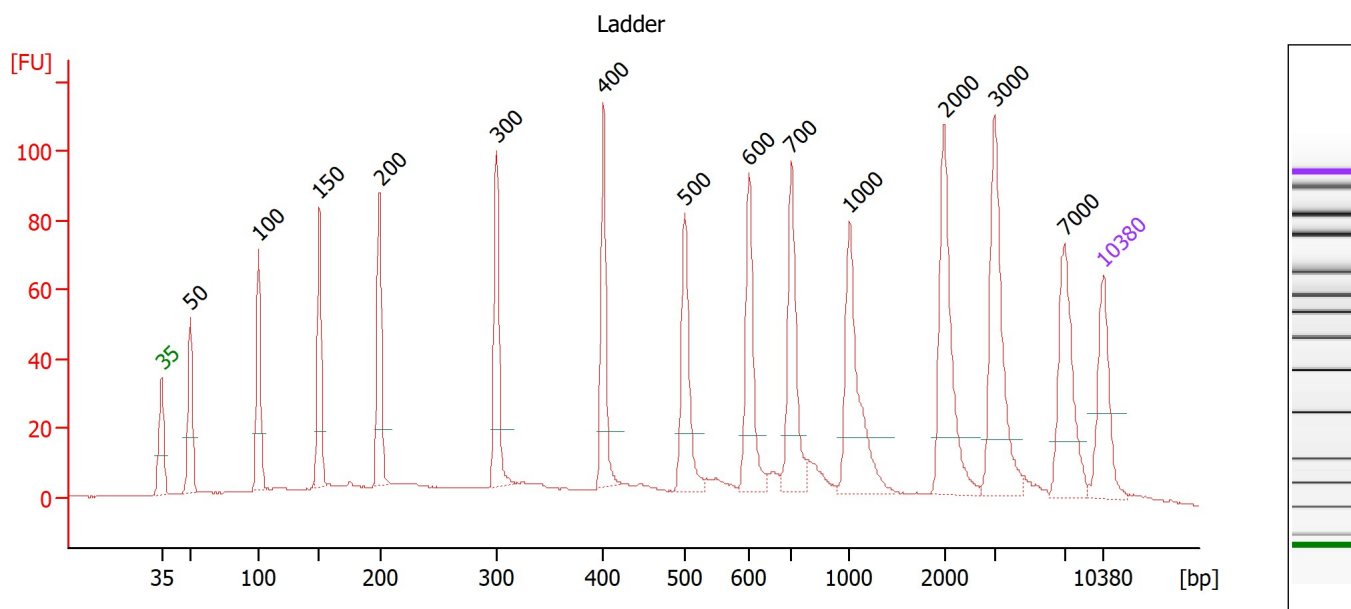

Overall Results for Ladder

Noise: 0.2


Peak table for Ladder

Peak	Size [bp]	Conc. [pg/μl]	Molarity [pmol/l]	Observations
1	35	125.00	5,411.3	Lower Marker
2	50	150.00	4,545.5	Ladder Peak
3	100	150.00	2,272.7	Ladder Peak
4	150	150.00	1,515.2	Ladder Peak
5	200	150.00	1,136.4	Ladder Peak
6	300	150.00	757.6	Ladder Peak
7	400	150.00	568.2	Ladder Peak
8	500	150.00	454.5	Ladder Peak
9	600	150.00	378.8	Ladder Peak
10	700	150.00	324.7	Ladder Peak
11	1,000	150.00	227.3	Ladder Peak
12	2,000	150.00	113.6	Ladder Peak
13	3,000	150.00	75.8	Ladder Peak
14	7,000	150.00	32.5	Ladder Peak
15	10,380	75.00	10.9	Upper Marker

Assay Class: High Sensitivity DNA Assay
 Data Path: C:\...pert\data\2025-06-02\10XSC0429_10XSC0430_ATAC_Libraries.xad

Created: 6/2/2025 2:37:26 PM
 Modified: 6/2/2025 3:21:59 PM

Electropherogram Summary Continued ...

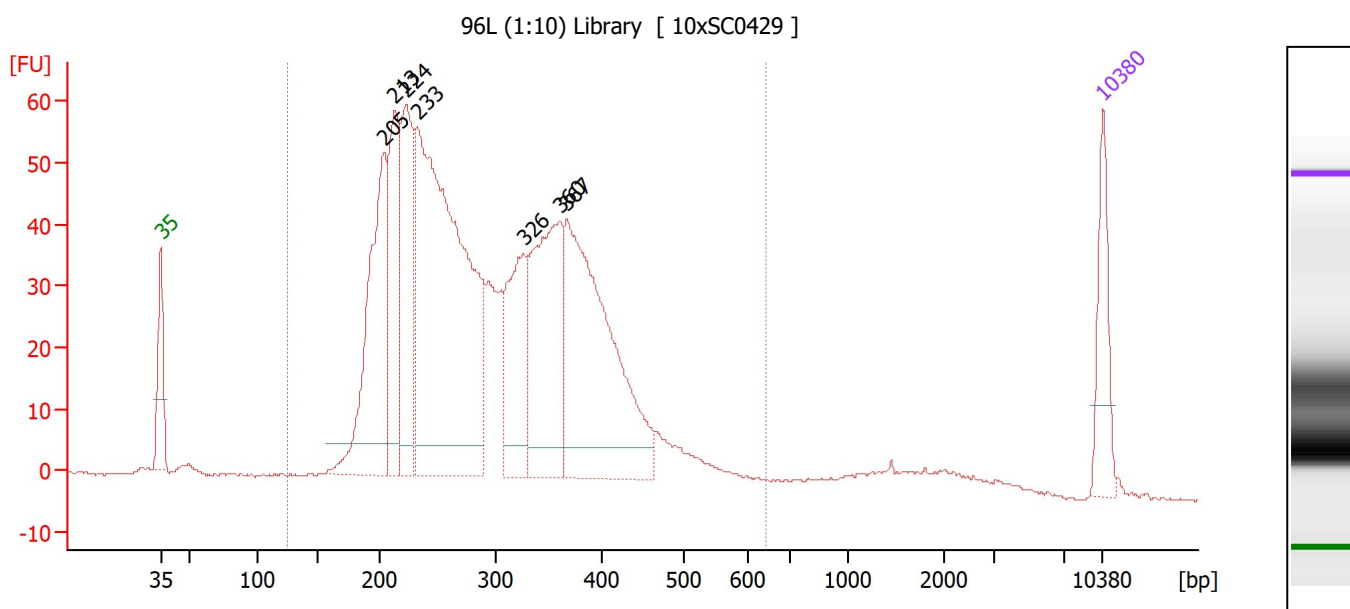

Overall Results for sample 2 : 96L (1:10) Library

Number of peaks found: 7 Corr. Area 1: 1,206.6
 Noise: 0.2

Peak table for sample 2 : 96L (1:10) Library

Peak	Size [bp]	Conc. [pg/μl]	Molarity [pmol/l]	Observations
1	35	125.00	5,411.3	Lower Marker
2	205	362.92	2,685.5	
3	213	217.24	1,543.6	
4	224	256.39	1,737.8	
5	233	907.19	5,893.2	
6	326	192.79	896.8	
7	360	312.41	1,314.4	
8	367	496.24	2,050.1	
9	10,380	75.00	10.9	Upper Marker


Region table for sample 2 : 96L (1:10) Library

From [bp]	To [bp]	Corr. Area	% of Total	Average Size [bp]	Size distribution in CV [%]	Conc. [pg/μl]	Molarity [pmol/l]	Color
124	641	1,206.6	95	303	27.7	3,058.68	16,930.0	Blue

Assay Class: High Sensitivity DNA Assay
 Data Path: C:\...pert\data\2025-06-02\10XSC0429_10XSC0430_ATAC_Libraries.xad

Created: 6/2/2025 2:37:26 PM
 Modified: 6/2/2025 3:21:59 PM

Electropherogram Summary Continued ...

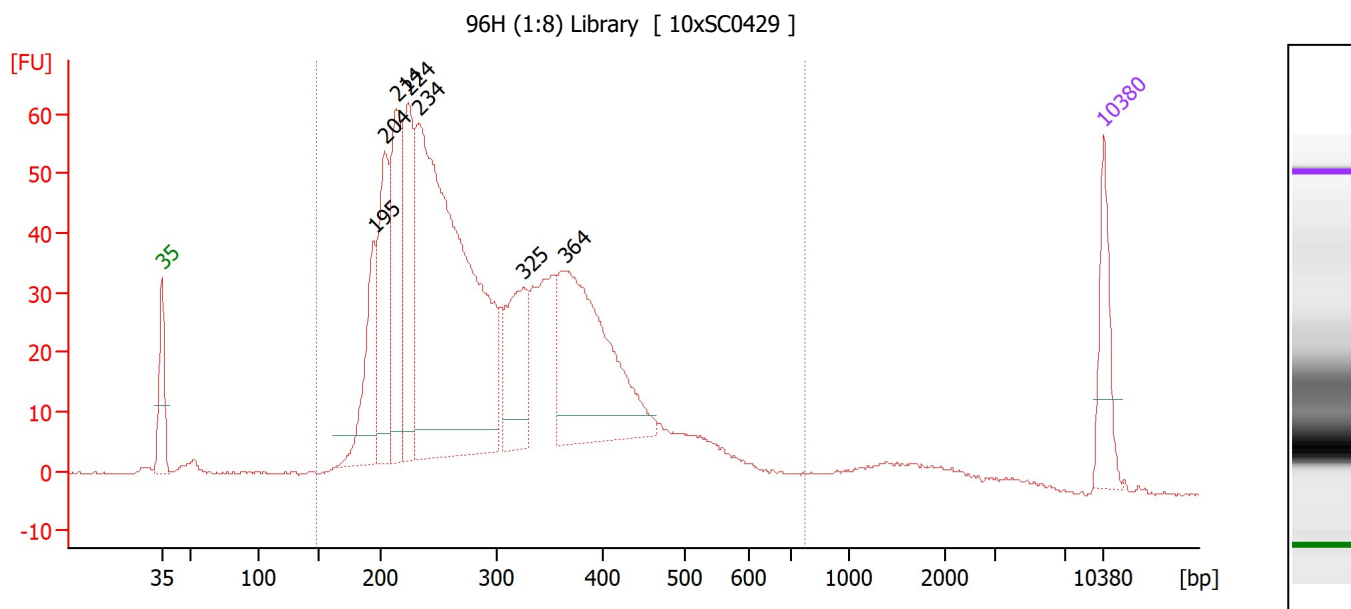

Overall Results for sample 3 : 96H (1:8) Library

Number of peaks found: 7 Corr. Area 1: 1,208.5
 Noise: 0.2

Peak table for sample 3 : 96H (1:8) Library

Peak	Size [bp]	Conc. [pg/μl]	Molarity [pmol/l]	Observations
1	35	125.00	5,411.3	Lower Marker
2	195	171.04	1,327.1	
3	204	215.32	1,596.2	
4	214	235.80	1,667.3	
5	224	251.79	1,701.8	
6	234	1,058.56	6,851.8	
7	325	171.88	800.6	
8	364	407.45	1,693.8	
9	10,380	75.00	10.9	Upper Marker


Region table for sample 3 : 96H (1:8) Library

From [bp]	To [bp]	Corr. Area	% of Total	Average Size [bp]	Size distribution in CV [%]	Conc. [pg/μl]	Molarity [pmol/l]	Color
147	767	1,208.5	92	311	32.1	3,306.75	18,223.0	Blue

Assay Class: High Sensitivity DNA Assay
 Data Path: C:\...pert\data\2025-06-02\10XSC0429_10XSC0430_ATAC_Libraries.xad

Created: 6/2/2025 2:37:26 PM
 Modified: 6/2/2025 3:21:59 PM

Electropherogram Summary Continued ...

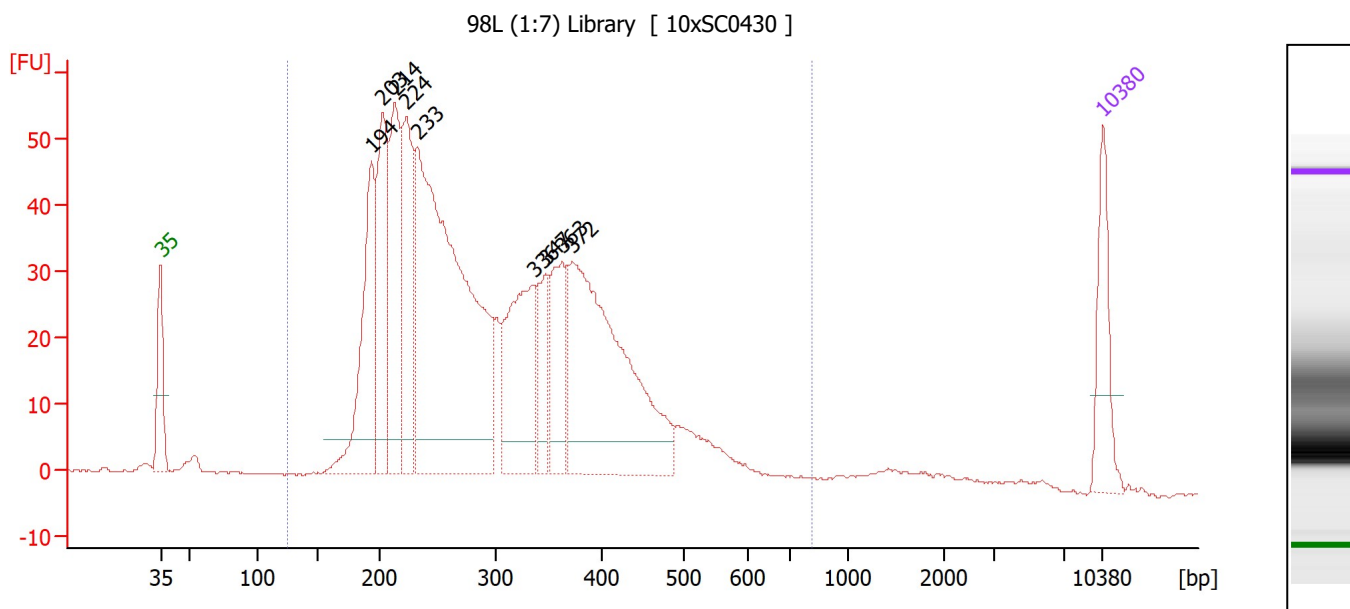

Overall Results for sample 4 : 98L (1:7) Library

Number of peaks found: 9 Corr. Area 1: 1,090.4
 Noise: 0.2

Peak table for sample 4 : 98L (1:7) Library

Peak	Size [bp]	Conc. [pg/μl]	Molarity [pmol/l]	Observations
1	35	125.00	5,411.3	Lower Marker
2	194	254.90	1,990.3	
3	203	216.10	1,610.2	
4	214	222.16	1,572.1	
5	224	210.71	1,427.7	
6	233	841.05	5,466.4	
7	336	223.07	1,007.1	
8	347	78.97	344.5	
9	363	124.59	520.6	
10	372	480.06	1,957.5	
11	10,380	75.00	10.9	Upper Marker


Region table for sample 4 : 98L (1:7) Library

From [bp]	To [bp]	Corr. Area	% of Total	Average Size [bp]	Size distribution in CV [%]	Conc. [pg/μl]	Molarity [pmol/l]	Color
124	806	1,090.4	95	310	32.4	2,909.08	16,199.6	Blue

Assay Class: High Sensitivity DNA Assay
 Data Path: C:\...pert\data\2025-06-02\10XSC0429_10XSC0430_ATAC_Libraries.xad

Created: 6/2/2025 2:37:26 PM
 Modified: 6/2/2025 3:21:59 PM

Electropherogram Summary Continued ...

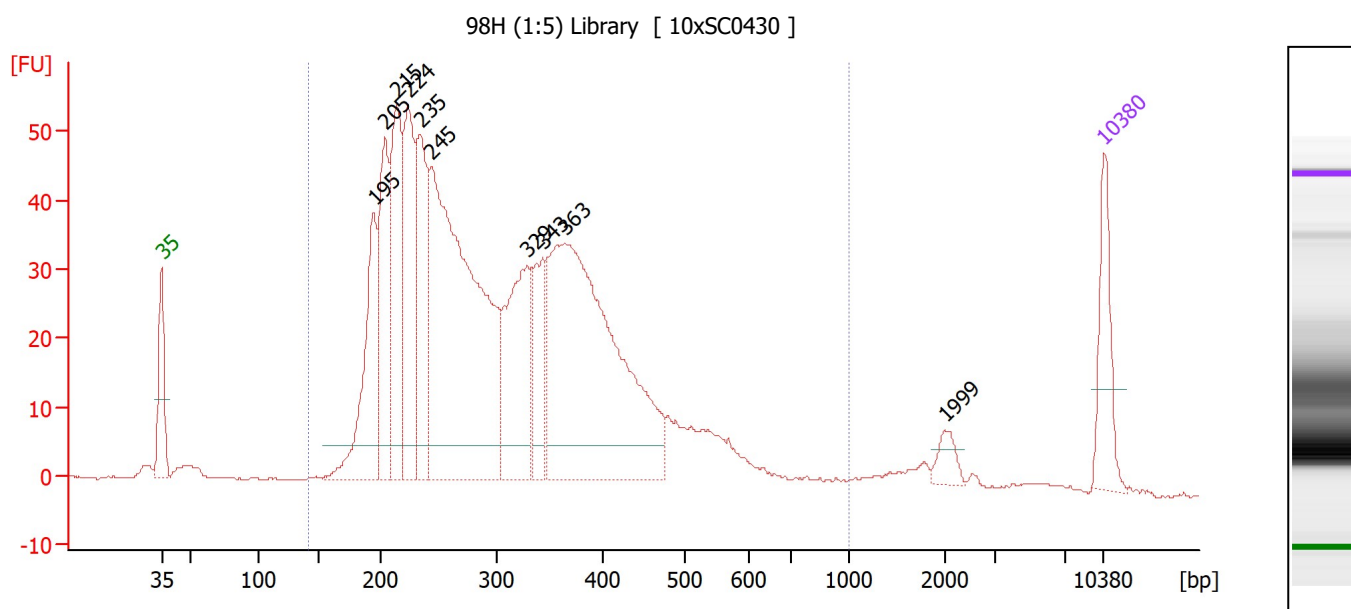

Overall Results for sample 5 : 98H (1:5) Library

Number of peaks found: 10 Corr. Area 1: 1,090.0
 Noise: 0.2

Peak table for sample 5 : 98H (1:5) Library

Peak	Size [bp]	Conc. [pg/μl]	Molarity [pmol/l]	Observations
1	35	125.00	5,411.3	Lower Marker
2	195	227.21	1,763.6	
3	205	204.75	1,514.8	
4	215	240.64	1,699.2	
5	224	237.96	1,607.1	
6	235	203.13	1,311.0	
7	245	791.20	4,902.2	
8	329	238.20	1,096.7	
9	343	115.80	511.1	
10	363	694.53	2,897.0	
11	1,999	21.63	16.4	Upper Marker
12	10,380	75.00	10.9	

Region table for sample 5 : 98H (1:5) Library

From [bp]	To [bp]	Corr. Area	% of Total	Average Size [bp]	Size distribution in CV [%]	Conc. [pg/μl]	Molarity [pmol/l]	Color
141	1,000	1,090.0	94	320	35.1	3,209.39	17,457.4	Blue

Assay Class: High Sensitivity DNA Assay
Data Path: C:\...pert\data\2025-06-02\10XSC0429_10XSC0430_ATAC_Libraries.xad

Created: 6/2/2025 2:37:26 PM
Modified: 6/2/2025 3:21:59 PM

Gel Image

Assay Class: High Sensitivity DNA Assay
 Data Path: C:\...pert\data\2025-06-02\10XSC0429_10XSC0430_ATAC_Libraries.xad

Created: 6/2/2025 2:37:26 PM
 Modified: 6/2/2025 3:21:59 PM

Run Logbook

Description	Number	Source	Category	Sub Category	Time	Time Zone	User	Host
Run ended on port 1 (Number of wells acquired: 12)		Instrument	Run		6/2/2025 3:18:52 PM	(GMT --07:00) Pacific Standard Time	sbsuser	DESKTOP-4UNVVOB
Run started on port 1 (File: C:\Program Files (x86)\Agilent\2100 bioanalyzer\2100 expert\data\2025-06-02\Bioanalyzer1_High Sensitivity DNA Assay_DE34903152_2025-06-02_001.xad)		Instrument	Run		6/2/2025 2:37:33 PM	(GMT --07:00) Pacific Standard Time	sbsuser	DESKTOP-4UNVVOB
Product Number : G2938C		Instrument	Run		6/2/2025 2:37:33 PM	(GMT --07:00) Pacific Standard Time	sbsuser	DESKTOP-4UNVVOB
Name :		Instrument	Run		6/2/2025 2:37:33 PM	(GMT --07:00) Pacific Standard Time	sbsuser	DESKTOP-4UNVVOB
Vendor : Agilent Technologies		Instrument	Run		6/2/2025 2:37:33 PM	(GMT --07:00) Pacific Standard Time	sbsuser	DESKTOP-4UNVVOB
Serial# : DE34903152		Instrument	Run		6/2/2025 2:37:33 PM	(GMT --07:00) Pacific Standard Time	sbsuser	DESKTOP-4UNVVOB
Firmware : C.01.069		Instrument	Run		6/2/2025 2:37:33 PM	(GMT --07:00) Pacific Standard Time	sbsuser	DESKTOP-4UNVVOB
Cartridge : Electrode		Instrument	Run		6/2/2025 2:37:33 PM	(GMT --07:00) Pacific Standard Time	sbsuser	DESKTOP-4UNVVOB