


Assay Class: High Sensitivity DNA Assay
Data Path: C:\... bioanalyzer\2100 expert\data\2025-07-15\10XSC0434_cDNA.xad

Created: 7/15/2025 3:04:30 PM
Modified: 7/15/2025 3:52:15 PM

Electrophoresis File Run Summary

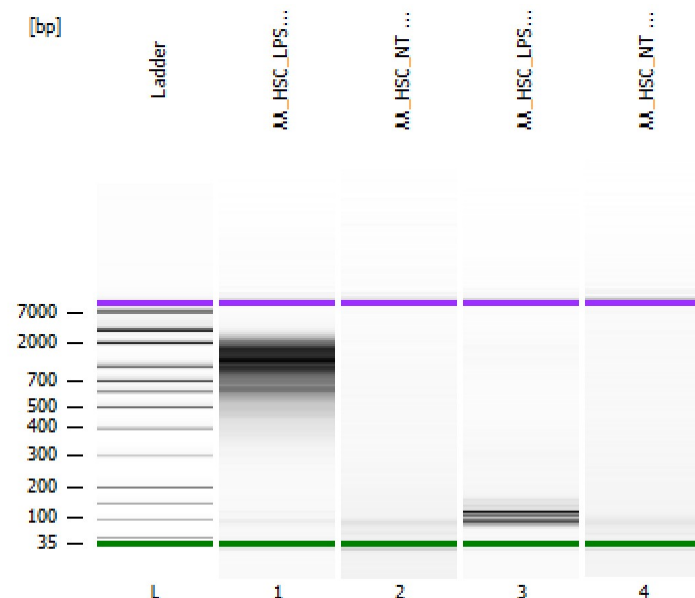

Instrument Information:





Instrument Name: DE34903152 Firmware: C.01.069
Serial#: DE34903152 Type: G2938C

Assay Information:

Assay Origin Path: C:\Program Files (x86)\Agilent\2100 bioanalyzer\2100 expert\assays\dsDNA\High Sensitivity DNA.xsy
Assay Class: High Sensitivity DNA Assay
Version: 1.03
Assay Comments: Copyright © 2003-2010 Agilent Technologies

Chip Information:

Chip Lot #:
Reagent Kit Lot #:
Chip Comments:

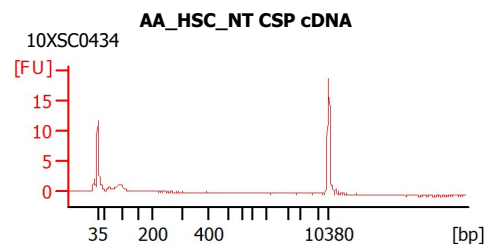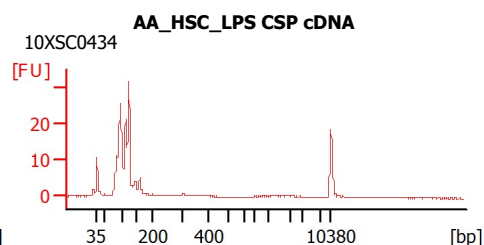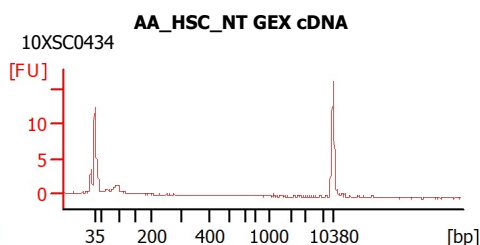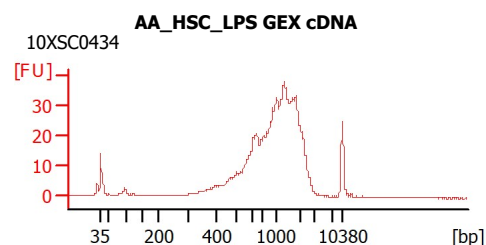

Assay Class: High Sensitivity DNA Assay
Data Path: C:\... bioanalyzer\2100 expert\data\2025-07-15\10XSC0434_cDNA.xad

Created: 7/15/2025 3:04:30 PM
Modified: 7/15/2025 3:52:15 PM

Electrophoresis File Run Summary (Chip Summary)

Sample Name	Sample Comment	Rest. Digest	Status	Observation	Result Label	Result Color
AA_HSC_LPS GEX cDNA	10XSC0434	☐		✓		
AA_HSC_NT GEX cDNA	10XSC0434	☐		✓		
AA_HSC_LPS CSP cDNA	10XSC0434	☐		✓		
AA_HSC_NT CSP cDNA	10XSC0434	☐		✓		
Ladder		☐		✓		

Chip Lot #**Reagent Kit Lot #****Chip Comments :**

Assay Class: High Sensitivity DNA Assay
Data Path: C:\... bioanalyzer\2100 expert\data\2025-07-15\10XSC0434_cDNA.xad

Created: 7/15/2025 3:04:30 PM
Modified: 7/15/2025 3:52:15 PM

Electrophoresis Assay Details

General Analysis Settings

Number of Available Sample and Ladder Wells (Max.) : 12
Minimum Visible Range [s] : 32
Maximum Visible Range [s] : 138
Start Analysis Time Range [s] : 33
End Analysis Time Range [s] : 137.5
Ladder Concentration [pg/μl] : 1950
Uses Standard Area for Ladder Fragments
Lower Marker Concentration [pg/μl] : 125
Upper Marker Concentration [pg/μl] : 75
Used Upper Marker for Quantitation
Standard Curve Fit is Point to Point
Show Data Aligned to Lower and Upper Marker

Integrator Settings

Integration Start Time [s] : 33.05
Integration End Time [s] : 137
Slope Threshold : 0.8
Height Threshold [FU] : 5
Area Threshold : 0.1
Width Threshold [s] : 0.6
Baseline Plateau [s] : 0.5

Filter Settings

Filter Width [s] : 0.5
Polynomial Order : 4


Ladder

Ladder Peak	Size	Area
1	35	160
2	50	210
3	100	208
4	150	221
5	200	242
6	300	270
7	400	305
8	500	306
9	600	336
10	700	321
11	1000	366
12	2000	413
13	3000	411
14	7000	400
15	10380	214

Assay Class: High Sensitivity DNA Assay
Data Path: C:\... bioanalyzer\2100 expert\data\2025-07-15\10XSC0434_cDNA.xad

Created: 7/15/2025 3:04:30 PM
Modified: 7/15/2025 3:52:15 PM

Electropherogram Summary

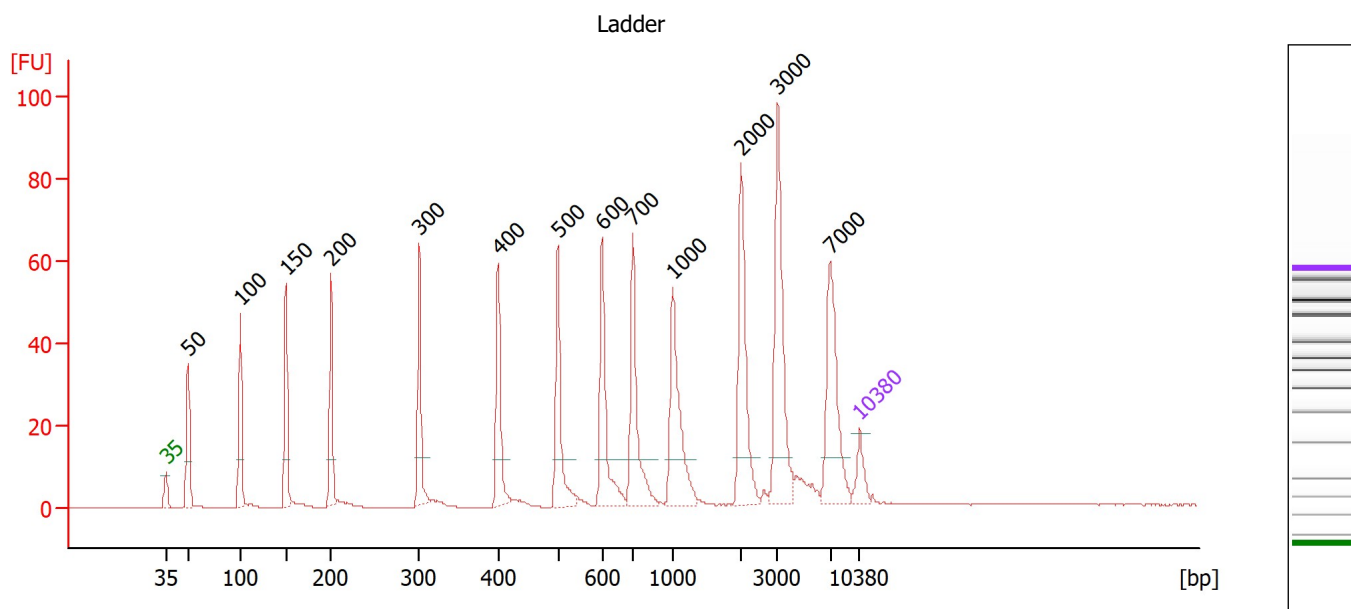

Overall Results for Ladder

Noise: 0.1


Peak table for Ladder

Peak	Size [bp]	Conc. [pg/μl]	Molarity [pmol/l]	Observations
1	35	125.00	5,411.3	Lower Marker
2	50	150.00	4,545.5	Ladder Peak
3	100	150.00	2,272.7	Ladder Peak
4	150	150.00	1,515.2	Ladder Peak
5	200	150.00	1,136.4	Ladder Peak
6	300	150.00	757.6	Ladder Peak
7	400	150.00	568.2	Ladder Peak
8	500	150.00	454.5	Ladder Peak
9	600	150.00	378.8	Ladder Peak
10	700	150.00	324.7	Ladder Peak
11	1,000	150.00	227.3	Ladder Peak
12	2,000	150.00	113.6	Ladder Peak
13	3,000	150.00	75.8	Ladder Peak
14	7,000	150.00	32.5	Ladder Peak
15	10,380	75.00	10.9	Upper Marker

Assay Class: High Sensitivity DNA Assay
 Data Path: C:\... bioanalyzer\2100 expert\data\2025-07-15\10XSC0434_cDNA.xad

Created: 7/15/2025 3:04:30 PM
 Modified: 7/15/2025 3:52:15 PM

Electropherogram Summary Continued ...

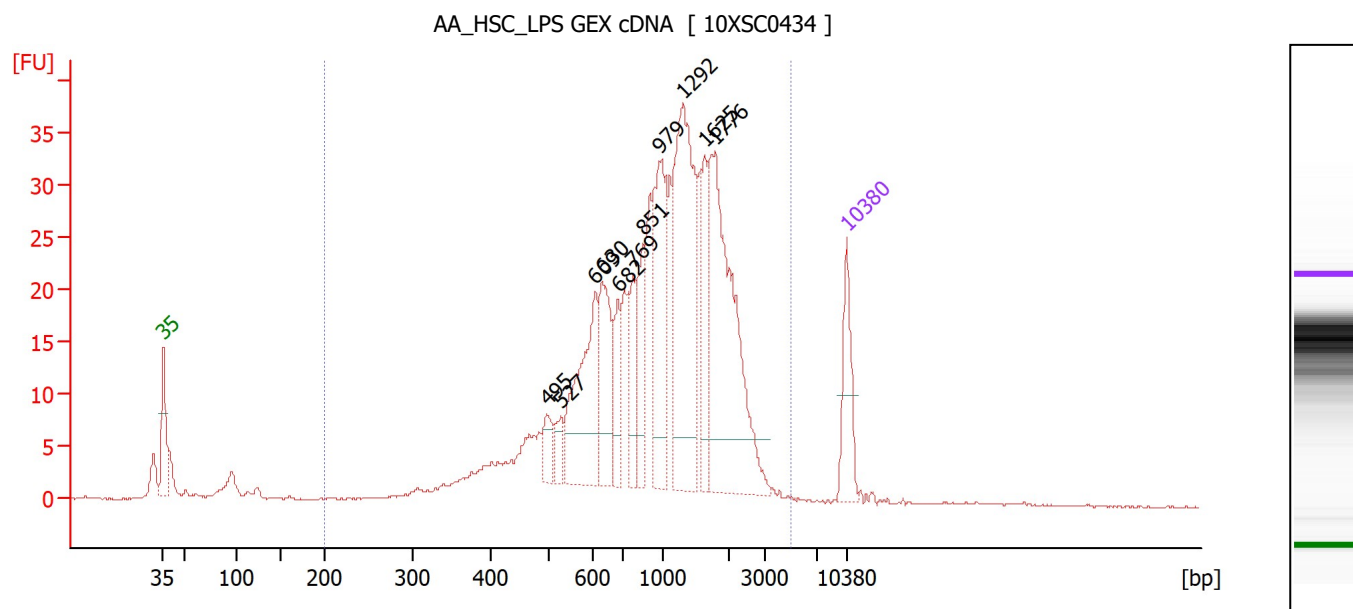

Overall Results for sample 1 : AA_HSC_LPS GEX cDNA

Number of peaks found: 11 Corr. Area 1: 557.4
 Noise: 0.1

Peak table for sample 1 : AA_HSC_LPS GEX cDNA

Peak	Size [bp]	Conc. [pg/μl]	Molarity [pmol/l]	Observations
1	35	125.00	5,411.3	Lower Marker
2	495	40.55	124.1	
3	527	34.09	98.1	
4	609	217.13	540.5	
5	630	134.32	323.2	
6	682	63.11	140.2	
7	769	75.11	148.0	
8	851	90.62	161.3	
9	979	211.16	326.8	
10	1,292	348.18	408.2	
11	1,625	87.16	81.2	
12	1,776	368.49	314.4	
13	10,380	75.00	10.9	Upper Marker


Region table for sample 1 : AA_HSC_LPS GEX cDNA

From [bp]	To [bp]	Corr. Area	% of Total	Average Size [bp]	Size distribution in CV [%]	Conc. [pg/μl]	Molarity [pmol/l]	Color
200	4,871	557.4	94	1,158	51.9	2,422.95	4,526.7	Blue

Assay Class: High Sensitivity DNA Assay
 Data Path: C:\... bioanalyzer\2100 expert\data\2025-07-15\10XSC0434_cDNA.xad

Created: 7/15/2025 3:04:30 PM
 Modified: 7/15/2025 3:52:15 PM

Electropherogram Summary Continued ...

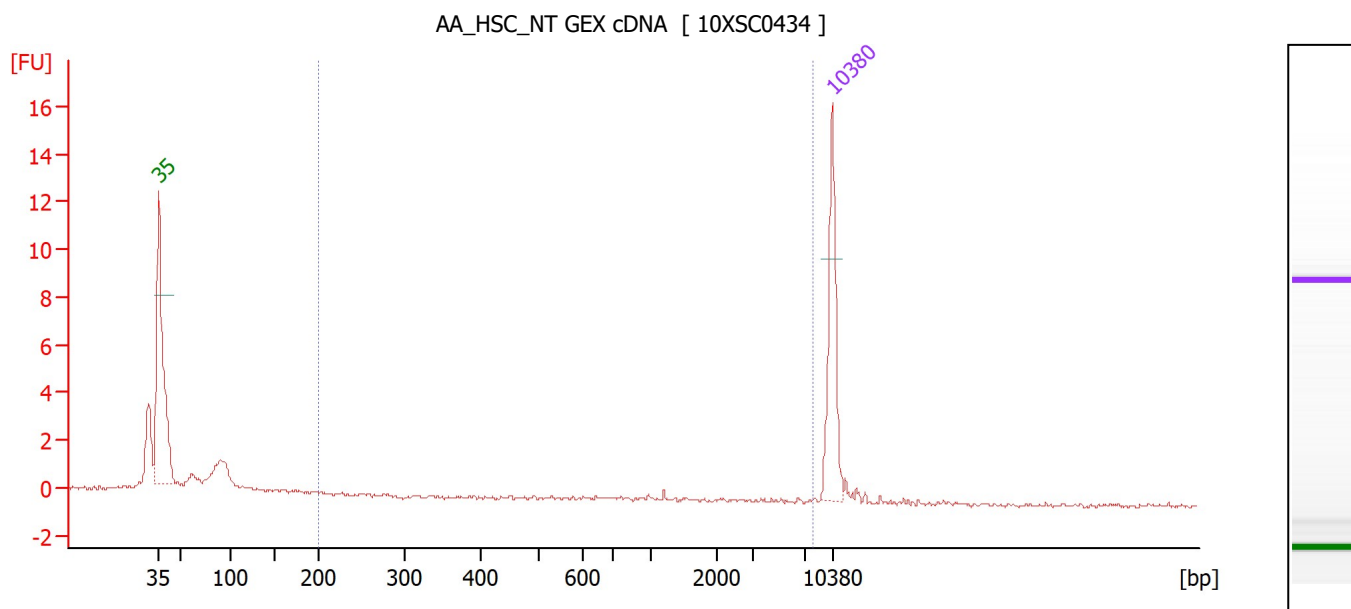

Overall Results for sample 2 : AA_HSC_NT GEX cDNA

Number of peaks found: 0 Corr. Area 1: 0.2
 Noise: 0.1

Peak table for sample 2 : AA_HSC_NT GEX cDNA

Peak	Size [bp]	Conc. [pg/μl]	Molarity [pmol/l]	Observations
1	35	125.00	5,411.3	Lower Marker
2	10,380	75.00	10.9	Upper Marker


Region table for sample 2 : AA_HSC_NT GEX cDNA

From [bp]	To [bp]	Corr. Area	% of Total	Average Size [bp]	Size distribution in CV [%]	Conc. [pg/μl]	Molarity [pmol/l]	Color
200	8,116	0.2	1	2,345	97.0	1.80	3.5	Blue

Assay Class: High Sensitivity DNA Assay
 Data Path: C:\... bioanalyzer\2100 expert\data\2025-07-15\10XSC0434_cDNA.xad

Created: 7/15/2025 3:04:30 PM
 Modified: 7/15/2025 3:52:15 PM

Electropherogram Summary Continued ...

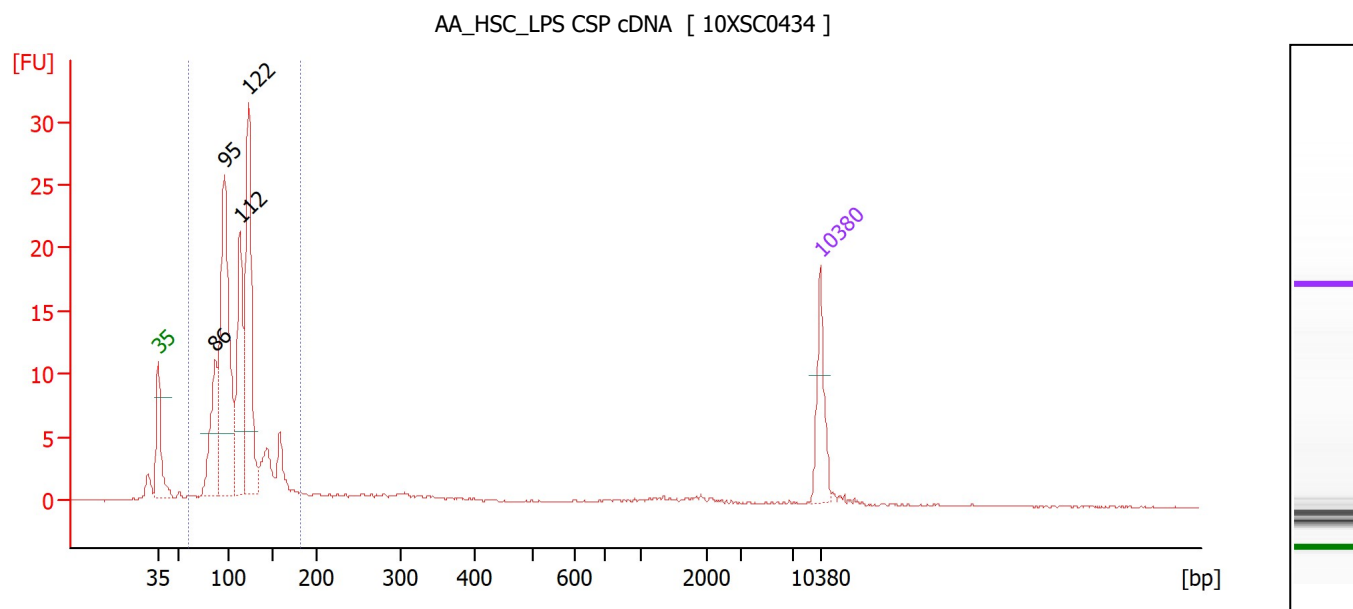

Overall Results for sample 3 : AA_HSC_LPS CSP cDNA

Number of peaks found: 4 Corr. Area 1: 162.9
 Noise: 0.1

Peak table for sample 3 : AA_HSC_LPS CSP cDNA

Peak	Size [bp]	Conc. [pg/μl]	Molarity [pmol/l]	Observations
1	35	125.00	5,411.3	Lower Marker
2	86	194.55	3,415.2	
3	95	586.78	9,318.8	
4	112	281.49	3,793.1	
5	122	448.03	5,555.4	
6	10,380	75.00	10.9	Upper Marker


Region table for sample 3 : AA_HSC_LPS CSP cDNA

From [bp]	To [bp]	Corr. Area	% of Total	Average Size [bp]	Size distribution in CV [%]	Conc. [pg/μl]	Molarity [pmol/l]	Color
60	183	162.9	86	111	18.8	1,762.43	24,820.6	Blue

Assay Class: High Sensitivity DNA Assay
 Data Path: C:\... bioanalyzer\2100 expert\data\2025-07-15\10XSC0434_cDNA.xad

Created: 7/15/2025 3:04:30 PM
 Modified: 7/15/2025 3:52:15 PM

Electropherogram Summary Continued ...

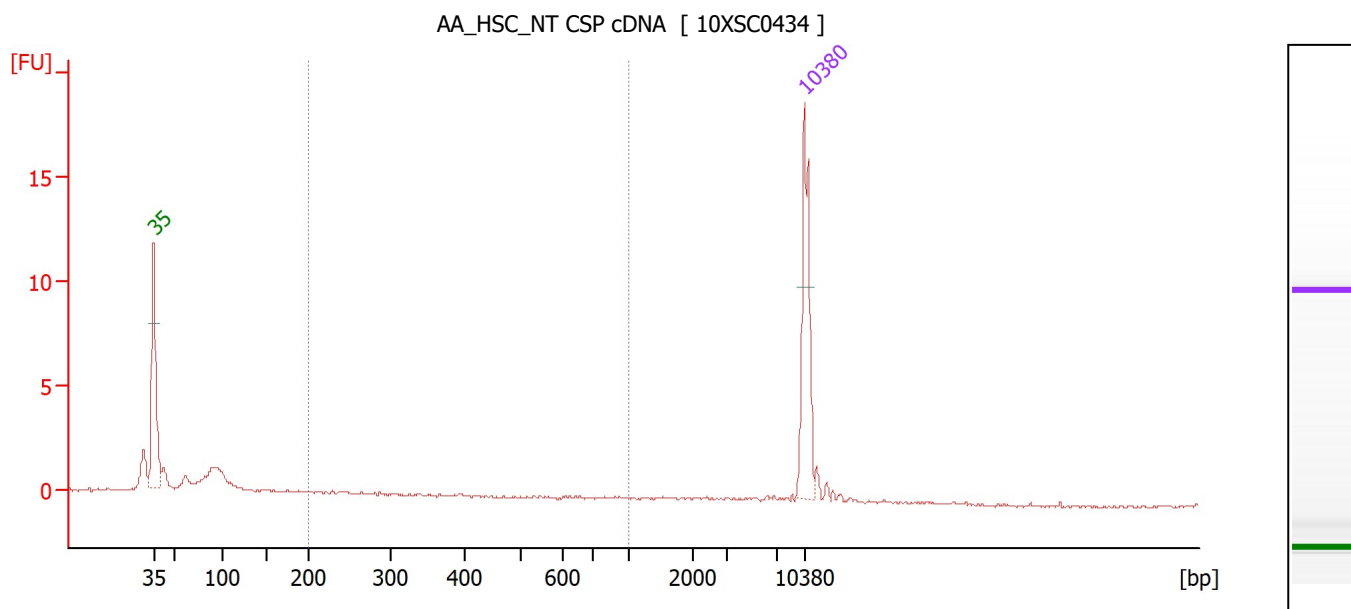

Overall Results for sample 4 : AA_HSC_NT CSP cDNA

Number of peaks found: 0 Corr. Area 1: 1.4
 Noise: 0.1

Peak table for sample 4 : AA_HSC_NT CSP cDNA


Peak	Size [bp]	Conc. [pg/μl]	Molarity [pmol/l]	Observations
1	35	125.00	5,411.3	Lower Marker
2	10,380	75.00	10.9	Upper Marker

Region table for sample 4 : AA_HSC_NT CSP cDNA

From [bp]	To [bp]	Corr. Area	% of Total	Average Size [bp]	Size distribution in CV [%]	Conc. [pg/μl]	Molarity [pmol/l]	Color
200	1,000	1.4	7	354	53.7	12.07	65.7	Blue

Assay Class: High Sensitivity DNA Assay
Data Path: C:\... bioanalyzer\2100 expert\data\2025-07-15\10XSC0434_cDNA.xad

Created: 7/15/2025 3:04:30 PM
Modified: 7/15/2025 3:52:15 PM

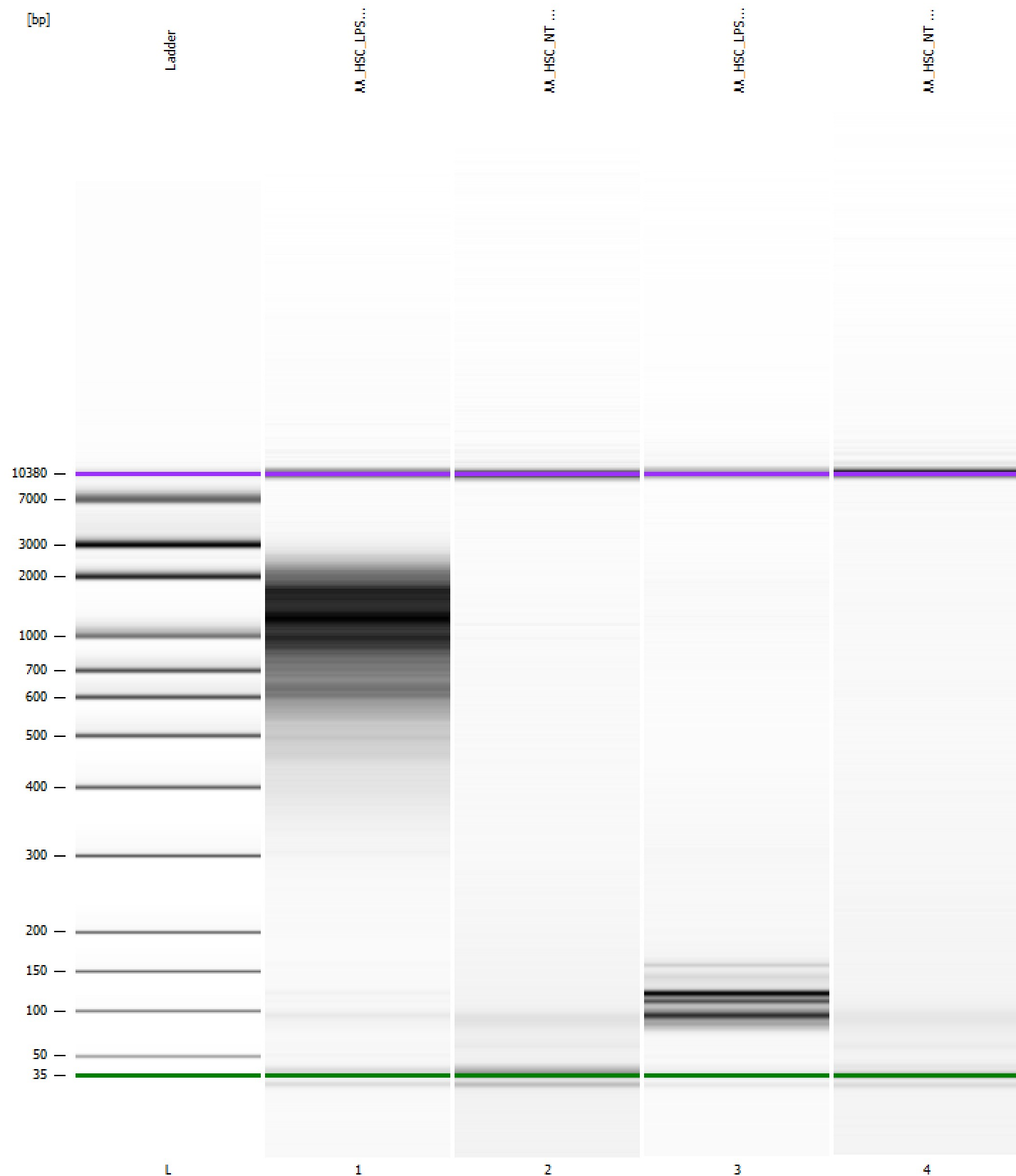Gel Image

Assay Class: High Sensitivity DNA Assay
 Data Path: C:\... bioanalyzer\2100 expert\data\2025-07-15\10XSC0434_cDNA.xad

Created: 7/15/2025 3:04:30 PM
 Modified: 7/15/2025 3:52:15 PM

Run Logbook

Description	Number	Source	Category	Sub Category	Time	Time Zone	User	Host
Run ended on port 1 (Number of wells acquired: 12)		Instrument	Run		7/15/2025 3:45:48 PM	(GMT --07:00) Pacific Standard Time	sbsuser	DESKTOP-4UNV VOB
Run started on port 1 (File: C:\Program Files (x86)\Agilent\2100 bioanalyzer\2100 expert\data\2025-07-15\Bioanalyzer1_High Sensitivity DNA Assay_DE34903152_2025-07-15_001.xad)		Instrument	Run		7/15/2025 3:04:36 PM	(GMT --07:00) Pacific Standard Time	sbsuser	DESKTOP-4UNV VOB
Product Number : G2938C		Instrument	Run		7/15/2025 3:04:36 PM	(GMT --07:00) Pacific Standard Time	sbsuser	DESKTOP-4UNV VOB
Name :		Instrument	Run		7/15/2025 3:04:36 PM	(GMT --07:00) Pacific Standard Time	sbsuser	DESKTOP-4UNV VOB
Vendor : Agilent Technologies		Instrument	Run		7/15/2025 3:04:36 PM	(GMT --07:00) Pacific Standard Time	sbsuser	DESKTOP-4UNV VOB
Serial# : DE34903152		Instrument	Run		7/15/2025 3:04:36 PM	(GMT --07:00) Pacific Standard Time	sbsuser	DESKTOP-4UNV VOB
Firmware : C.01.069		Instrument	Run		7/15/2025 3:04:36 PM	(GMT --07:00) Pacific Standard Time	sbsuser	DESKTOP-4UNV VOB
Cartridge : Electrode		Instrument	Run		7/15/2025 3:04:36 PM	(GMT --07:00) Pacific Standard Time	sbsuser	DESKTOP-4UNV VOB