


Assay Class: High Sensitivity DNA Assay
Data Path: C:\...gh Sensitivity DNA Assay_DE13701287_2025-09-09_16-06-55.xad

Created: 9/9/2025 4:06:55 PM
Modified: 9/10/2025 5:43:51 AM

Electrophoresis File Run Summary

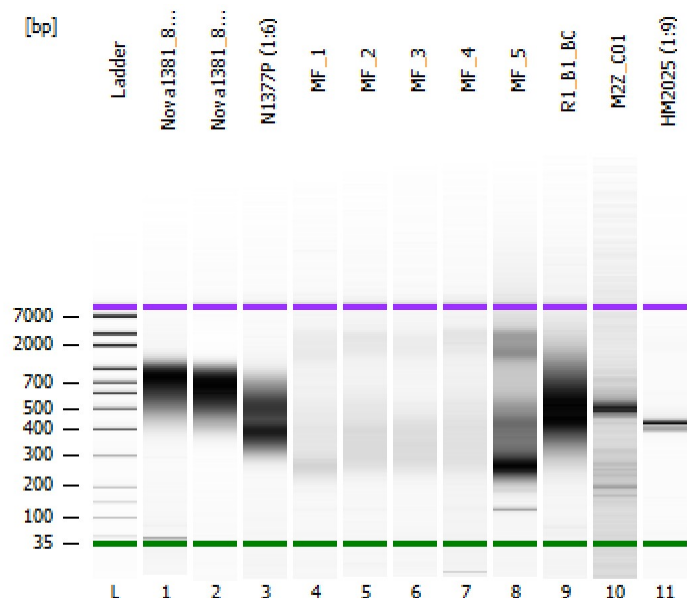

Instrument Information:












Instrument Name: BioA#2
Serial#: DE13701287
Firmware: C.01.069
Type: G2938B

Assay Information:

Assay Origin Path: C:\Program Files (x86)\Agilent\2100 bioanalyzer\2100 expert\assays\dsDNA\High Sensitivity DNA.xsy
Assay Class: High Sensitivity DNA Assay
Version: 1.03
Assay Comments: Copyright © 2003-2010 Agilent Technologies

Chip Information:

Chip Lot #:
Reagent Kit Lot #:
Chip Comments:

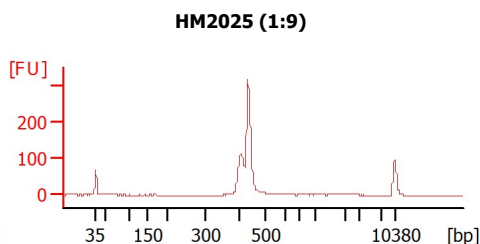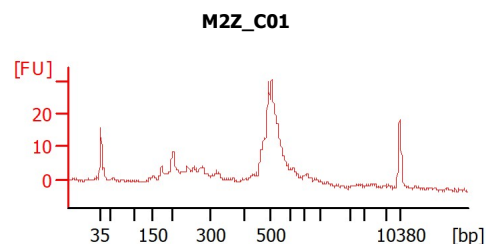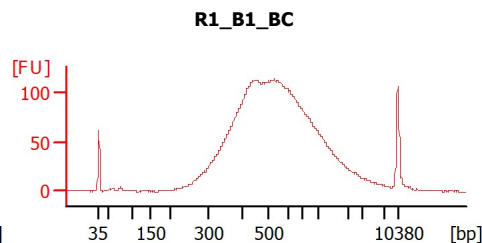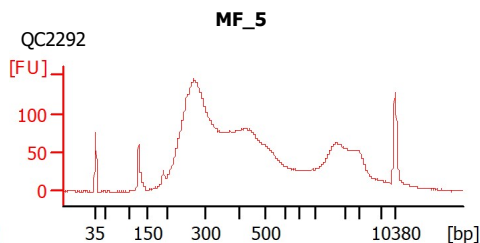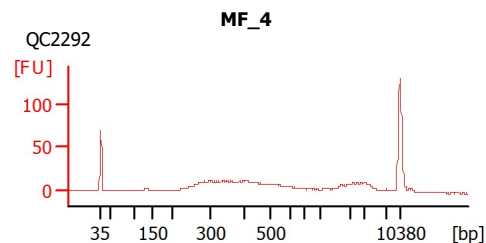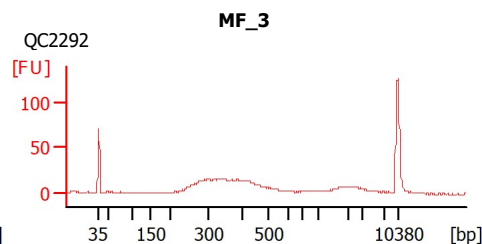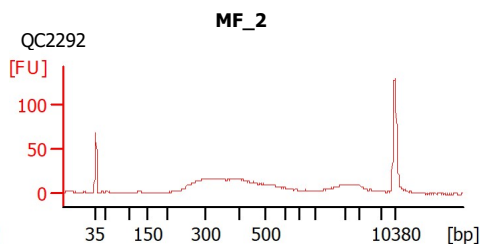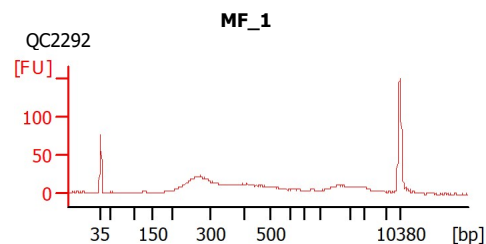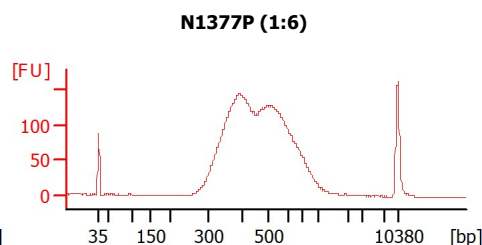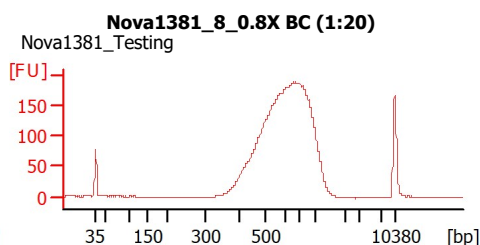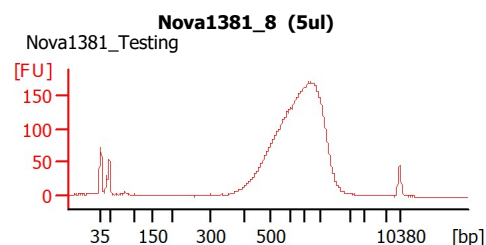

Assay Class: High Sensitivity DNA Assay
Data Path: C:\...gh Sensitivity DNA Assay_DE13701287_2025-09-09_16-06-55.xad

Created: 9/9/2025 4:06:55 PM
Modified: 9/10/2025 5:43:51 AM

Electrophoresis File Run Summary (Chip Summary)

Sample Name	Sample Comment	Rest. Digest	Status	Observation	Result Label	Result Color
Nova1381_8 (5ul)	Nova1381_Testing	☐	✓			
Nova1381_8_0.8X BC (1:20)	Nova1381_Testing	☐	✓			
N1377P (1:6)		☐	✓			
MF_1	QC2292	☐	✓			
MF_2	QC2292	☐	✓			
MF_3	QC2292	☐	✓			
MF_4	QC2292	☐	✓			
MF_5	QC2292	☐	✓			
R1_B1_BC		☐	✓			
M2Z_C01		☐	✓			
HM2025 (1:9)		☐	✓			
Ladder		☐	✓			

Chip Lot #**Reagent Kit Lot #****Chip Comments :**

Assay Class: High Sensitivity DNA Assay
Data Path: C:\...gh Sensitivity DNA Assay_DE13701287_2025-09-09_16-06-55.xad

Created: 9/9/2025 4:06:55 PM
Modified: 9/10/2025 5:43:51 AM

Electrophoresis Assay Details

General Analysis Settings

Number of Available Sample and Ladder Wells (Max.) : 12
Minimum Visible Range [s] : 32
Maximum Visible Range [s] : 138
Start Analysis Time Range [s] : 33
End Analysis Time Range [s] : 137.5
Ladder Concentration [pg/ μ l] : 1950
Uses Standard Area for Ladder Fragments
Lower Marker Concentration [pg/ μ l] : 125
Upper Marker Concentration [pg/ μ l] : 75
Used Upper Marker for Quantitation
Standard Curve Fit is Point to Point
Show Data Aligned to Lower and Upper Marker

Integrator Settings

Integration Start Time [s] : 33.05
Integration End Time [s] : 137
Slope Threshold : 0.8
Height Threshold [FU] : 5
Area Threshold : 0.1
Width Threshold [s] : 0.6
Baseline Plateau [s] : 0.5

Filter Settings

Filter Width [s] : 0.5
Polynomial Order : 4


Ladder

Ladder Peak	Size	Area
1	35	160
2	50	210
3	100	208
4	150	221
5	200	242
6	300	270
7	400	305
8	500	306
9	600	336
10	700	321
11	1000	366
12	2000	413
13	3000	411
14	7000	400
15	10380	214

Assay Class: High Sensitivity DNA Assay
 Data Path: C:\...gh Sensitivity DNA Assay_DE13701287_2025-09-09_16-06-55.xad

Created: 9/9/2025 4:06:55 PM
 Modified: 9/10/2025 5:43:51 AM

Electropherogram Summary

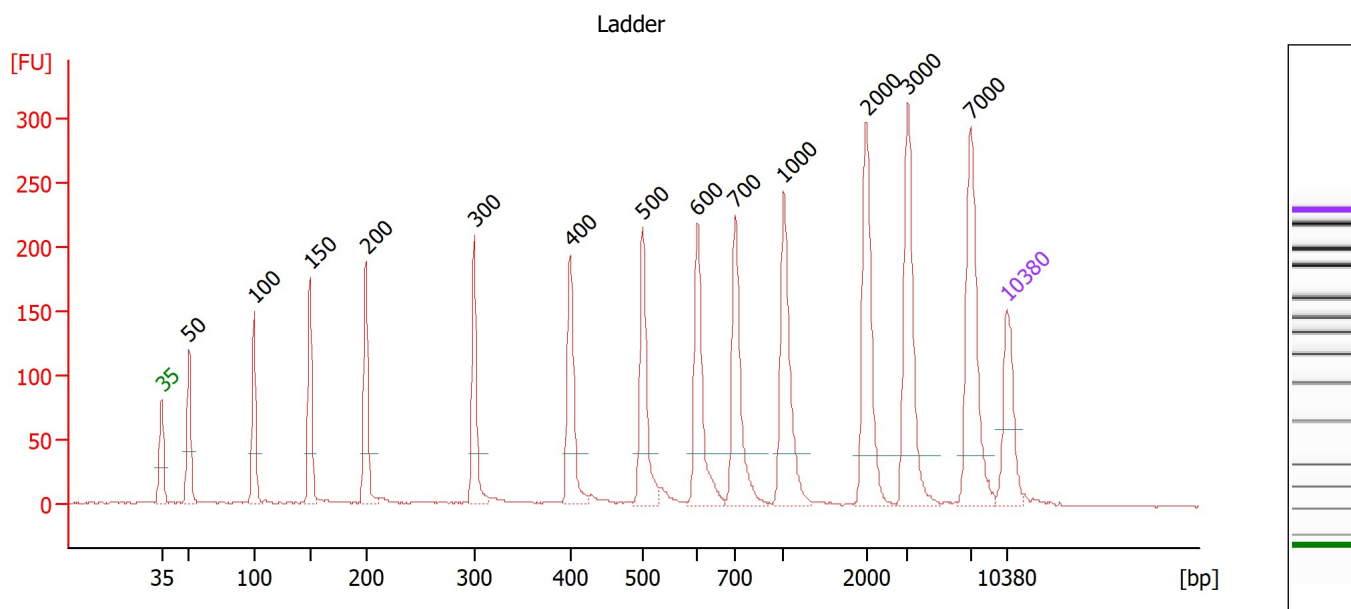

Overall Results for Ladder

Noise: 0.6


Peak table for Ladder

Peak	Size [bp]	Conc. [pg/μl]	Molarity [pmol/l]	Observations
1	35	125.00	5,411.3	Lower Marker
2	50	150.00	4,545.5	Ladder Peak
3	100	150.00	2,272.7	Ladder Peak
4	150	150.00	1,515.2	Ladder Peak
5	200	150.00	1,136.4	Ladder Peak
6	300	150.00	757.6	Ladder Peak
7	400	150.00	568.2	Ladder Peak
8	500	150.00	454.5	Ladder Peak
9	600	150.00	378.8	Ladder Peak
10	700	150.00	324.7	Ladder Peak
11	1,000	150.00	227.3	Ladder Peak
12	2,000	150.00	113.6	Ladder Peak
13	3,000	150.00	75.8	Ladder Peak
14	7,000	150.00	32.5	Ladder Peak
15	10,380	75.00	10.9	Upper Marker

Assay Class: High Sensitivity DNA Assay
 Data Path: C:\...gh Sensitivity DNA Assay_DE13701287_2025-09-09_16-06-55.xad

Created: 9/9/2025 4:06:55 PM
 Modified: 9/10/2025 5:43:51 AM

Electropherogram Summary Continued ...

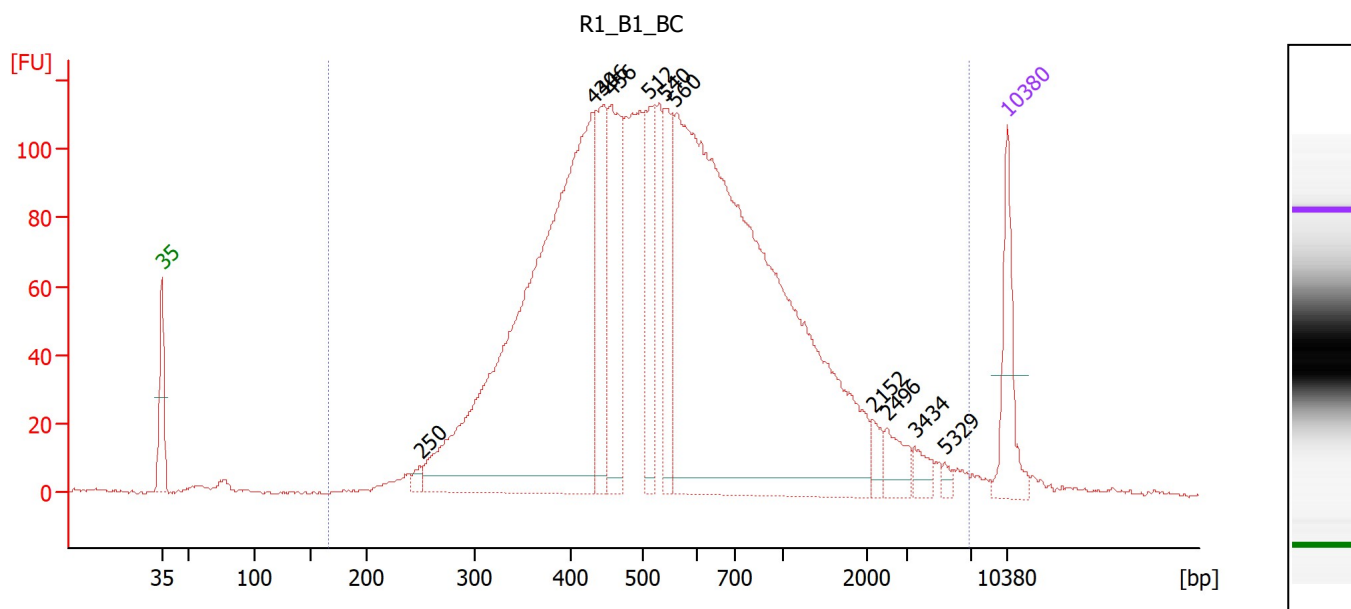

Overall Results for sample 9 : R1_B1_BC

Number of peaks found: 11 Corr. Area 1: 3,076.8
 Noise: 0.5

Peak table for sample 9 : R1_B1_BC

Peak	Size [bp]	Conc. [pg/μl]	Molarity [pmol/l]	Observations
1	35	125.00	5,411.3	Lower Marker
2	250	12.44	75.4	
3	430	1,073.36	3,780.0	
4	446	160.82	546.6	
5	456	195.27	648.2	
6	512	116.00	343.4	
7	540	110.56	310.0	
8	560	1,364.53	3,692.5	
9	2,152	17.54	12.4	
10	2,496	31.91	19.4	
11	3,434	17.53	7.7	
12	5,329	6.90	2.0	
13	10,380	75.00	10.9	Upper Marker


Region table for sample 9 : R1_B1_BC

From [bp]	To [bp]	Corr. Area	% of Total	Average Size [bp]	Size distribution in CV [%]	Conc. [pg/μl]	Molarity [pmol/l]	Color
166	6,971	3,076.8	97	752	99.4	3,579.87	10,956.2	Blue

Assay Class: High Sensitivity DNA Assay
 Data Path: C:\...gh Sensitivity DNA Assay_DE13701287_2025-09-09_16-06-55.xad

Created: 9/9/2025 4:06:55 PM
 Modified: 9/10/2025 5:43:51 AM

Electropherogram Summary Continued ...

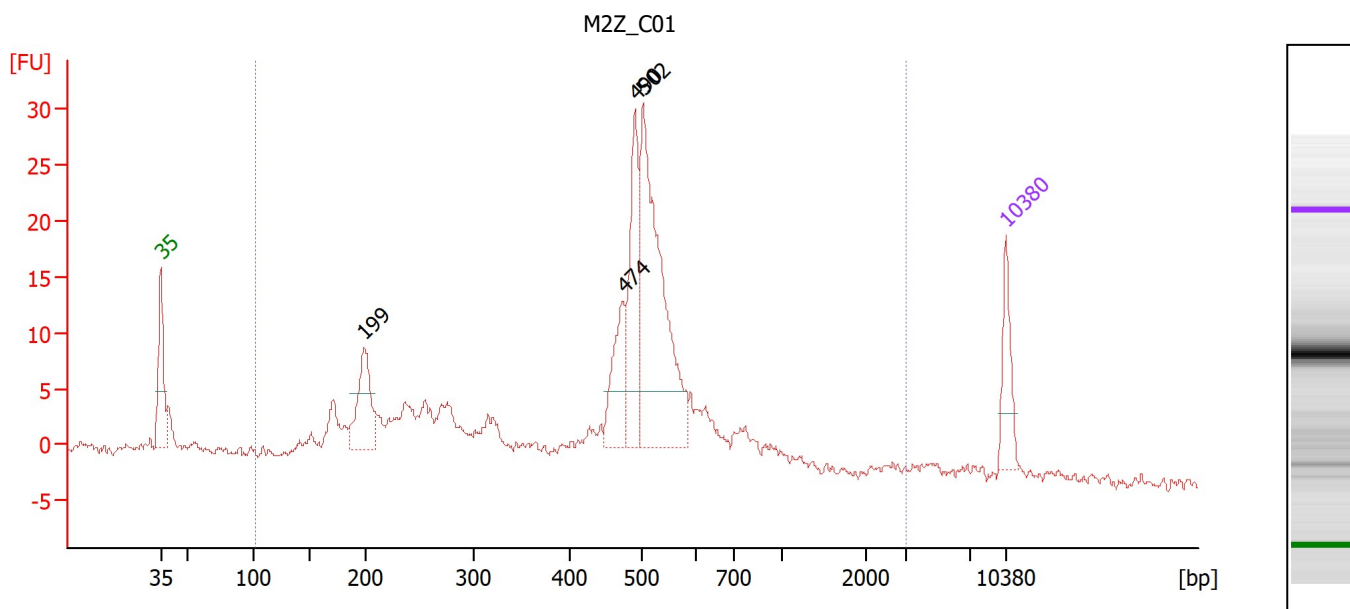

Overall Results for sample 10 : M2Z_C01

Number of peaks found: 4 Corr. Area 1: 283.9
 Noise: 0.5

Peak table for sample 10 : M2Z_C01

Peak	Size [bp]	Conc. [pg/μl]	Molarity [pmol/l]	Observations
1	35	125.00	5,411.3	Lower Marker
2	199	172.28	1,308.8	
3	474	134.06	428.8	
4	490	221.12	683.6	
5	502	555.54	1,676.0	
6	10,380	75.00	10.9	Upper Marker


Region table for sample 10 : M2Z_C01

From [bp]	To [bp]	Corr. Area	% of Total	Average Size [bp]	Size distribution in CV [%]	Conc. [pg/μl]	Molarity [pmol/l]	Color
103	2,986	283.9	90	502	64.9	2,326.11	10,151.5	Blue

Assay Class: High Sensitivity DNA Assay
 Data Path: C:\...gh Sensitivity DNA Assay_DE13701287_2025-09-09_16-06-55.xad

Created: 9/9/2025 4:06:55 PM
 Modified: 9/10/2025 5:43:51 AM

Electropherogram Summary Continued ...

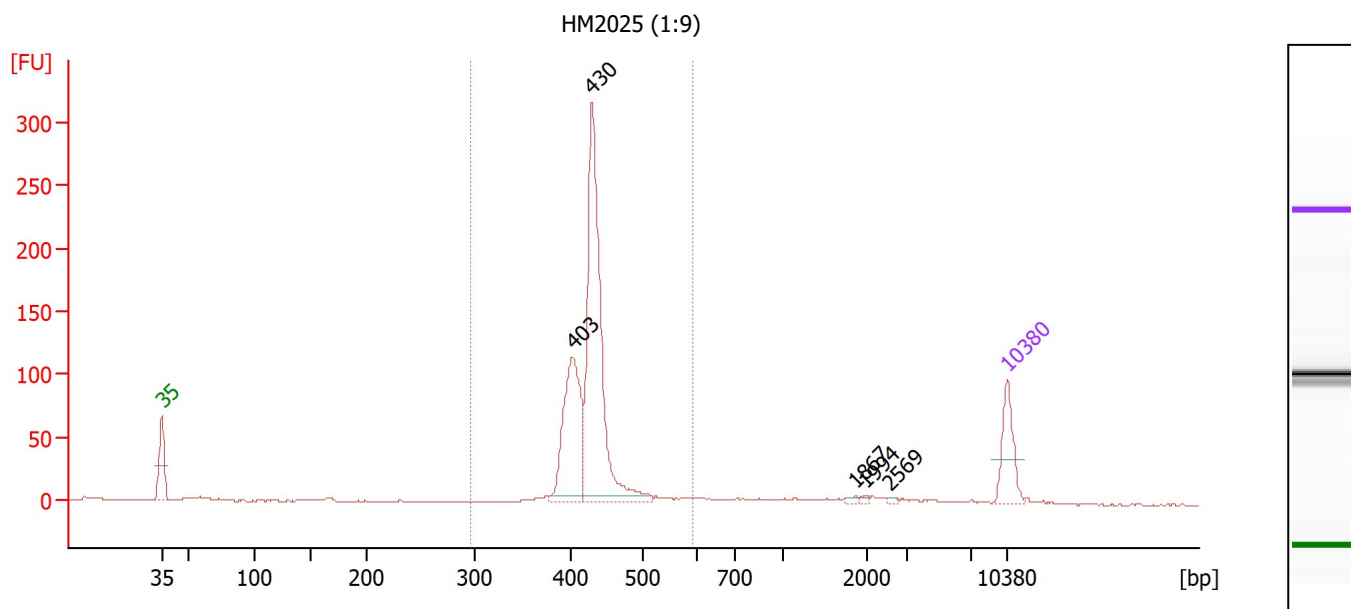

Overall Results for sample 11 : HM2025 (1:9)

Number of peaks found: 5 Corr. Area 1: 750.6
 Noise: 0.5

Peak table for sample 11 : HM2025 (1:9)


Peak	Size [bp]	Conc. [pg/μl]	Molarity [pmol/l]	Observations
1	35	125.00	5,411.3	Lower Marker
2	403	275.52	1,036.6	
3	430	566.15	1,996.4	
4	1,867	4.54	3.7	
5	1,994	3.48	2.6	
6	2,569	3.84	2.3	
7	10,380	75.00	10.9	Upper Marker

Region table for sample 11 : HM2025 (1:9)

From [bp]	To [bp]	Corr. Area	% of Total	Average Size [bp]	Size distribution in CV [%]	Conc. [pg/μl]	Molarity [pmol/l]	Color
296	592	750.6	85	428	6.9	874.47	3,105.7	Blue

Assay Class: High Sensitivity DNA Assay
Data Path: C:\...gh Sensitivity DNA Assay_DE13701287_2025-09-09_16-06-55.xad

Created: 9/9/2025 4:06:55 PM
Modified: 9/10/2025 5:43:51 AM

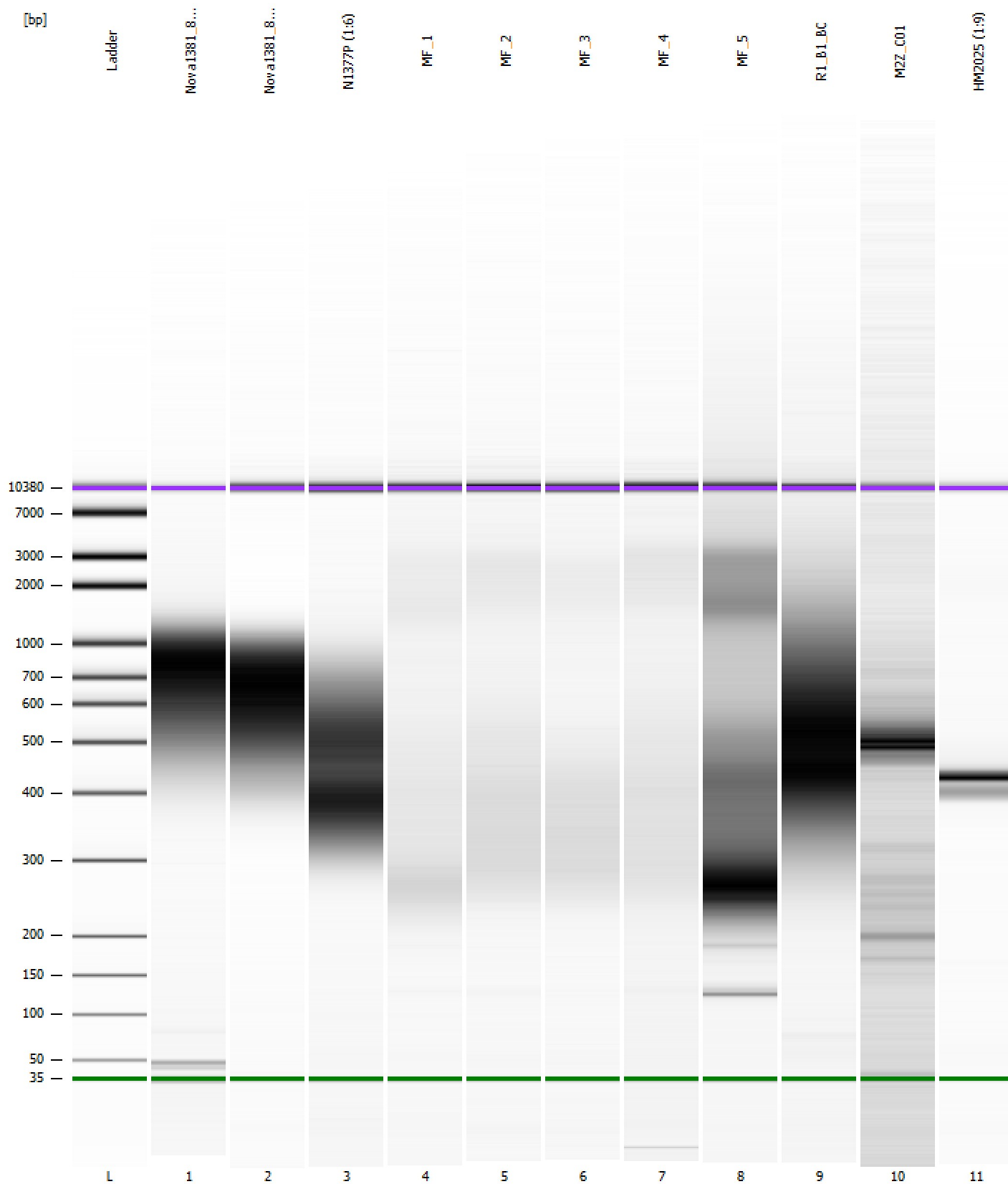Gel Image

Assay Class: High Sensitivity DNA Assay
 Data Path: C:\...gh Sensitivity DNA Assay_DE13701287_2025-09-09_16-06-55.xad

Created: 9/9/2025 4:06:55 PM
 Modified: 9/10/2025 5:43:51 AM

Run Logbook

Description	Number	Source	Category	Sub Category	Time	Time Zone	User	Host
Run ended on port 3 (Number of wells acquired: 12)		Instrument	Run		9/9/2025 4:48:18 PM	(GMT --07:00) Pacific Standard Time	dnat	CBS-GC1402RT ALP
Run started on port 3 (File: C:\Program Files (x86)\Agilent\2100 bioanalyzer\2100 expert\Data\2025-09-09\2100 expert_High Sensitivity DNA Assay_DE13701287_2025-09-09_16-06-55.xad)		Instrument	Run		9/9/2025 4:07:01 PM	(GMT --07:00) Pacific Standard Time	dnat	CBS-GC1402RT ALP
Product Number : G2938B		Instrument	Run		9/9/2025 4:07:01 PM	(GMT --07:00) Pacific Standard Time	dnat	CBS-GC1402RT ALP
Name : BioA#2		Instrument	Run		9/9/2025 4:07:01 PM	(GMT --07:00) Pacific Standard Time	dnat	CBS-GC1402RT ALP
Vendor : Agilent Technologies		Instrument	Run		9/9/2025 4:07:01 PM	(GMT --07:00) Pacific Standard Time	dnat	CBS-GC1402RT ALP
Serial# : DE13701287		Instrument	Run		9/9/2025 4:07:01 PM	(GMT --07:00) Pacific Standard Time	dnat	CBS-GC1402RT ALP
Firmware : C.01.069		Instrument	Run		9/9/2025 4:07:01 PM	(GMT --07:00) Pacific Standard Time	dnat	CBS-GC1402RT ALP
Cartridge : Electrode		Instrument	Run		9/9/2025 4:07:01 PM	(GMT --07:00) Pacific Standard Time	dnat	CBS-GC1402RT ALP