


Assay Class: High Sensitivity DNA Assay
Data Path: C:\...100 bioanalyzer\2100 expert\data\2025-04-03\36_NARMS_BC.xad

Created: 4/3/2025 11:16:18 AM
Modified: 4/3/2025 12:24:01 PM

Electrophoresis File Run Summary

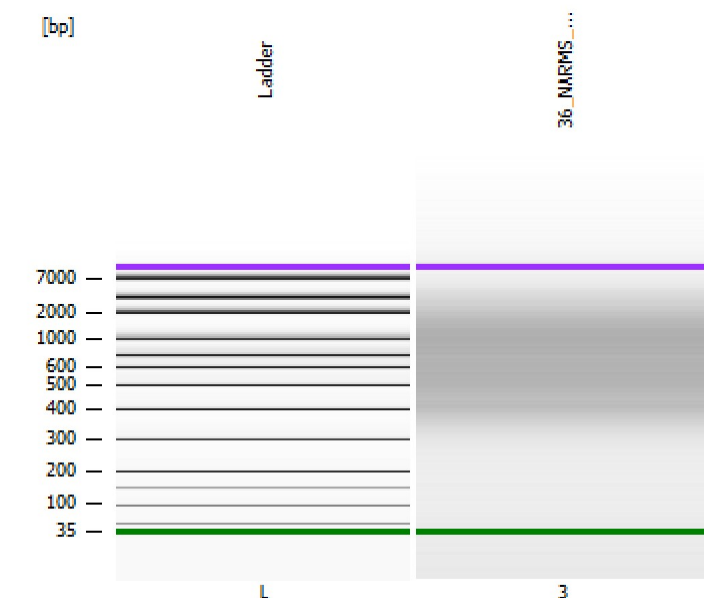

Instrument Information:


Instrument Name: DE34903152 Firmware: C.01.069
Serial#: DE34903152 Type: G2938C

Assay Information:

Assay Origin Path: C:\Program Files (x86)\Agilent\2100 bioanalyzer\2100 expert\assays\dsDNA\High Sensitivity DNA.xsy
Assay Class: High Sensitivity DNA Assay
Version: 1.03
Assay Comments: Copyright © 2003-2010 Agilent Technologies

Chip Information:

Chip Lot #:
Reagent Kit Lot #:
Chip Comments:

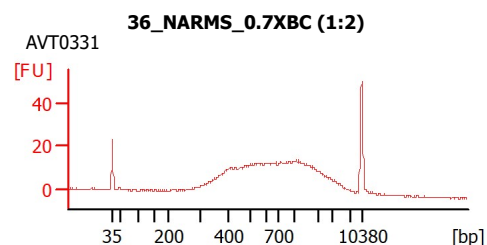

Assay Class: High Sensitivity DNA Assay
Data Path: C:\...100 bioanalyzer\2100 expert\data\2025-04-03\36_NARMS_BC.xad

Created: 4/3/2025 11:16:18 AM
Modified: 4/3/2025 12:24:01 PM

Electrophoresis File Run Summary (Chip Summary)

Sample Name	Sample Comment	Rest. Digest	Status	Observation	Result Label	Result Color
36_NARMS_0.7XBC (1:2)	AVT0331	☐	✓			
Ladder		☐	✓			

Chip Lot #**Reagent Kit Lot #****Chip Comments :**

Assay Class: High Sensitivity DNA Assay
Data Path: C:\...100 bioanalyzer\2100 expert\data\2025-04-03\36_NARMS_BC.xad

Created: 4/3/2025 11:16:18 AM
Modified: 4/3/2025 12:24:01 PM

Electrophoresis Assay Details

General Analysis Settings

Number of Available Sample and Ladder Wells (Max.) : 12
Minimum Visible Range [s] : 32
Maximum Visible Range [s] : 138
Start Analysis Time Range [s] : 33
End Analysis Time Range [s] : 137.5
Ladder Concentration [pg/μl] : 1950
Uses Standard Area for Ladder Fragments
Lower Marker Concentration [pg/μl] : 125
Upper Marker Concentration [pg/μl] : 75
Used Upper Marker for Quantitation
Standard Curve Fit is Point to Point
Show Data Aligned to Lower and Upper Marker

Integrator Settings

Integration Start Time [s] : 33.05
Integration End Time [s] : 137
Slope Threshold : 0.8
Height Threshold [FU] : 5
Area Threshold : 0.1
Width Threshold [s] : 0.6
Baseline Plateau [s] : 0.5

Filter Settings

Filter Width [s] : 0.5
Polynomial Order : 4


Ladder

Ladder Peak	Size	Area
1	35	160
2	50	210
3	100	208
4	150	221
5	200	242
6	300	270
7	400	305
8	500	306
9	600	336
10	700	321
11	1000	366
12	2000	413
13	3000	411
14	7000	400
15	10380	214

Assay Class: High Sensitivity DNA Assay
 Data Path: C:\...100 bioanalyzer\2100 expert\data\2025-04-03\36_NARMS_BC.xad

Created: 4/3/2025 11:16:18 AM
 Modified: 4/3/2025 12:24:01 PM

Electropherogram Summary

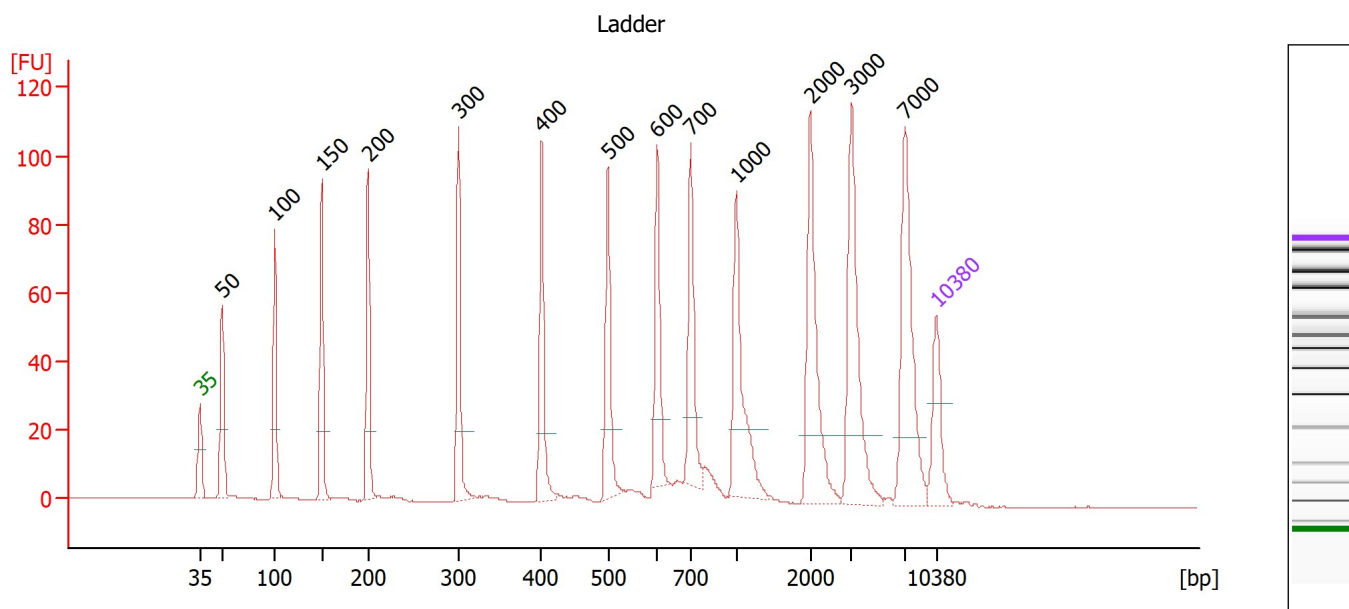

Overall Results for Ladder


Noise: 0.1

Peak table for Ladder

Peak		Size [bp]	Conc. [pg/μl]	Molarity [pmol/l]	Observations
1	◀	35	125.00	5,411.3	Lower Marker
2	█	50	150.00	4,545.5	Ladder Peak
3	█	100	150.00	2,272.7	Ladder Peak
4	█	150	150.00	1,515.2	Ladder Peak
5	█	200	150.00	1,136.4	Ladder Peak
6	█	300	150.00	757.6	Ladder Peak
7	█	400	150.00	568.2	Ladder Peak
8	█	500	150.00	454.5	Ladder Peak
9	█	600	150.00	378.8	Ladder Peak
10	█	700	150.00	324.7	Ladder Peak
11	█	1,000	150.00	227.3	Ladder Peak
12	█	2,000	150.00	113.6	Ladder Peak
13	█	3,000	150.00	75.8	Ladder Peak
14	█	7,000	150.00	32.5	Ladder Peak
15	▶	10,380	75.00	10.9	Upper Marker

Assay Class: High Sensitivity DNA Assay
Data Path: C:\...100 bioanalyzer\2100 expert\data\2025-04-03\36_NARMS_BC.xad

Created: 4/3/2025 11:16:18 AM
Modified: 4/3/2025 12:24:01 PM

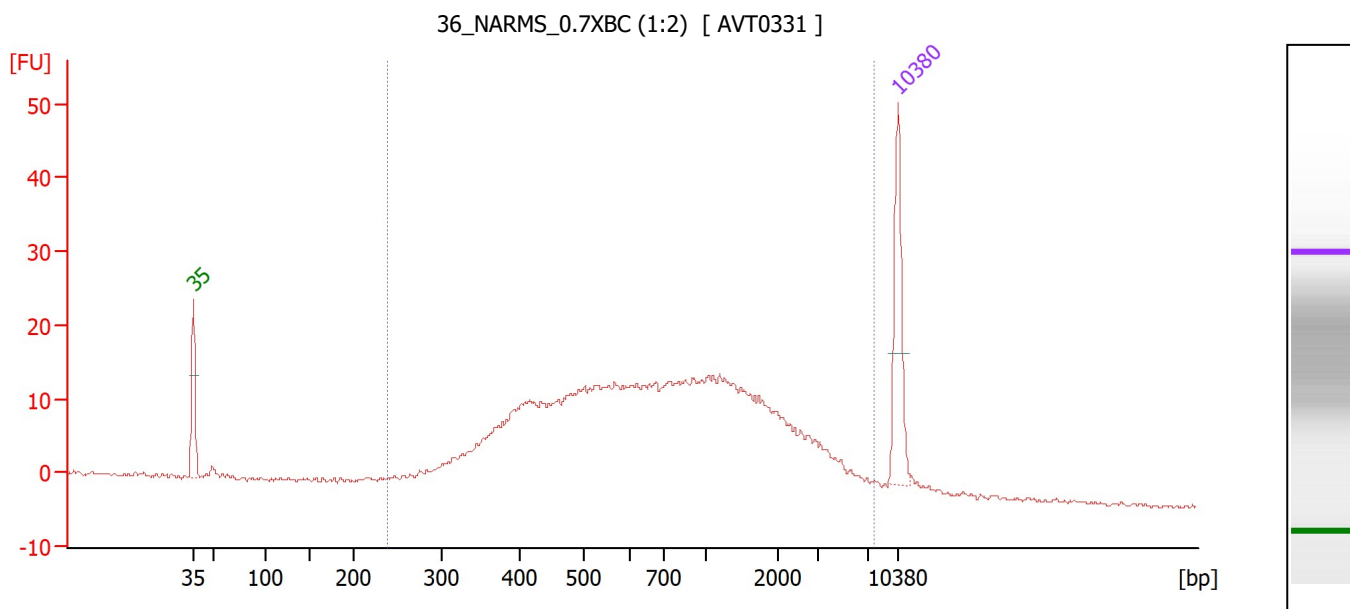Electropherogram Summary Continued ...**Overall Results for sample 3 : 36_NARMS_0.7XBC (1:2)**

Number of peaks found: 0 Corr. Area 1: 505.2
Noise: 0.3

Peak table for sample 3 : 36_NARMS_0.7XBC (1:2)


Peak	Size [bp]	Conc. [pg/μl]	Molarity [pmol/l]	Observations
1	35	125.00	5,411.3	Lower Marker
2	10,380	75.00	10.9	Upper Marker

Region table for sample 3 : 36_NARMS_0.7XBC (1:2)

From [bp]	To [bp]	Corr. Area	% of Total	Average Size [bp]	Size distribution in CV [%]	Conc. [pg/μl]	Molarity [pmol/l]	Color
239	7,656	505.2	96	1,118	96.4	1,281.61	3,196.6	Blue

Assay Class: High Sensitivity DNA Assay
Data Path: C:\...100 bioanalyzer\2100 expert\data\2025-04-03\36_NARMS_BC.xad

Created: 4/3/2025 11:16:18 AM
Modified: 4/3/2025 12:24:01 PM

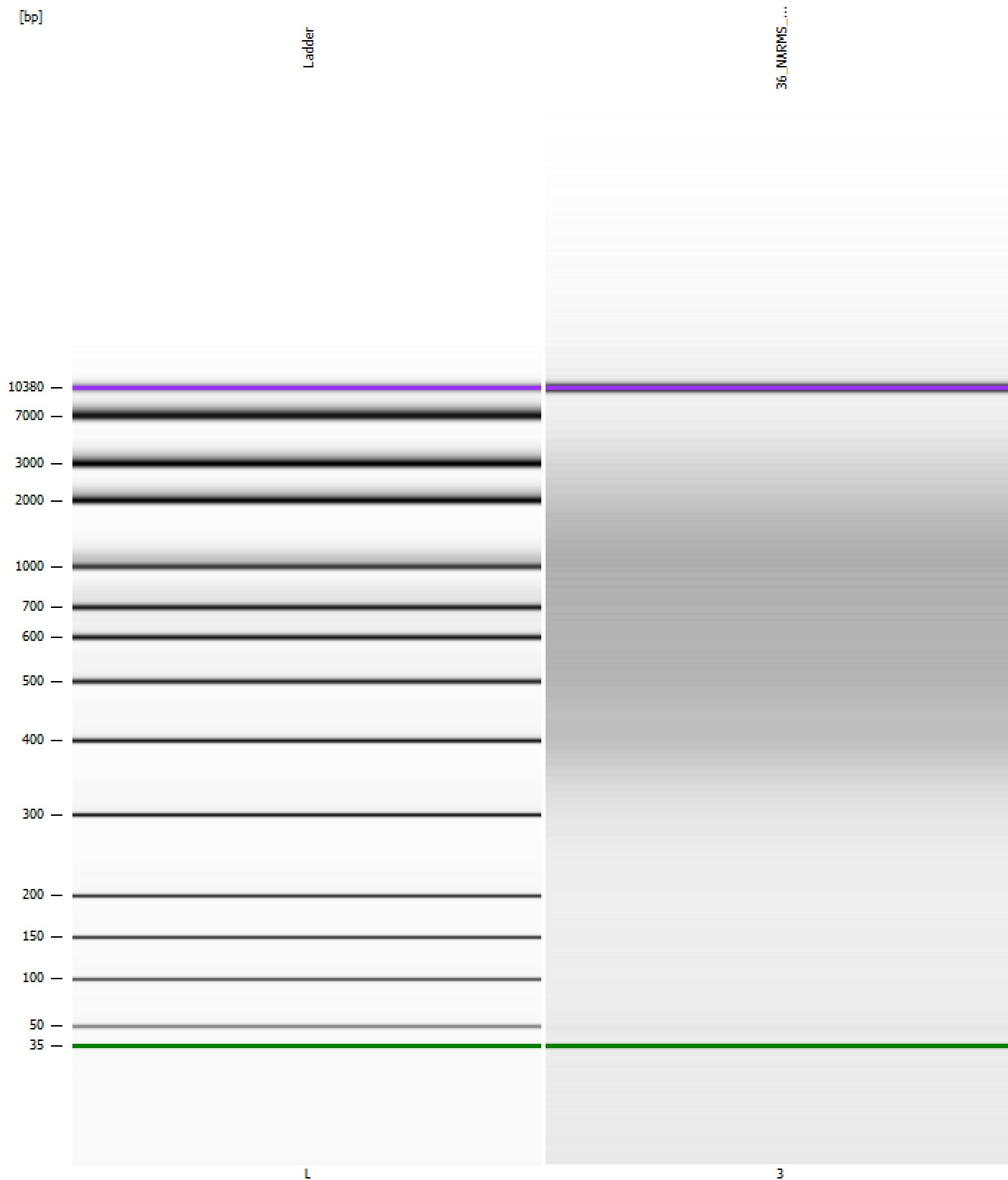Gel Image

Assay Class: High Sensitivity DNA Assay
 Data Path: C:\...100 bioanalyzer\2100 expert\data\2025-04-03\36_NARMS_BC.xad

Created: 4/3/2025 11:16:18 AM
 Modified: 4/3/2025 12:24:01 PM

Run Logbook

Description	Number	Source	Category	Sub Category	Time	Time Zone	User	Host
Run ended on port 1 (Number of wells acquired: 12)		Instrument	Run		4/3/2025 11:57:39 AM	(GMT --07:00) Pacific Standard Time	sbsuser	DESKTOP-4UNV VOB
Run started on port 1 (File: C:\Program Files (x86)\Agilent\2100 bioanalyzer\2100 expert\data\2025-04-03\Bioanalyzer1_High Sensitivity DNA Assay_DE34903152_2025-04-03_001.xad)		Instrument	Run		4/3/2025 11:16:25 AM	(GMT --07:00) Pacific Standard Time	sbsuser	DESKTOP-4UNV VOB
Product Number : G2938C		Instrument	Run		4/3/2025 11:16:25 AM	(GMT --07:00) Pacific Standard Time	sbsuser	DESKTOP-4UNV VOB
Name :		Instrument	Run		4/3/2025 11:16:25 AM	(GMT --07:00) Pacific Standard Time	sbsuser	DESKTOP-4UNV VOB
Vendor : Agilent Technologies		Instrument	Run		4/3/2025 11:16:25 AM	(GMT --07:00) Pacific Standard Time	sbsuser	DESKTOP-4UNV VOB
Serial# : DE34903152		Instrument	Run		4/3/2025 11:16:25 AM	(GMT --07:00) Pacific Standard Time	sbsuser	DESKTOP-4UNV VOB
Firmware : C.01.069		Instrument	Run		4/3/2025 11:16:25 AM	(GMT --07:00) Pacific Standard Time	sbsuser	DESKTOP-4UNV VOB
Cartridge : Electrode		Instrument	Run		4/3/2025 11:16:25 AM	(GMT --07:00) Pacific Standard Time	sbsuser	DESKTOP-4UNV VOB