


Assay Class: High Sensitivity DNA Assay
Data Path: C:\...2100 expert\data\2025-05-19\A0682_A0683_A0684_Rapid_PCR.xad

Created: 5/19/2025 12:50:37 PM
Modified: 5/19/2025 1:38:19 PM

Electrophoresis File Run Summary

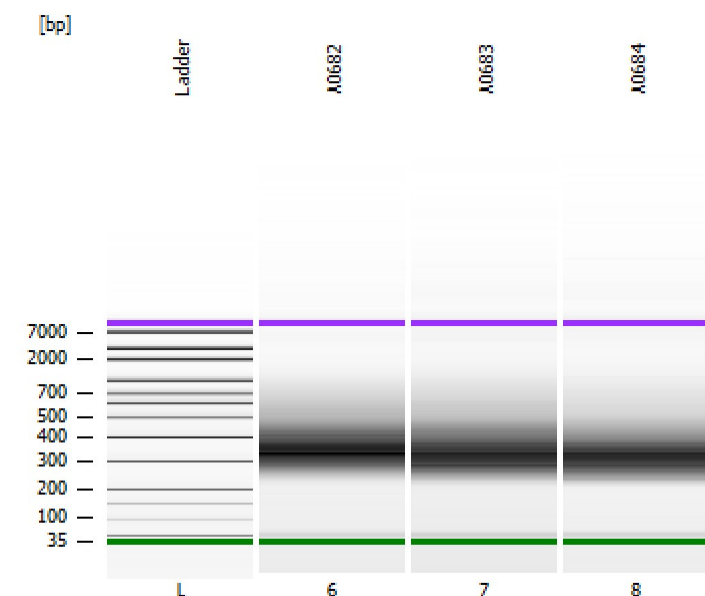

Instrument Information:




Instrument Name: DE34903152 Firmware: C.01.069
Serial#: DE34903152 Type: G2938C

Assay Information:

Assay Origin Path: C:\Program Files (x86)\Agilent\2100 bioanalyzer\2100 expert\assays\dsDNA\High Sensitivity DNA.xsy
Assay Class: High Sensitivity DNA Assay
Version: 1.03
Assay Comments: Copyright © 2003-2010 Agilent Technologies

Chip Information:

Chip Lot #:
Reagent Kit Lot #:
Chip Comments:

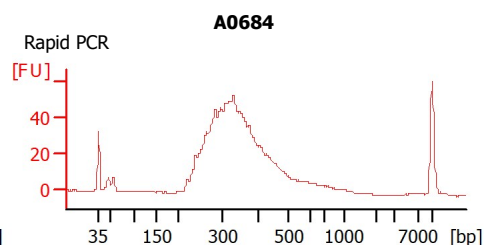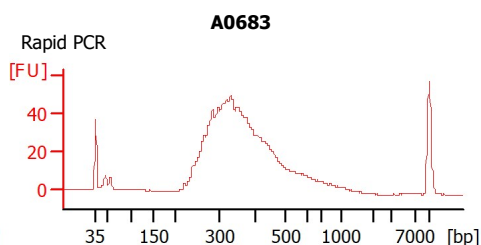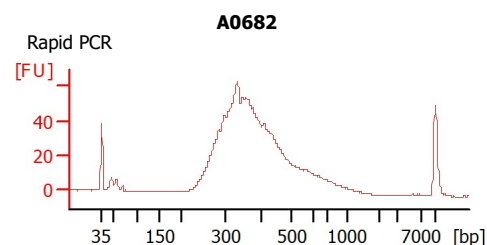

Assay Class: High Sensitivity DNA Assay
Data Path: C:\...2100 expert\data\2025-05-19\A0682_A0683_A0684_Rapid_PCR.xad

Created: 5/19/2025 12:50:37 PM
Modified: 5/19/2025 1:38:19 PM

Electrophoresis File Run Summary (Chip Summary)

Sample Name	Sample Comment	Rest. Digest	Status	Observation	Result Label	Result Color
A0682	Rapid PCR	☐	✓			
A0683	Rapid PCR	☐	✓			
A0684	Rapid PCR	☐	✓			
Ladder		☐	✓			

Chip Lot #**Reagent Kit Lot #****Chip Comments :**

Assay Class: High Sensitivity DNA Assay
Data Path: C:\...2100 expert\data\2025-05-19\A0682_A0683_A0684_Rapid_PCR.xad

Created: 5/19/2025 12:50:37 PM
Modified: 5/19/2025 1:38:19 PM

Electrophoresis Assay Details

General Analysis Settings

Number of Available Sample and Ladder Wells (Max.) : 12
Minimum Visible Range [s] : 32
Maximum Visible Range [s] : 138
Start Analysis Time Range [s] : 33
End Analysis Time Range [s] : 137.5
Ladder Concentration [pg/ μ l] : 1950
Uses Standard Area for Ladder Fragments
Lower Marker Concentration [pg/ μ l] : 125
Upper Marker Concentration [pg/ μ l] : 75
Used Upper Marker for Quantitation
Standard Curve Fit is Point to Point
Show Data Aligned to Lower and Upper Marker

Integrator Settings

Integration Start Time [s] : 33.05
Integration End Time [s] : 137
Slope Threshold : 0.8
Height Threshold [FU] : 5
Area Threshold : 0.1
Width Threshold [s] : 0.6
Baseline Plateau [s] : 0.5

Filter Settings

Filter Width [s] : 0.5
Polynomial Order : 4


Ladder

Ladder Peak	Size	Area
1	35	160
2	50	210
3	100	208
4	150	221
5	200	242
6	300	270
7	400	305
8	500	306
9	600	336
10	700	321
11	1000	366
12	2000	413
13	3000	411
14	7000	400
15	10380	214

Assay Class: High Sensitivity DNA Assay
 Data Path: C:\...2100 expert\data\2025-05-19\A0682_A0683_A0684_Rapid_PCR.xad

Created: 5/19/2025 12:50:37 PM
 Modified: 5/19/2025 1:38:19 PM

Electropherogram Summary

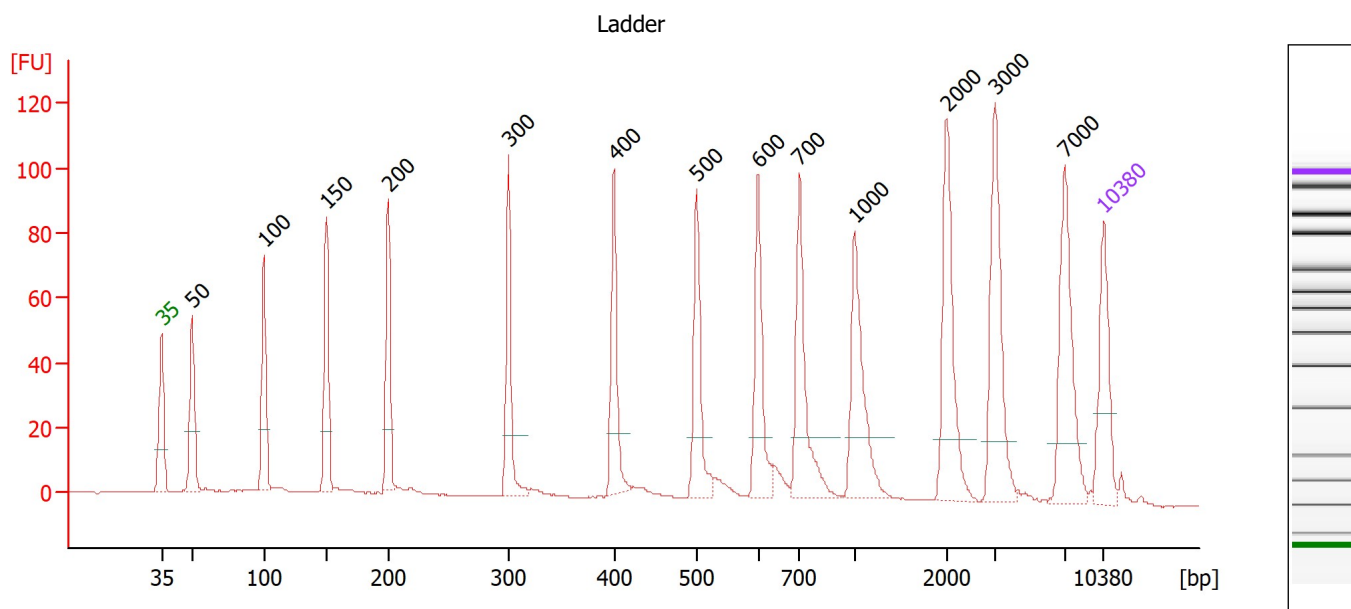

Overall Results for Ladder

Noise: 0.2


Peak table for Ladder

Peak	Size [bp]	Conc. [pg/μl]	Molarity [pmol/l]	Observations
1	35	125.00	5,411.3	Lower Marker
2	50	150.00	4,545.5	Ladder Peak
3	100	150.00	2,272.7	Ladder Peak
4	150	150.00	1,515.2	Ladder Peak
5	200	150.00	1,136.4	Ladder Peak
6	300	150.00	757.6	Ladder Peak
7	400	150.00	568.2	Ladder Peak
8	500	150.00	454.5	Ladder Peak
9	600	150.00	378.8	Ladder Peak
10	700	150.00	324.7	Ladder Peak
11	1,000	150.00	227.3	Ladder Peak
12	2,000	150.00	113.6	Ladder Peak
13	3,000	150.00	75.8	Ladder Peak
14	7,000	150.00	32.5	Ladder Peak
15	10,380	75.00	10.9	Upper Marker

Assay Class: High Sensitivity DNA Assay
 Data Path: C:\...2100 expert\data\2025-05-19\A0682_A0683_A0684_Rapid_PCR.xad

Created: 5/19/2025 12:50:37 PM
 Modified: 5/19/2025 1:38:19 PM

Electropherogram Summary Continued ...

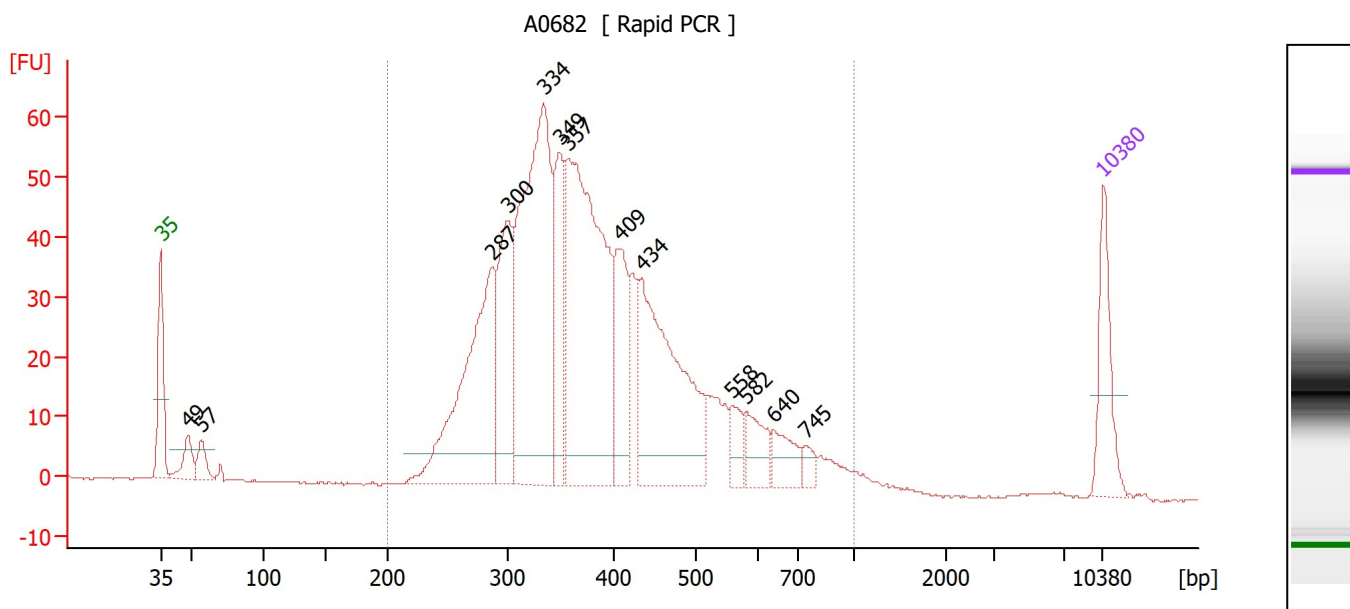

Overall Results for sample 6 : A0682

Number of peaks found: 13 Corr. Area 1: 1,002.8
 Noise: 0.2

Peak table for sample 6 : A0682

Peak	Size [bp]	Conc. [pg/μl]	Molarity [pmol/l]	Observations
1	35	125.00	5,411.3	Lower Marker
2	49	43.74	1,366.1	
3	57	32.94	874.2	
4	287	366.38	1,932.2	
5	300	196.67	994.4	
6	334	548.70	2,487.9	
7	349	143.76	624.0	
8	357	547.00	2,321.9	
9	409	131.41	486.6	
10	434	336.37	1,174.1	
11	558	32.44	88.1	
12	582	47.31	123.1	
13	640	40.53	96.0	
14	745	17.66	35.9	
15	10,380	75.00	10.9	Upper Marker


Region table for sample 6 : A0682

From [bp]	To [bp]	Corr. Area	% of Total	Average Size [bp]	Size distribution in CV [%]	Conc. [pg/μl]	Molarity [pmol/l]	Color
200	1,000	1,002.8	97	396	29.4	2,596.19	10,882.8	Blue

Assay Class: High Sensitivity DNA Assay
 Data Path: C:\...2100 expert\data\2025-05-19\A0682_A0683_A0684_Rapid_PCR.xad

Created: 5/19/2025 12:50:37 PM
 Modified: 5/19/2025 1:38:19 PM

Electropherogram Summary Continued ...

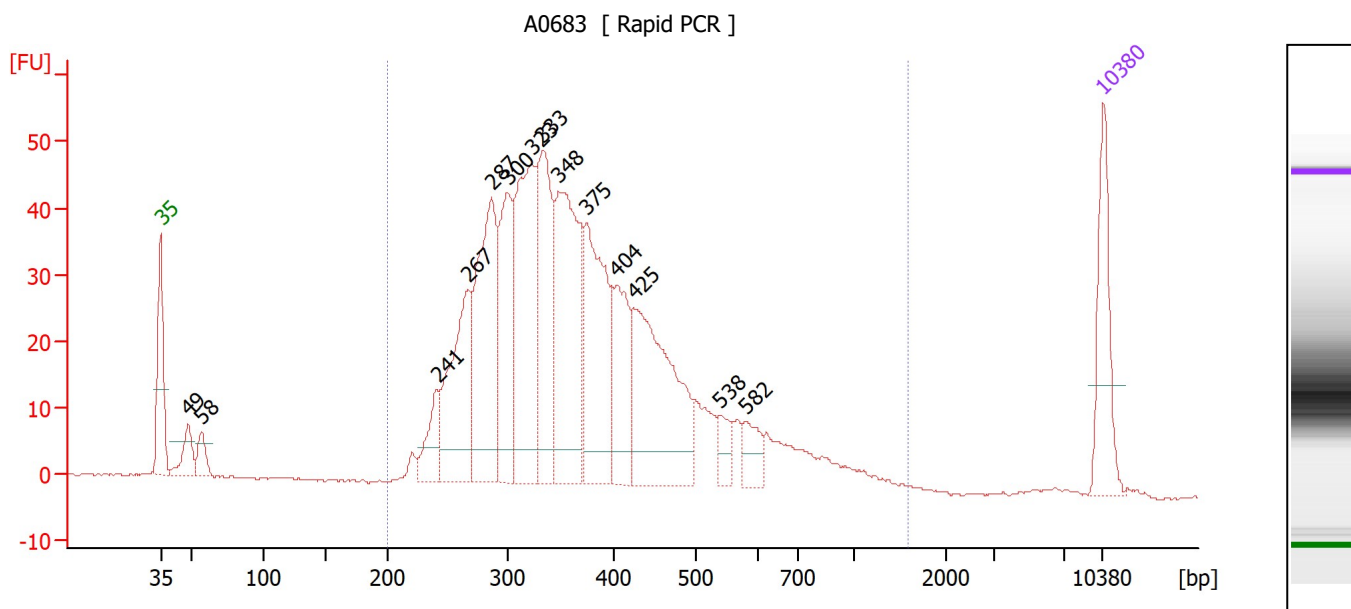

Overall Results for sample 7 : A0683

Number of peaks found: 14 Corr. Area 1: 898.0
 Noise: 0.2

Peak table for sample 7 : A0683

Peak	Size [bp]	Conc. [pg/μl]	Molarity [pmol/l]	Observations
1	35	125.00	5,411.3	Lower Marker
2	49	38.88	1,209.1	
3	58	27.97	736.9	
4	241	51.31	322.2	
5	267	176.18	999.5	
6	287	248.81	1,315.4	
7	300	155.41	785.3	
8	323	252.75	1,185.8	
9	333	181.63	827.7	
10	348	255.37	1,110.4	
11	375	205.16	829.5	
12	404	108.78	407.6	
13	425	236.65	844.4	
14	538	26.74	75.4	
15	582	30.66	79.9	
16	10,380	75.00	10.9	Upper Marker


Region table for sample 7 : A0683

From [bp]	To [bp]	Corr. Area	% of Total	Average Size [bp]	Size distribution in CV [%]	Conc. [pg/μl]	Molarity [pmol/l]	Color
200	1,593	898.0	97	389	36.5	2,099.63	9,157.0	Blue

Assay Class: High Sensitivity DNA Assay
 Data Path: C:\...2100 expert\data\2025-05-19\A0682_A0683_A0684_Rapid_PCR.xad

Created: 5/19/2025 12:50:37 PM
 Modified: 5/19/2025 1:38:19 PM

Electropherogram Summary Continued ...

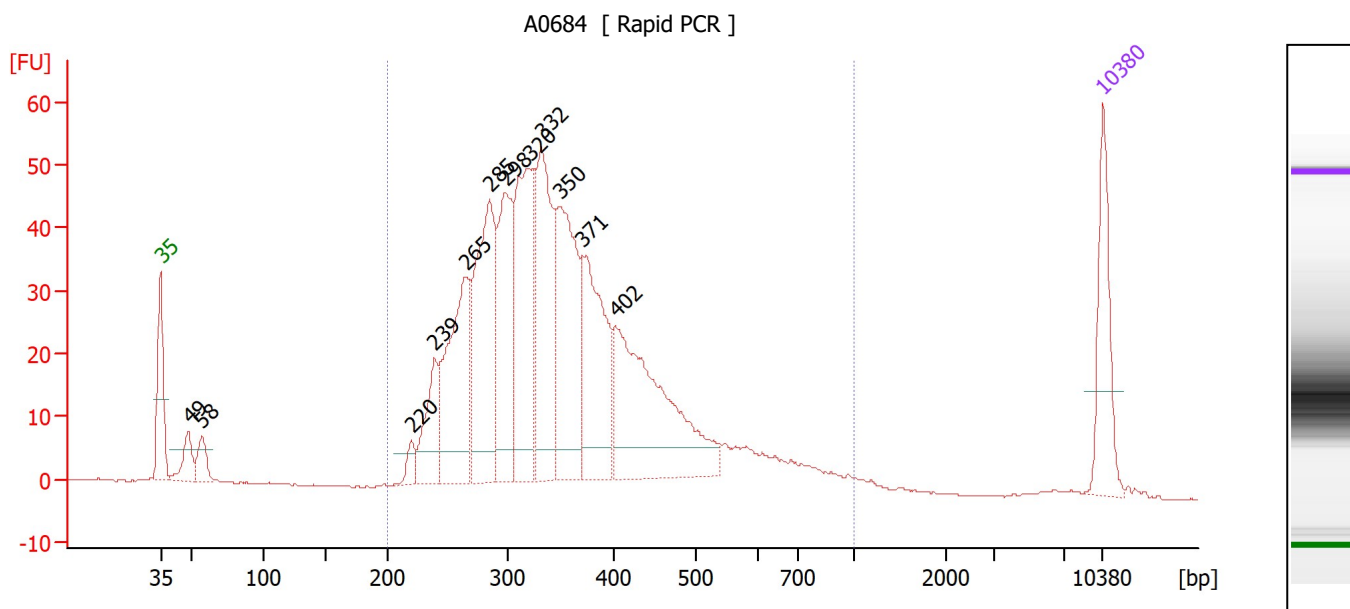

Overall Results for sample 8 : A0684

Number of peaks found: 12 Corr. Area 1: 875.1
 Noise: 0.1

Peak table for sample 8 : A0684


Peak	Size [bp]	Conc. [pg/μl]	Molarity [pmol/l]	Observations
1	35	125.00	5,411.3	Lower Marker
2	49	39.57	1,226.3	
3	58	31.08	813.7	
4	220	19.78	136.1	
5	239	83.39	528.1	
6	265	204.85	1,172.8	
7	285	244.47	1,299.1	
8	298	169.67	863.2	
9	320	231.99	1,098.0	
10	332	215.57	983.1	
11	350	217.33	940.8	
12	371	180.93	739.0	
13	402	272.58	1,027.1	
14	10,380	75.00	10.9	Upper Marker

Region table for sample 8 : A0684

From [bp]	To [bp]	Corr. Area	% of Total	Average Size [bp]	Size distribution in CV [%]	Conc. [pg/μl]	Molarity [pmol/l]	Color
200	1,000	875.1	97	363	29.8	2,017.86	9,207.1	Blue

Assay Class: High Sensitivity DNA Assay
Data Path: C:\...2100 expert\data\2025-05-19\A0682_A0683_A0684_Rapid_PCR.xad

Created: 5/19/2025 12:50:37 PM
Modified: 5/19/2025 1:38:19 PM

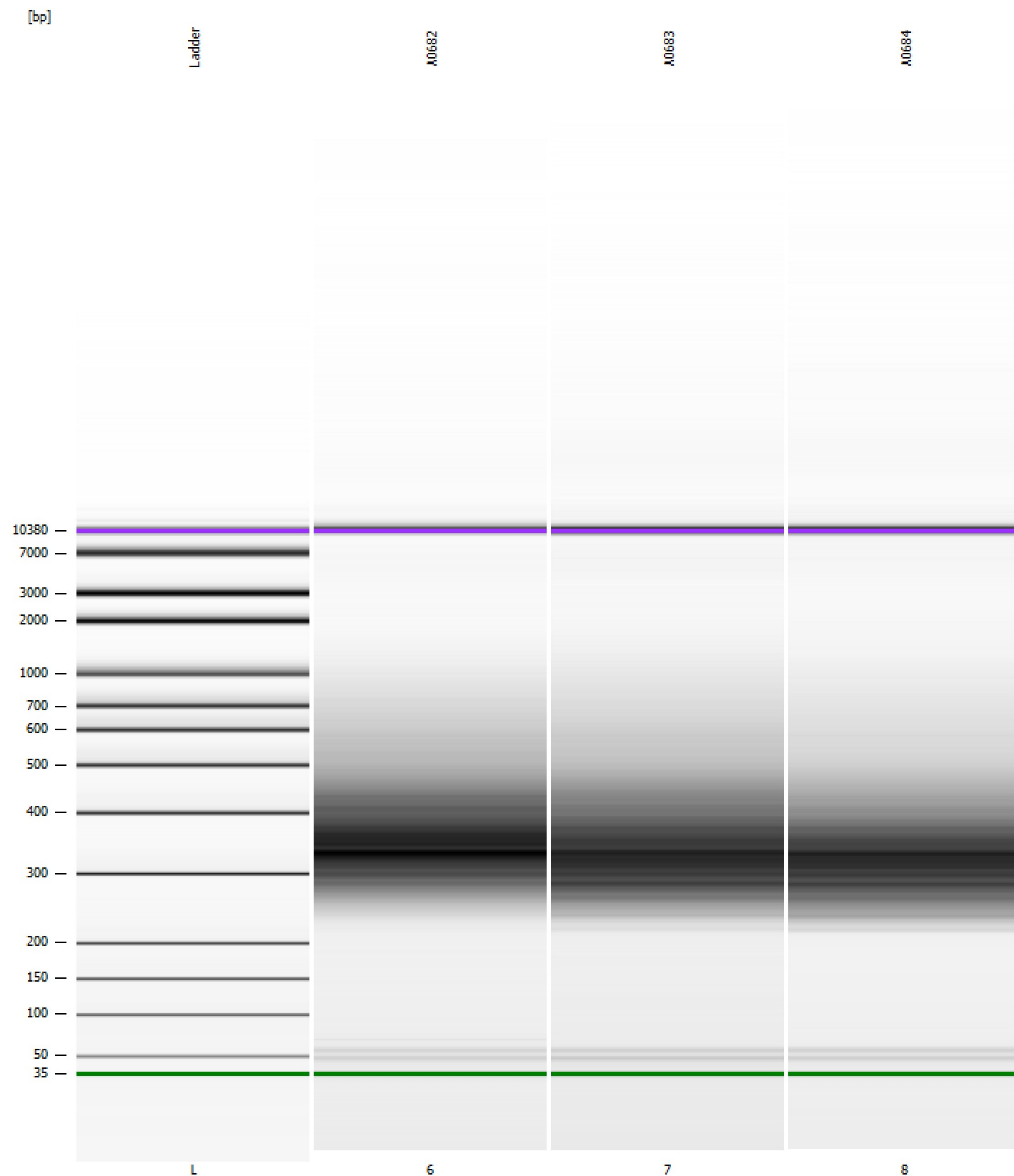Gel Image

Assay Class: High Sensitivity DNA Assay
 Data Path: C:\...2100 expert\data\2025-05-19\A0682_A0683_A0684_Rapid_PCR.xad

Created: 5/19/2025 12:50:37 PM
 Modified: 5/19/2025 1:38:19 PM

Run Logbook

Description	Number	Source	Category	Sub Category	Time	Time Zone	User	Host
Run ended on port 1 (Number of wells acquired: 12)		Instrument	Run		5/19/2025 1:32:03 PM	(GMT --07:00) Pacific Standard Time	sbsuser	DESKTOP-4UNV VOB
Run started on port 1 (File: C:\Program Files (x86)\Agilent\2100 bioanalyzer\2100 expert\data\2025-05-19\Bioanalyzer1_High Sensitivity DNA Assay_DE34903152_2025-05-19_001.xad)		Instrument	Run		5/19/2025 12:50:47 PM	(GMT --07:00) Pacific Standard Time	sbsuser	DESKTOP-4UNV VOB
Product Number : G2938C		Instrument	Run		5/19/2025 12:50:47 PM	(GMT --07:00) Pacific Standard Time	sbsuser	DESKTOP-4UNV VOB
Name :		Instrument	Run		5/19/2025 12:50:47 PM	(GMT --07:00) Pacific Standard Time	sbsuser	DESKTOP-4UNV VOB
Vendor : Agilent Technologies		Instrument	Run		5/19/2025 12:50:47 PM	(GMT --07:00) Pacific Standard Time	sbsuser	DESKTOP-4UNV VOB
Serial# : DE34903152		Instrument	Run		5/19/2025 12:50:47 PM	(GMT --07:00) Pacific Standard Time	sbsuser	DESKTOP-4UNV VOB
Firmware : C.01.069		Instrument	Run		5/19/2025 12:50:47 PM	(GMT --07:00) Pacific Standard Time	sbsuser	DESKTOP-4UNV VOB
Cartridge : Electrode		Instrument	Run		5/19/2025 12:50:46 PM	(GMT --07:00) Pacific Standard Time	sbsuser	DESKTOP-4UNV VOB