


Assay Class: High Sensitivity DNA Assay
 Data Path: C:\Users\sbsuser\Desktop\2023-03-23\BatchSeq114.xad

Created: 3/23/2023 11:11:01 AM
 Modified: 3/23/2023 12:09:58 PM

Electrophoresis File Run Summary

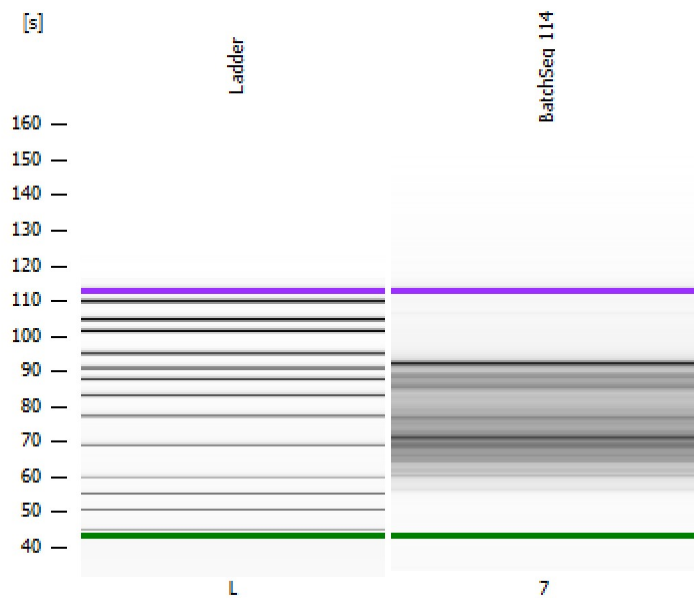

Instrument Information:

Instrument Name: DE34903152 Firmware: C.01.069
 Serial#: DE34903152 Type: G2938C


Assay Information:

Assay Origin Path: C:\Program Files (x86)\Agilent\2100 bioanalyzer\2100 expert\assays\dsDNA\High Sensitivity DNA.xsy
 Assay Class: High Sensitivity DNA Assay
 Version: 1.03
 Assay Comments: Copyright © 2003-2010 Agilent Technologies

Chip Information:

Chip Lot #:
 Reagent Kit Lot #:
 Chip Comments:

BatchSeq 114

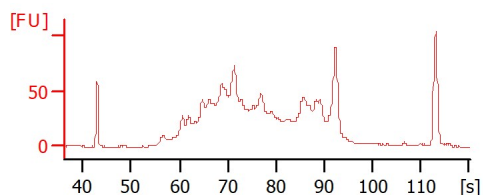

Assay Class: High Sensitivity DNA Assay
Data Path: C:\Users\sbsuser\Desktop\2023-03-23\BatchSeq114.xad

Created: 3/23/2023 11:11:01 AM
Modified: 3/23/2023 12:09:58 PM

Electrophoresis File Run Summary (Chip Summary)

Sample Name	Sample Comment	Rest. Digest	Status	Observation	Result Label	Result Color
BatchSeq 114		<input type="checkbox"/>	✓			
Ladder		<input type="checkbox"/>	✓			

Chip Lot #

Reagent Kit Lot #

Chip Comments :

Assay Class: High Sensitivity DNA Assay
Data Path: C:\Users\sbsuser\Desktop\2023-03-23\BatchSeq114.xad

Created: 3/23/2023 11:11:01 AM
Modified: 3/23/2023 12:09:58 PM

Electrophoresis Assay Details

General Analysis Settings

Number of Available Sample and Ladder Wells (Max.) : 12
Minimum Visible Range [s] : 32
Maximum Visible Range [s] : 138
Start Analysis Time Range [s] : 33
End Analysis Time Range [s] : 137.5
Ladder Concentration [pg/μl] : 1950
Uses Standard Area for Ladder Fragments
Lower Marker Concentration [pg/μl] : 125
Upper Marker Concentration [pg/μl] : 75
Used Upper Marker for Quantitation
Standard Curve Fit is Point to Point
Show Data Aligned to Lower and Upper Marker

Integrator Settings

Integration Start Time [s] : 33.05
Integration End Time [s] : 137
Slope Threshold : 0.8
Height Threshold [FU] : 5
Area Threshold : 0.1
Width Threshold [s] : 0.6
Baseline Plateau [s] : 0.5

Filter Settings

Filter Width [s] : 0.5
Polynomial Order : 4


Ladder

Ladder Peak	Size	Area
1	35	160
2	50	210
3	100	208
4	150	221
5	200	242
6	300	270
7	400	305
8	500	306
9	600	336
10	700	321
11	1000	366
12	2000	413
13	3000	411
14	7000	400
15	10380	214

Assay Class: High Sensitivity DNA Assay
 Data Path: C:\Users\sbsuser\Desktop\2023-03-23\BatchSeq114.xad

Created: 3/23/2023 11:11:01 AM
 Modified: 3/23/2023 12:09:58 PM

Electropherogram Summary

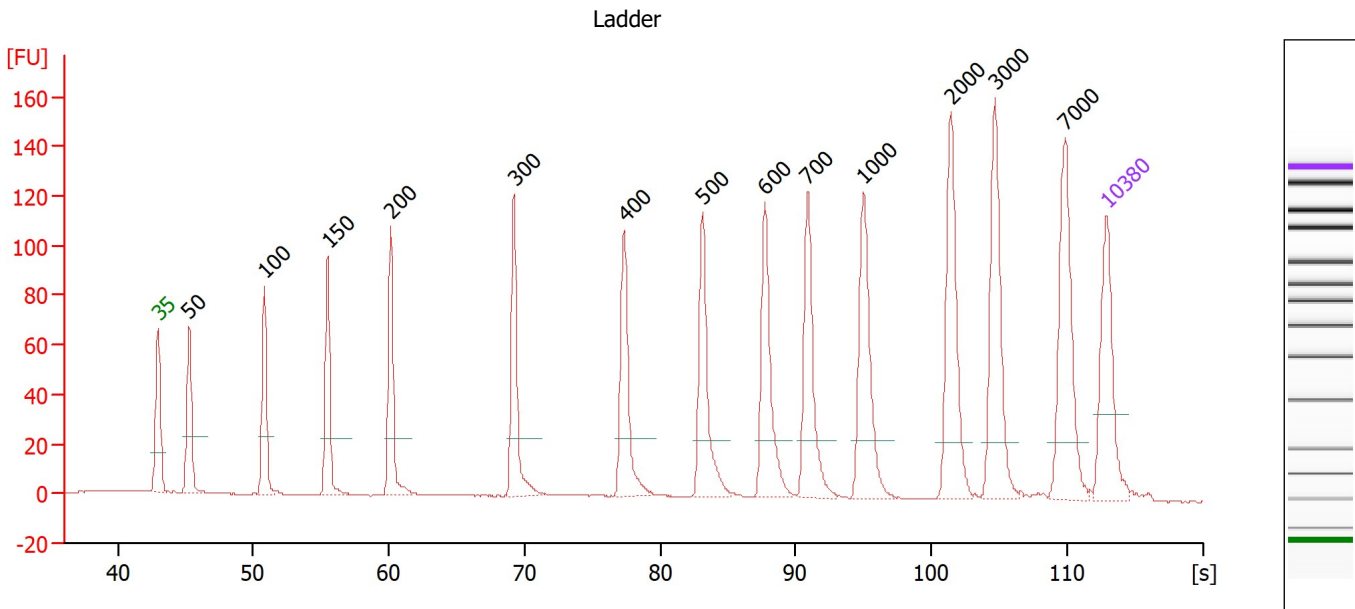

Overall Results for Ladder

Noise: 0.2


Peak table for Ladder

Peak	Size [bp]	Conc. [pg/μl]	Molarity [pmol/l]	Observations
1	35	125.00	5,411.3	Lower Marker
2	50	150.00	4,545.5	Ladder Peak
3	100	150.00	2,272.7	Ladder Peak
4	150	150.00	1,515.2	Ladder Peak
5	200	150.00	1,136.4	Ladder Peak
6	300	150.00	757.6	Ladder Peak
7	400	150.00	568.2	Ladder Peak
8	500	150.00	454.5	Ladder Peak
9	600	150.00	378.8	Ladder Peak
10	700	150.00	324.7	Ladder Peak
11	1,000	150.00	227.3	Ladder Peak
12	2,000	150.00	113.6	Ladder Peak
13	3,000	150.00	75.8	Ladder Peak
14	7,000	150.00	32.5	Ladder Peak
15	10,380	75.00	10.9	Upper Marker

Assay Class: High Sensitivity DNA Assay
 Data Path: C:\Users\sbsuser\Desktop\2023-03-23\BatchSeq114.xad

Created: 3/23/2023 11:11:01 AM
 Modified: 3/23/2023 12:09:58 PM

Electropherogram Summary Continued ...

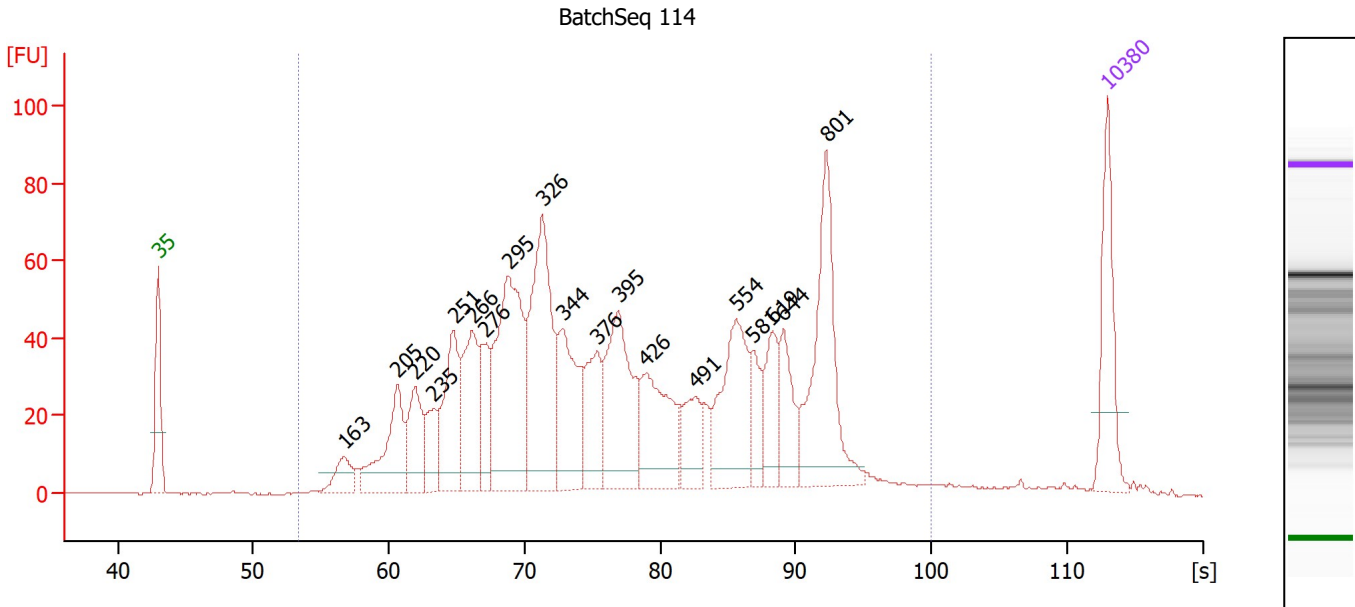

Overall Results for sample 7 : BatchSeq 114

Number of peaks found: 19 Corr. Area 1: 1,638.0
 Noise: 0.1

Peak table for sample 7 : BatchSeq 114

Peak	Size [bp]	Conc. [pg/μl]	Molarity [pmol/l]	Observations
1	35	125.00	5,411.3	Lower Marker
2	163	39.65	368.7	
3	205	114.66	845.6	
4	220	76.81	529.8	
5	235	44.95	290.4	
6	251	120.19	724.9	
7	266	118.71	676.1	
8	276	64.17	352.4	
9	295	253.25	1,300.3	
10	326	235.96	1,097.7	
11	344	126.39	556.1	
12	376	80.88	326.3	
13	395	158.02	606.2	
14	426	119.14	423.9	
15	491	51.93	160.4	
16	554	135.21	369.7	
17	581	34.62	90.3	
18	619	56.35	137.9	
19	644	57.63	135.6	
20	801	177.98	336.7	
21	10,380	75.00	10.9	Upper Marker


Region table for sample 7 : BatchSeq 114

From [s]	To [s]	Corr. Area	% of Total	Average Size [bp]	Size distribution in CV [%]	Conc. [pg/μl]	Molarity [pmol/l]	Color
53.40	100.03	1,638.0	97	448	47.9	2,246.91	9,919.2	Blue

Assay Class: High Sensitivity DNA Assay
Data Path: C:\Users\sbsuser\Desktop\2023-03-23\BatchSeq114.xad

Created: 3/23/2023 11:11:01 AM
Modified: 3/23/2023 12:09:58 PM

Gel Image

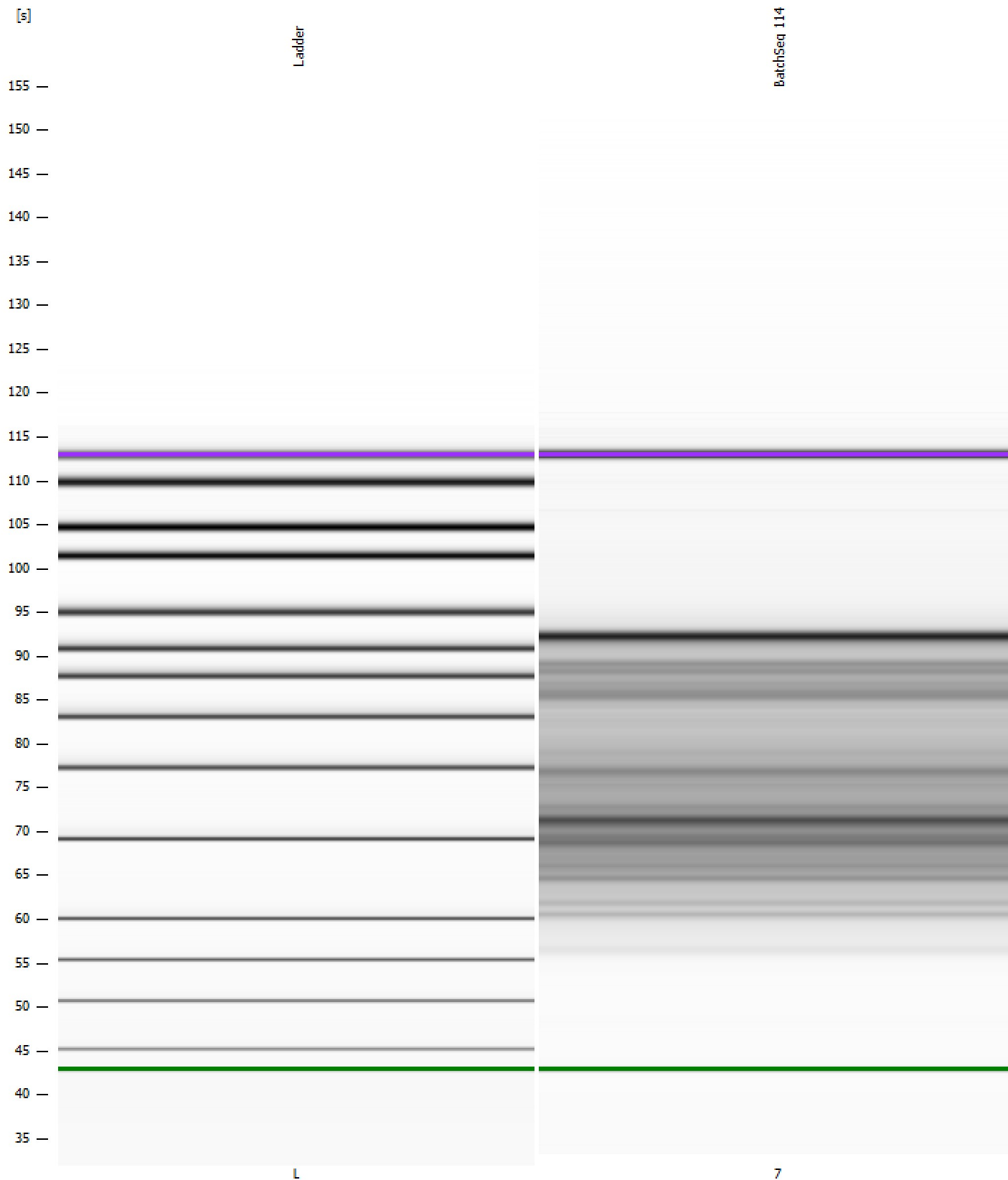

Assay Class: High Sensitivity DNA Assay
 Data Path: C:\Users\sbsuser\Desktop\2023-03-23\BatchSeq114.xad

Created: 3/23/2023 11:11:01 AM
 Modified: 3/23/2023 12:09:58 PM

Run Logbook

Description	Number	Source	Category	Sub Category	Time	Time Zone	User	Host
Run ended on port 1 (Number of wells acquired: 12)		Instrument	Run		3/23/2023 11:52:21 AM	(GMT --07:00) Pacific Standard Time	sbsuser	DESKTOP-4UNV VOB
Run started on port 1 (File: C:\Users\sbsuser\Desktop\2023-03-23\Bioanalyzer1_High Sensitivity DNA Assay_2023-03-23_001.xad)		Instrument	Run		3/23/2023 11:11:06 AM	(GMT --07:00) Pacific Standard Time	sbsuser	DESKTOP-4UNV VOB
Product Number : G2938C		Instrument	Run		3/23/2023 11:11:06 AM	(GMT --07:00) Pacific Standard Time	sbsuser	DESKTOP-4UNV VOB
Name :		Instrument	Run		3/23/2023 11:11:06 AM	(GMT --07:00) Pacific Standard Time	sbsuser	DESKTOP-4UNV VOB
Vendor : Agilent Technologies		Instrument	Run		3/23/2023 11:11:06 AM	(GMT --07:00) Pacific Standard Time	sbsuser	DESKTOP-4UNV VOB
Serial# : DE34903152		Instrument	Run		3/23/2023 11:11:06 AM	(GMT --07:00) Pacific Standard Time	sbsuser	DESKTOP-4UNV VOB
Firmware : C.01.069		Instrument	Run		3/23/2023 11:11:06 AM	(GMT --07:00) Pacific Standard Time	sbsuser	DESKTOP-4UNV VOB
Cartridge : Electrode		Instrument	Run		3/23/2023 11:11:06 AM	(GMT --07:00) Pacific Standard Time	sbsuser	DESKTOP-4UNV VOB