


Assay Class: High Sensitivity DNA Assay
Data Path: C:\...0 expert\data\2023-03-16\BioA_DNA_HS_N778P_N808P_UC390S.xad

Created: 3/16/2023 2:14:26 PM
Modified: 3/16/2023 3:19:06 PM

Electrophoresis File Run Summary

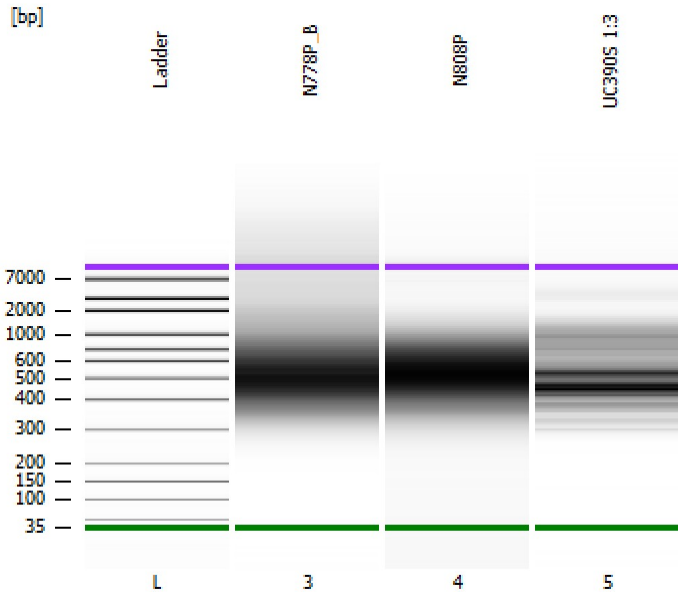

Instrument Information:

Instrument Name: DE34903152 Firmware: C.01.069
Serial#: DE34903152 Type: G2938C

Assay Information:

Assay Origin Path: C:\Program Files (x86)\Agilent\2100 bioanalyzer\2100 expert\assays\dsDNA\High Sensitivity DNA.xsy
Assay Class: High Sensitivity DNA Assay
Version: 1.03
Assay Comments: Copyright © 2003-2010 Agilent Technologies




Chip Information:

Chip Lot #:
Reagent Kit Lot #:
Chip Comments:

N778P_B

N808P

UC390S 1:3

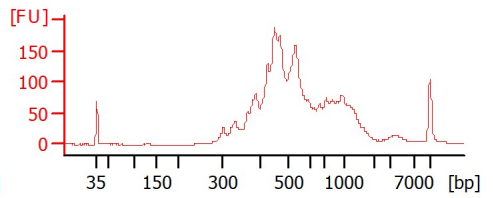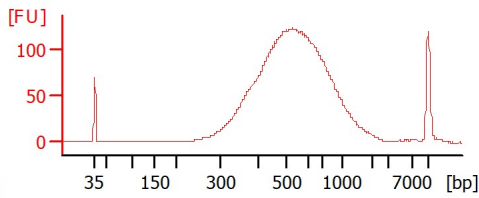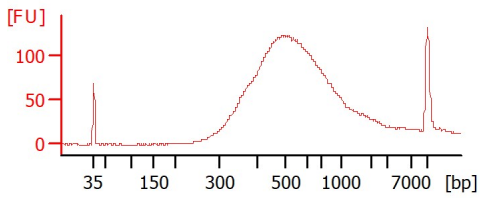

Assay Class: High Sensitivity DNA Assay
 Data Path: C:\...0 expert\data\2023-03-16\BioA_DNA_HS_N778P_N808P_UC390S.xad

Created: 3/16/2023 2:14:26 PM
 Modified: 3/16/2023 3:19:06 PM

Electrophoresis File Run Summary (Chip Summary)

Sample Name	Sample Comment	Rest. Digest	Status	Observation	Result Label	Result Color
N778P_B		<input type="checkbox"/>				
N808P		<input type="checkbox"/>				
UC390S 1:3		<input type="checkbox"/>				
Ladder		<input type="checkbox"/>				

Chip Lot #

Reagent Kit Lot #

Chip Comments :

Assay Class: High Sensitivity DNA Assay
Data Path: C:\...0 expert\data\2023-03-16\BioA_DNA_HS_N778P_N808P_UC390S.xad

Created: 3/16/2023 2:14:26 PM
Modified: 3/16/2023 3:19:06 PM

Electrophoresis Assay Details

General Analysis Settings

Number of Available Sample and Ladder Wells (Max.) : 12
Minimum Visible Range [s] : 32
Maximum Visible Range [s] : 138
Start Analysis Time Range [s] : 33
End Analysis Time Range [s] : 137.5
Ladder Concentration [pg/ μ l] : 1950
Uses Standard Area for Ladder Fragments
Lower Marker Concentration [pg/ μ l] : 125
Upper Marker Concentration [pg/ μ l] : 75
Used Upper Marker for Quantitation
Standard Curve Fit is Point to Point
Show Data Aligned to Lower and Upper Marker

Integrator Settings

Integration Start Time [s] : 33.05
Integration End Time [s] : 137
Slope Threshold : 0.8
Height Threshold [FU] : 5
Area Threshold : 0.1
Width Threshold [s] : 0.6
Baseline Plateau [s] : 0.5

Filter Settings

Filter Width [s] : 0.5
Polynomial Order : 4


Ladder

Ladder Peak	Size	Area
1	35	160
2	50	210
3	100	208
4	150	221
5	200	242
6	300	270
7	400	305
8	500	306
9	600	336
10	700	321
11	1000	366
12	2000	413
13	3000	411
14	7000	400
15	10380	214

Assay Class: High Sensitivity DNA Assay
 Data Path: C:\...0 expert\data\2023-03-16\BioA_DNA_HS_N778P_N808P_UC390S.xad

Created: 3/16/2023 2:14:26 PM
 Modified: 3/16/2023 3:19:06 PM

Electropherogram Summary

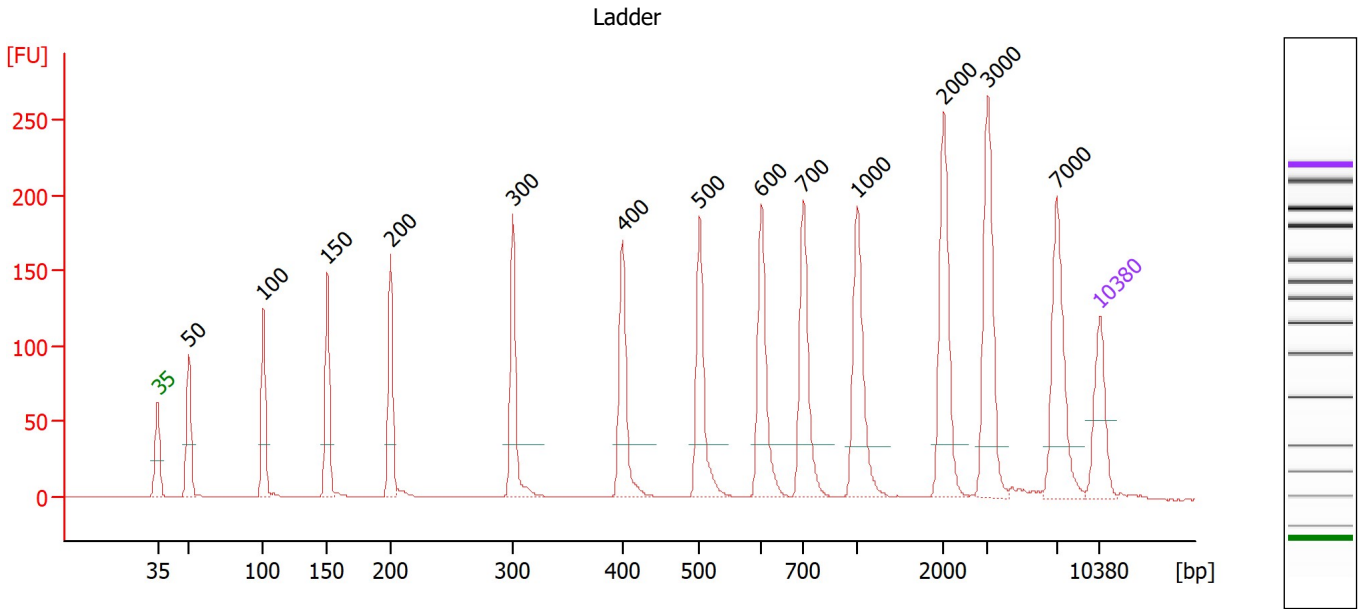

Overall Results for Ladder

Noise: 0.2


Peak table for Ladder

Peak	Size [bp]	Conc. [pg/μl]	Molarity [pmol/l]	Observations
1	35	125.00	5,411.3	Lower Marker
2	50	150.00	4,545.5	Ladder Peak
3	100	150.00	2,272.7	Ladder Peak
4	150	150.00	1,515.2	Ladder Peak
5	200	150.00	1,136.4	Ladder Peak
6	300	150.00	757.6	Ladder Peak
7	400	150.00	568.2	Ladder Peak
8	500	150.00	454.5	Ladder Peak
9	600	150.00	378.8	Ladder Peak
10	700	150.00	324.7	Ladder Peak
11	1,000	150.00	227.3	Ladder Peak
12	2,000	150.00	113.6	Ladder Peak
13	3,000	150.00	75.8	Ladder Peak
14	7,000	150.00	32.5	Ladder Peak
15	10,380	75.00	10.9	Upper Marker

Assay Class: High Sensitivity DNA Assay
 Data Path: C:\...0 expert\data\2023-03-16\BioA_DNA_HS_N778P_N808P_UC390S.xad

Created: 3/16/2023 2:14:26 PM
 Modified: 3/16/2023 3:19:06 PM

Electropherogram Summary Continued ...

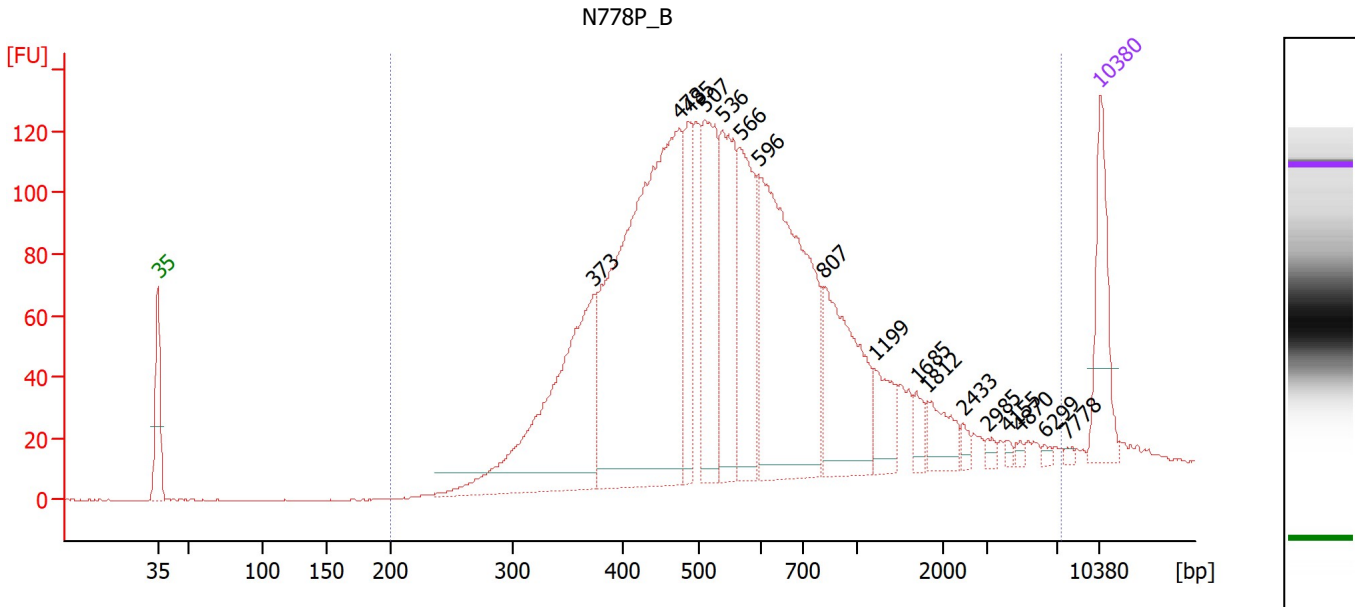

Overall Results for sample 3 : N778P_B

Number of peaks found: 17 Corr. Area 1: 2,958.8
 Noise: 0.2

Peak table for sample 3 : N778P_B

Peak	Size [bp]	Conc. [pg/μl]	Molarity [pmol/l]	Observations
1	35	125.00	5,411.3	Lower Marker
2	373	349.82	1,419.5	
3	473	705.74	2,259.0	
4	485	92.47	288.9	
5	507	179.99	537.8	
6	536	172.23	487.0	
7	566	158.97	425.7	
8	596	389.67	991.0	
9	807	178.28	334.7	
10	1,199	45.07	56.9	
11	1,685	17.04	15.3	
12	1,812	35.04	29.3	
13	2,433	7.59	4.7	
14	2,985	5.65	2.9	
15	4,155	3.74	1.4	
16	4,870	3.83	1.2	
17	6,299	3.27	0.8	
18	7,778	2.65	0.5	
19	10,380	75.00	10.9	Upper Marker


Region table for sample 3 : N778P_B

From [bp]	To [bp]	Corr. Area	% of Total	Average Size [bp]	Size distribution in CV [%]	Conc. [pg/μl]	Molarity [pmol/l]	Color
200	7,289	2,958.8	95	843	100.0	2,788.57	8,139.4	Blue

Assay Class: High Sensitivity DNA Assay
 Data Path: C:\...0 expert\data\2023-03-16\BioA_DNA_HS_N778P_N808P_UC390S.xad

Created: 3/16/2023 2:14:26 PM
 Modified: 3/16/2023 3:19:06 PM

Electropherogram Summary Continued ...

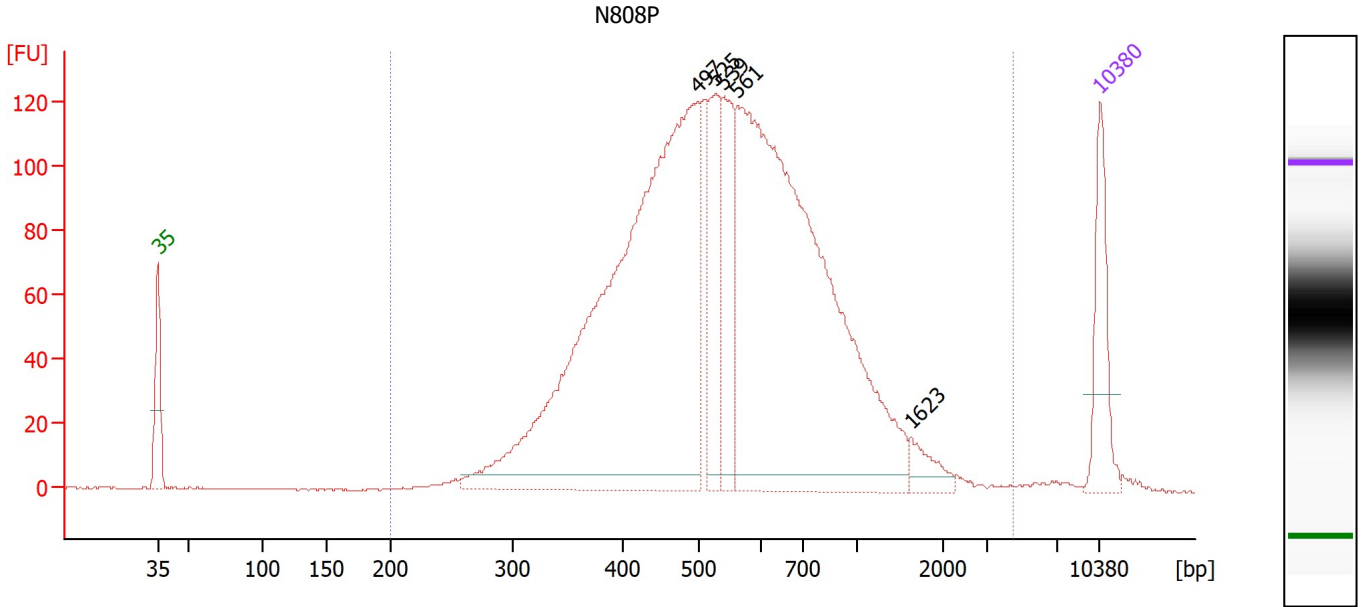

Overall Results for sample 4 : N808P

Number of peaks found: 5 Corr. Area 1: 2,647.6
 Noise: 0.3

Peak table for sample 4 : N808P

Peak	Size [bp]	Conc. [pg/µl]	Molarity [pmol/l]	Observations
1	35	125.00	5,411.3	Lower Marker
2	497	1,099.92	3,350.5	
3	525	143.09	413.0	
4	539	131.80	370.8	
5	561	962.36	2,598.5	
6	1,623	27.78	25.9	
7	10,380	75.00	10.9	Upper Marker


Region table for sample 4 : N808P

From [bp]	To [bp]	Corr. Area	% of Total	Average Size [bp]	Size distribution in CV [%]	Conc. [pg/µl]	Molarity [pmol/l]	Color
200	4,417	2,647.6	98	612	54.6	2,511.30	7,567.9	Blue

Assay Class: High Sensitivity DNA Assay
 Data Path: C:\...0 expert\data\2023-03-16\BioA_DNA_HS_N778P_N808P_UC390S.xad

Created: 3/16/2023 2:14:26 PM
 Modified: 3/16/2023 3:19:06 PM

Electropherogram Summary Continued ...

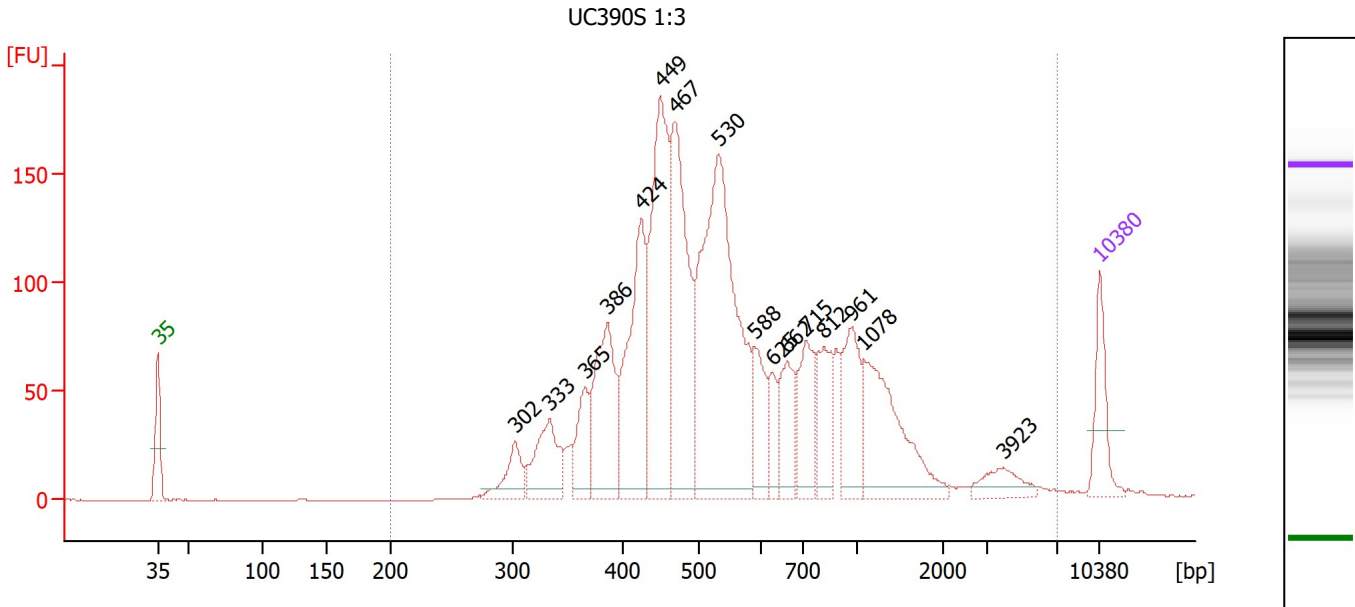

Overall Results for sample 5 : UC390S 1:3

Number of peaks found: 16 Corr. Area 1: 2,728.6
 Noise: 0.2

Peak table for sample 5 : UC390S 1:3

Peak	Size [bp]	Conc. [pg/μl]	Molarity [pmol/l]	Observations
1	35	125.00	5,411.3	Lower Marker
2	302	78.20	392.1	
3	333	138.00	627.9	
4	365	95.89	397.9	
5	386	235.30	923.4	
6	424	315.52	1,127.0	
7	449	430.22	1,453.2	
8	467	392.99	1,274.2	
9	530	693.74	1,982.3	
10	588	95.29	245.5	
11	625	55.68	134.9	
12	662	100.91	231.1	
13	715	122.78	260.0	
14	812	99.48	185.6	
15	961	135.03	212.9	
16	1,078	226.91	319.0	
17	3,923	42.97	16.6	
18	10,380	75.00	10.9	Upper Marker


Region table for sample 5 : UC390S 1:3

From [bp]	To [bp]	Corr. Area	% of Total	Average Size [bp]	Size distribution in CV [%]	Conc. [pg/μl]	Molarity [pmol/l]	Color
200	7,010	2,728.6	99	718	95.2	3,398.30	10,023.6	Blue

Assay Class: High Sensitivity DNA Assay
Data Path: C:\...0 expert\data\2023-03-16\BioA_DNA_HS_N778P_N808P_UC390S.xad

Created: 3/16/2023 2:14:26 PM
Modified: 3/16/2023 3:19:06 PM

Gel Image

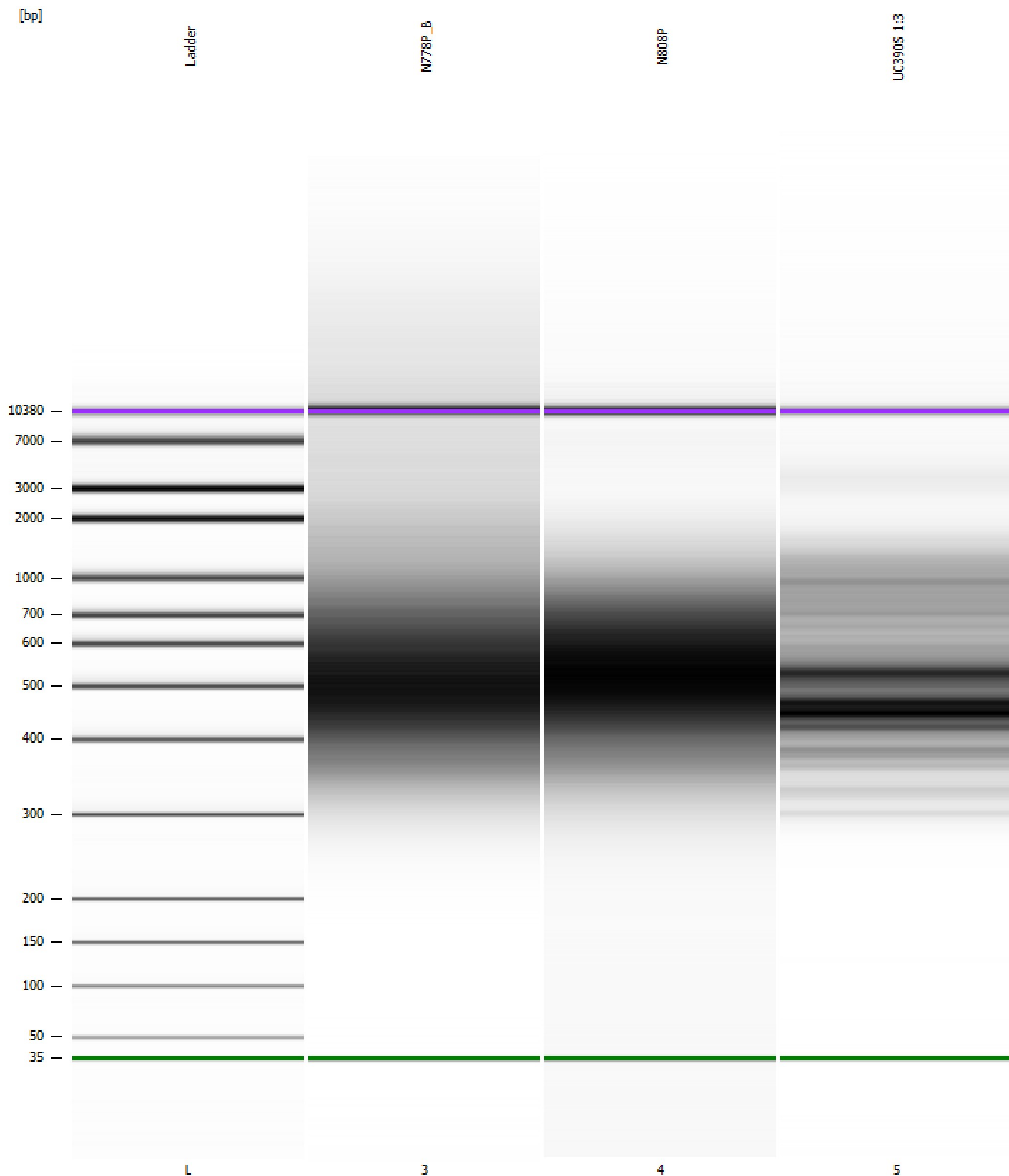

Assay Class: High Sensitivity DNA Assay
 Data Path: C:\...0 expert\data\2023-03-16\BioA_DNA_HS_N778P_N808P_UC390S.xad

Created: 3/16/2023 2:14:26 PM
 Modified: 3/16/2023 3:19:06 PM

Run Logbook

Description	Number	Source	Category	Sub Category	Time	Time Zone	User	Host
Run ended on port 1 (Number of wells acquired: 12)		Instrument	Run		3/16/2023 2:55:44 PM	(GMT --07:00) Pacific Standard Time	sbsuser	DESKTOP-4UNV VOB
Run started on port 1 (File: C:\Users\sbsuser\Desktop\2023-03-16\Bioanalyzer1_High Sensitivity DNA Assay_2023-03-16_001.xad)		Instrument	Run		3/16/2023 2:14:32 PM	(GMT --07:00) Pacific Standard Time	sbsuser	DESKTOP-4UNV VOB
Product Number : G2938C		Instrument	Run		3/16/2023 2:14:32 PM	(GMT --07:00) Pacific Standard Time	sbsuser	DESKTOP-4UNV VOB
Name :		Instrument	Run		3/16/2023 2:14:32 PM	(GMT --07:00) Pacific Standard Time	sbsuser	DESKTOP-4UNV VOB
Vendor : Agilent Technologies		Instrument	Run		3/16/2023 2:14:32 PM	(GMT --07:00) Pacific Standard Time	sbsuser	DESKTOP-4UNV VOB
Serial# : DE34903152		Instrument	Run		3/16/2023 2:14:32 PM	(GMT --07:00) Pacific Standard Time	sbsuser	DESKTOP-4UNV VOB
Firmware : C.01.069		Instrument	Run		3/16/2023 2:14:32 PM	(GMT --07:00) Pacific Standard Time	sbsuser	DESKTOP-4UNV VOB
Cartridge : Electrode		Instrument	Run		3/16/2023 2:14:32 PM	(GMT --07:00) Pacific Standard Time	sbsuser	DESKTOP-4UNV VOB