


Assay Class: High Sensitivity DNA Assay  
 Data Path: C:\...lyzer1\_High Sensitivity DNA Assay\_2023-05-17\_001\_QC1482.xad

Created: 5/17/2023 3:54:12 PM  
 Modified: 5/17/2023 4:57:51 PM

**Electrophoresis File Run Summary**

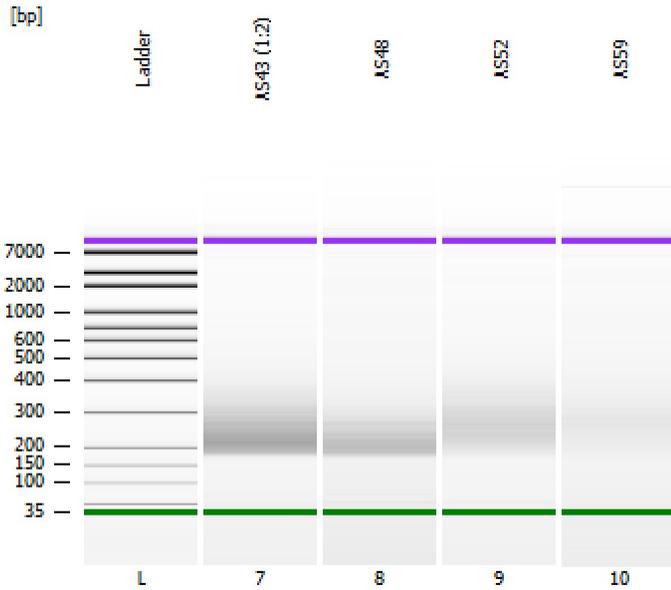

Instrument Information:





Instrument Name: DE34903152      Firmware: C.01.069  
 Serial#: DE34903152      Type: G2938C

Assay Information:

Assay Origin Path: C:\Program Files (x86)\Agilent\2100 bioanalyzer\2100 expert\assays\dsDNA\High Sensitivity DNA.xsy  
 Assay Class: High Sensitivity DNA Assay  
 Version: 1.03  
 Assay Comments: Copyright © 2003-2010 Agilent Technologies

Chip Information:

Chip Lot #:   
 Reagent Kit Lot #:   
 Chip Comments:

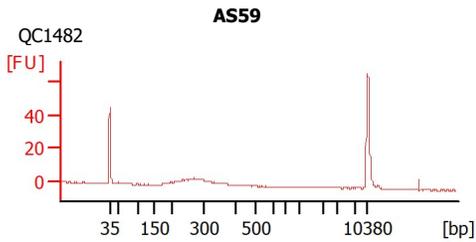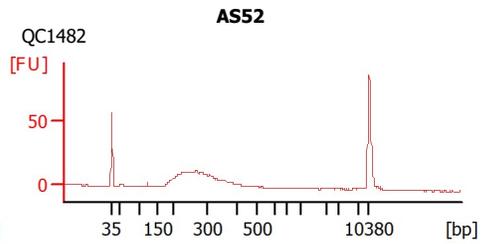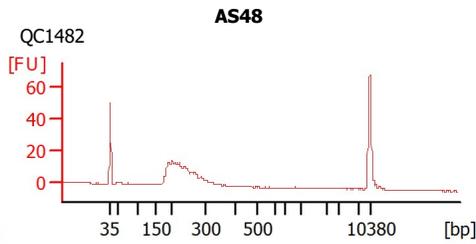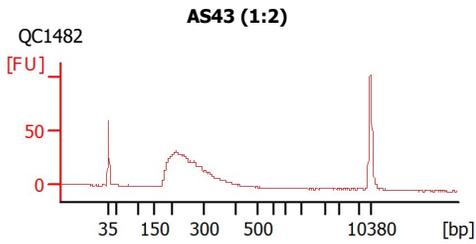

Assay Class: High Sensitivity DNA Assay  
 Data Path: C:\...lyzer1\_High Sensitivity DNA Assay\_2023-05-17\_001\_QC1482.xad

Created: 5/17/2023 3:54:12 PM  
 Modified: 5/17/2023 4:57:51 PM

**Electrophoresis File Run Summary (Chip Summary)**

Sample Name	Sample Comment	Rest. Digest	Status	Observation	Result Label	Result Color
AS43 (1:2)	QC1482	<input type="checkbox"/>		✓		
AS48	QC1482	<input type="checkbox"/>		✓		
AS52	QC1482	<input type="checkbox"/>		✓		
AS59	QC1482	<input type="checkbox"/>		✓		
Ladder		<input type="checkbox"/>		✓		

Chip Lot #

Reagent Kit Lot #

Chip Comments :

Assay Class: High Sensitivity DNA Assay  
 Data Path: C:\...lyzer1\_High Sensitivity DNA Assay\_2023-05-17\_001\_QC1482.xad

Created: 5/17/2023 3:54:12 PM  
 Modified: 5/17/2023 4:57:51 PM

**Electrophoresis Assay Details**

**General Analysis Settings**

Number of Available Sample and Ladder Wells (Max.) : 12  
 Minimum Visible Range [s] : 32  
 Maximum Visible Range [s] : 138  
 Start Analysis Time Range [s] : 33  
 End Analysis Time Range [s] : 137.5  
 Ladder Concentration [pg/μl] : 1950  
 Uses Standard Area for Ladder Fragments  
 Lower Marker Concentration [pg/μl] : 125  
 Upper Marker Concentration [pg/μl] : 75  
 Used Upper Marker for Quantitation  
 Standard Curve Fit is Point to Point  
 Show Data Aligned to Lower and Upper Marker

**Integrator Settings**

Integration Start Time [s] : 33.05  
 Integration End Time [s] : 137  
 Slope Threshold : 0.8  
 Height Threshold [FU] : 5  
 Area Threshold : 0.1  
 Width Threshold [s] : 0.6  
 Baseline Plateau [s] : 0.5

**Filter Settings**

Filter Width [s] : 0.5  
 Polynomial Order : 4


**Ladder**

Ladder Peak	Size	Area
1	35	160
2	50	210
3	100	208
4	150	221
5	200	242
6	300	270
7	400	305
8	500	306
9	600	336
10	700	321
11	1000	366
12	2000	413
13	3000	411
14	7000	400
15	10380	214

Assay Class: High Sensitivity DNA Assay  
 Data Path: C:\...lyzer1\_High Sensitivity DNA Assay\_2023-05-17\_001\_QC1482.xad

Created: 5/17/2023 3:54:12 PM  
 Modified: 5/17/2023 4:57:51 PM

**Electropherogram Summary**

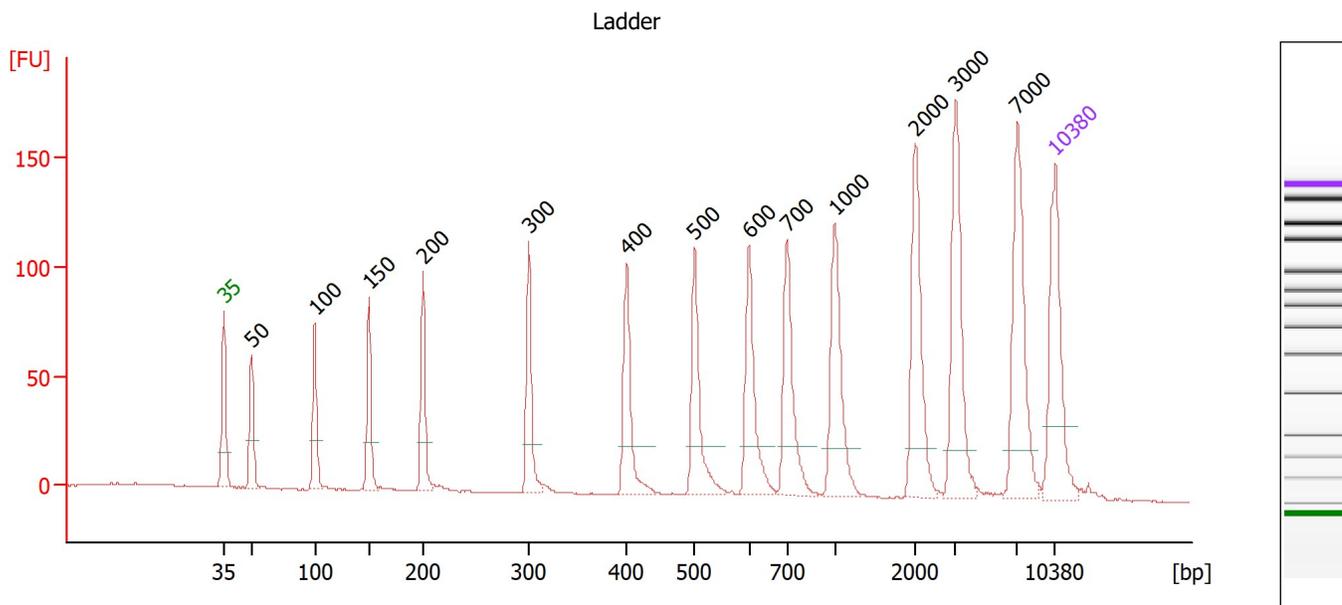

**Overall Results for Ladder**

Noise: 0.3


**Peak table for Ladder**

Peak	Size [bp]	Conc. [pg/μl]	Molarity [pmol/l]	Observations
1	35	125.00	5,411.3	Lower Marker
2	50	150.00	4,545.5	Ladder Peak
3	100	150.00	2,272.7	Ladder Peak
4	150	150.00	1,515.2	Ladder Peak
5	200	150.00	1,136.4	Ladder Peak
6	300	150.00	757.6	Ladder Peak
7	400	150.00	568.2	Ladder Peak
8	500	150.00	454.5	Ladder Peak
9	600	150.00	378.8	Ladder Peak
10	700	150.00	324.7	Ladder Peak
11	1,000	150.00	227.3	Ladder Peak
12	2,000	150.00	113.6	Ladder Peak
13	3,000	150.00	75.8	Ladder Peak
14	7,000	150.00	32.5	Ladder Peak
15	10,380	75.00	10.9	Upper Marker

Assay Class: High Sensitivity DNA Assay  
 Data Path: C:\...lyzer1\_High Sensitivity DNA Assay\_2023-05-17\_001\_QC1482.xad

Created: 5/17/2023 3:54:12 PM  
 Modified: 5/17/2023 4:57:51 PM

**Electropherogram Summary Continued ...**

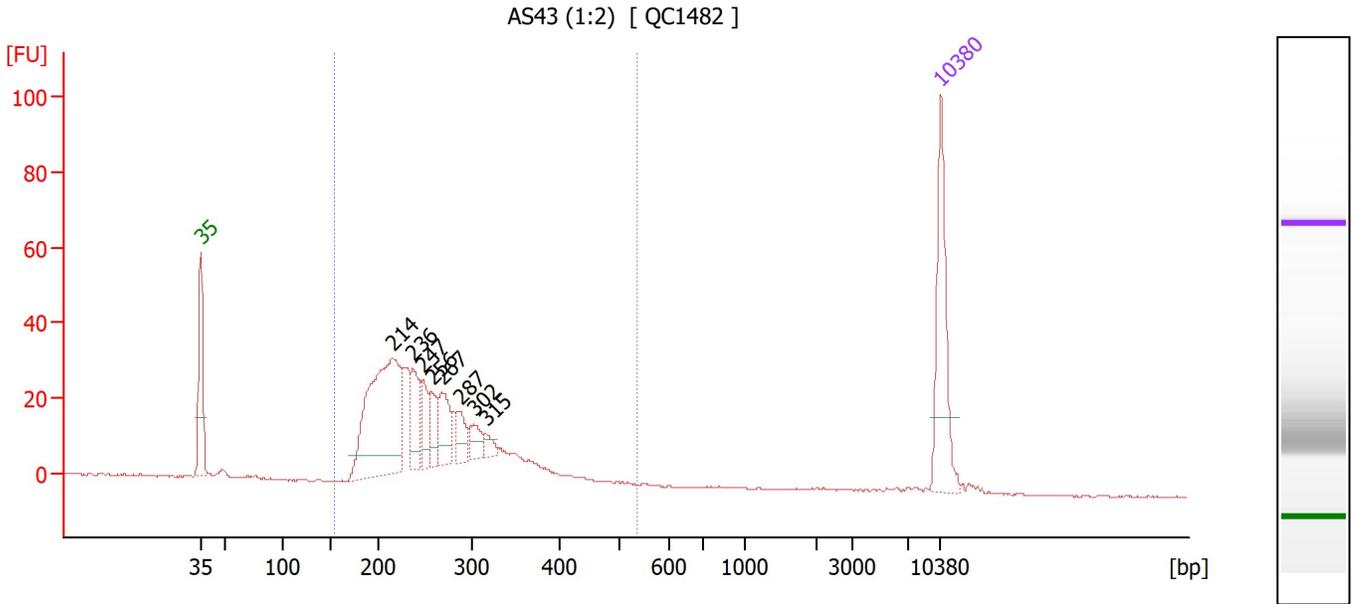

**Overall Results for sample 7 : AS43 (1:2)**

Number of peaks found: 8                      Corr. Area 1: 506.0  
 Noise: 0.2

**Peak table for sample 7 : AS43 (1:2)**

Peak	Size [bp]	Conc. [pg/μl]	Molarity [pmol/l]	Observations
1	35	125.00	5,411.3	Lower Marker
2	214	249.79	1,765.2	
3	236	48.93	313.6	
4	247	33.95	208.6	
5	256	29.42	174.0	
6	267	46.34	263.0	
7	287	24.38	128.8	
8	302	19.65	98.6	
9	315	10.44	50.2	
10	10,380	75.00	10.9	Upper Marker


**Region table for sample 7 : AS43 (1:2)**

From [bp]	To [bp]	Corr. Area	% of Total	Average Size [bp]	Size distribution in CV [%]	Conc. [pg/μl]	Molarity [pmol/l]	Color
155	536	506.0	98	259	22.8	690.64	4,266.2	Blue

Assay Class: High Sensitivity DNA Assay  
 Data Path: C:\...lyzer1\_High Sensitivity DNA Assay\_2023-05-17\_001\_QC1482.xad

Created: 5/17/2023 3:54:12 PM  
 Modified: 5/17/2023 4:57:51 PM

**Electropherogram Summary Continued ...**

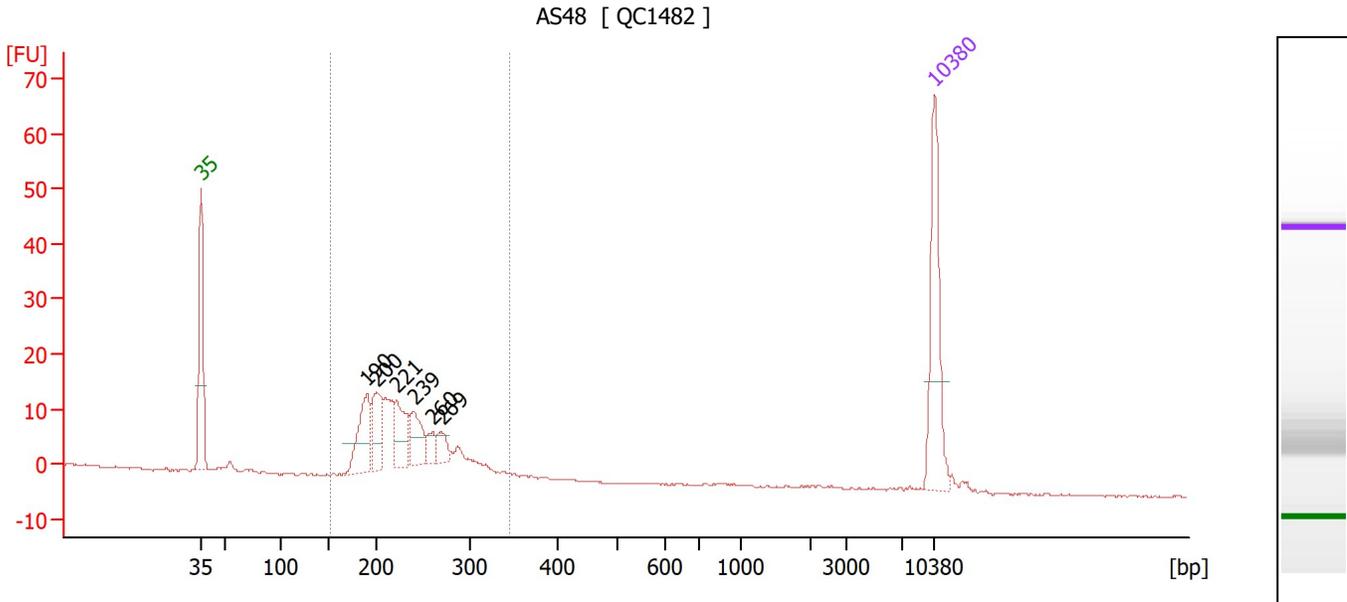

**Overall Results for sample 8 : AS48**

Number of peaks found: 6                      Corr. Area 1: 171.4  
 Noise: 0.2

**Peak table for sample 8 : AS48**

Peak	Size [bp]	Conc. [pg/μl]	Molarity [pmol/l]	Observations
1	35	125.00	5,411.3	Lower Marker
2	190	66.87	534.6	
3	200	52.67	400.0	
4	221	48.12	330.5	
5	239	40.78	258.8	
6	260	14.34	83.6	
7	269	17.94	101.0	
8	10,380	75.00	10.9	Upper Marker


**Region table for sample 8 : AS48**

From [bp]	To [bp]	Corr. Area	% of Total	Average Size [bp]	Size distribution in CV [%]	Conc. [pg/μl]	Molarity [pmol/l]	Color
150	345	171.4	97	231	16.0	362.63	2,439.1	Blue

Assay Class: High Sensitivity DNA Assay  
 Data Path: C:\...lyzer1\_High Sensitivity DNA Assay\_2023-05-17\_001\_QC1482.xad

Created: 5/17/2023 3:54:12 PM  
 Modified: 5/17/2023 4:57:51 PM

**Electropherogram Summary Continued ...**

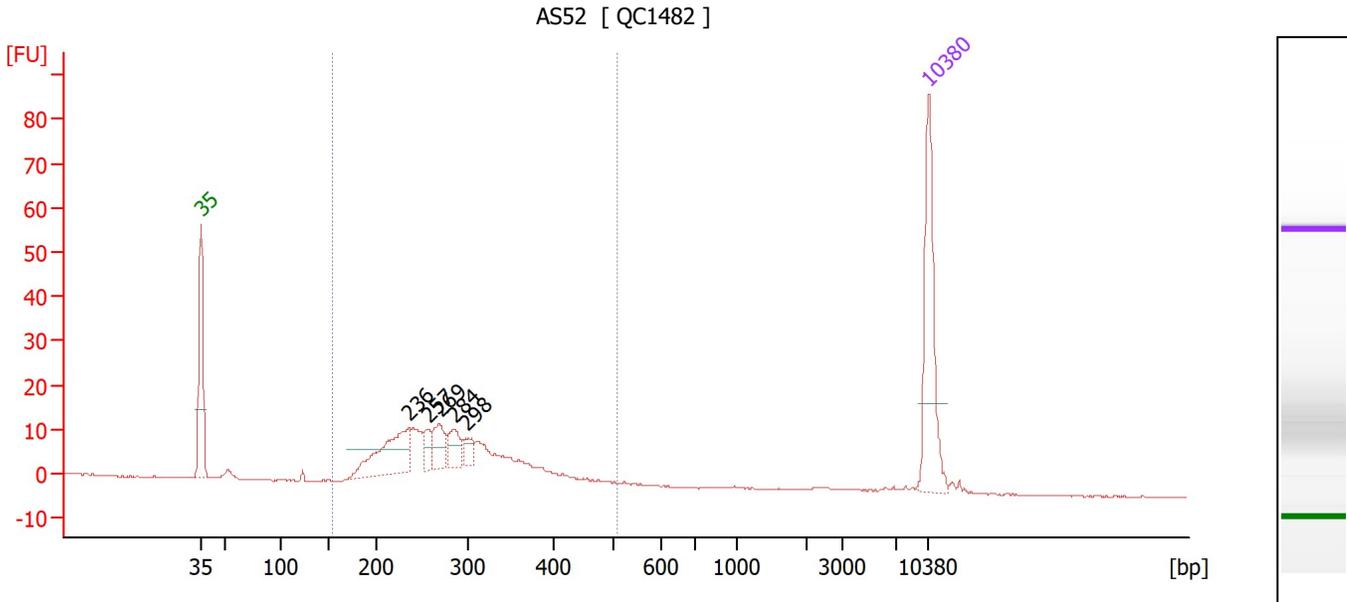

**Overall Results for sample 9 : AS52**

Number of peaks found: 5                      Corr. Area 1: 224.2  
 Noise: 0.2

**Peak table for sample 9 : AS52**

Peak	Size [bp]	Conc. [pg/μl]	Molarity [pmol/l]	Observations
1	35	125.00	5,411.3	Lower Marker
2	236	84.71	544.7	
3	257	16.26	95.7	
4	269	27.91	157.3	
5	284	24.93	133.0	
6	298	12.23	62.1	
7	10,380	75.00	10.9	Upper Marker


**Region table for sample 9 : AS52**

From [bp]	To [bp]	Corr. Area	% of Total	Average Size [bp]	Size distribution in CV [%]	Conc. [pg/μl]	Molarity [pmol/l]	Color
155	510	224.2	96	284	22.1	346.89	1,958.7	Blue

Assay Class: High Sensitivity DNA Assay  
 Data Path: C:\...lyzer1\_High Sensitivity DNA Assay\_2023-05-17\_001\_QC1482.xad

Created: 5/17/2023 3:54:12 PM  
 Modified: 5/17/2023 4:57:51 PM

**Electropherogram Summary Continued ...**

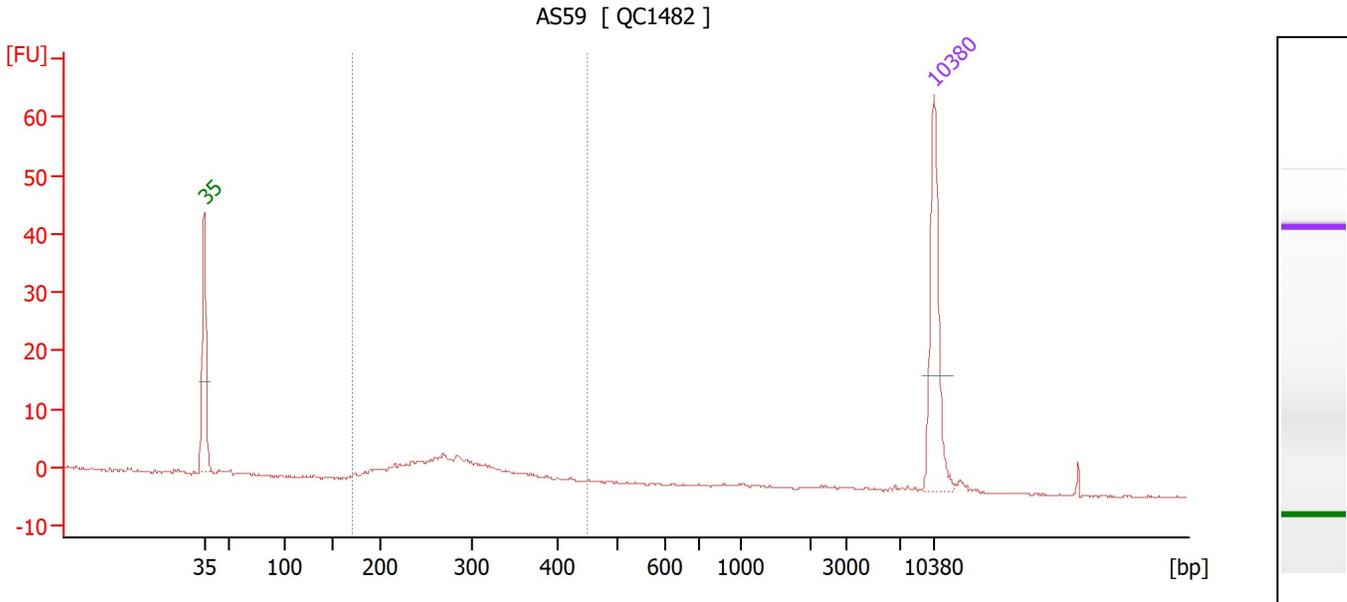

**Overall Results for sample 10 : AS59**

Number of peaks found: 0                      Corr. Area 1: 54.3  
 Noise: 0.2

**Peak table for sample 10 : AS59**

Peak	Size [bp]	Conc. [pg/μl]	Molarity [pmol/l]	Observations
1	35	125.00	5,411.3	Lower Marker
2	10,380	75.00	10.9	Upper Marker


**Region table for sample 10 : AS59**

From [bp]	To [bp]	Corr. Area	% of Total	Average Size [bp]	Size distribution in CV [%]	Conc. [pg/μl]	Molarity [pmol/l]	Color
172	450	54.3	86	282	18.4	116.71	653.4	Blue

Assay Class: High Sensitivity DNA Assay  
Data Path: C:\...lyzer1\_High Sensitivity DNA Assay\_2023-05-17\_001\_QC1482.xad

Created: 5/17/2023 3:54:12 PM  
Modified: 5/17/2023 4:57:51 PM

**Gel Image**

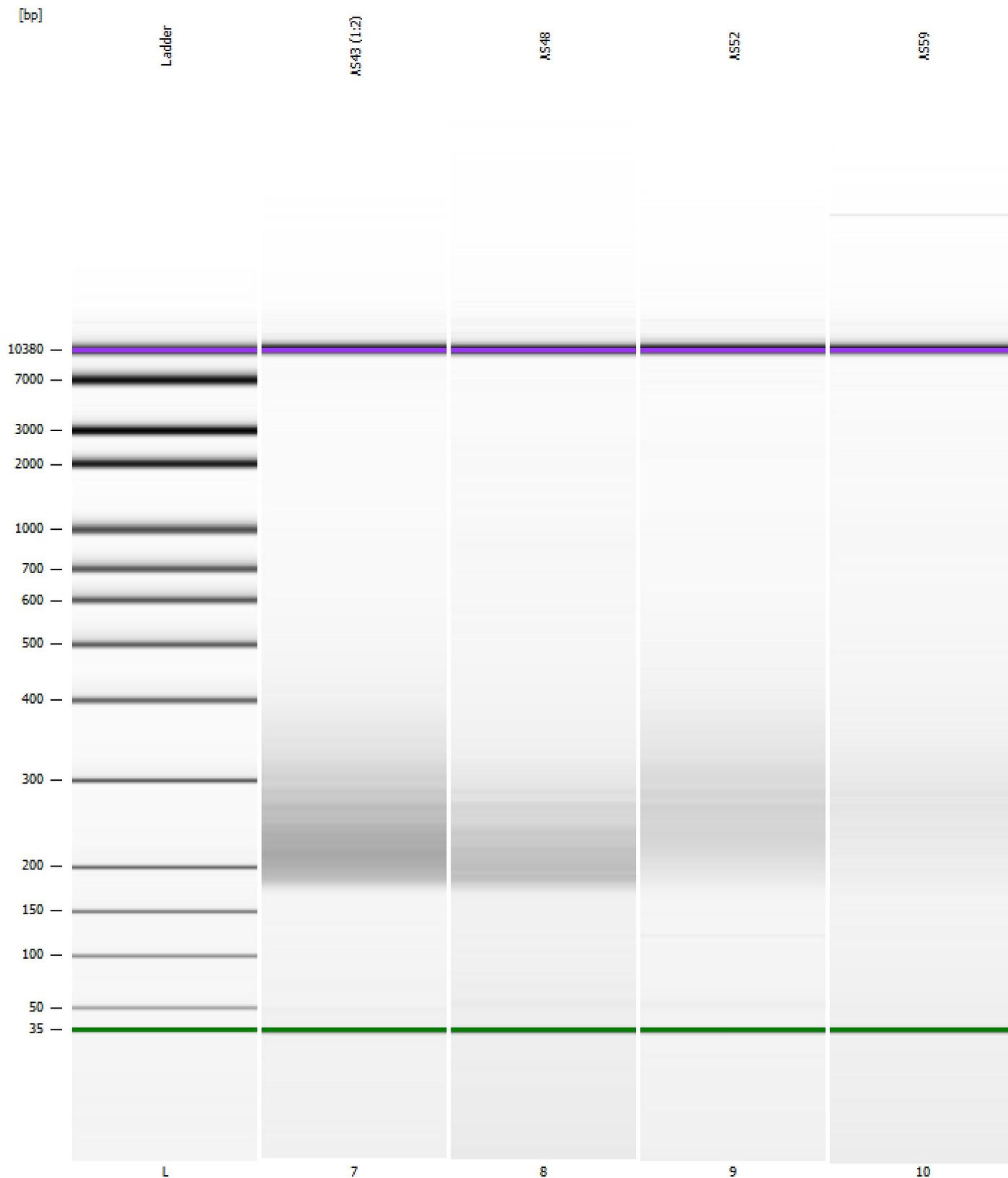

Assay Class: High Sensitivity DNA Assay Created: 5/17/2023 3:54:12 PM  
 Data Path: C:\...lyzer1\_High Sensitivity DNA Assay\_2023-05-17\_001\_QC1482.xad Modified: 5/17/2023 4:57:51 PM

**Run Logbook**

Description	Number	Source	Category	Sub Category	Time	Time Zone	User	Host
Run ended on port 1 (Number of wells acquired: 12)		Instrument	Run		5/17/2023 4:35:35 PM	(GMT --07:00) Pacific Standard Time	sbsuser	DESKTOP-4UNV VOB
Run started on port 1 (File: C:\Users\sbsuser\Desktop\2023-05-17\Bioanalyzer1_High Sensitivity DNA Assay_2023-05-17_001.xad)		Instrument	Run		5/17/2023 3:54:21 PM	(GMT --07:00) Pacific Standard Time	sbsuser	DESKTOP-4UNV VOB
Product Number : G2938C		Instrument	Run		5/17/2023 3:54:21 PM	(GMT --07:00) Pacific Standard Time	sbsuser	DESKTOP-4UNV VOB
Name :		Instrument	Run		5/17/2023 3:54:21 PM	(GMT --07:00) Pacific Standard Time	sbsuser	DESKTOP-4UNV VOB
Vendor : Agilent Technologies		Instrument	Run		5/17/2023 3:54:21 PM	(GMT --07:00) Pacific Standard Time	sbsuser	DESKTOP-4UNV VOB
Serial# : DE34903152		Instrument	Run		5/17/2023 3:54:21 PM	(GMT --07:00) Pacific Standard Time	sbsuser	DESKTOP-4UNV VOB
Firmware : C.01.069		Instrument	Run		5/17/2023 3:54:21 PM	(GMT --07:00) Pacific Standard Time	sbsuser	DESKTOP-4UNV VOB
Cartridge : Electrode		Instrument	Run		5/17/2023 3:54:21 PM	(GMT --07:00) Pacific Standard Time	sbsuser	DESKTOP-4UNV VOB