


Assay Class: High Sensitivity DNA Assay
Data Path: C:\...pert\data\2025-05-28\N1293P_N1297P_N1295P_N1302P_N1294P.xad

Created: 5/28/2025 3:25:46 PM
Modified: 5/28/2025 4:46:05 PM

Electrophoresis File Run Summary

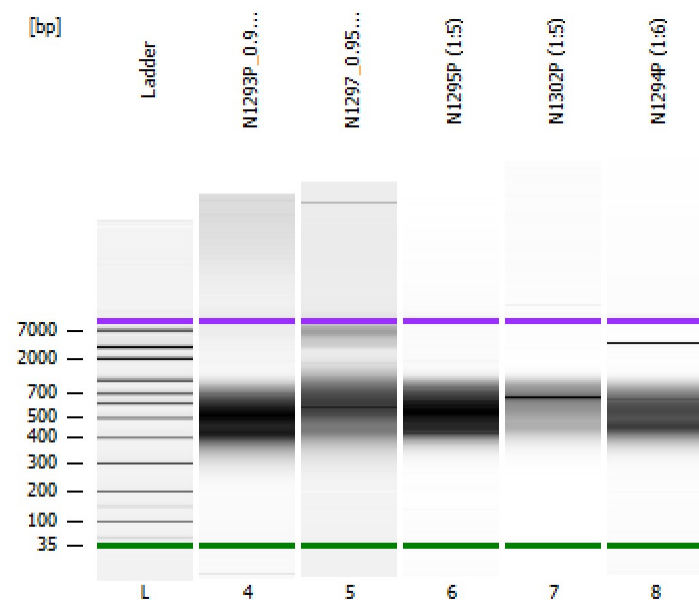

Instrument Information:

Instrument Name: DE34903152

Firmware: C.01.069

Serial#: DE34903152

Type: G2938C

Assay Information:

Assay Origin Path: C:\Program Files (x86)\Agilent\2100 bioanalyzer\2100 expert\assays\dsDNA\High Sensitivity DNA.xsy

Assay Class: High Sensitivity DNA Assay

Version: 1.03

Assay Comments: Copyright © 2003-2010 Agilent Technologies


Chip Information:

Chip Lot #:


Reagent Kit Lot #:

Chip Comments:


N1293P_0.95xBC (1:10)

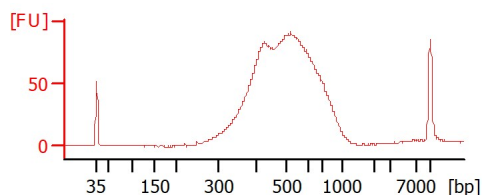


N1297_0.95xBC (1:4)

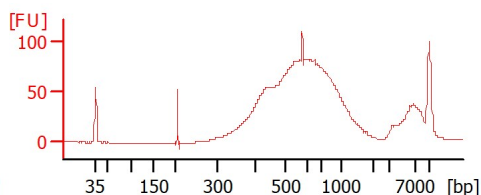


N1295P (1:5)

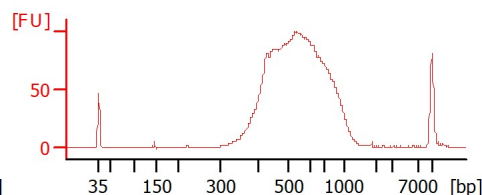

N1302P (1:5)

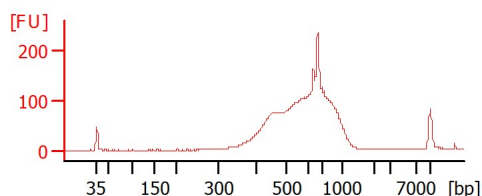

N1294P (1:6)

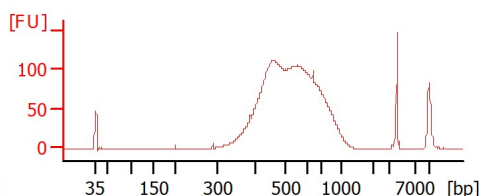

Assay Class: High Sensitivity DNA Assay
Data Path: C:\...pert\data\2025-05-28\N1293P_N1297P_N1295P_N1302P_N1294P.xad

Created: 5/28/2025 3:25:46 PM
Modified: 5/28/2025 4:46:05 PM

Electrophoresis File Run Summary (Chip Summary)

Sample Name	Sample Comment	Rest. Digest	Status	Observation	Result Label	Result Color
N1293P_0.95xBC (1:10)		☐	✓			
N1297_0.95xBC (1:4)		☐	✓			
N1295P (1:5)		☐	✓			
N1302P (1:5)		☐	✓			
N1294P (1:6)		☐	✓			
Ladder		☐	✓			

Chip Lot #**Reagent Kit Lot #****Chip Comments :**

Assay Class: High Sensitivity DNA Assay
Data Path: C:\...pert\data\2025-05-28\N1293P_N1297P_N1295P_N1302P_N1294P.xad

Created: 5/28/2025 3:25:46 PM
Modified: 5/28/2025 4:46:05 PM

Electrophoresis Assay Details

General Analysis Settings

Number of Available Sample and Ladder Wells (Max.) : 12

Minimum Visible Range [s] : 32

Maximum Visible Range [s] : 138

Start Analysis Time Range [s] : 33

End Analysis Time Range [s] : 137.5

Ladder Concentration [pg/μl] : 1950

Uses Standard Area for Ladder Fragments

Lower Marker Concentration [pg/μl] : 125

Upper Marker Concentration [pg/μl] : 75

Used Upper Marker for Quantitation

Standard Curve Fit is Point to Point

Show Data Aligned to Lower and Upper Marker

Integrator Settings

Integration Start Time [s] : 33.05

Integration End Time [s] : 137

Slope Threshold : 0.8

Height Threshold [FU] : 5

Area Threshold : 0.1

Width Threshold [s] : 0.6

Baseline Plateau [s] : 0.5

Filter Settings

Filter Width [s] : 0.5

Polynomial Order : 4


Ladder

Ladder Peak	Size	Area
1	35	160
2	50	210
3	100	208
4	150	221
5	200	242
6	300	270
7	400	305
8	500	306
9	600	336
10	700	321
11	1000	366
12	2000	413
13	3000	411
14	7000	400
15	10380	214

Assay Class: High Sensitivity DNA Assay
 Data Path: C:\...pert\data\2025-05-28\N1293P_N1297P_N1295P_N1302P_N1294P.xad

Created: 5/28/2025 3:25:46 PM
 Modified: 5/28/2025 4:46:05 PM

Electropherogram Summary

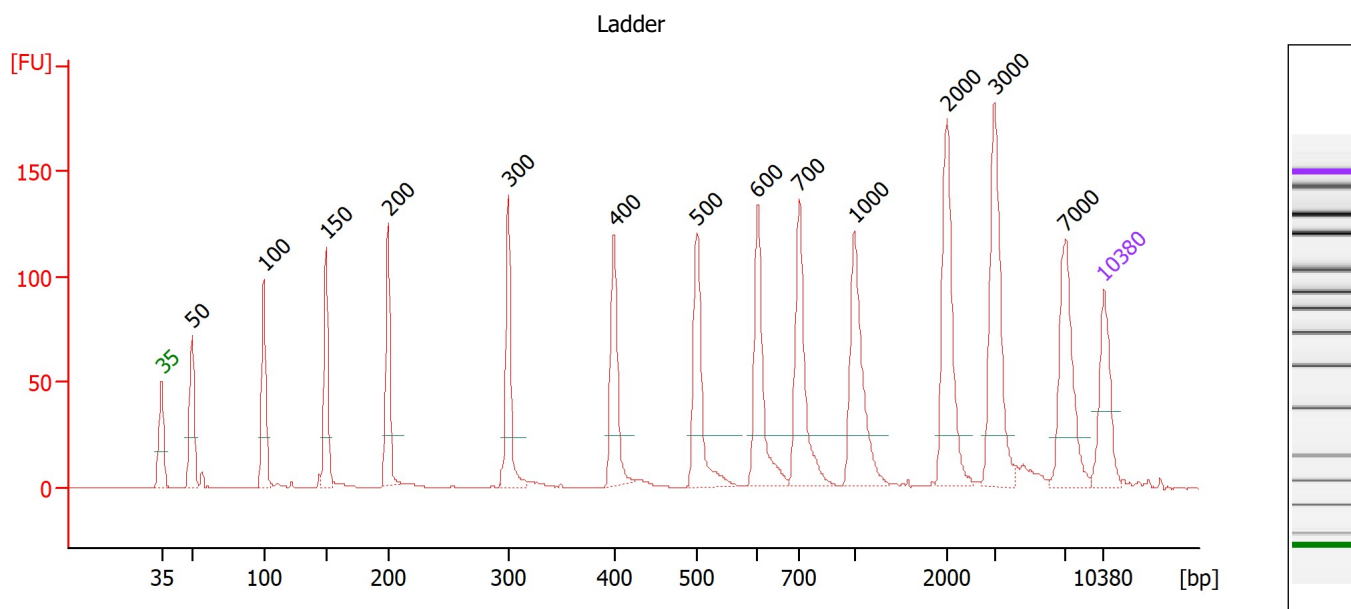

Overall Results for Ladder

Noise: 0.1


Peak table for Ladder

Peak	Size [bp]	Conc. [pg/μl]	Molarity [pmol/l]	Observations
1	35	125.00	5,411.3	Lower Marker
2	50	150.00	4,545.5	Ladder Peak
3	100	150.00	2,272.7	Ladder Peak
4	150	150.00	1,515.2	Ladder Peak
5	200	150.00	1,136.4	Ladder Peak
6	300	150.00	757.6	Ladder Peak
7	400	150.00	568.2	Ladder Peak
8	500	150.00	454.5	Ladder Peak
9	600	150.00	378.8	Ladder Peak
10	700	150.00	324.7	Ladder Peak
11	1,000	150.00	227.3	Ladder Peak
12	2,000	150.00	113.6	Ladder Peak
13	3,000	150.00	75.8	Ladder Peak
14	7,000	150.00	32.5	Ladder Peak
15	10,380	75.00	10.9	Upper Marker

Assay Class: High Sensitivity DNA Assay
 Data Path: C:\...pert\data\2025-05-28\N1293P_N1297P_N1295P_N1302P_N1294P.xad

Created: 5/28/2025 3:25:46 PM
 Modified: 5/28/2025 4:46:05 PM

Electropherogram Summary Continued ...

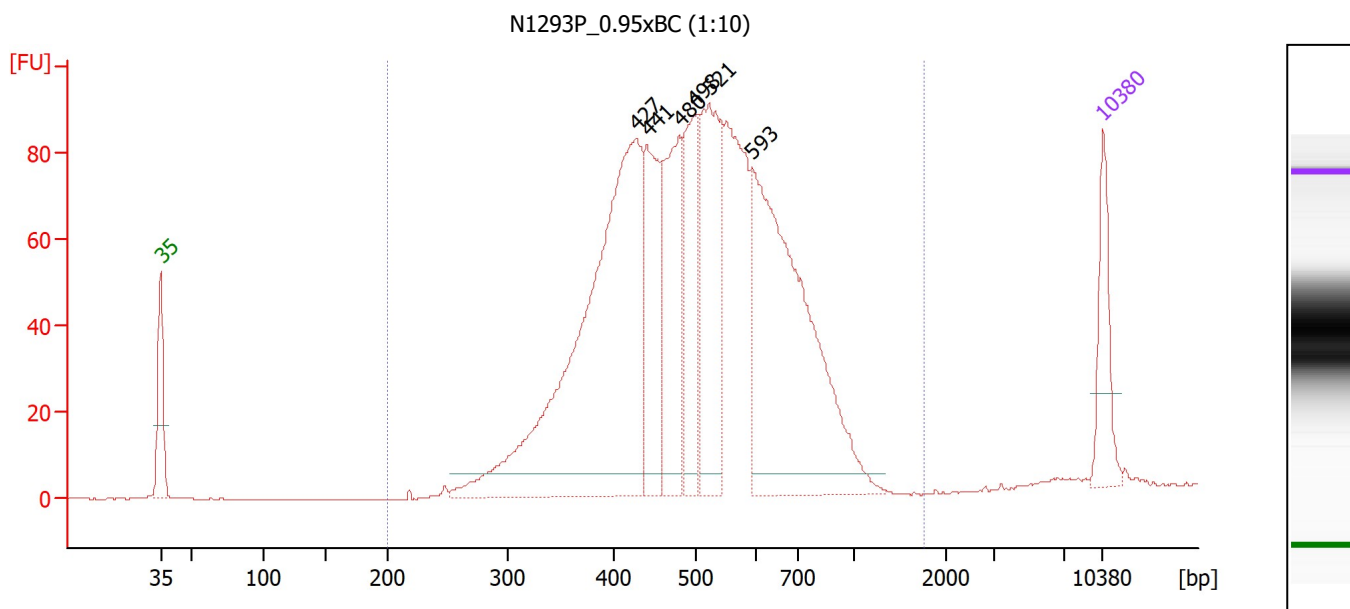

Overall Results for sample 4 : N1293P 0.95xBC (1:10)

Number of peaks found: 6 Corr. Area 1: 1,572.2
 Noise: 0.1

Peak table for sample 4 : N1293P 0.95xBC (1:10)

Peak	Size [bp]	Conc. [pg/μl]	Molarity [pmol/l]	Observations
1	35	125.00	5,411.3	Lower Marker
2	427	860.94	3,053.2	
3	441	201.92	693.7	
4	480	225.30	711.7	
5	498	184.32	560.8	
6	521	274.27	797.1	
7	593	574.50	1,466.8	
8	10,380	75.00	10.9	Upper Marker


Region table for sample 4 : N1293P 0.95xBC (1:10)

From [bp]	To [bp]	Corr. Area	% of Total	Average Size [bp]	Size distribution in CV [%]	Conc. [pg/μl]	Molarity [pmol/l]	Color
200	1,767	1,572.2	99	518	24.8	2,509.67	7,915.2	Blue

Assay Class: High Sensitivity DNA Assay
 Data Path: C:\...pert\data\2025-05-28\N1293P_N1297P_N1295P_N1302P_N1294P.xad

Created: 5/28/2025 3:25:46 PM
 Modified: 5/28/2025 4:46:05 PM

Electropherogram Summary Continued ...

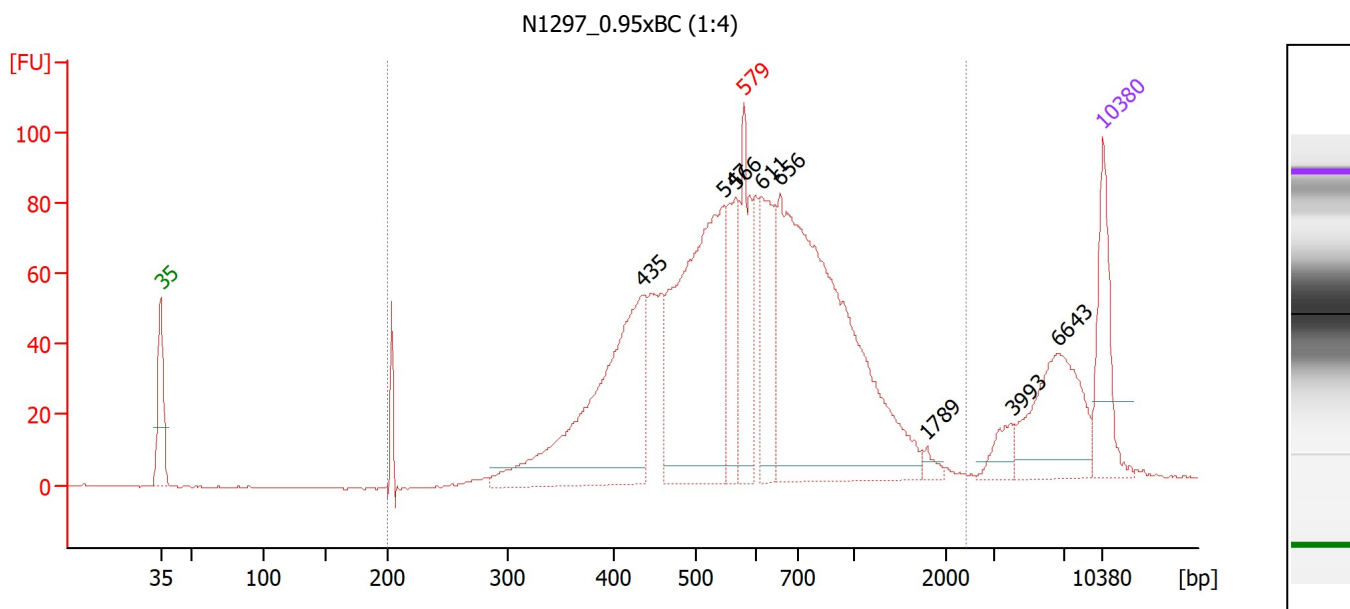

Overall Results for sample 5 : N1297 0.95xBC (1:4)

Number of peaks found: 8 Corr. Area 1: 1,609.5
 Noise: 0.2

Peak table for sample 5 : N1297 0.95xBC (1:4)

Peak	Size [bp]	Conc. [pg/μl]	Molarity [pmol/l]	Observations
1	35	125.00	5,411.3	Lower Marker
2	435	370.47	1,290.6	
3	547	405.94	1,124.8	
4	566	97.47	260.8	
5	579	123.09	321.9	excluded peak
6	611	116.66	289.1	
7	656	614.73	1,419.3	
8	1,789	8.89	7.5	
9	3,993	23.28	8.8	
10	6,643	125.40	28.6	
11	10,380	75.00	10.9	Upper Marker


Region table for sample 5 : N1297 0.95xBC (1:4)

From [bp]	To [bp]	Corr. Area	% of Total	Average Size [bp]	Size distribution in CV [%]	Conc. [pg/μl]	Molarity [pmol/l]	Color
200	2,412	1,609.5	89	650	44.2	1,927.46	5,408.5	Blue

Assay Class: High Sensitivity DNA Assay
 Data Path: C:\...pert\data\2025-05-28\N1293P_N1297P_N1295P_N1302P_N1294P.xad

Created: 5/28/2025 3:25:46 PM
 Modified: 5/28/2025 4:46:05 PM

Electropherogram Summary Continued ...

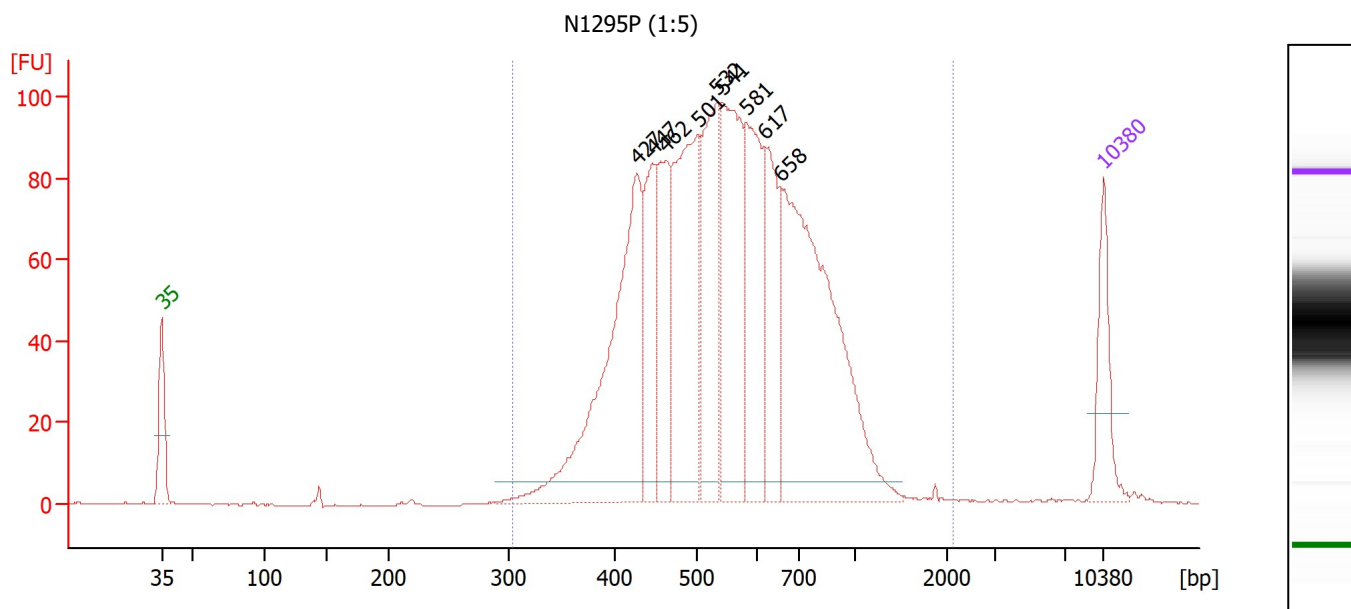

Overall Results for sample 6 : N1295P (1:5)

Number of peaks found: 9 Corr. Area 1: 1,715.0
 Noise: 0.2

Peak table for sample 6 : N1295P (1:5)

Peak	Size [bp]	Conc. [pg/μl]	Molarity [pmol/l]	Observations
1	35	125.00	5,411.3	Lower Marker
2	427	490.01	1,736.7	
3	447	164.16	556.6	
4	462	166.99	547.2	
5	501	346.05	1,045.6	
6	532	247.34	704.3	
7	541	298.58	835.6	
8	581	211.14	550.2	
9	617	165.08	405.2	
10	658	586.34	1,350.9	
11	10,380	75.00	10.9	Upper Marker


Region table for sample 6 : N1295P (1:5)

From [bp]	To [bp]	Corr. Area	% of Total	Average Size [bp]	Size distribution in CV [%]	Conc. [pg/μl]	Molarity [pmol/l]	Color
304	2,144	1,715.0	98	585	32.0	2,755.75	7,858.9	Blue

Assay Class: High Sensitivity DNA Assay
 Data Path: C:\...pert\data\2025-05-28\N1293P_N1297P_N1295P_N1302P_N1294P.xad

Created: 5/28/2025 3:25:46 PM
 Modified: 5/28/2025 4:46:05 PM

Electropherogram Summary Continued ...

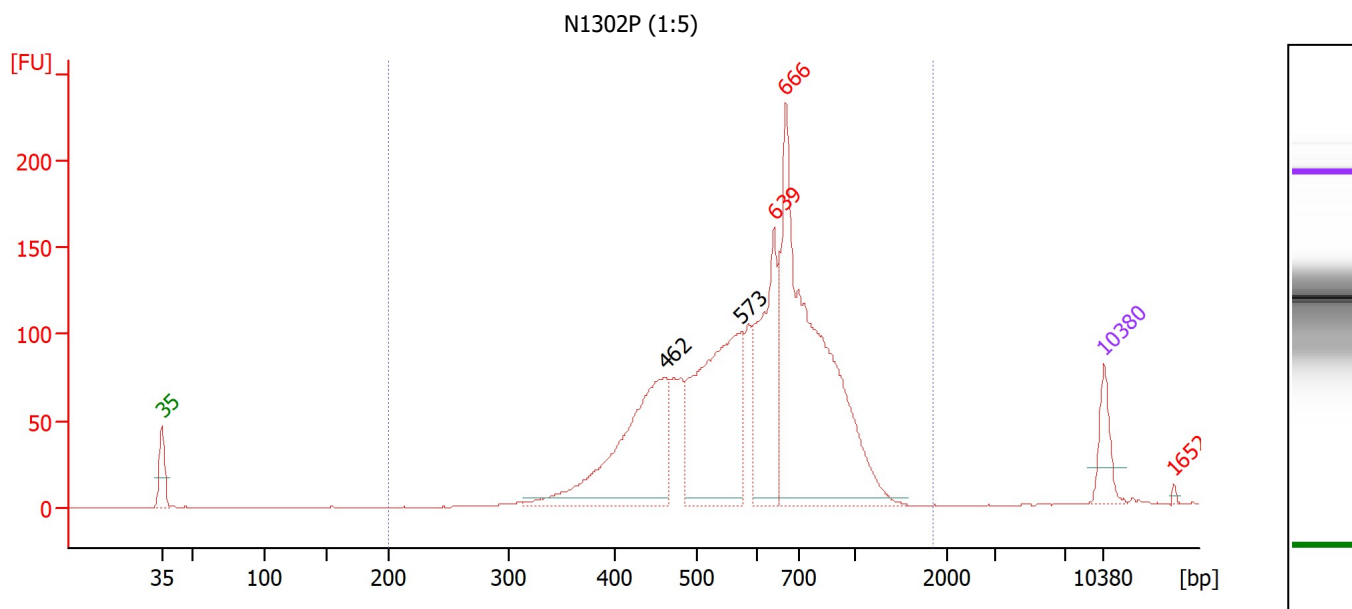

Overall Results for sample 7 : N1302P (1:5)

Number of peaks found: 2 Corr. Area 1: 1,922.9
 Noise: 0.3

Peak table for sample 7 : N1302P (1:5)

Peak	Size [bp]	Conc. [pg/μl]	Molarity [pmol/l]	Observations
1	35	125.00	5,411.3	Lower Marker
2	462	627.50	2,057.4	
3	573	636.07	1,682.6	
4	639	409.36	970.6	excluded peak
5	666	1,101.59	2,504.6	excluded peak
6	10,380	75.00	10.9	Upper Marker
7	16,521	0.00	0.0	excluded peak
8	65,646	0.00	0.0	excluded peak


Region table for sample 7 : N1302P (1:5)

From [bp]	To [bp]	Corr. Area	% of Total	Average Size [bp]	Size distribution in CV [%]	Conc. [pg/μl]	Molarity [pmol/l]	Color
200	1,849	1,922.9	98	625	28.3	3,073.54	8,179.8	Blue

Assay Class: High Sensitivity DNA Assay
 Data Path: C:\...pert\data\2025-05-28\N1293P_N1297P_N1295P_N1302P_N1294P.xad

Created: 5/28/2025 3:25:46 PM
 Modified: 5/28/2025 4:46:05 PM

Electropherogram Summary Continued ...

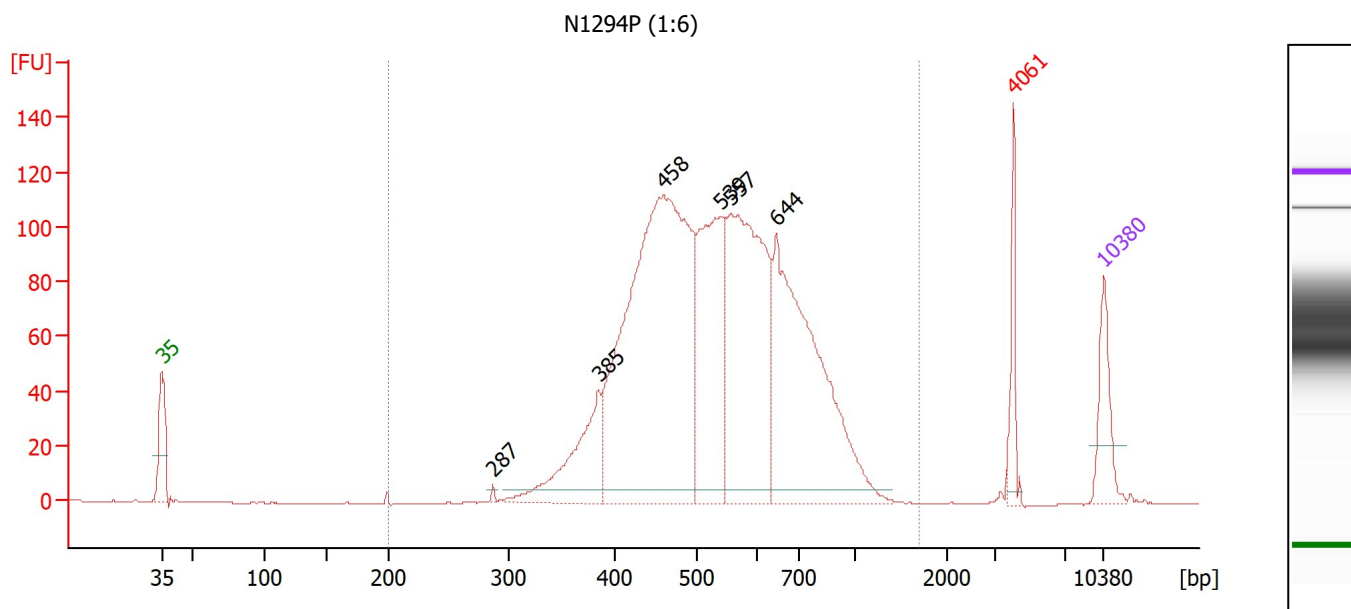

Overall Results for sample 8 : N1294P (1:6)

Number of peaks found: 6 Corr. Area 1: 1,863.5
 Noise: 0.2

Peak table for sample 8 : N1294P (1:6)


Peak	Size [bp]	Conc. [pg/μl]	Molarity [pmol/l]	Observations
1	35	125.00	5,411.3	Lower Marker
2	287	5.04	26.6	
3	385	206.15	810.9	
4	458	1,146.84	3,792.9	
5	539	394.13	1,108.6	
6	557	574.25	1,563.3	
7	644	608.08	1,429.6	
8	4,061	56.22	21.0	excluded peak
9	10,380	75.00	10.9	Upper Marker

Region table for sample 8 : N1294P (1:6)

From [bp]	To [bp]	Corr. Area	% of Total	Average Size [bp]	Size distribution in CV [%]	Conc. [pg/μl]	Molarity [pmol/l]	Color
200	1,702	1,863.5	97	554	26.4	2,925.00	8,631.3	Blue

Assay Class: High Sensitivity DNA Assay
Data Path: C:\...pert\data\2025-05-28\N1293P_N1297P_N1295P_N1302P_N1294P.xad

Created: 5/28/2025 3:25:46 PM
Modified: 5/28/2025 4:46:05 PM

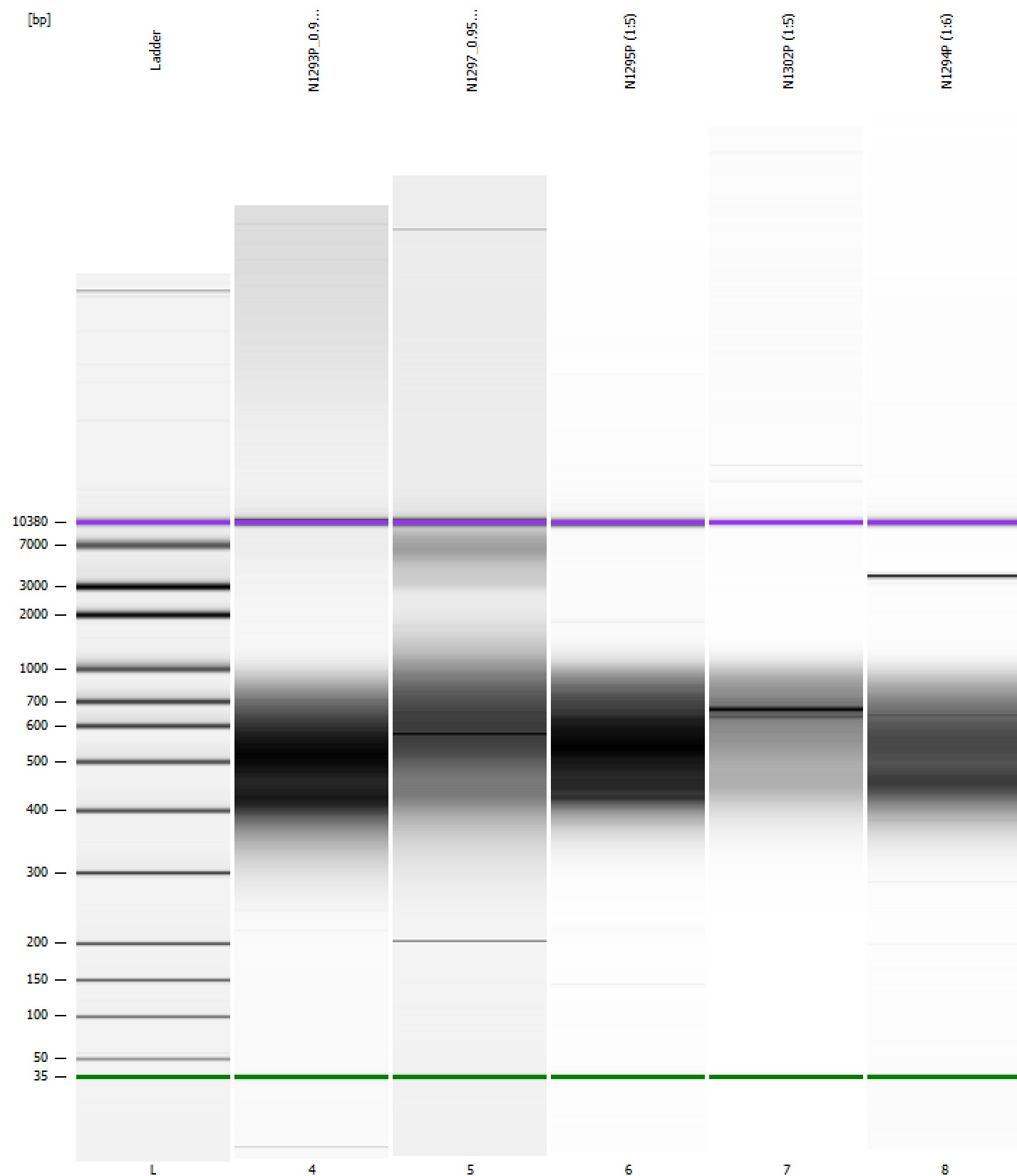Gel Image

Assay Class: High Sensitivity DNA Assay
 Data Path: C:\...pert\data\2025-05-28\N1293P_N1297P_N1295P_N1302P_N1294P.xad

Created: 5/28/2025 3:25:46 PM
 Modified: 5/28/2025 4:46:05 PM

Run Logbook

Description	Number	Source	Category	Sub Category	Time	Time Zone	User	Host
Run ended on port 1 (Number of wells acquired: 12)		Instrument	Run		5/28/2025 4:07:13 PM	(GMT --07:00) Pacific Standard Time	sbsuser	DESKTOP-4UNV VOB
Run started on port 1 (File: C:\Program Files (x86)\Agilent\2100 bioanalyzer\2100 expert\data\2025-05-28\Bioanalyzer1_High Sensitivity DNA Assay_DE34903152_2025-05-28_001.xad)		Instrument	Run		5/28/2025 3:25:53 PM	(GMT --07:00) Pacific Standard Time	sbsuser	DESKTOP-4UNV VOB
Product Number : G2938C		Instrument	Run		5/28/2025 3:25:53 PM	(GMT --07:00) Pacific Standard Time	sbsuser	DESKTOP-4UNV VOB
Name :		Instrument	Run		5/28/2025 3:25:53 PM	(GMT --07:00) Pacific Standard Time	sbsuser	DESKTOP-4UNV VOB
Vendor : Agilent Technologies		Instrument	Run		5/28/2025 3:25:53 PM	(GMT --07:00) Pacific Standard Time	sbsuser	DESKTOP-4UNV VOB
Serial# : DE34903152		Instrument	Run		5/28/2025 3:25:53 PM	(GMT --07:00) Pacific Standard Time	sbsuser	DESKTOP-4UNV VOB
Firmware : C.01.069		Instrument	Run		5/28/2025 3:25:53 PM	(GMT --07:00) Pacific Standard Time	sbsuser	DESKTOP-4UNV VOB
Cartridge : Electrode		Instrument	Run		5/28/2025 3:25:53 PM	(GMT --07:00) Pacific Standard Time	sbsuser	DESKTOP-4UNV VOB