

Assay Class: High Sensitivity DNA Assay
 Data Path: C:\...100 expert\data\2023-09-19\N927P_N949P_N965P_N966P_BioA.xad

Created: 9/19/2023 11:57:54 AM
 Modified: 9/19/2023 12:56:46 PM

Electrophoresis File Run Summary

Instrument Information:


Instrument Name: DE34903152 Firmware: C.01.069
 Serial#: DE34903152 Type: G2938C

Assay Information:

Assay Origin Path: C:\Program Files (x86)\Agilent\2100 bioanalyzer\2100 expert\assays\dsDNA\High Sensitivity DNA.xsy
 Assay Class: High Sensitivity DNA Assay
 Version: 1.03
 Assay Comments: Copyright © 2003-2010 Agilent Technologies

Chip Information:




Chip Lot #:
 Reagent Kit Lot #:
 Chip Comments:

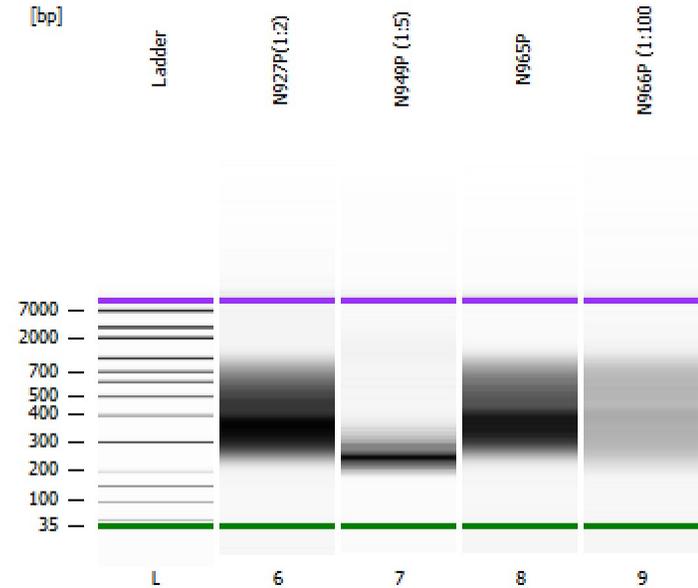


N927P (1:2)

N949P (1:5)

N965P

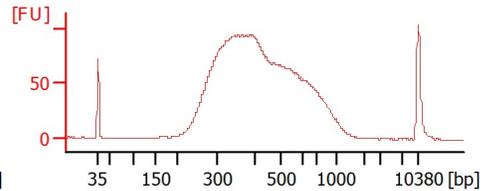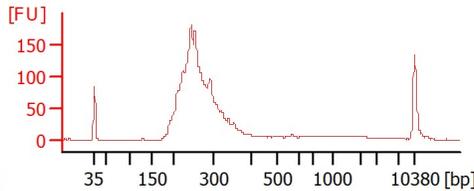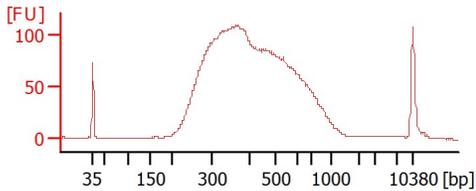

N966P (1:100)

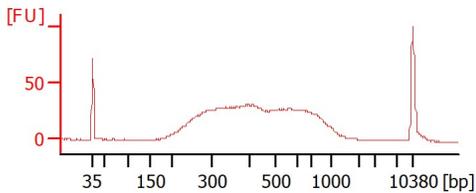

Assay Class: High Sensitivity DNA Assay
 Data Path: C:\...100 expert\data\2023-09-19\N927P_N949P_N965P_N966P_BioA.xad

Created: 9/19/2023 11:57:54 AM
 Modified: 9/19/2023 12:56:46 PM

Electrophoresis File Run Summary (Chip Summary)

Sample Name	Sample Comment	Rest. Digest	Status	Observation	Result Label	Result Color
N927P(1:2)		<input type="checkbox"/>		✓		
N949P (1:5)		<input type="checkbox"/>		✓		
N965P		<input type="checkbox"/>		✓		
N966P (1:100)		<input type="checkbox"/>		✓		
Ladder		<input type="checkbox"/>		✓		

Chip Lot #

Reagent Kit Lot #

Chip Comments :

Assay Class: High Sensitivity DNA Assay
Data Path: C:\...100 expert\data\2023-09-19\N927P_N949P_N965P_N966P_BioA.xad

Created: 9/19/2023 11:57:54 AM
Modified: 9/19/2023 12:56:46 PM

Electrophoresis Assay Details

General Analysis Settings

Number of Available Sample and Ladder Wells (Max.) : 12
Minimum Visible Range [s] : 32
Maximum Visible Range [s] : 138
Start Analysis Time Range [s] : 33
End Analysis Time Range [s] : 137.5
Ladder Concentration [pg/ μ l] : 1950
Uses Standard Area for Ladder Fragments
Lower Marker Concentration [pg/ μ l] : 125
Upper Marker Concentration [pg/ μ l] : 75
Used Upper Marker for Quantitation
Standard Curve Fit is Point to Point
Show Data Aligned to Lower and Upper Marker

Integrator Settings

Integration Start Time [s] : 33.05
Integration End Time [s] : 137
Slope Threshold : 0.8
Height Threshold [FU] : 5
Area Threshold : 0.1
Width Threshold [s] : 0.6
Baseline Plateau [s] : 0.5

Filter Settings

Filter Width [s] : 0.5
Polynomial Order : 4


Ladder

Ladder Peak	Size	Area
1	35	160
2	50	210
3	100	208
4	150	221
5	200	242
6	300	270
7	400	305
8	500	306
9	600	336
10	700	321
11	1000	366
12	2000	413
13	3000	411
14	7000	400
15	10380	214

Assay Class: High Sensitivity DNA Assay
 Data Path: C:\...100 expert\data\2023-09-19\N927P_N949P_N965P_N966P_BioA.xad

Created: 9/19/2023 11:57:54 AM
 Modified: 9/19/2023 12:56:46 PM

Electropherogram Summary

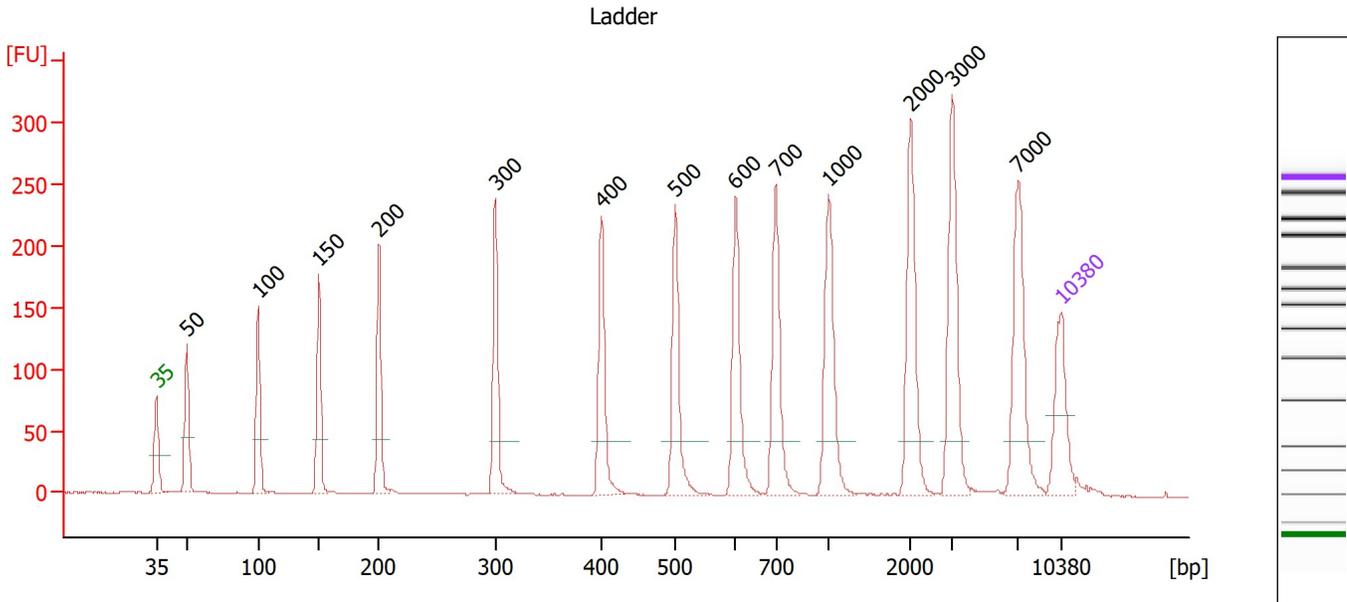

Overall Results for Ladder

Noise: 0.2


Peak table for Ladder

Peak	Size [bp]	Conc. [pg/μl]	Molarity [pmol/l]	Observations
1	35	125.00	5,411.3	Lower Marker
2	50	150.00	4,545.5	Ladder Peak
3	100	150.00	2,272.7	Ladder Peak
4	150	150.00	1,515.2	Ladder Peak
5	200	150.00	1,136.4	Ladder Peak
6	300	150.00	757.6	Ladder Peak
7	400	150.00	568.2	Ladder Peak
8	500	150.00	454.5	Ladder Peak
9	600	150.00	378.8	Ladder Peak
10	700	150.00	324.7	Ladder Peak
11	1,000	150.00	227.3	Ladder Peak
12	2,000	150.00	113.6	Ladder Peak
13	3,000	150.00	75.8	Ladder Peak
14	7,000	150.00	32.5	Ladder Peak
15	10,380	75.00	10.9	Upper Marker

Assay Class: High Sensitivity DNA Assay
 Data Path: C:\...100 expert\data\2023-09-19\N927P_N949P_N965P_N966P_BioA.xad

Created: 9/19/2023 11:57:54 AM
 Modified: 9/19/2023 12:56:46 PM

Electropherogram Summary Continued ...

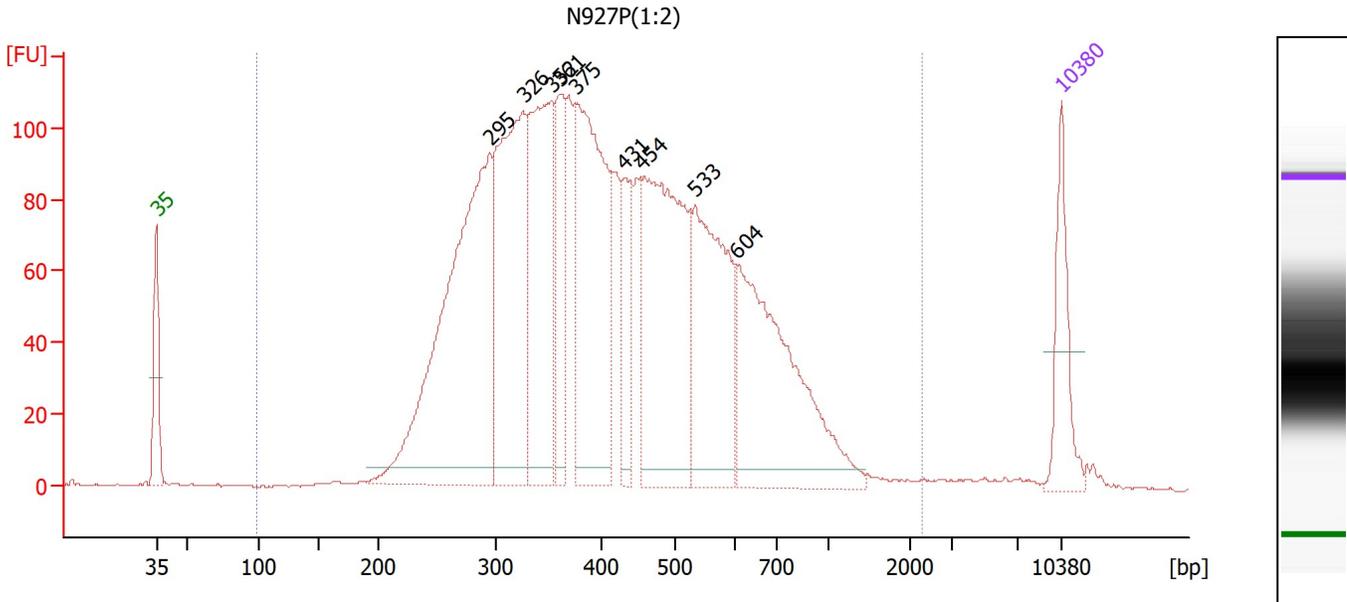

Overall Results for sample 6 : N927P(1:2)

Number of peaks found: 9 Corr. Area 1: 3,115.1
 Noise: 0.2

Peak table for sample 6 : N927P(1:2)

Peak	Size [bp]	Conc. [pg/μl]	Molarity [pmol/l]	Observations
1	35	125.00	5,411.3	Lower Marker
2	295	664.31	3,413.2	
3	326	415.00	1,930.1	
4	352	335.39	1,442.5	
5	361	118.98	499.8	
6	375	399.94	1,614.6	
7	431	83.91	295.3	
8	454	414.90	1,385.9	
9	533	283.83	807.3	
10	604	340.63	854.4	
11	10,380	75.00	10.9	Upper Marker


Region table for sample 6 : N927P(1:2)

From [bp]	To [bp]	Corr. Area	% of Total	Average Size [bp]	Size distribution in CV [%]	Conc. [pg/μl]	Molarity [pmol/l]	Color
99	2,288	3,115.1	98	448	44.2	3,442.67	13,889.7	Blue

Assay Class: High Sensitivity DNA Assay
 Data Path: C:\...100 expert\data\2023-09-19\N927P_N949P_N965P_N966P_BioA.xad

Created: 9/19/2023 11:57:54 AM
 Modified: 9/19/2023 12:56:46 PM

Electropherogram Summary Continued ...

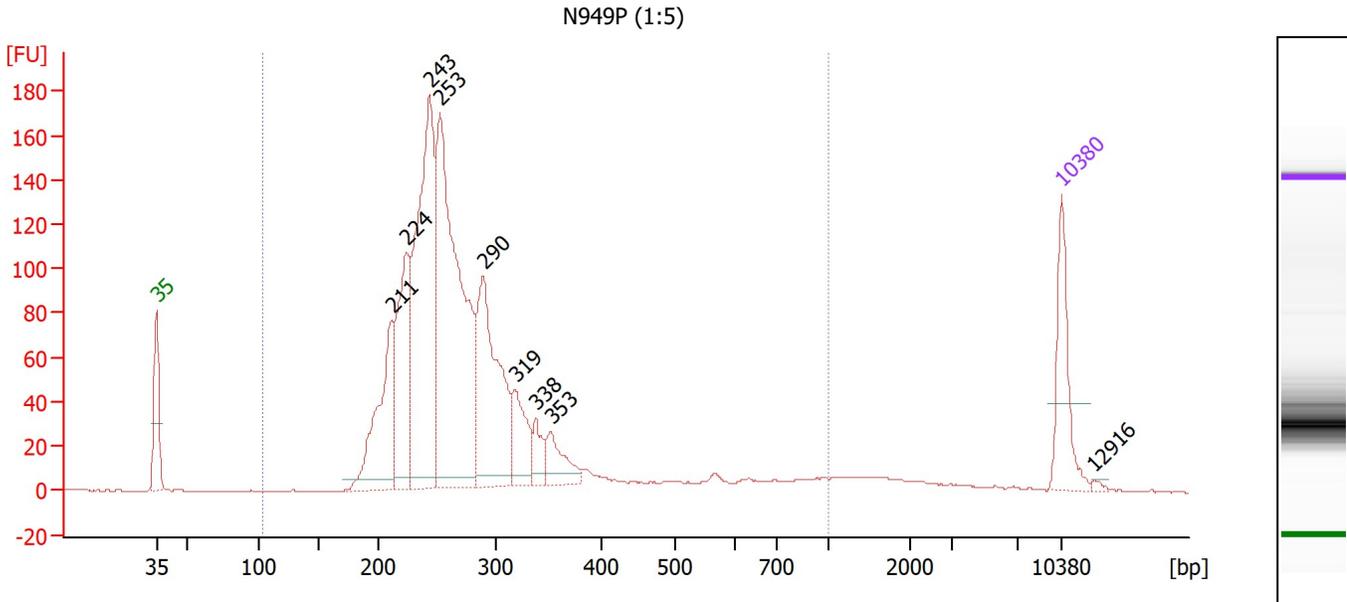

Overall Results for sample 7 : N949P (1:5)

Number of peaks found: 9 Corr. Area 1: 1,871.6
 Noise: 0.2

Peak table for sample 7 : N949P (1:5)

Peak	Size [bp]	Conc. [pg/μl]	Molarity [pmol/l]	Observations
1	35	125.00	5,411.3	Lower Marker
2	211	204.59	1,471.1	
3	224	207.92	1,407.8	
4	243	467.46	2,910.2	
5	253	589.40	3,536.5	
6	290	277.37	1,450.4	
7	319	77.23	367.0	
8	338	34.98	156.7	
9	353	53.08	228.1	
10	10,380	75.00	10.9	Upper Marker
11	12,916	0.00	0.0	


Region table for sample 7 : N949P (1:5)

From [bp]	To [bp]	Corr. Area	% of Total	Average Size [bp]	Size distribution in CV [%]	Conc. [pg/μl]	Molarity [pmol/l]	Color
104	1,000	1,871.6	94	290	38.6	2,071.88	11,898.1	Blue

Assay Class: High Sensitivity DNA Assay
 Data Path: C:\...100 expert\data\2023-09-19\N927P_N949P_N965P_N966P_BioA.xad

Created: 9/19/2023 11:57:54 AM
 Modified: 9/19/2023 12:56:46 PM

Electropherogram Summary Continued ...

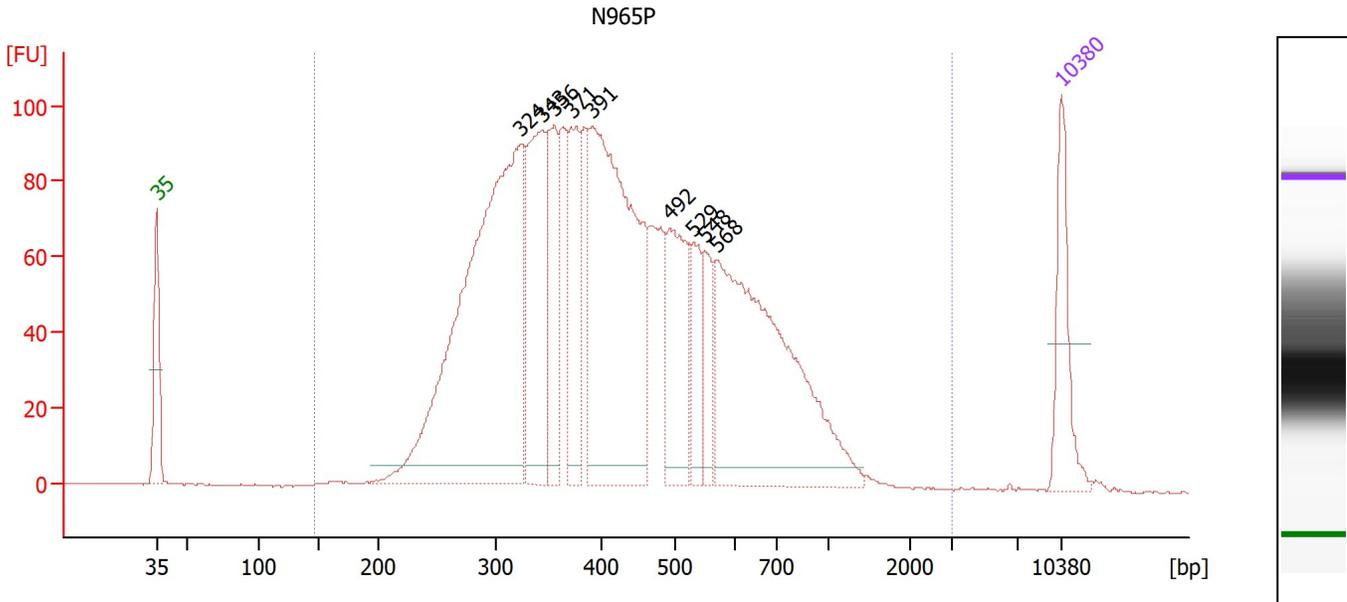

Overall Results for sample 8 : N965P

Number of peaks found: 9 Corr. Area 1: 2,589.4
 Noise: 0.3

Peak table for sample 8 : N965P

Peak	Size [bp]	Conc. [pg/μl]	Molarity [pmol/l]	Observations
1	35	125.00	5,411.3	Lower Marker
2	324	785.07	3,674.0	
3	343	234.79	1,036.2	
4	356	129.37	551.2	
5	371	133.46	545.3	
6	391	539.36	2,089.5	
7	492	159.35	490.6	
8	529	79.68	228.4	
9	548	57.09	157.9	
10	568	465.18	1,240.6	
11	10,380	75.00	10.9	Upper Marker


Region table for sample 8 : N965P

From [bp]	To [bp]	Corr. Area	% of Total	Average Size [bp]	Size distribution in CV [%]	Conc. [pg/μl]	Molarity [pmol/l]	Color
147	2,999	2,589.4	100	454	41.9	2,993.53	11,773.4	Blue

Assay Class: High Sensitivity DNA Assay
 Data Path: C:\...100 expert\data\2023-09-19\N927P_N949P_N965P_N966P_BioA.xad

Created: 9/19/2023 11:57:54 AM
 Modified: 9/19/2023 12:56:46 PM

Electropherogram Summary Continued ...

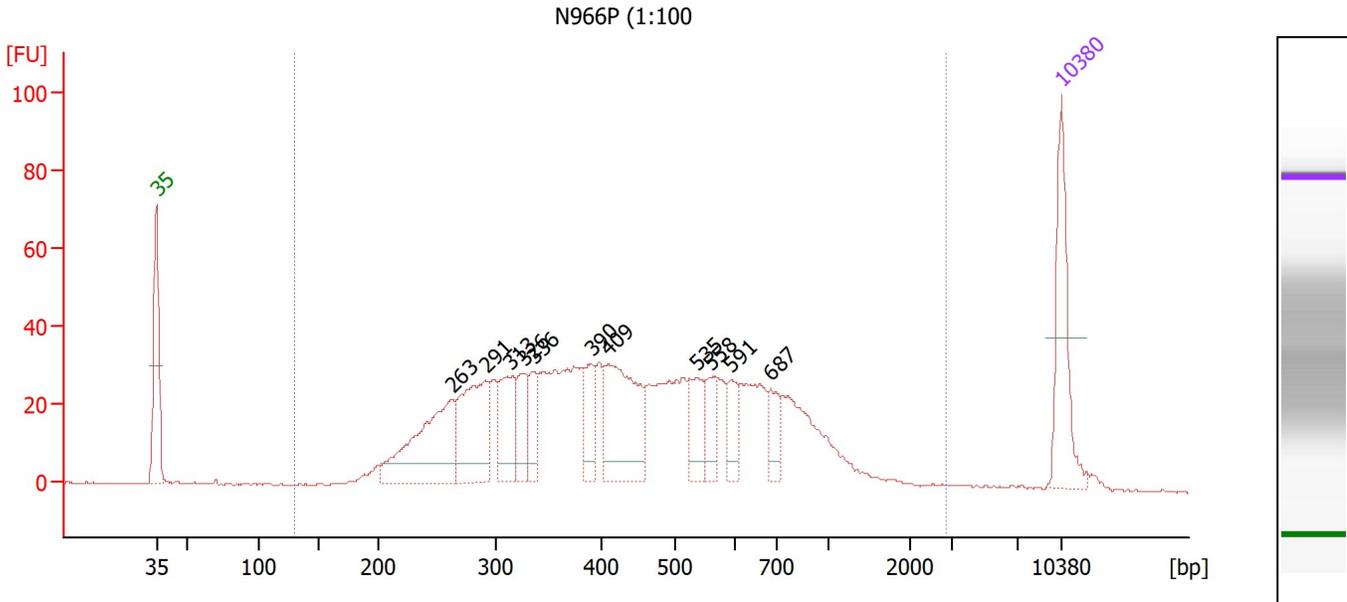

Overall Results for sample 9 : N966P (1:100)

Number of peaks found: 11 Corr. Area 1: 1,106.6
 Noise: 0.2

Peak table for sample 9 : N966P (1:100)

Peak	Size [bp]	Conc. [pg/μl]	Molarity [pmol/l]	Observations
1	35	125.00	5,411.3	Lower Marker
2	263	155.93	899.3	
3	291	122.11	635.6	
4	313	68.59	331.5	
5	326	41.87	194.7	
6	336	38.43	173.3	
7	390	39.10	151.8	
8	409	129.93	481.3	
9	535	46.43	131.4	
10	558	28.90	78.4	
11	591	26.45	67.7	
12	687	26.94	59.4	
13	10,380	75.00	10.9	Upper Marker


Region table for sample 9 : N966P (1:100)

From [bp]	To [bp]	Corr. Area	% of Total	Average Size [bp]	Size distribution in CV [%]	Conc. [pg/μl]	Molarity [pmol/l]	Color
130	2,875	1,106.6	98	494	55.1	1,408.89	5,604.0	Blue

Assay Class: High Sensitivity DNA Assay
Data Path: C:\...100 expert\data\2023-09-19\N927P_N949P_N965P_N966P_BioA.xad

Created: 9/19/2023 11:57:54 AM
Modified: 9/19/2023 12:56:46 PM

Gel Image

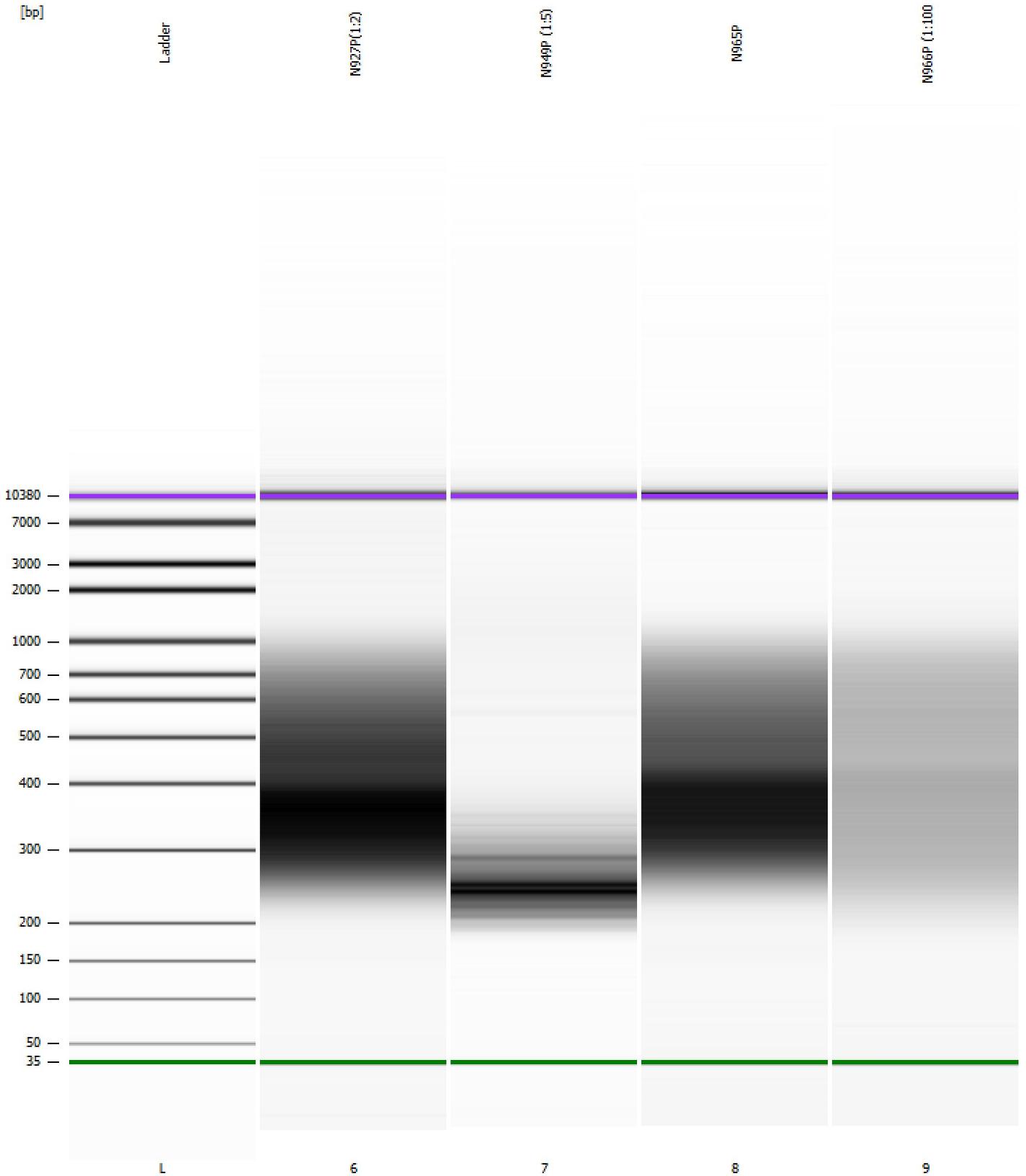


Assay Class: High Sensitivity DNA Assay
Data Path: C:\...100 expert\data\2023-09-19\N927P_N949P_N965P_N966P_BioA.xad

Created: 9/19/2023 11:57:54 AM
Modified: 9/19/2023 12:56:46 PM

Curves

Standard Curve

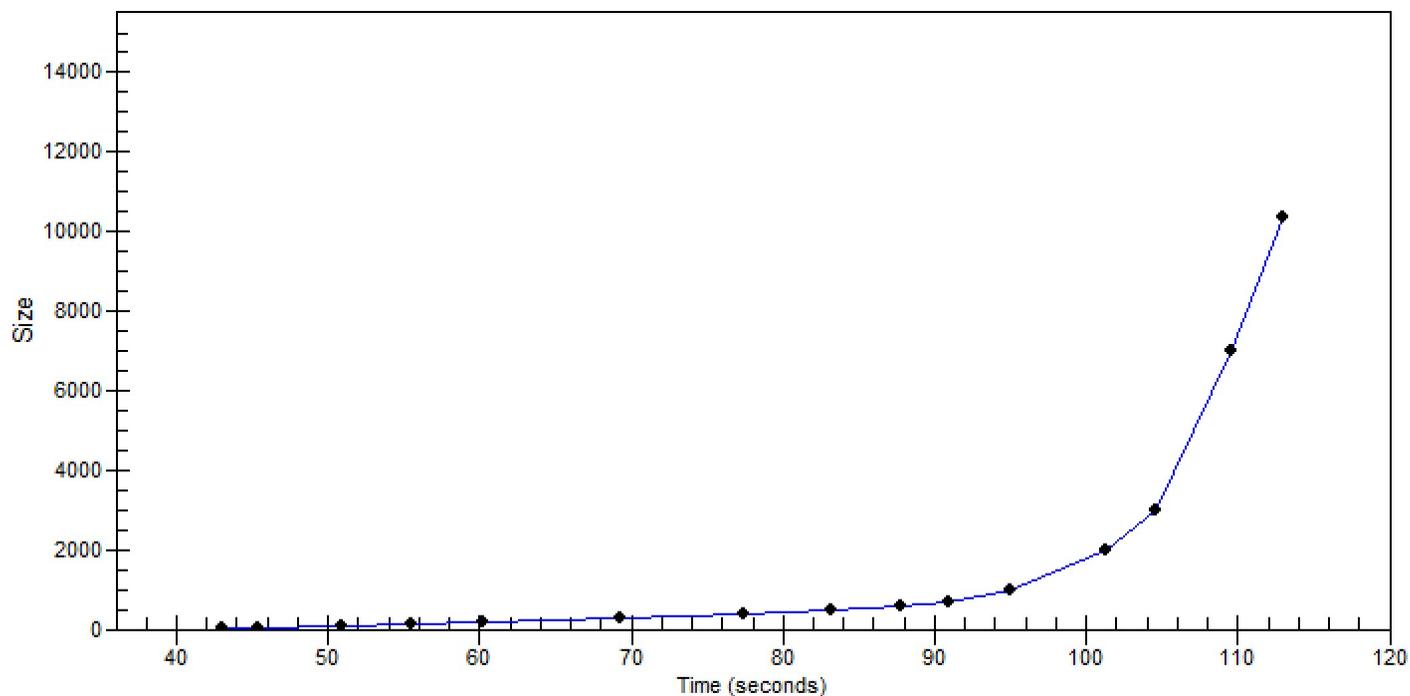

Assay Class: High Sensitivity DNA Assay
 Data Path: C:\...100 expert\data\2023-09-19\N927P_N949P_N965P_N966P_BioA.xad

Created: 9/19/2023 11:57:54 AM
 Modified: 9/19/2023 12:56:46 PM

Run Logbook

Description	Number	Source	Category	Sub Category	Time	Time Zone	User	Host
Run ended on port 1 (Number of wells acquired: 12)		Instrument	Run		9/19/2023 12:39:13 PM	(GMT --07:00) Pacific Standard Time	sbsuser	DESKTOP-4UNV VOB
Run started on port 1 (File: C:\Program Files (x86)\Agilent\2100 bioanalyzer\2100 expert\data\2023-09-19\Bioanalyzer1_High Sensitivity DNA Assay_DE34903152_2023-09-19_001.xad)		Instrument	Run		9/19/2023 11:58:01 AM	(GMT --07:00) Pacific Standard Time	sbsuser	DESKTOP-4UNV VOB
Product Number : G2938C		Instrument	Run		9/19/2023 11:58:01 AM	(GMT --07:00) Pacific Standard Time	sbsuser	DESKTOP-4UNV VOB
Name :		Instrument	Run		9/19/2023 11:58:01 AM	(GMT --07:00) Pacific Standard Time	sbsuser	DESKTOP-4UNV VOB
Vendor : Agilent Technologies		Instrument	Run		9/19/2023 11:58:01 AM	(GMT --07:00) Pacific Standard Time	sbsuser	DESKTOP-4UNV VOB
Serial# : DE34903152		Instrument	Run		9/19/2023 11:58:01 AM	(GMT --07:00) Pacific Standard Time	sbsuser	DESKTOP-4UNV VOB
Firmware : C.01.069		Instrument	Run		9/19/2023 11:58:01 AM	(GMT --07:00) Pacific Standard Time	sbsuser	DESKTOP-4UNV VOB
Cartridge : Electrode		Instrument	Run		9/19/2023 11:58:01 AM	(GMT --07:00) Pacific Standard Time	sbsuser	DESKTOP-4UNV VOB