


Assay Class: High Sensitivity DNA Assay
Data Path: C:\...zer\2100 expert\data\2024-01-03\New Projects_01-13-2024.xad

Created: 1/3/2024 4:39:30 PM
Modified: 1/3/2024 5:29:16 PM

Electrophoresis File Run Summary

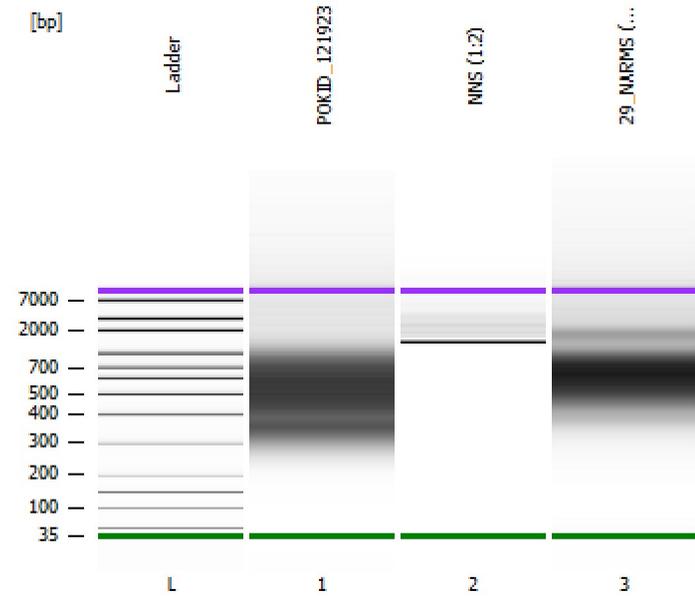

Instrument Information:

Instrument Name: DE34903152 Firmware: C.01.069
Serial#: DE34903152 Type: G2938C

Assay Information:

Assay Origin Path: C:\Program Files (x86)\Agilent\2100 bioanalyzer\2100 expert\assays\dsDNA\High Sensitivity DNA.xsy
Assay Class: High Sensitivity DNA Assay
Version: 1.03
Assay Comments: Copyright © 2003-2010 Agilent Technologies




Chip Information:

Chip Lot #:
Reagent Kit Lot #:
Chip Comments:

POKID_121923

NNS (1:2)

29_NARMS (1:20)

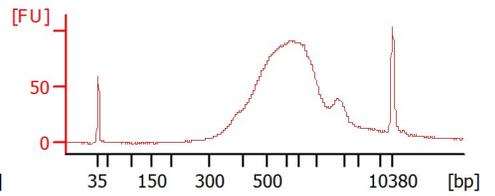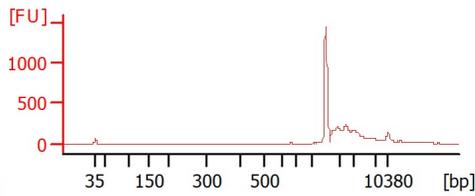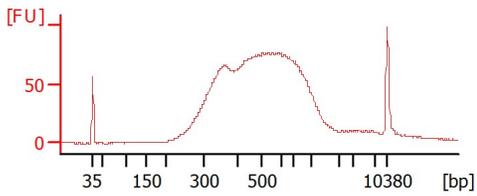

Assay Class: High Sensitivity DNA Assay
 Data Path: C:\...zer\2100 expert\data\2024-01-03\New Projects_01-13-2024.xad

Created: 1/3/2024 4:39:30 PM
 Modified: 1/3/2024 5:29:16 PM

Electrophoresis File Run Summary (Chip Summary)

Sample Name	Sample Comment	Rest. Digest	Status	Observation	Result Label	Result Color
POKID_121923		<input type="checkbox"/>	✓			
NNS (1:2)		<input type="checkbox"/>	✓			
29_NARMS (1:20)		<input type="checkbox"/>	✓			
Ladder		<input type="checkbox"/>	✓			

Chip Lot #

Reagent Kit Lot #

Chip Comments :

Assay Class: High Sensitivity DNA Assay
Data Path: C:\...zer\2100 expert\data\2024-01-03\New Projects_01-13-2024.xad

Created: 1/3/2024 4:39:30 PM
Modified: 1/3/2024 5:29:16 PM

Electrophoresis Assay Details

General Analysis Settings

Number of Available Sample and Ladder Wells (Max.) : 12
Minimum Visible Range [s] : 32
Maximum Visible Range [s] : 138
Start Analysis Time Range [s] : 33
End Analysis Time Range [s] : 137.5
Ladder Concentration [pg/μl] : 1950
Uses Standard Area for Ladder Fragments
Lower Marker Concentration [pg/μl] : 125
Upper Marker Concentration [pg/μl] : 75
Used Upper Marker for Quantitation
Standard Curve Fit is Point to Point
Show Data Aligned to Lower and Upper Marker

Integrator Settings

Integration Start Time [s] : 33.05
Integration End Time [s] : 137
Slope Threshold : 0.8
Height Threshold [FU] : 5
Area Threshold : 0.1
Width Threshold [s] : 0.6
Baseline Plateau [s] : 0.5

Filter Settings

Filter Width [s] : 0.5
Polynomial Order : 4


Ladder

Ladder Peak	Size	Area
1	35	160
2	50	210
3	100	208
4	150	221
5	200	242
6	300	270
7	400	305
8	500	306
9	600	336
10	700	321
11	1000	366
12	2000	413
13	3000	411
14	7000	400
15	10380	214

Assay Class: High Sensitivity DNA Assay
 Data Path: C:\...zer\2100 expert\data\2024-01-03\New Projects_01-13-2024.xad

Created: 1/3/2024 4:39:30 PM
 Modified: 1/3/2024 5:29:16 PM

Electropherogram Summary

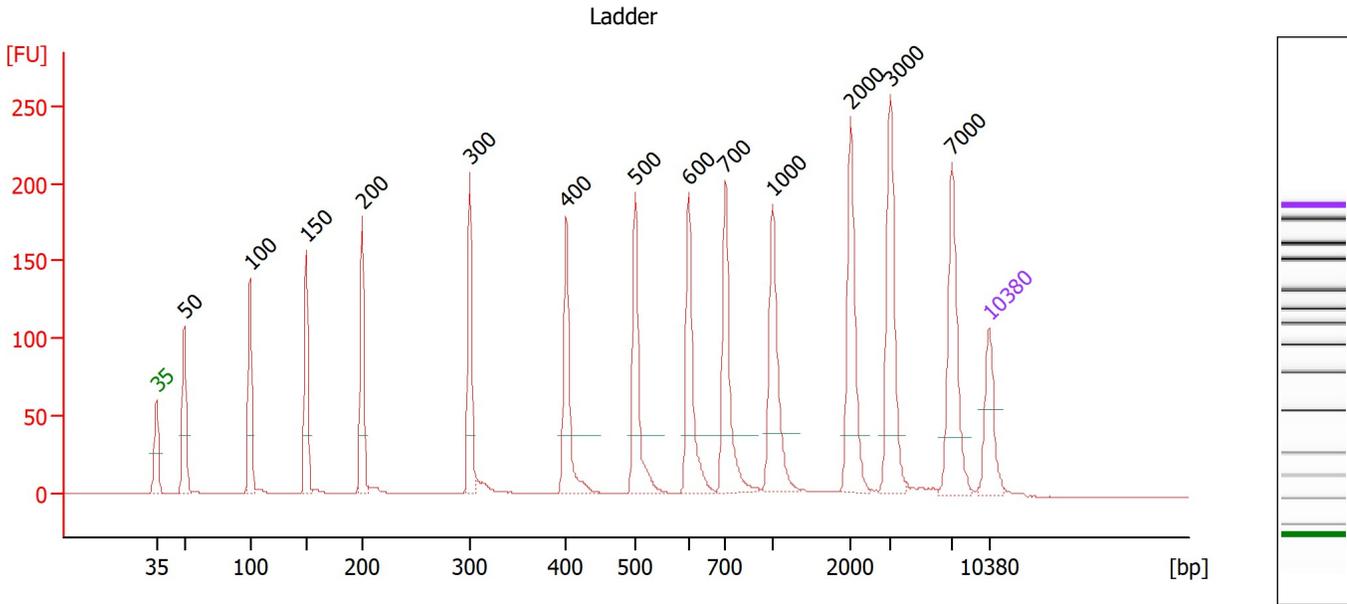

Overall Results for Ladder

Noise: 0.2


Peak table for Ladder

Peak	Size [bp]	Conc. [pg/μl]	Molarity [pmol/l]	Observations
1	35	125.00	5,411.3	Lower Marker
2	50	150.00	4,545.5	Ladder Peak
3	100	150.00	2,272.7	Ladder Peak
4	150	150.00	1,515.2	Ladder Peak
5	200	150.00	1,136.4	Ladder Peak
6	300	150.00	757.6	Ladder Peak
7	400	150.00	568.2	Ladder Peak
8	500	150.00	454.5	Ladder Peak
9	600	150.00	378.8	Ladder Peak
10	700	150.00	324.7	Ladder Peak
11	1,000	150.00	227.3	Ladder Peak
12	2,000	150.00	113.6	Ladder Peak
13	3,000	150.00	75.8	Ladder Peak
14	7,000	150.00	32.5	Ladder Peak
15	10,380	75.00	10.9	Upper Marker

Assay Class: High Sensitivity DNA Assay
 Data Path: C:\...zer\2100 expert\data\2024-01-03\New Projects_01-13-2024.xad

Created: 1/3/2024 4:39:30 PM
 Modified: 1/3/2024 5:29:16 PM

Electropherogram Summary Continued ...

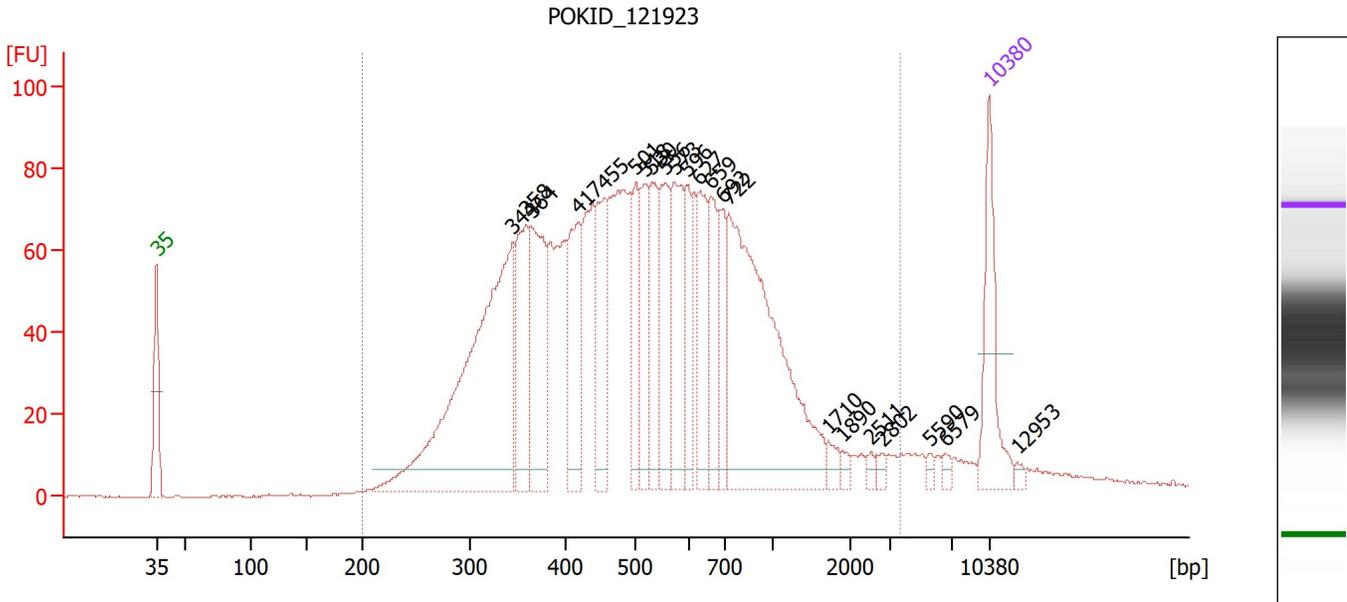

Overall Results for sample 1 : POKID_121923

Number of peaks found: 22 Corr. Area 1: 2,319.5
 Noise: 0.1

Peak table for sample 1 : POKID_121923

Peak	Size [bp]	Conc. [pg/μl]	Molarity [pmol/l]	Observations
1	35	125.00	5,411.3	Lower Marker
2	344	493.93	2,176.4	
3	358	129.34	546.7	
4	364	149.28	621.8	
5	417	117.80	428.0	
6	455	109.29	363.9	
7	501	75.45	228.4	
8	518	73.99	216.2	
9	530	95.17	271.9	
10	556	91.43	249.0	
11	573	110.43	292.0	
12	596	65.47	166.6	
13	627	88.43	213.7	
14	659	76.96	177.1	
15	692	65.35	143.1	
16	722	403.88	847.3	
17	1,710	11.25	10.0	
18	1,890	7.07	5.7	
19	2,511	5.63	3.4	
20	2,802	6.37	3.4	
21	5,590	5.03	1.4	
22	6,579	6.40	1.5	
23	10,380	75.00	10.9	Upper Marker
24	12,953	0.00	0.0	

Assay Class: High Sensitivity DNA Assay
Data Path: C:\...zer\2100 expert\data\2024-01-03\New Projects_01-13-2024.xad

Created: 1/3/2024 4:39:30 PM
Modified: 1/3/2024 5:29:16 PM

Electropherogram Summary Continued ...

... Region table for sample 1 :


POKID_121923

From [bp]	To [bp]	Corr. Area	% of Total	Average Size [bp]	Size distribution in CV [%]	Conc. [pg/ μ l]	Molarity [pmol/l]	Color
200	3,655	2,319.5	95	613	67.2	2,830.07	9,355.4	

Assay Class: High Sensitivity DNA Assay
 Data Path: C:\...zer\2100 expert\data\2024-01-03\New Projects_01-13-2024.xad

Created: 1/3/2024 4:39:30 PM
 Modified: 1/3/2024 5:29:16 PM

Electropherogram Summary Continued ...

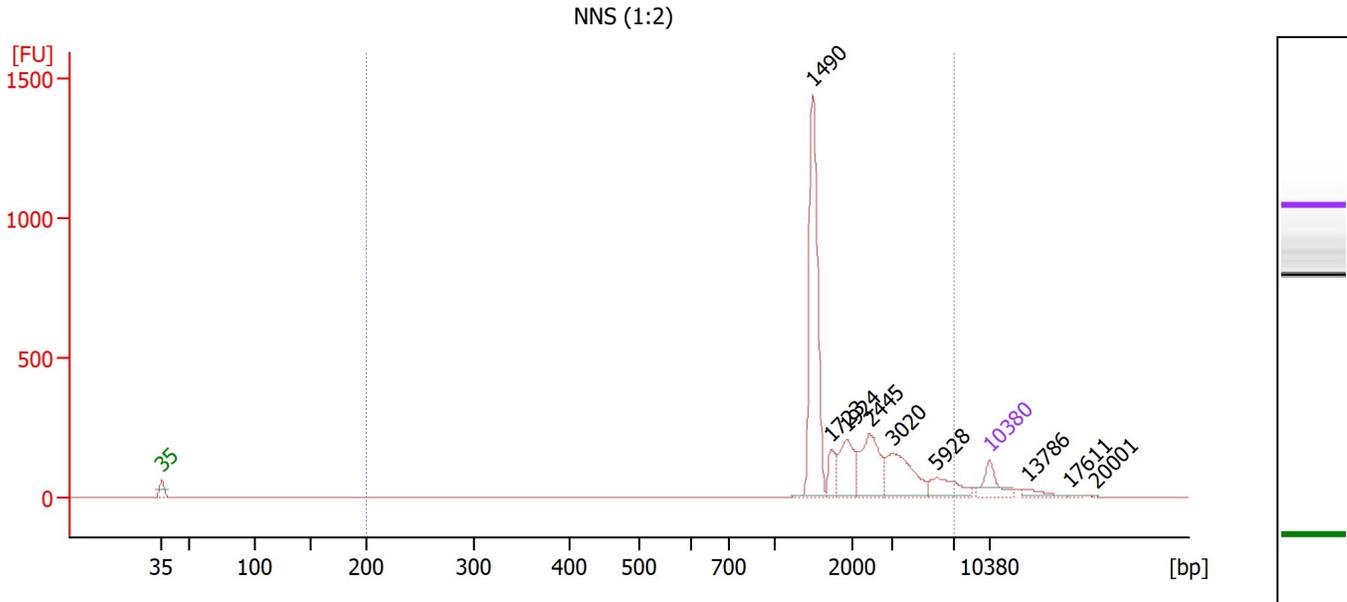

Overall Results for sample 2 : NNS (1:2)

Number of peaks found: 9 Corr. Area 1: 2,596.1
 Noise: 0.3

Peak table for sample 2 : NNS (1:2)

Peak	Size [bp]	Conc. [pg/μl]	Molarity [pmol/l]	Observations
1	35	125.00	5,411.3	Lower Marker
2	1,490	620.36	630.9	
3	1,723	62.83	55.2	
4	1,924	143.81	113.3	
5	2,445	195.63	121.2	
6	3,020	194.23	97.5	
7	5,928	90.90	23.2	
8	10,380	75.00	10.9	Upper Marker
9	13,786	0.00	0.0	
10	17,611	0.00	0.0	
11	20,001	0.00	0.0	


Region table for sample 2 : NNS (1:2)

From [bp]	To [bp]	Corr. Area	% of Total	Average Size [bp]	Size distribution in CV [%]	Conc. [pg/μl]	Molarity [pmol/l]	Color
200	7,140	2,596.1	93	2,321	54.1	1,288.98	1,024.8	Blue

Assay Class: High Sensitivity DNA Assay
 Data Path: C:\...zer\2100 expert\data\2024-01-03\New Projects_01-13-2024.xad

Created: 1/3/2024 4:39:30 PM
 Modified: 1/3/2024 5:29:16 PM

Electropherogram Summary Continued ...

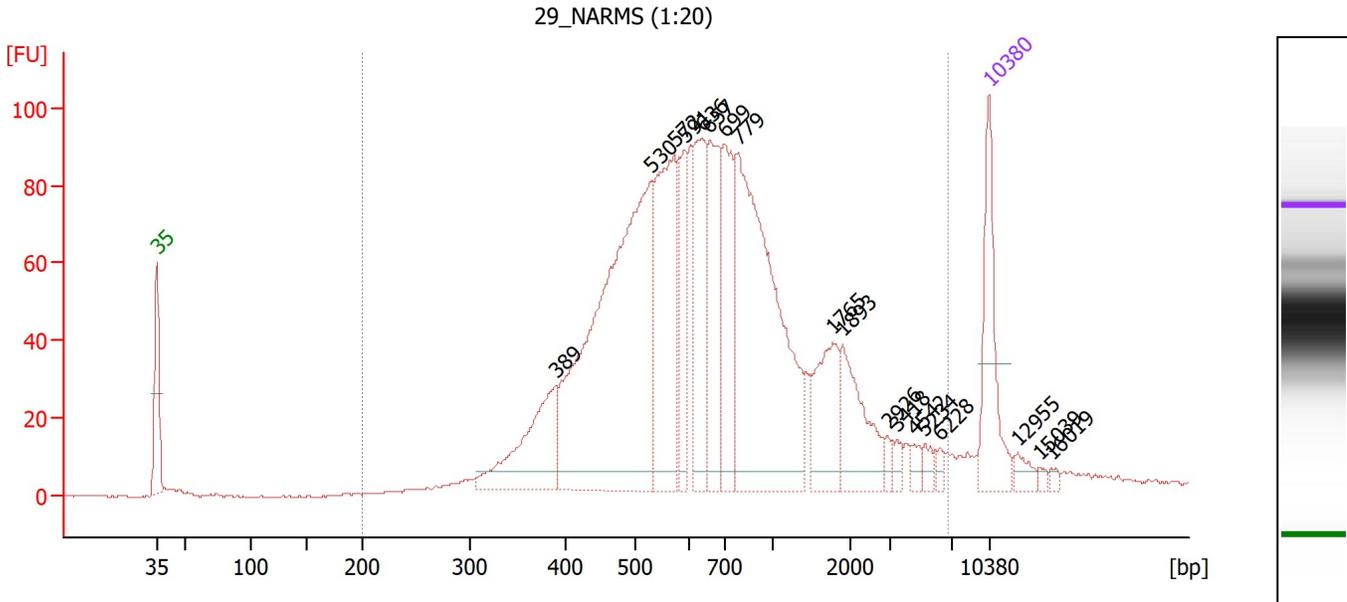

Overall Results for sample 3 : 29_NARMS (1:20)

Number of peaks found: 18 Corr. Area 1: 1,968.2
 Noise: 0.3

Peak table for sample 3 : 29_NARMS (1:20)

Peak	Size [bp]	Conc. [pg/μl]	Molarity [pmol/l]	Observations
1	35	125.00	5,411.3	Lower Marker
2	389	123.51	481.3	
3	530	510.35	1,459.2	
4	572	193.04	511.1	
5	591	70.87	181.7	
6	636	112.20	267.2	
7	657	122.43	282.5	
8	699	104.56	226.6	
9	779	364.56	709.0	
10	1,765	68.45	58.8	
11	1,893	70.90	56.7	
12	2,926	7.55	3.9	
13	3,418	7.04	3.1	
14	4,542	7.97	2.7	
15	5,234	8.49	2.5	
16	6,228	6.31	1.5	
17	10,380	75.00	10.9	Upper Marker
18	12,955	0.00	0.0	
19	15,039	0.00	0.0	
20	16,019	0.00	0.0	


Region table for sample 3 : 29_NARMS (1:20)

From [bp]	To [bp]	Corr. Area	% of Total	Average Size [bp]	Size distribution in CV [%]	Conc. [pg/μl]	Molarity [pmol/l]	Color
200	6,751	1,968.2	94	937	96.7	1,985.28	5,005.1	Blue

Assay Class: High Sensitivity DNA Assay
Data Path: C:\...zer\2100 expert\data\2024-01-03\New Projects_01-13-2024.xad

Created: 1/3/2024 4:39:30 PM
Modified: 1/3/2024 5:29:16 PM

Gel Image

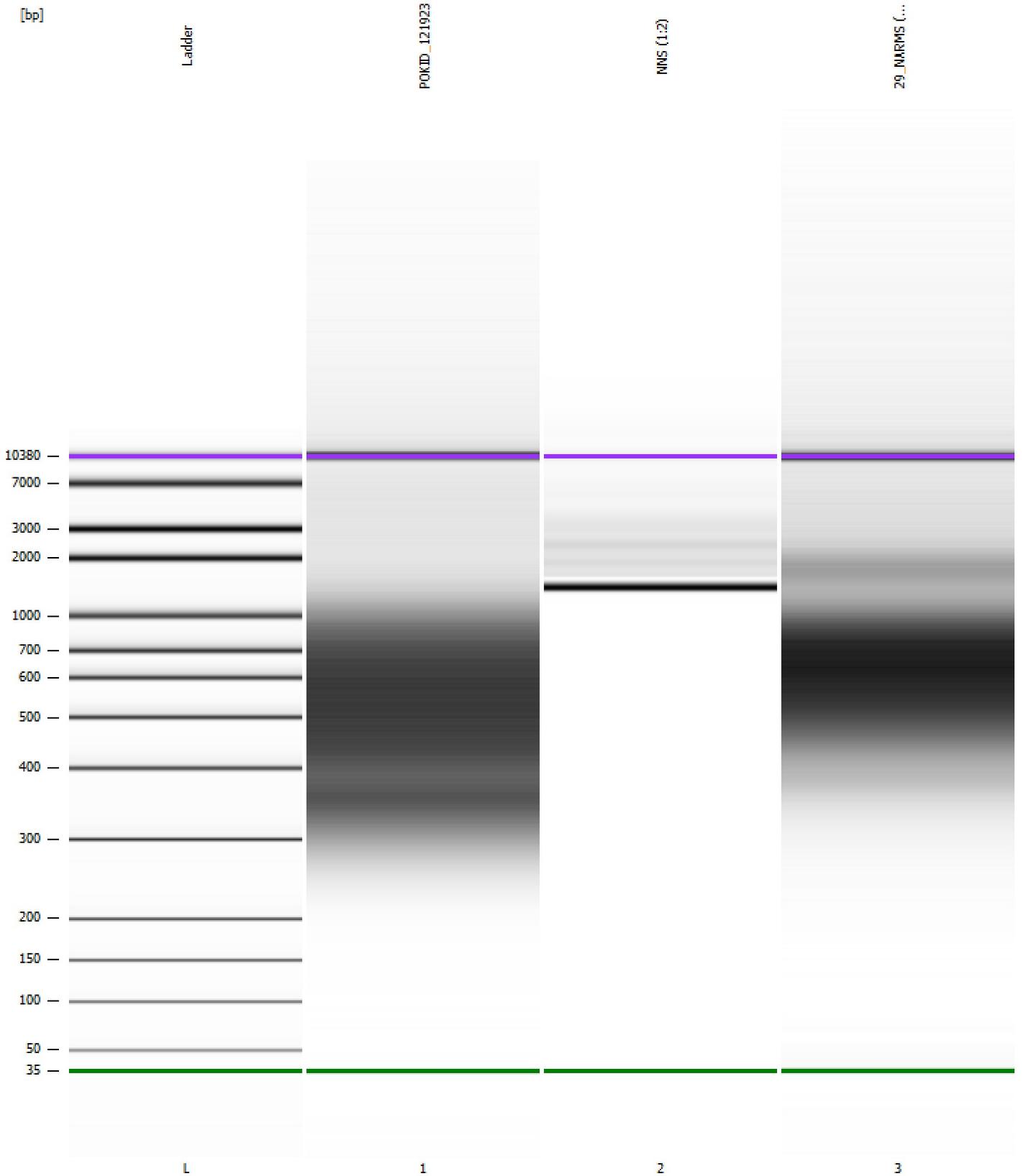


Assay Class: High Sensitivity DNA Assay
Data Path: C:\...zer\2100 expert\data\2024-01-03\New Projects_01-13-2024.xad

Created: 1/3/2024 4:39:30 PM
Modified: 1/3/2024 5:29:16 PM

Curves

Standard Curve

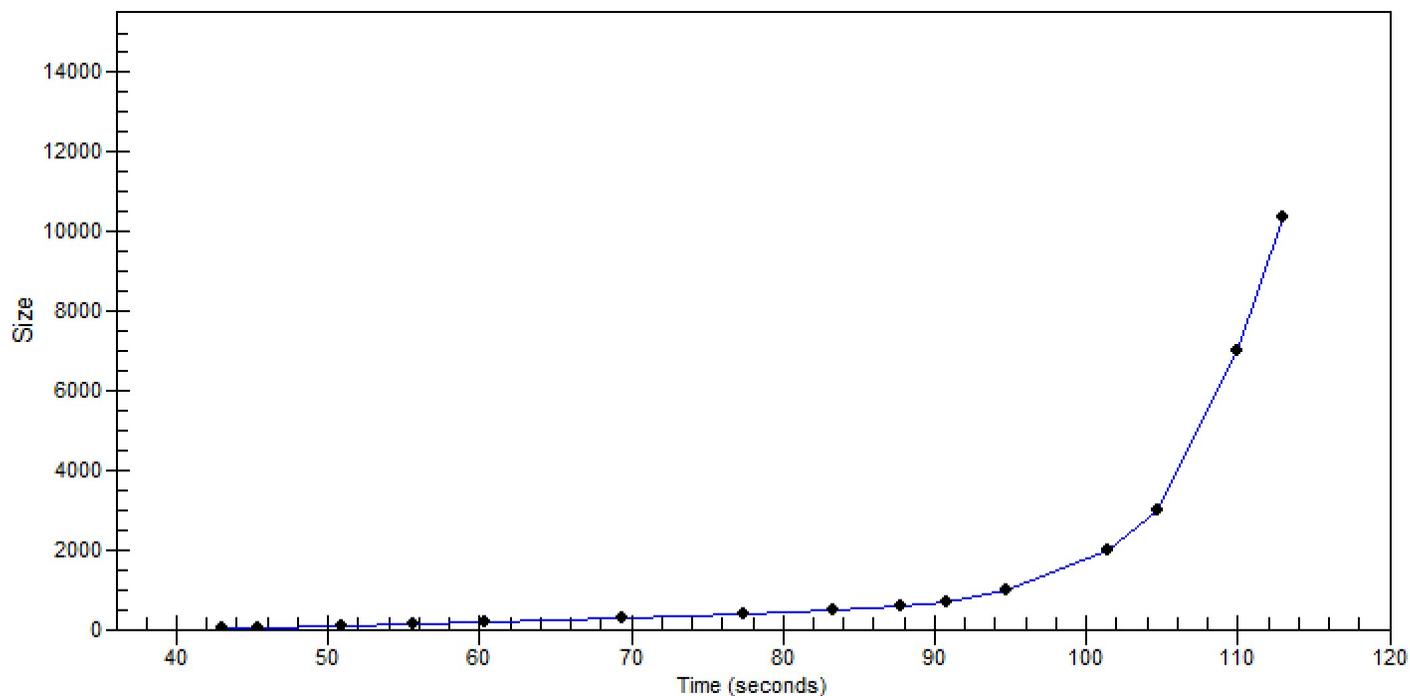

Assay Class: High Sensitivity DNA Assay
 Data Path: C:\...zer\2100 expert\data\2024-01-03\New Projects_01-13-2024.xad

Created: 1/3/2024 4:39:30 PM
 Modified: 1/3/2024 5:29:16 PM

Run Logbook

Description	Number	Source	Category	Sub Category	Time	Time Zone	User	Host
Run ended on port 1 (Number of wells acquired: 12)		Instrument	Run		1/3/2024 5:20:49 PM	(GMT --08:00) Pacific Standard Time	sbsuser	DESKTOP-4UNV VOB
Run started on port 1 (File: C:\Program Files (x86)\Agilent\2100 bioanalyzer\2100 expert\data\2024-01-03\Bioanalyzer1_High Sensitivity DNA Assay_DE34903152_2024-01-03_001.xad)		Instrument	Run		1/3/2024 4:39:37 PM	(GMT --08:00) Pacific Standard Time	sbsuser	DESKTOP-4UNV VOB
Product Number : G2938C		Instrument	Run		1/3/2024 4:39:37 PM	(GMT --08:00) Pacific Standard Time	sbsuser	DESKTOP-4UNV VOB
Name :		Instrument	Run		1/3/2024 4:39:37 PM	(GMT --08:00) Pacific Standard Time	sbsuser	DESKTOP-4UNV VOB
Vendor : Agilent Technologies		Instrument	Run		1/3/2024 4:39:37 PM	(GMT --08:00) Pacific Standard Time	sbsuser	DESKTOP-4UNV VOB
Serial# : DE34903152		Instrument	Run		1/3/2024 4:39:37 PM	(GMT --08:00) Pacific Standard Time	sbsuser	DESKTOP-4UNV VOB
Firmware : C.01.069		Instrument	Run		1/3/2024 4:39:37 PM	(GMT --08:00) Pacific Standard Time	sbsuser	DESKTOP-4UNV VOB
Cartridge : Electrode		Instrument	Run		1/3/2024 4:39:36 PM	(GMT --08:00) Pacific Standard Time	sbsuser	DESKTOP-4UNV VOB