


Assay Class: High Sensitivity DNA Assay
Data Path: C:\...2100 expert\data\2023-03-17\PB896_SPKv3_Amplicon_Lib_n2.xad

Created: 3/17/2023 4:37:39 PM
Modified: 3/20/2023 8:59:45 AM

Electrophoresis File Run Summary

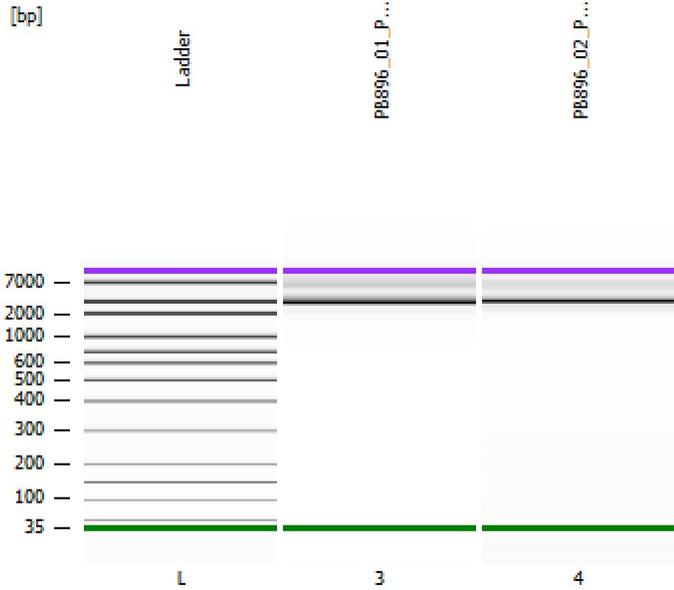

Instrument Information:

Instrument Name: DE34903152 Firmware: C.01.069
Serial#: DE34903152 Type: G2938C

Assay Information:



Assay Origin Path: C:\Program Files (x86)\Agilent\2100 bioanalyzer\2100 expert\assays\dsDNA\High Sensitivity DNA.xsy
Assay Class: High Sensitivity DNA Assay
Version: 1.03
Assay Comments: Copyright © 2003-2010 Agilent Technologies

Chip Information:

Chip Lot #:
Reagent Kit Lot #:
Chip Comments:

PB896_01_PBEH1_1_SPKv3_Amplicon_Lib

PB896_02_PBEH1_2_SPKv3_Amplicon_Lib

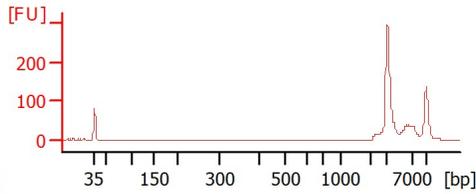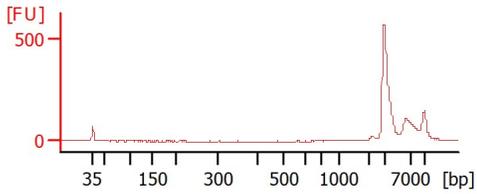

Assay Class: High Sensitivity DNA Assay
 Data Path: C:\...2100 expert\data\2023-03-17\PB896_SPKv3_Amplicon_Lib_n2.xad

Created: 3/17/2023 4:37:39 PM
 Modified: 3/20/2023 8:59:45 AM

Electrophoresis File Run Summary (Chip Summary)

Sample Name	Sample Comment	Rest. Digest	Status	Observation	Result Label	Result Color
PB896_01_PBEH1_1_SP Kv3_Amplicon_Lib		<input type="checkbox"/>	✓			
PB896_02_PBEH1_2_SP Kv3_Amplicon_Lib		<input type="checkbox"/>	✓			
Ladder		<input type="checkbox"/>	✓			

Chip Lot #

Reagent Kit Lot #

Chip Comments :

Assay Class: High Sensitivity DNA Assay
Data Path: C:\...2100 expert\data\2023-03-17\PB896_SPKv3_Amplicon_Lib_n2.xad

Created: 3/17/2023 4:37:39 PM
Modified: 3/20/2023 8:59:45 AM

Electrophoresis Assay Details

General Analysis Settings

Number of Available Sample and Ladder Wells (Max.) : 12
Minimum Visible Range [s] : 32
Maximum Visible Range [s] : 138
Start Analysis Time Range [s] : 33
End Analysis Time Range [s] : 137.5
Ladder Concentration [pg/ μ l] : 1950
Uses Standard Area for Ladder Fragments
Lower Marker Concentration [pg/ μ l] : 125
Upper Marker Concentration [pg/ μ l] : 75
Used Upper Marker for Quantitation
Standard Curve Fit is Point to Point
Show Data Aligned to Lower and Upper Marker

Integrator Settings

Integration Start Time [s] : 33.05
Integration End Time [s] : 137
Slope Threshold : 0.8
Height Threshold [FU] : 5
Area Threshold : 0.1
Width Threshold [s] : 0.6
Baseline Plateau [s] : 0.5

Filter Settings

Filter Width [s] : 0.5
Polynomial Order : 4


Ladder

Ladder Peak	Size	Area
1	35	160
2	50	210
3	100	208
4	150	221
5	200	242
6	300	270
7	400	305
8	500	306
9	600	336
10	700	321
11	1000	366
12	2000	413
13	3000	411
14	7000	400
15	10380	214

Assay Class: High Sensitivity DNA Assay
 Data Path: C:\...2100 expert\data\2023-03-17\PB896_SPKv3_Amplicon_Lib_n2.xad

Created: 3/17/2023 4:37:39 PM
 Modified: 3/20/2023 8:59:45 AM

Electropherogram Summary

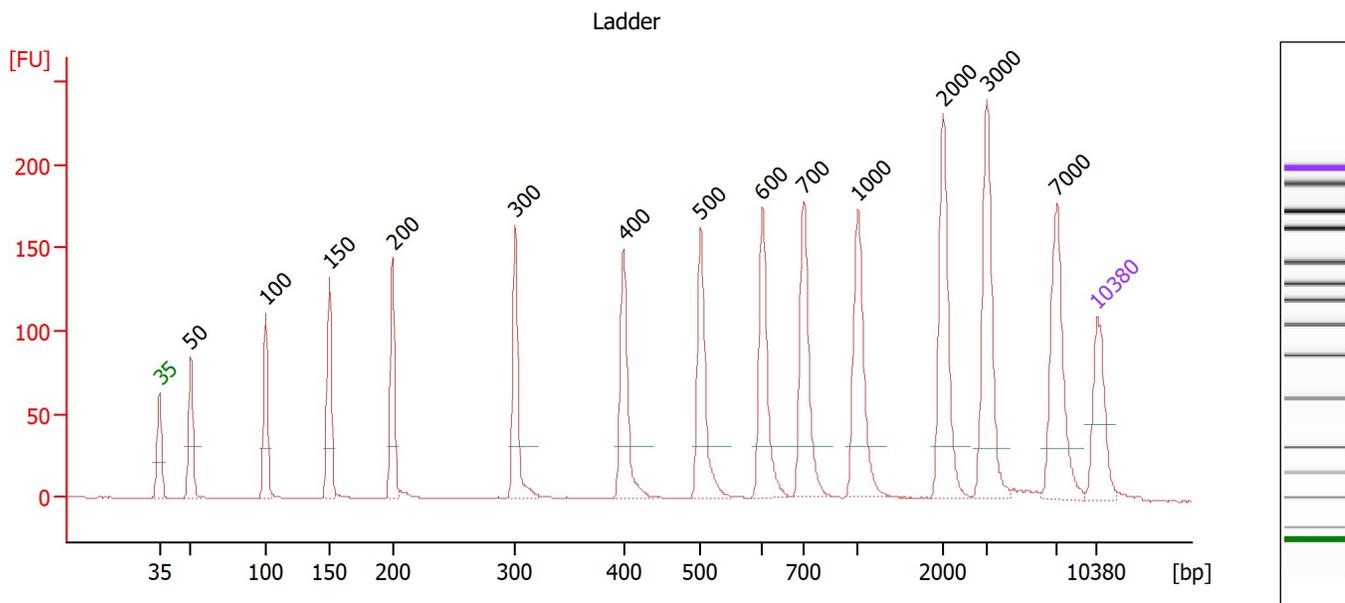

Overall Results for Ladder

Noise: 0.2


Peak table for Ladder

Peak	Size [bp]	Conc. [pg/μl]	Molarity [pmol/l]	Observations
1	35	125.00	5,411.3	Lower Marker
2	50	150.00	4,545.5	Ladder Peak
3	100	150.00	2,272.7	Ladder Peak
4	150	150.00	1,515.2	Ladder Peak
5	200	150.00	1,136.4	Ladder Peak
6	300	150.00	757.6	Ladder Peak
7	400	150.00	568.2	Ladder Peak
8	500	150.00	454.5	Ladder Peak
9	600	150.00	378.8	Ladder Peak
10	700	150.00	324.7	Ladder Peak
11	1,000	150.00	227.3	Ladder Peak
12	2,000	150.00	113.6	Ladder Peak
13	3,000	150.00	75.8	Ladder Peak
14	7,000	150.00	32.5	Ladder Peak
15	10,380	75.00	10.9	Upper Marker

Assay Class: High Sensitivity DNA Assay
 Data Path: C:\...2100 expert\data\2023-03-17\PB896_SPKv3_Amplicon_Lib_n2.xad

Created: 3/17/2023 4:37:39 PM
 Modified: 3/20/2023 8:59:45 AM

Electropherogram Summary Continued ...

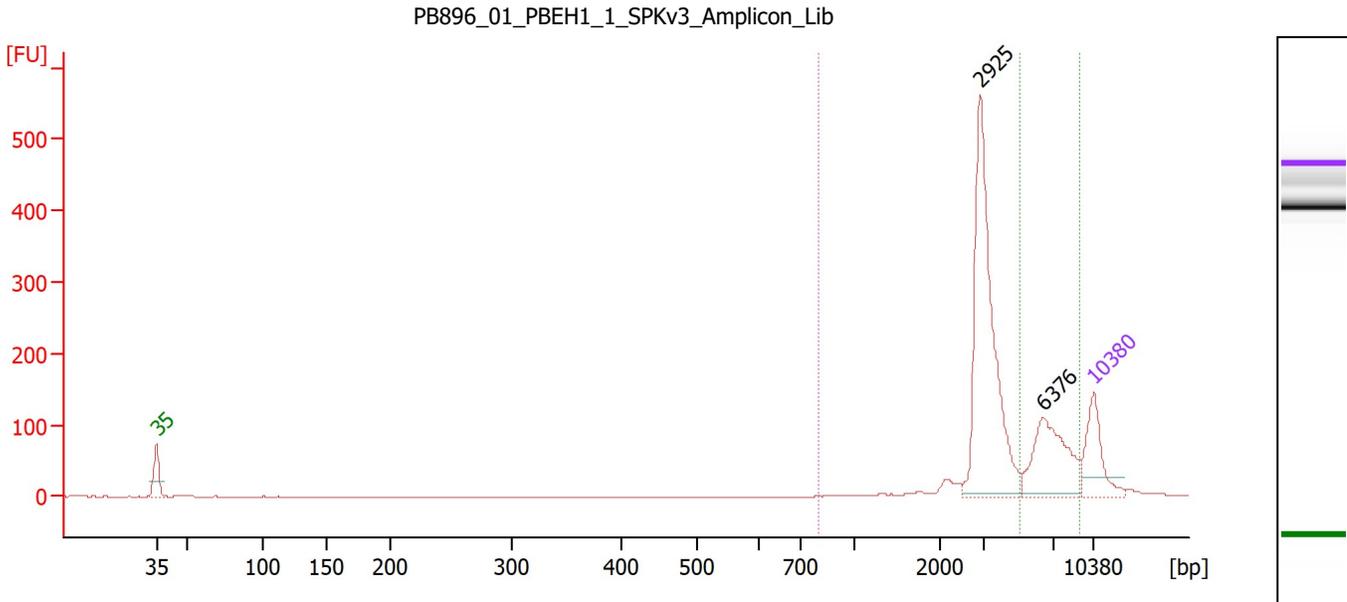

Overall Results for sample 3 : PB896_01_PBEH1_1_SPKv3_Amplicon_Lib

Number of peaks found:	2	Corr. Area 2:	859.6
Noise:	0.3	Corr. Area 3:	297.7
Corr. Area 1:	1,157.3		

Peak table for sample 3 : PB896_01_PBEH1_1_SPKv3_Amplicon_Lib

Peak	Size [bp]	Conc. [pg/μl]	Molarity [pmol/l]	Observations
1	35	125.00	5,411.3	Lower Marker
2	2,925	366.15	189.7	
3	6,376	138.45	32.9	
4	10,380	75.00	10.9	Upper Marker

Region table for sample 3 : PB896_01_PBEH1_1_SPKv3_Amplicon_Lib


From [bp]	To [bp]	Corr. Area	% of Total	Average Size [bp]	Size distribution in CV [%]	Conc. [pg/μl]	Molarity [pmol/l]	Color
800	9,200	1,157.3	97	4,160	43.8	527.34	225.6	Blue
800	5,000	859.6	72	3,157	19.3	390.23	195.0	Pink
5,000	9,200	297.7	25	6,911	15.0	137.11	30.6	Green

Assay Class: High Sensitivity DNA Assay
 Data Path: C:\...2100 expert\data\2023-03-17\PB896_SPKv3_Amplicon_Lib_n2.xad

Created: 3/17/2023 4:37:39 PM
 Modified: 3/20/2023 8:59:45 AM

Electropherogram Summary Continued ...

PB896_02_PBEH1_2_SPKv3_Amplicon_Lib

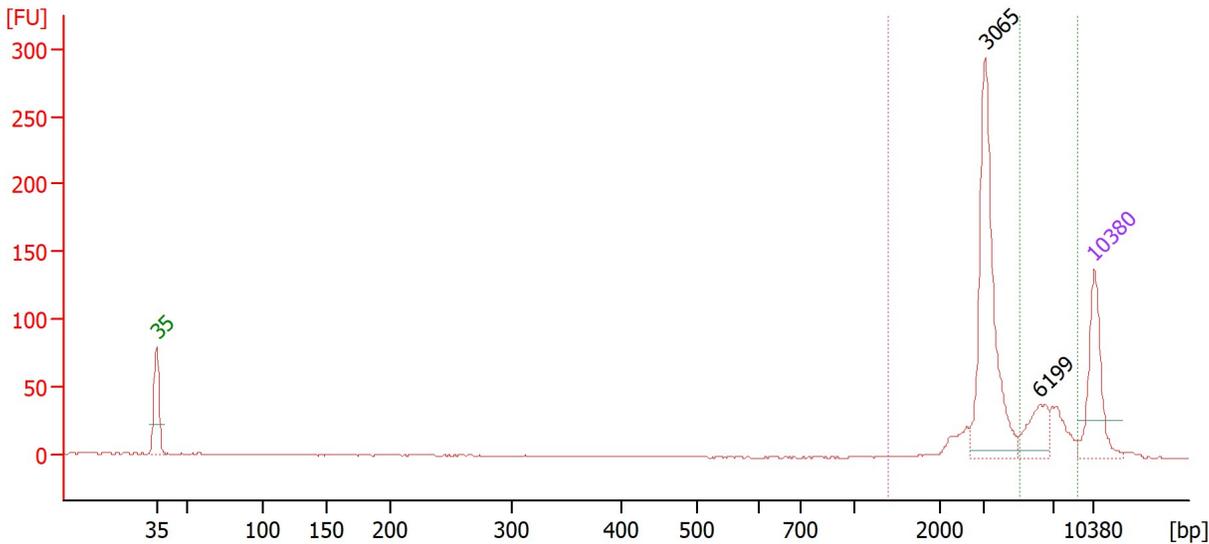

Overall Results for sample 4 : PB896_02_PBEH1_2_SPKv3_Amplicon_Lib

Number of peaks found: 2 Corr. Area 2: 373.6
 Noise: 0.3 Corr. Area 3: 113.3
 Corr. Area 1: 486.9

Peak table for sample 4 : PB896_02_PBEH1_2_SPKv3_Amplicon_Lib

Peak	Size [bp]	Conc. [pg/μl]	Molarity [pmol/l]	Observations
1	35	125.00	5,411.3	Lower Marker
2	3,065	200.56	99.1	
3	6,199	39.28	9.6	
4	10,380	75.00	10.9	Upper Marker


Region table for sample 4 : PB896_02_PBEH1_2_SPKv3_Amplicon_Lib

From [bp]	To [bp]	Corr. Area	% of Total	Average Size [bp]	Size distribution in CV [%]	Conc. [pg/μl]	Molarity [pmol/l]	Color
1,400	9,000	486.9	96	4,048	39.9	292.02	122.7	Blue
1,400	5,000	373.6	73	3,219	16.5	223.13	106.8	Pink
5,000	9,000	113.3	22	6,661	14.3	68.89	15.9	Green

Assay Class: High Sensitivity DNA Assay
Data Path: C:\...2100 expert\data\2023-03-17\PB896_SPKv3_Amplicon_Lib_n2.xad

Created: 3/17/2023 4:37:39 PM
Modified: 3/20/2023 8:59:45 AM

Gel Image

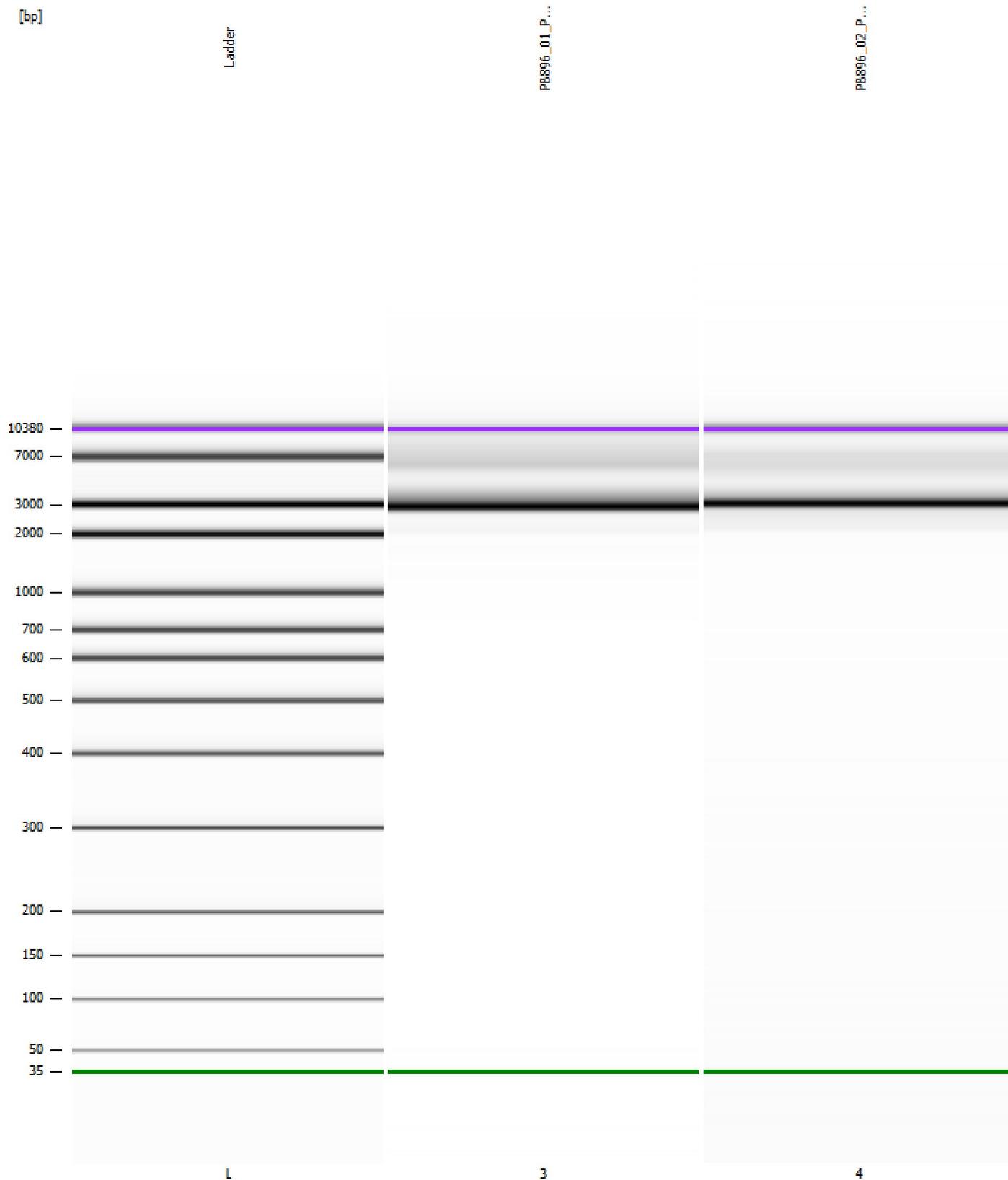

Assay Class: High Sensitivity DNA Assay
 Data Path: C:\...2100 expert\data\2023-03-17\PB896_SPKv3_Amplicon_Lib_n2.xad

Created: 3/17/2023 4:37:39 PM
 Modified: 3/20/2023 8:59:45 AM

Run Logbook

Description	Number	Source	Category	Sub Category	Time	Time Zone	User	Host
Run ended on port 1 (Number of wells acquired: 12)		Instrument	Run		3/17/2023 5:18:58 PM	(GMT --07:00) Pacific Standard Time	sbsuser	DESKTOP-4UNV VOB
Run started on port 1 (File: C:\Users\sbsuser\Desktop\2023-03-17\Bioanalyzer1_High Sensitivity DNA Assay_2023-03-17_001.xad)		Instrument	Run		3/17/2023 4:37:46 PM	(GMT --07:00) Pacific Standard Time	sbsuser	DESKTOP-4UNV VOB
Product Number : G2938C		Instrument	Run		3/17/2023 4:37:46 PM	(GMT --07:00) Pacific Standard Time	sbsuser	DESKTOP-4UNV VOB
Name :		Instrument	Run		3/17/2023 4:37:46 PM	(GMT --07:00) Pacific Standard Time	sbsuser	DESKTOP-4UNV VOB
Vendor : Agilent Technologies		Instrument	Run		3/17/2023 4:37:46 PM	(GMT --07:00) Pacific Standard Time	sbsuser	DESKTOP-4UNV VOB
Serial# : DE34903152		Instrument	Run		3/17/2023 4:37:46 PM	(GMT --07:00) Pacific Standard Time	sbsuser	DESKTOP-4UNV VOB
Firmware : C.01.069		Instrument	Run		3/17/2023 4:37:46 PM	(GMT --07:00) Pacific Standard Time	sbsuser	DESKTOP-4UNV VOB
Cartridge : Electrode		Instrument	Run		3/17/2023 4:37:46 PM	(GMT --07:00) Pacific Standard Time	sbsuser	DESKTOP-4UNV VOB