


Assay Class: High Sensitivity DNA Assay
Data Path: C:\...bioanalyzer\2100 expert\data\2025-06-18\QC2229_DV8-DV12.xad

Created: 6/18/2025 1:06:43 PM
Modified: 6/18/2025 2:40:02 PM

Electrophoresis File Run Summary

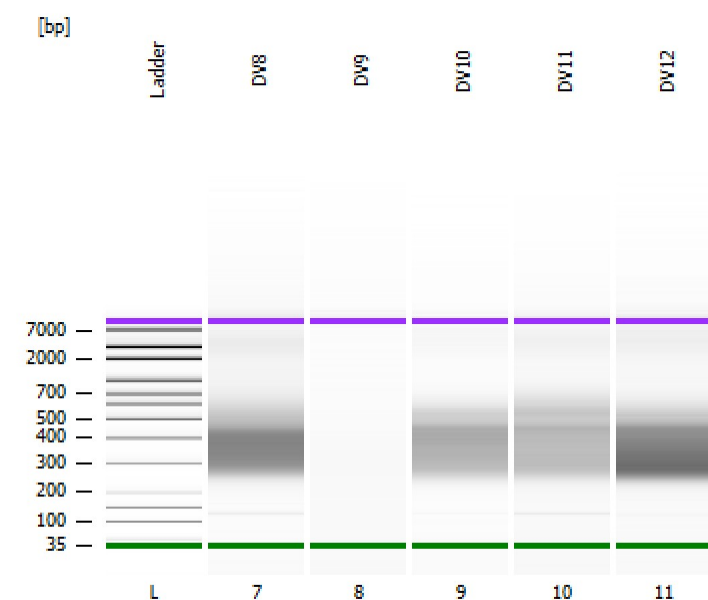

Instrument Information:






Instrument Name: DE34903152 Firmware: C.01.069
Serial#: DE34903152 Type: G2938C

Assay Information:

Assay Origin Path: C:\Program Files (x86)\Agilent\2100 bioanalyzer\2100 expert\assays\dsDNA\High Sensitivity DNA.xsy
Assay Class: High Sensitivity DNA Assay
Version: 1.03
Assay Comments: Copyright © 2003-2010 Agilent Technologies

Chip Information:

Chip Lot #:
Reagent Kit Lot #:
Chip Comments:

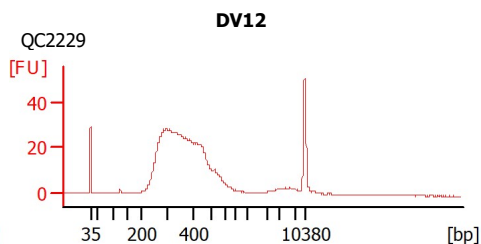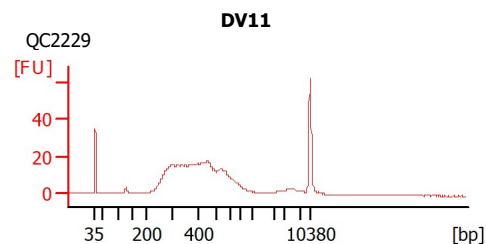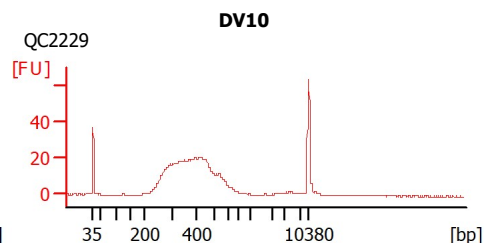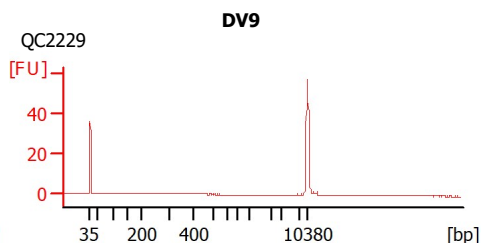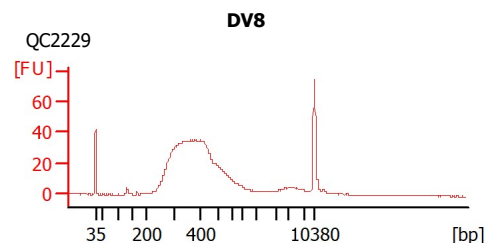

Assay Class: High Sensitivity DNA Assay
Data Path: C:\...bioanalyzer\2100 expert\data\2025-06-18\QC2229_DV8-DV12.xad

Created: 6/18/2025 1:06:43 PM
Modified: 6/18/2025 2:40:02 PM

Electrophoresis File Run Summary (Chip Summary)

Sample Name	Sample Comment	Rest. Digest	Status	Observation	Result Label	Result Color
DV8	QC2229	☐	✓			
DV9	QC2229	☐	✓			
DV10	QC2229	☐	✓			
DV11	QC2229	☐	✓			
DV12	QC2229	☐	✓			
Ladder		☐	✓			

Chip Lot #**Reagent Kit Lot #****Chip Comments :**

Assay Class: High Sensitivity DNA Assay
Data Path: C:\...bioanalyzer\2100 expert\data\2025-06-18\QC2229_DV8-DV12.xad

Created: 6/18/2025 1:06:43 PM
Modified: 6/18/2025 2:40:02 PM

Electrophoresis Assay Details

General Analysis Settings

Number of Available Sample and Ladder Wells (Max.) : 12
Minimum Visible Range [s] : 32
Maximum Visible Range [s] : 138
Start Analysis Time Range [s] : 33
End Analysis Time Range [s] : 137.5
Ladder Concentration [pg/μl] : 1950
Uses Standard Area for Ladder Fragments
Lower Marker Concentration [pg/μl] : 125
Upper Marker Concentration [pg/μl] : 75
Used Upper Marker for Quantitation
Standard Curve Fit is Point to Point
Show Data Aligned to Lower and Upper Marker

Integrator Settings

Integration Start Time [s] : 33.05
Integration End Time [s] : 137
Slope Threshold : 0.8
Height Threshold [FU] : 5
Area Threshold : 0.1
Width Threshold [s] : 0.6
Baseline Plateau [s] : 0.5

Filter Settings

Filter Width [s] : 0.5
Polynomial Order : 4


Ladder

Ladder Peak	Size	Area
1	35	160
2	50	210
3	100	208
4	150	221
5	200	242
6	300	270
7	400	305
8	500	306
9	600	336
10	700	321
11	1000	366
12	2000	413
13	3000	411
14	7000	400
15	10380	214

Assay Class: High Sensitivity DNA Assay
 Data Path: C:\...bioanalyzer\2100 expert\data\2025-06-18\QC2229_DV8-DV12.xad

Created: 6/18/2025 1:06:43 PM
 Modified: 6/18/2025 2:40:02 PM

Electropherogram Summary

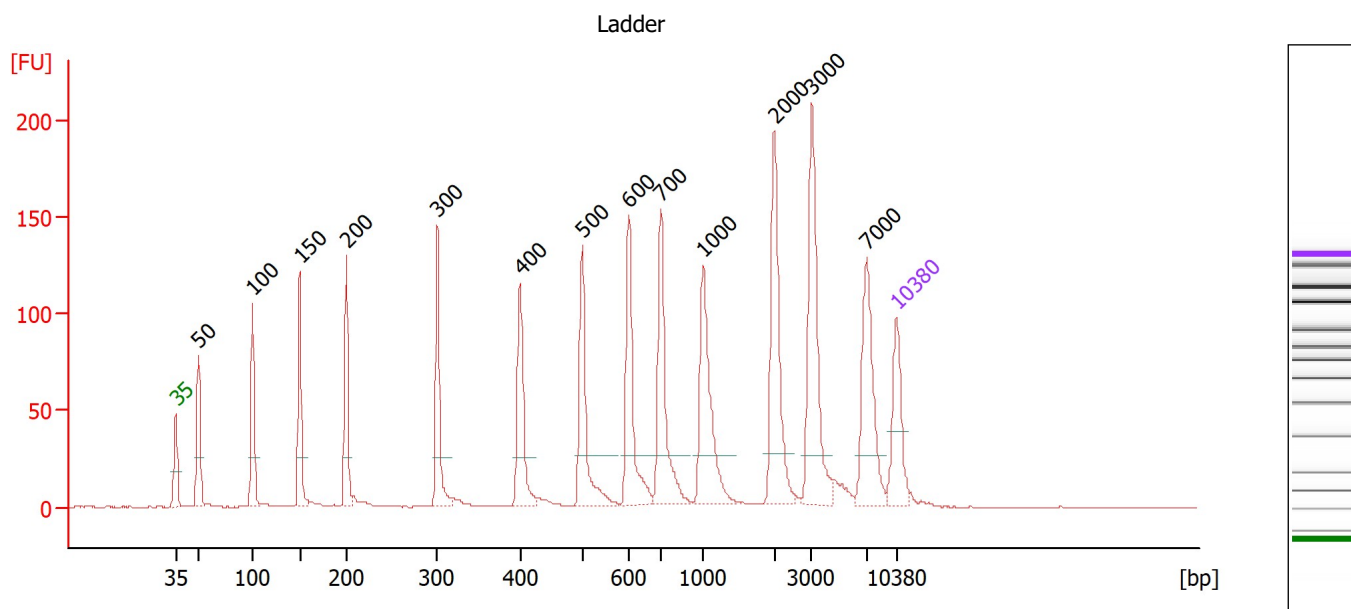

Overall Results for Ladder


Noise: 0.1

Peak table for Ladder

Peak	Size [bp]	Conc. [pg/μl]	Molarity [pmol/l]	Observations
1	35	125.00	5,411.3	Lower Marker
2	50	150.00	4,545.5	Ladder Peak
3	100	150.00	2,272.7	Ladder Peak
4	150	150.00	1,515.2	Ladder Peak
5	200	150.00	1,136.4	Ladder Peak
6	300	150.00	757.6	Ladder Peak
7	400	150.00	568.2	Ladder Peak
8	500	150.00	454.5	Ladder Peak
9	600	150.00	378.8	Ladder Peak
10	700	150.00	324.7	Ladder Peak
11	1,000	150.00	227.3	Ladder Peak
12	2,000	150.00	113.6	Ladder Peak
13	3,000	150.00	75.8	Ladder Peak
14	7,000	150.00	32.5	Ladder Peak
15	10,380	75.00	10.9	Upper Marker

Assay Class: High Sensitivity DNA Assay
Data Path: C:\...bioanalyzer\2100 expert\data\2025-06-18\QC2229_DV8-DV12.xad

Created: 6/18/2025 1:06:43 PM
Modified: 6/18/2025 2:40:02 PM

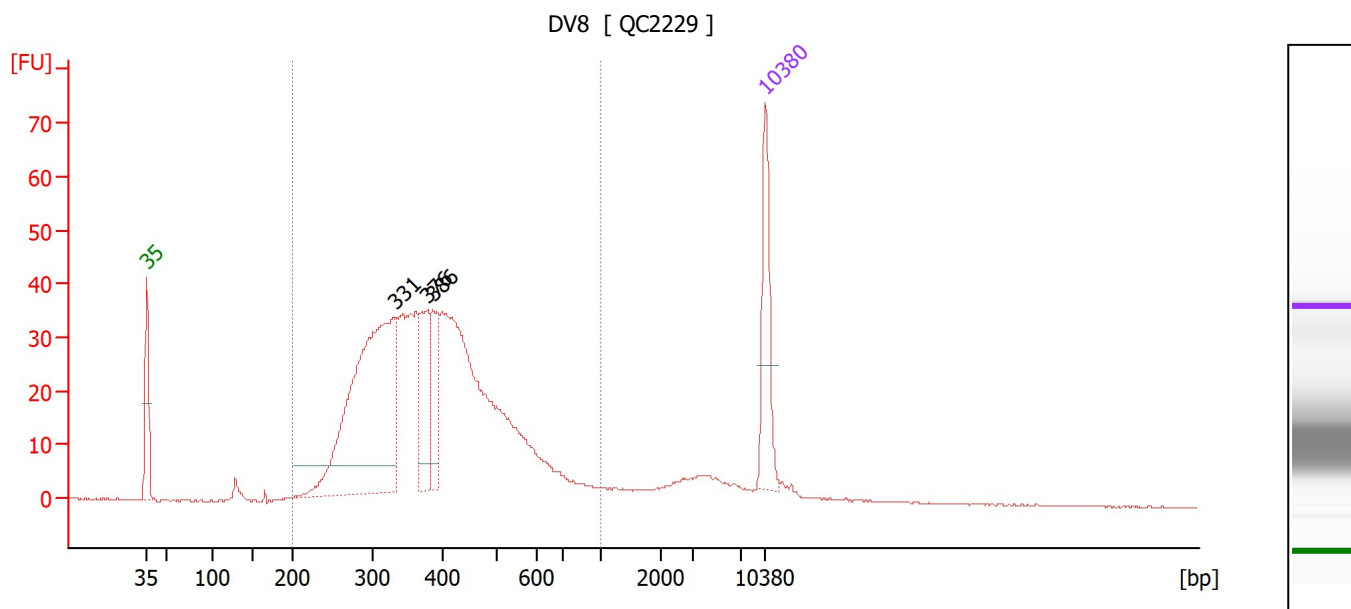Electropherogram Summary Continued ...**Overall Results for sample 7 : DV8**

Number of peaks found: 3 Corr. Area 1: 807.8
Noise: 0.2

Peak table for sample 7 : DV8


Peak	Size [bp]	Conc. [pg/μl]	Molarity [pmol/l]	Observations
1	35	125.00	5,411.3	Lower Marker
2	331	464.70	2,126.7	
3	376	95.33	384.0	
4	386	66.37	260.6	
5	10,380	75.00	10.9	Upper Marker

Region table for sample 7 : DV8

From [bp]	To [bp]	Corr. Area	% of Total	Average Size [bp]	Size distribution in CV [%]	Conc. [pg/μl]	Molarity [pmol/l]	Color
200	1,000	807.8	91	401	28.8	1,471.50	6,109.2	Blue

Assay Class: High Sensitivity DNA Assay
Data Path: C:\...bioanalyzer\2100 expert\data\2025-06-18\QC2229_DV8-DV12.xad

Created: 6/18/2025 1:06:43 PM
Modified: 6/18/2025 2:40:02 PM

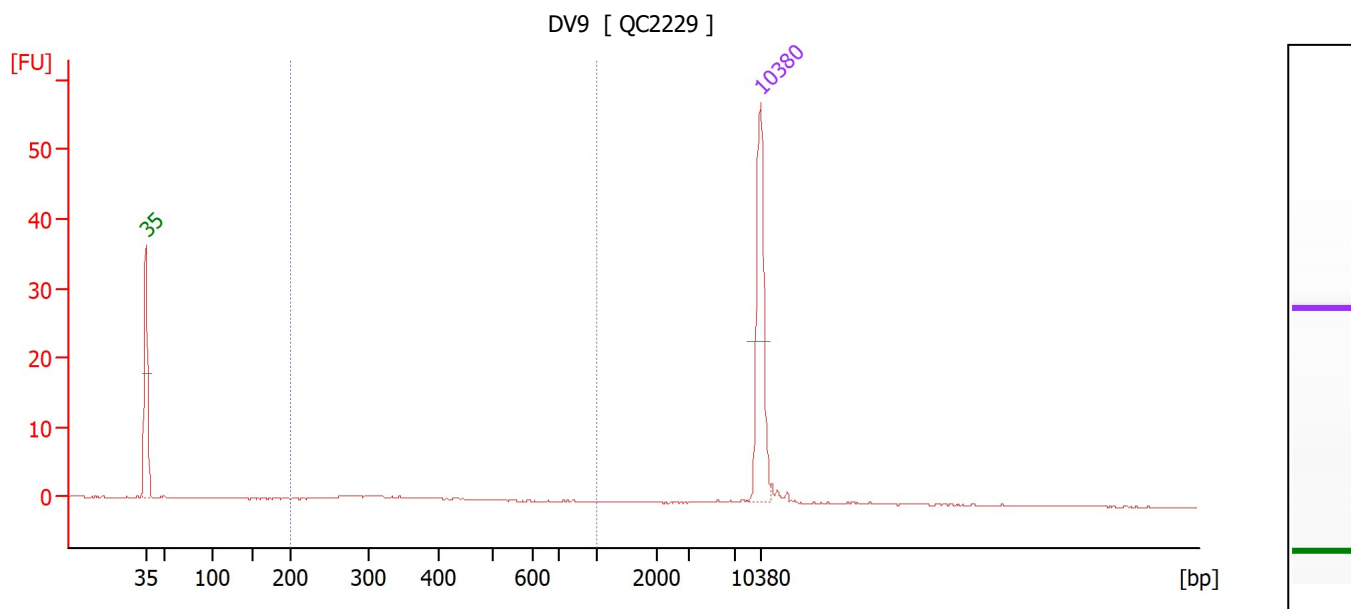Electropherogram Summary Continued ...**Overall Results for sample 8 : DV9**

Number of peaks found: 0 Corr. Area 1: 7.6
Noise: 0.1

Peak table for sample 8 : DV9

Peak	Size [bp]	Conc. [pg/μl]	Molarity [pmol/l]	Observations
1	35	125.00	5,411.3	Lower Marker
2	10,380	75.00	10.9	Upper Marker


Region table for sample 8 : DV9

From [bp]	To [bp]	Corr. Area	% of Total	Average Size [bp]	Size distribution in CV [%]	Conc. [pg/μl]	Molarity [pmol/l]	Color
200	1,000	7.6	53	350	27.5	18.25	85.2	Blue

Assay Class: High Sensitivity DNA Assay
 Data Path: C:\...bioanalyzer\2100 expert\data\2025-06-18\QC2229_DV8-DV12.xad

Created: 6/18/2025 1:06:43 PM
 Modified: 6/18/2025 2:40:02 PM

Electropherogram Summary Continued ...

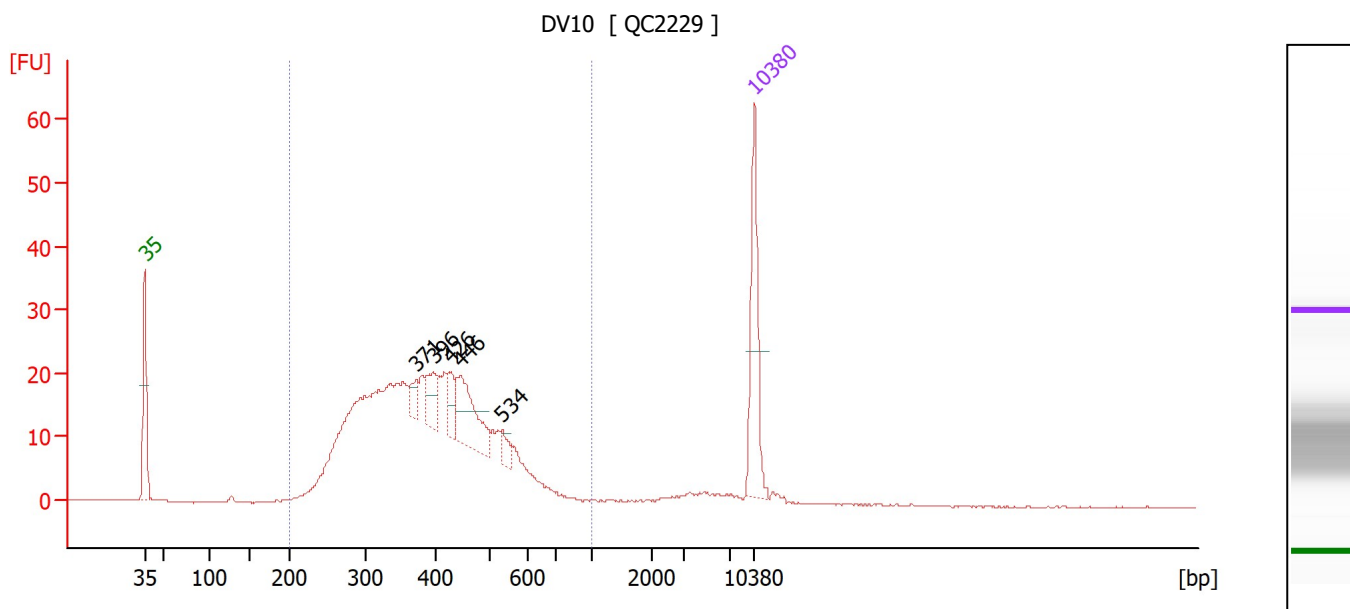

Overall Results for sample 9 : DV10

Number of peaks found: 5 Corr. Area 1: 472.4
 Noise: 0.1

Peak table for sample 9 : DV10


Peak	Size [bp]	Conc. [pg/μl]	Molarity [pmol/l]	Observations
1	35	125.00	5,411.3	Lower Marker
2	371	16.24	66.3	
3	396	31.63	121.1	
4	426	23.82	84.8	
5	446	70.78	240.5	
6	534	13.59	38.5	
7	10,380	75.00	10.9	Upper Marker

Region table for sample 9 : DV10

From [bp]	To [bp]	Corr. Area	% of Total	Average Size [bp]	Size distribution in CV [%]	Conc. [pg/μl]	Molarity [pmol/l]	Color
200	1,000	472.4	95	397	25.9	1,070.47	4,459.2	Blue

Assay Class: High Sensitivity DNA Assay
Data Path: C:\...bioanalyzer\2100 expert\data\2025-06-18\QC2229_DV8-DV12.xad

Created: 6/18/2025 1:06:43 PM
Modified: 6/18/2025 2:40:02 PM

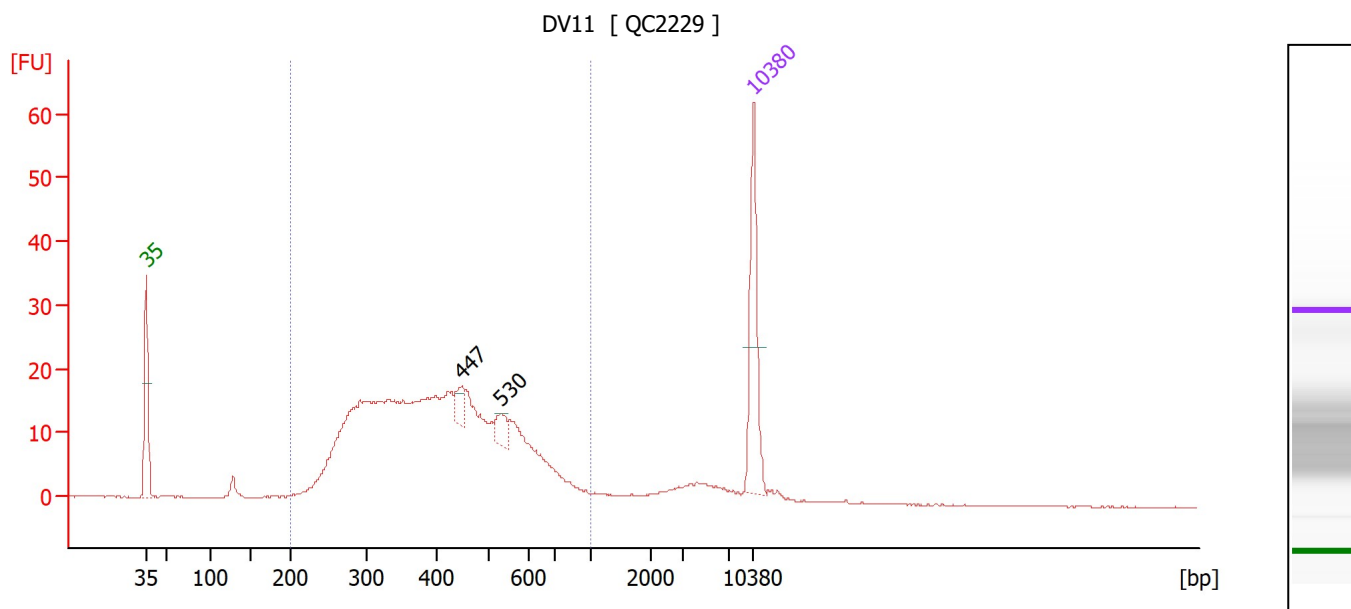Electropherogram Summary Continued ...**Overall Results for sample 10 : DV11**

Number of peaks found: 2 Corr. Area 1: 455.1
Noise: 0.1

Peak table for sample 10 : DV11

Peak	Size [bp]	Conc. [pg/μl]	Molarity [pmol/l]	Observations
1	35	125.00	5,411.3	Lower Marker
2	447	16.36	55.4	
3	530	18.39	52.5	
4	10,380	75.00	10.9	Upper Marker


Region table for sample 10 : DV11

From [bp]	To [bp]	Corr. Area	% of Total	Average Size [bp]	Size distribution in CV [%]	Conc. [pg/μl]	Molarity [pmol/l]	Color
200	1,000	455.1	91	422	30.6	1,033.49	4,184.5	Blue

Assay Class: High Sensitivity DNA Assay
 Data Path: C:\...bioanalyzer\2100 expert\data\2025-06-18\QC2229_DV8-DV12.xad

Created: 6/18/2025 1:06:43 PM
 Modified: 6/18/2025 2:40:02 PM

Electropherogram Summary Continued ...

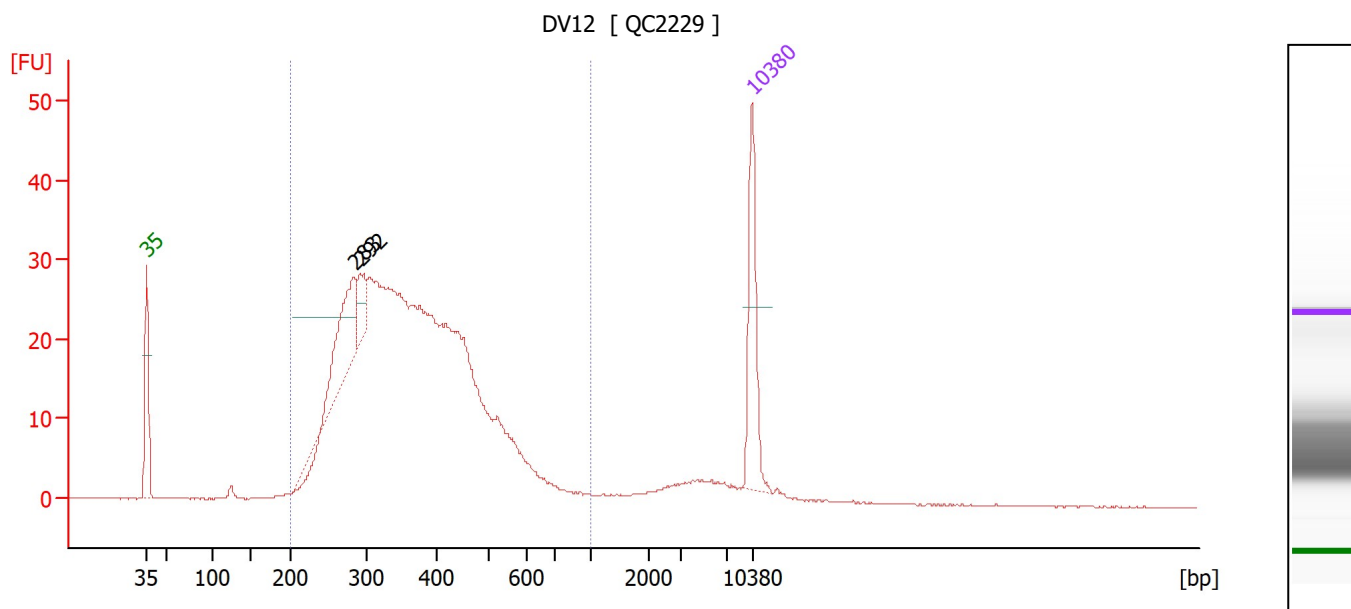

Overall Results for sample 11 : DV12

Number of peaks found: 2 Corr. Area 1: 637.5
 Noise: 0.1

Peak table for sample 11 : DV12


Peak	Size [bp]	Conc. [pg/μl]	Molarity [pmol/l]	Observations
1	35	125.00	5,411.3	Lower Marker
2	283	118.45	633.9	
3	292	34.74	180.1	
4	10,380	75.00	10.9	Upper Marker

Region table for sample 11 : DV12

From [bp]	To [bp]	Corr. Area	% of Total	Average Size [bp]	Size distribution in CV [%]	Conc. [pg/μl]	Molarity [pmol/l]	Color
200	1,000	637.5	94	376	27.9	1,725.78	7,629.1	Blue

Assay Class: High Sensitivity DNA Assay
Data Path: C:\...bioanalyzer\2100 expert\data\2025-06-18\QC2229_DV8-DV12.xad

Created: 6/18/2025 1:06:43 PM
Modified: 6/18/2025 2:40:02 PM

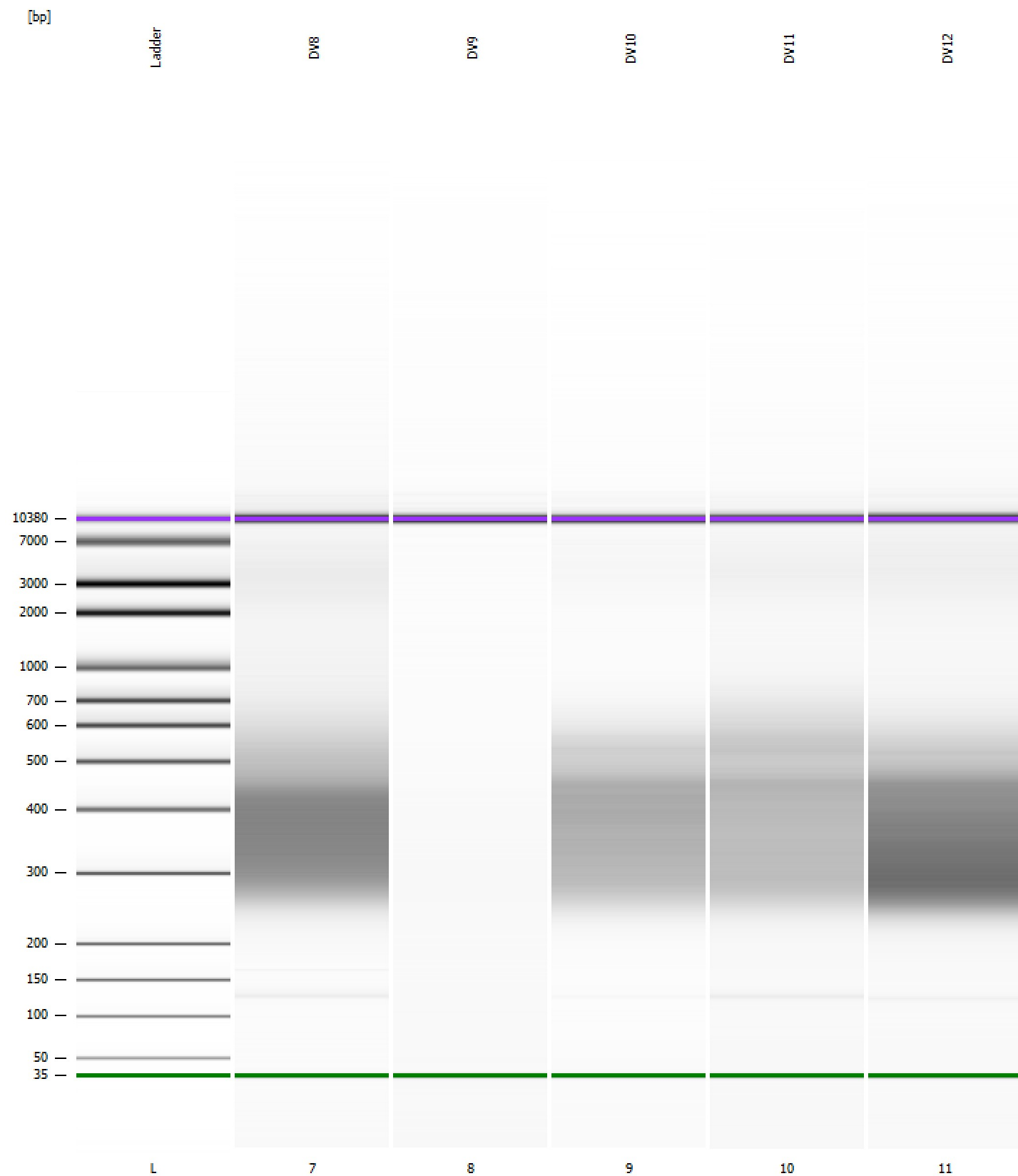Gel Image

Assay Class: High Sensitivity DNA Assay
 Data Path: C:\...bioanalyzer\2100 expert\data\2025-06-18\QC2229_DV8-DV12.xad

Created: 6/18/2025 1:06:43 PM
 Modified: 6/18/2025 2:40:02 PM

Run Logbook

Description	Number	Source	Category	Sub Category	Time	Time Zone	User	Host
Run ended on port 1 (Number of wells acquired: 12)		Instrument	Run		6/18/2025 1:48:00 PM	(GMT --07:00) Pacific Standard Time	sbsuser	DESKTOP-4UNV VOB
Run started on port 1 (File: C:\Program Files (x86)\Agilent\2100 bioanalyzer\2100 expert\data\2025-06-18\Bioanalyzer1_High Sensitivity DNA Assay_DE34903152_2025-06-18_001.xad)		Instrument	Run		6/18/2025 1:06:45 PM	(GMT --07:00) Pacific Standard Time	sbsuser	DESKTOP-4UNV VOB
Product Number : G2938C		Instrument	Run		6/18/2025 1:06:45 PM	(GMT --07:00) Pacific Standard Time	sbsuser	DESKTOP-4UNV VOB
Name :		Instrument	Run		6/18/2025 1:06:45 PM	(GMT --07:00) Pacific Standard Time	sbsuser	DESKTOP-4UNV VOB
Vendor : Agilent Technologies		Instrument	Run		6/18/2025 1:06:45 PM	(GMT --07:00) Pacific Standard Time	sbsuser	DESKTOP-4UNV VOB
Serial# : DE34903152		Instrument	Run		6/18/2025 1:06:45 PM	(GMT --07:00) Pacific Standard Time	sbsuser	DESKTOP-4UNV VOB
Firmware : C.01.069		Instrument	Run		6/18/2025 1:06:45 PM	(GMT --07:00) Pacific Standard Time	sbsuser	DESKTOP-4UNV VOB
Cartridge : Electrode		Instrument	Run		6/18/2025 1:06:45 PM	(GMT --07:00) Pacific Standard Time	sbsuser	DESKTOP-4UNV VOB