


Assay Class: High Sensitivity DNA Assay  
Data Path: C:\...Settings\Bioanalyzer\2015-08-06\2015-08-06\_002\_MiSeq155.xad

Created: 8/6/2015 3:53:09 PM  
Modified: 8/6/2015 5:11:05 PM

**Electrophoresis File Run Summary**

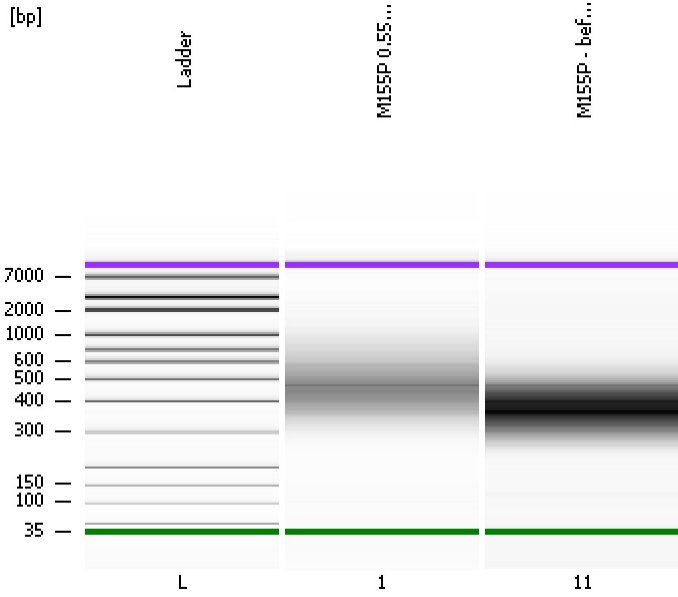

Instrument Information:

Instrument Name: DE13701086      Firmware: C.01.069  
Serial#: DE13701086      Type: G2938B

Assay Information:



Assay Origin Path: C:\Program Files\Agilent\2100 bioanalyzer\2100 expert\assays\dsDNA\High Sensitivity DNA.xsy  
Assay Class: High Sensitivity DNA Assay  
Version: 1.03  
Assay Comments: Copyright © 2003-2010 Agilent Technologies

Chip Information:

Chip Lot #:  
Reagent Kit Lot #:  
Chip Comments:

**M155P 0.55x size selection**

**M155P - before size selection**

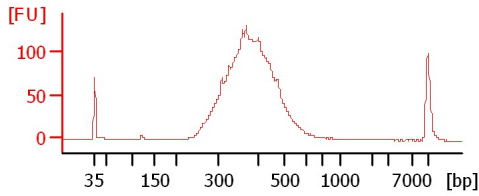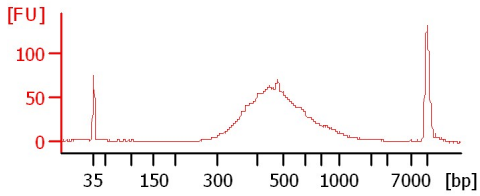

Assay Class: High Sensitivity DNA Assay  
Data Path: C:\...Settings\Bioanalyzer\2015-08-06\2015-08-06\_002\_MiSeq155.xad

Created: 8/6/2015 3:53:09 PM  
Modified: 8/6/2015 5:11:05 PM

**Electrophoresis File Run Summary (Chip Summary)**

Sample Name	Sample Comment	Rest. Digest	Status	Observation	Result Label	Result Color
M155P 0.55x size selection		<input type="checkbox"/>	✓			
M155P - before size selection		<input type="checkbox"/>	✓			
Ladder		<input type="checkbox"/>	✓			
<b>Chip Lot #</b>				<b>Reagent Kit Lot #</b>		

**Chip Comments :**

Assay Class: High Sensitivity DNA Assay  
Data Path: C:\...Settings\Bioanalyzer\2015-08-06\2015-08-06\_002\_MiSeq155.xad

Created: 8/6/2015 3:53:09 PM  
Modified: 8/6/2015 5:11:05 PM

**Electrophoresis Assay Details**

**General Analysis Settings**

Number of Available Sample and Ladder Wells (Max.) : 12  
Minimum Visible Range [s] : 32  
Maximum Visible Range [s] : 138  
Start Analysis Time Range [s] : 33  
End Analysis Time Range [s] : 137.5  
Ladder Concentration [pg/μl] : 1950  
Uses Standard Area for Ladder Fragments  
Lower Marker Concentration [pg/μl] : 125  
Upper Marker Concentration [pg/μl] : 75  
Used Upper Marker for Quantitation  
Standard Curve Fit is Point to Point  
Show Data Aligned to Lower and Upper Marker

**Integrator Settings**

Integration Start Time [s] : 33.05  
Integration End Time [s] : 137  
Slope Threshold : 0.8  
Height Threshold [FU] : 5  
Area Threshold : 0.1  
Width Threshold [s] : 0.6  
Baseline Plateau [s] : 0.5

**Filter Settings**

Filter Width [s] : 0.5  
Polynomial Order : 4


**Ladder**

Ladder Peak	Size	Area
1	35	160
2	50	210
3	100	208
4	150	221
5	200	242
6	300	270
7	400	305
8	500	306
9	600	336
10	700	321
11	1000	366
12	2000	413
13	3000	411
14	7000	400
15	10380	214

Assay Class: High Sensitivity DNA Assay  
 Data Path: C:\...Settings\Bioanalyzer\2015-08-06\2015-08-06\_002\_MiSeq155.xad

Created: 8/6/2015 3:53:09 PM  
 Modified: 8/6/2015 5:11:05 PM

**Electropherogram Summary**

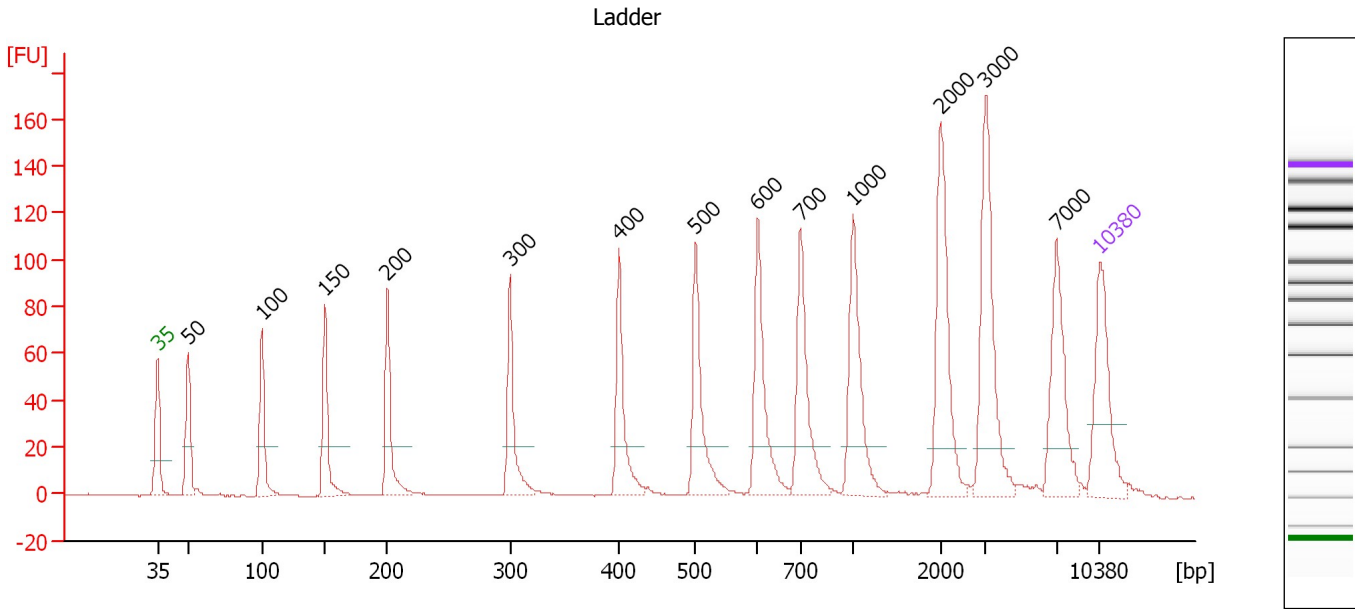

**Overall Results for Ladder**

Noise: 0.2


**Peak table for Ladder**

Peak	Size [bp]	Conc. [pg/μl]	Molarity [pmol/l]	Observations	Aligned Migration Time [s]
1	35	125.00	5,411.3	Lower Marker	43.00
2	50	150.00	4,545.5	Ladder Peak	45.29
3	100	150.00	2,272.7	Ladder Peak	50.74
4	150	150.00	1,515.2	Ladder Peak	55.45
5	200	150.00	1,136.4	Ladder Peak	60.07
6	300	150.00	757.6	Ladder Peak	69.16
7	400	150.00	568.2	Ladder Peak	77.28
8	500	150.00	454.5	Ladder Peak	82.97
9	600	150.00	378.8	Ladder Peak	87.58
10	700	150.00	324.7	Ladder Peak	90.75
11	1,000	150.00	227.3	Ladder Peak	94.67
12	2,000	150.00	113.6	Ladder Peak	101.15
13	3,000	150.00	75.8	Ladder Peak	104.51
14	7,000	150.00	32.5	Ladder Peak	109.78
15	10,380	75.00	10.9	Upper Marker	113.00

Assay Class: High Sensitivity DNA Assay  
 Data Path: C:\...Settings\Bioanalyzer\2015-08-06\2015-08-06\_002\_MiSeq155.xad

Created: 8/6/2015 3:53:09 PM  
 Modified: 8/6/2015 5:11:05 PM

**Electropherogram Summary Continued ...**

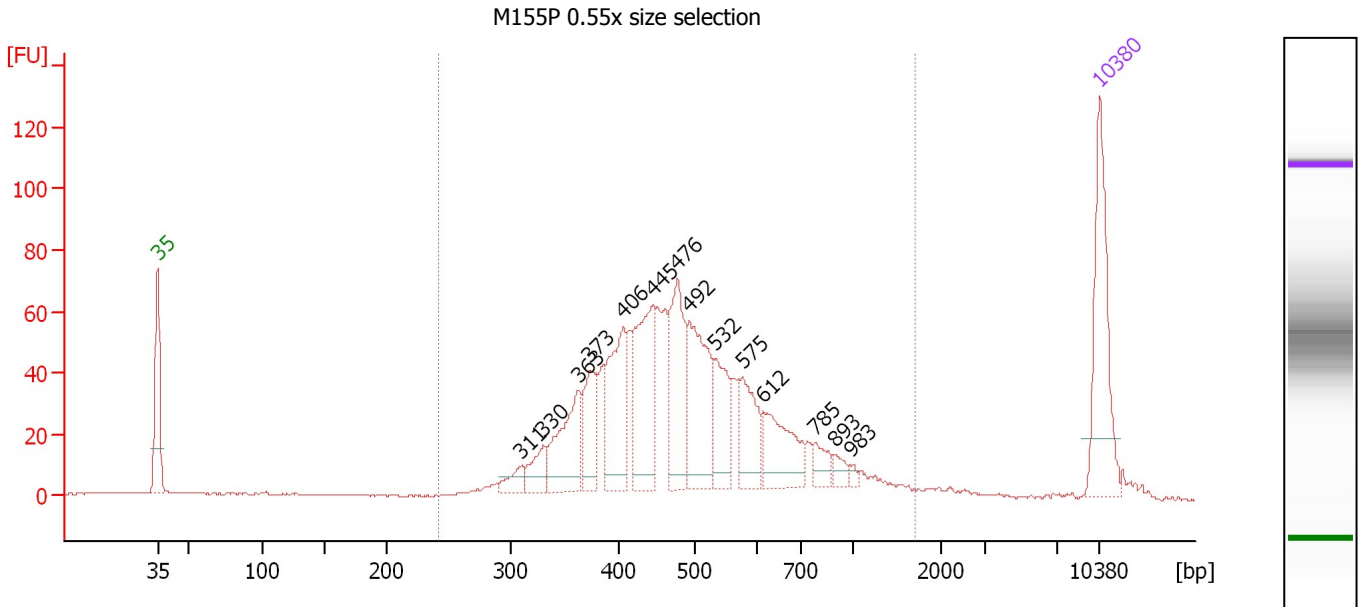

**Overall Results for sample 1 : M155P 0.55x size selection**

Number of peaks found: 14                      Corr. Area 1: 1,167.1  
 Noise: 0.2

**Peak table for sample 1 : M155P 0.55x size selection**

Peak	Size [bp]	Conc. [pg/μl]	Molarity [pmol/l]	Observations	Aligned Migration Time [s]
1	35	125.00	5,411.3	Lower Marker	43.00
2	311	15.58	76.0		70.03
3	330	23.84	109.4		71.63
4	363	72.36	302.3		74.25
5	373	47.61	193.3		75.10
6	406	86.20	321.6		77.63
7	445	105.00	357.6		79.83
8	476	87.28	277.7		81.61
9	492	98.11	302.3		82.50
10	532	51.16	145.6		84.47
11	575	47.17	124.3		86.43
12	612	54.61	135.1		87.98
13	785	14.79	28.5		91.87
14	893	8.96	15.2		93.27
15	983	3.03	4.7		94.45
16	10,380	75.00	10.9	Upper Marker	113.00


**Region table for sample 1 : M155P 0.55x size selection**

From [bp]	Average Size [bp]	To [bp]	Corr. Area	Molarity [pmol/l]	Co % of Ior Total	Size distribution in CV [%]	Conc. [pg/μl]
243	513	1,706	1,167.1	3,154.4	96	36.3	955.09

Assay Class: High Sensitivity DNA Assay  
 Data Path: C:\...Settings\Bioanalyzer\2015-08-06\2015-08-06\_002\_MiSeq155.xad

Created: 8/6/2015 3:53:09 PM  
 Modified: 8/6/2015 5:11:05 PM

**Electropherogram Summary Continued ...**

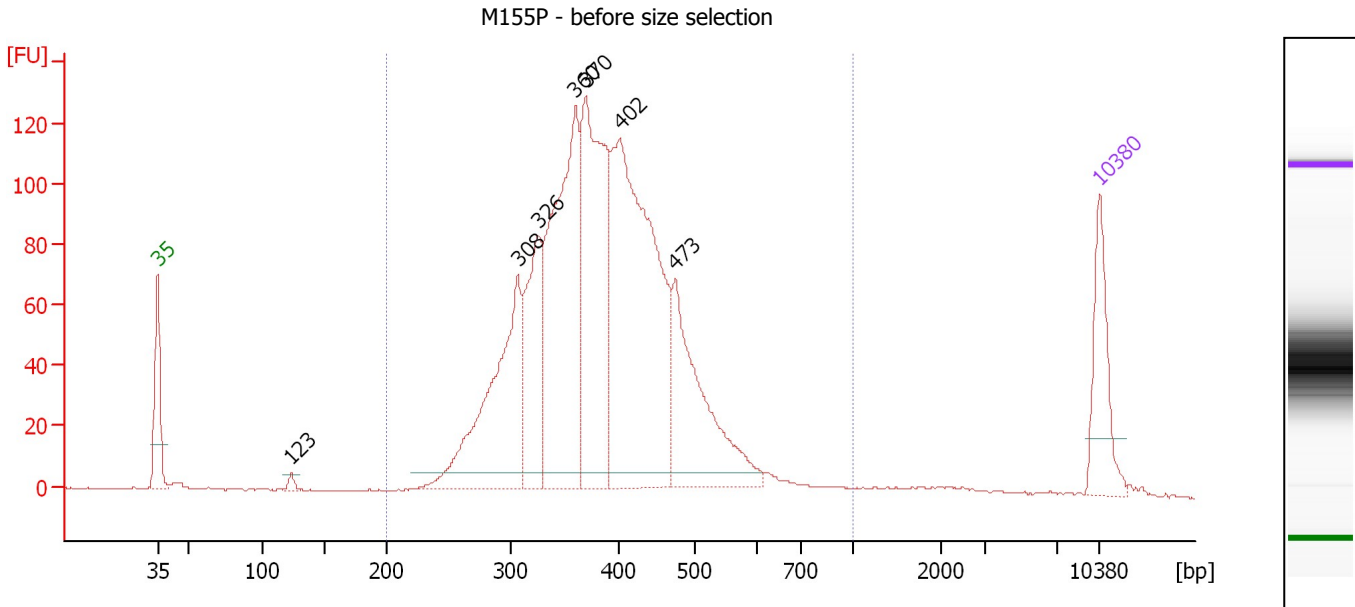

**Overall Results for sample 11 : M155P - before size selection**

Number of peaks found: 7                      Corr. Area 1: 1,872.1  
 Noise: 0.1

**Peak table for sample 11 : M155P - before size selection**

Peak	Size [bp]	Conc. [pg/μl]	Molarity [pmol/l]	Observations	Aligned Migration Time [s]
1	35	125.00	5,411.3	Lower Marker	43.00
2	123	7.62	93.6		52.94
3	308	307.44	1,513.9		69.79
4	326	172.70	803.3		71.25
5	360	408.75	1,720.8		74.02
6	370	367.24	1,505.6		74.81
7	402	561.35	2,117.5		77.37
8	473	228.89	732.6		81.45
9	10,380	75.00	10.9	Upper Marker	113.00

**Region table for sample 11 : M155P - before size selection**

From [bp]	Average Size [bp]	To [bp]	Corr. Area	Molarity [pmol/l]	Co % of Ior Total	Size distribution in CV [%]	Conc. [pg/μl]
200	389	1,000	1,872.1	8,108.0	100	17.9	2,000.00

Assay Class: High Sensitivity DNA Assay  
Data Path: C:\...Settings\Bioanalyzer\2015-08-06\2015-08-06\_002\_MiSeq155.xad

Created: 8/6/2015 3:53:09 PM  
Modified: 8/6/2015 5:11:05 PM

**Gel Image**

