


Assay Class: High Sensitivity DNA Assay  
Data Path: C:\...ings\Bioanalyzer\2015-11-19\2015-11-19\_002\_H241P\_4\_5\_17.xad

Created: 11/19/2015 10:44:58 AM  
Modified: 11/19/2015 11:27:44 AM

**Electrophoresis File Run Summary**

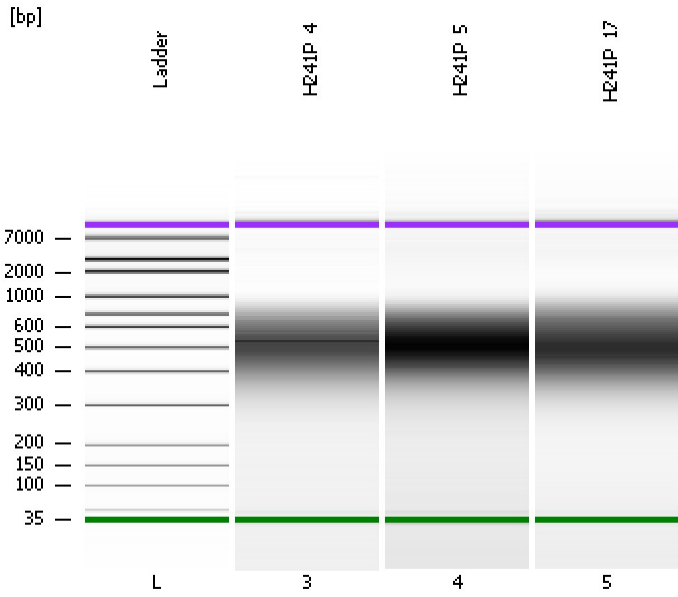

Instrument Information:

Instrument Name: DE13701086      Firmware: C.01.069  
Serial#: DE13701086      Type: G2938B

Assay Information:

Assay Origin Path: C:\Program Files\Agilent\2100 bioanalyzer\2100 expert\assays\dsDNA\High Sensitivity DNA.xsy  
Assay Class: High Sensitivity DNA Assay  
Version: 1.03  
Assay Comments: Copyright © 2003-2010 Agilent Technologies




Chip Information:

Chip Lot #:  
Reagent Kit Lot #:  
Chip Comments:

**H241P\_4**

**H241P\_5**

**H241P\_17**

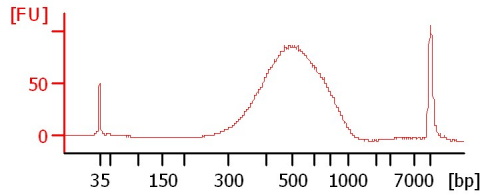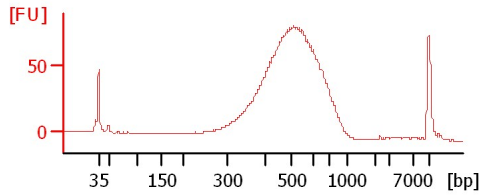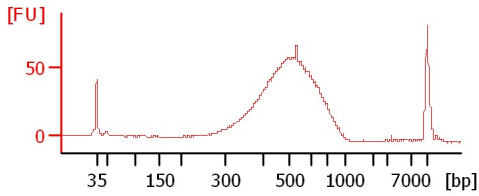

Assay Class: High Sensitivity DNA Assay  
 Data Path: C:\...ings\Bioanalyzer\2015-11-19\2015-11-19\_002\_H241P\_4\_5\_17.xad

Created: 11/19/2015 10:44:58 AM  
 Modified: 11/19/2015 11:27:44 AM

**Electrophoresis File Run Summary (Chip Summary)**

Sample Name	Sample Comment	Rest. Digest	Status	Observation	Result Label	Result Color
H241P_4		<input type="checkbox"/>	✓			
H241P_5		<input type="checkbox"/>	✓			
H241P_17		<input type="checkbox"/>	✓			
Ladder		<input type="checkbox"/>	✓			

**Chip Lot #** **Reagent Kit Lot #**

**Chip Comments :**

Assay Class: High Sensitivity DNA Assay  
Data Path: C:\...ings\Bioanalyzer\2015-11-19\2015-11-19\_002\_H241P\_4\_5\_17.xad

Created: 11/19/2015 10:44:58 AM  
Modified: 11/19/2015 11:27:44 AM

**Electrophoresis Assay Details**

**General Analysis Settings**

Number of Available Sample and Ladder Wells (Max.) : 12  
Minimum Visible Range [s] : 32  
Maximum Visible Range [s] : 138  
Start Analysis Time Range [s] : 33  
End Analysis Time Range [s] : 137.5  
Ladder Concentration [pg/μl] : 1950  
Uses Standard Area for Ladder Fragments  
Lower Marker Concentration [pg/μl] : 125  
Upper Marker Concentration [pg/μl] : 75  
Used Upper Marker for Quantitation  
Standard Curve Fit is Point to Point  
Show Data Aligned to Lower and Upper Marker

**Integrator Settings**

Integration Start Time [s] : 33.05  
Integration End Time [s] : 137  
Slope Threshold : 0.8  
Height Threshold [FU] : 5  
Area Threshold : 0.1  
Width Threshold [s] : 0.6  
Baseline Plateau [s] : 0.5

**Filter Settings**

Filter Width [s] : 0.5  
Polynomial Order : 4


**Ladder**

Ladder Peak	Size	Area
1	35	160
2	50	210
3	100	208
4	150	221
5	200	242
6	300	270
7	400	305
8	500	306
9	600	336
10	700	321
11	1000	366
12	2000	413
13	3000	411
14	7000	400
15	10380	214

Assay Class: High Sensitivity DNA Assay  
 Data Path: C:\...ings\Bioanalyzer\2015-11-19\2015-11-19\_002\_H241P\_4\_5\_17.xad

Created: 11/19/2015 10:44:58 AM  
 Modified: 11/19/2015 11:27:44 AM

**Electropherogram Summary**

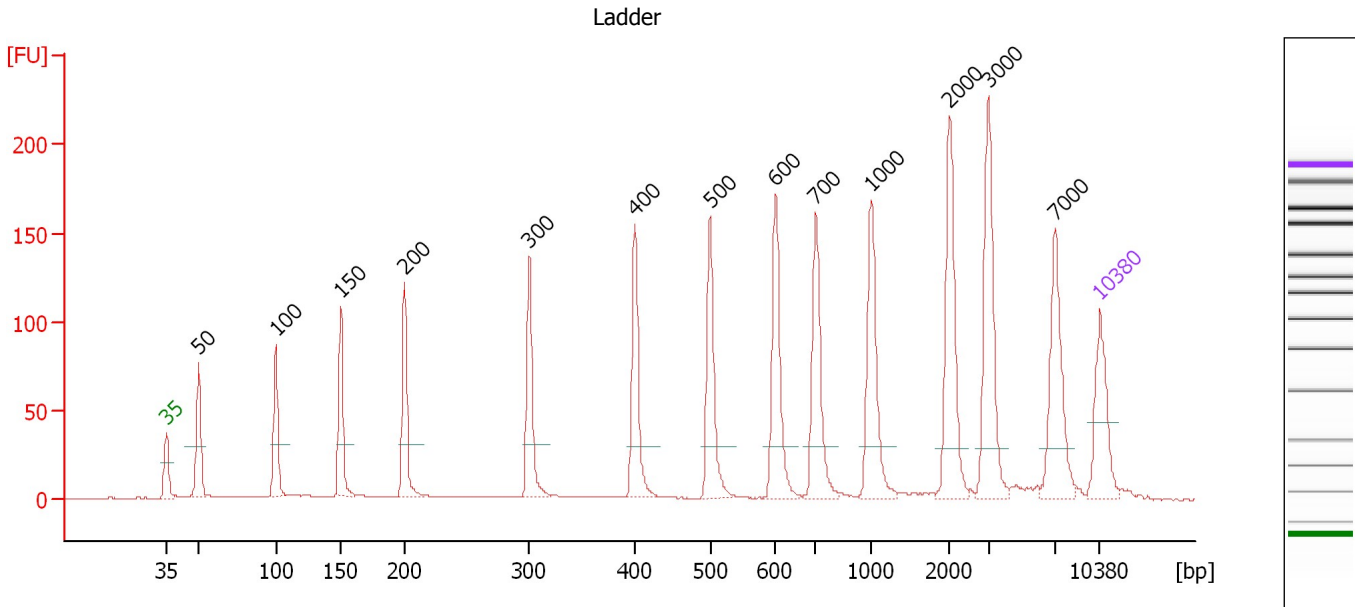

**Overall Results for Ladder**

Noise: 0.2


**Peak table for Ladder**

Peak	Size [bp]	Conc. [pg/μl]	Molarity [pmol/l]	Observations	Aligned Migration Time [s]
1	35	125.00	5,411.3	Lower Marker	43.00
2	50	150.00	4,545.5	Ladder Peak	45.43
3	100	150.00	2,272.7	Ladder Peak	51.20
4	150	150.00	1,515.2	Ladder Peak	56.11
5	200	150.00	1,136.4	Ladder Peak	60.88
6	300	150.00	757.6	Ladder Peak	70.21
7	400	150.00	568.2	Ladder Peak	78.14
8	500	150.00	454.5	Ladder Peak	83.77
9	600	150.00	378.8	Ladder Peak	88.68
10	700	150.00	324.7	Ladder Peak	91.69
11	1,000	150.00	227.3	Ladder Peak	95.84
12	2,000	150.00	113.6	Ladder Peak	101.74
13	3,000	150.00	75.8	Ladder Peak	104.62
14	7,000	150.00	32.5	Ladder Peak	109.62
15	10,380	75.00	10.9	Upper Marker	113.00

Assay Class: High Sensitivity DNA Assay  
 Data Path: C:\...ings\Bioanalyzer\2015-11-19\2015-11-19\_002\_H241P\_4\_5\_17.xad

Created: 11/19/2015 10:44:58 AM  
 Modified: 11/19/2015 11:27:44 AM

**Electropherogram Summary Continued ...**

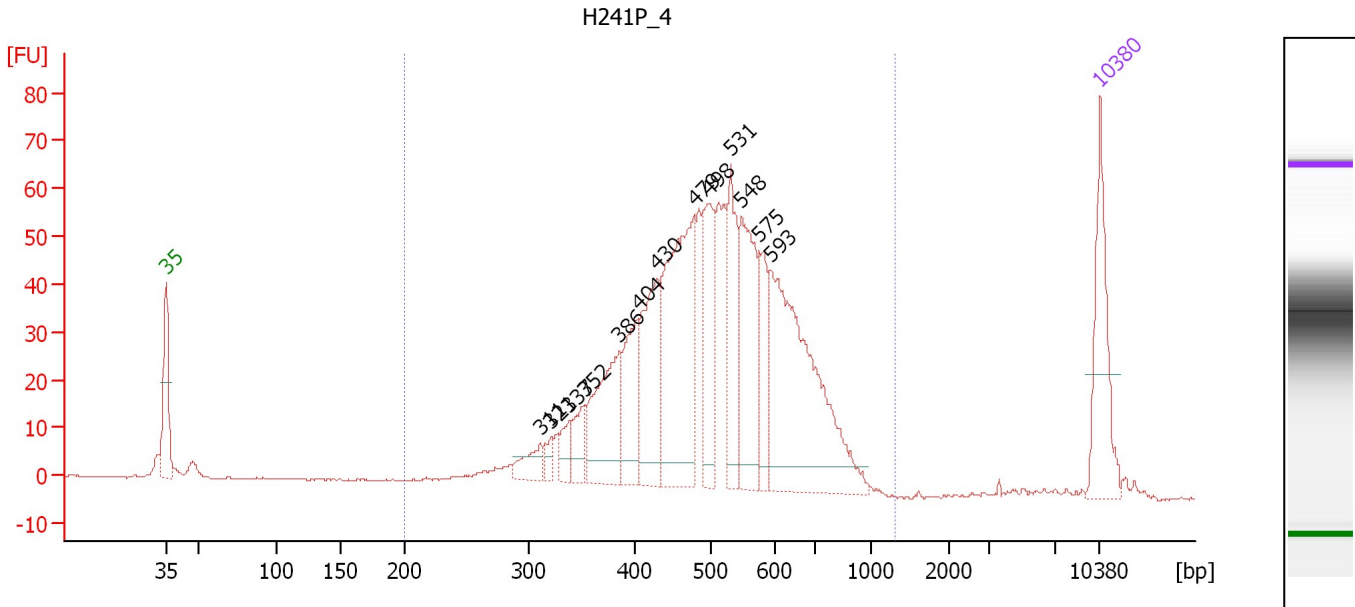

**Overall Results for sample 3 : H241P 4**

Number of peaks found: 13                      Corr. Area 1: 1,062.6  
 Noise: 0.2

**Peak table for sample 3 : H241P 4**

Peak	Size [bp]	Conc. [pg/μl]	Molarity [pmol/l]	Observations	Aligned Migration Time [s]
1	35	125.00	5,411.3	Lower Marker	43.00
2	311	27.78	135.4		71.06
3	321	12.99	61.3		71.89
4	337	23.88	107.3		73.16
5	352	28.92	124.6		74.30
6	386	111.59	438.5		77.00
7	404	77.94	292.2		78.37
8	430	112.22	395.4		79.83
9	479	232.05	734.5		82.56
10	498	98.08	298.3		83.66
11	531	90.63	258.5		85.30
12	548	125.71	347.6		86.12
13	575	58.52	154.2		87.45
14	593	267.26	683.3		88.31
15	10,380	75.00	10.9	Upper Marker	113.00


**Region table for sample 3 : H241P 4**

From [bp]	To [bp]	Average Size [bp]	Corr. Area	Molarity [pmol/l]	Co Conc. [pg/μl]	% of Total	Size distribution in CV [%]
200	1,307	507	1,062.6	4,850.2	1,493.85	96	23.7

Assay Class: High Sensitivity DNA Assay  
 Data Path: C:\...ings\Bioanalyzer\2015-11-19\2015-11-19\_002\_H241P\_4\_5\_17.xad

Created: 11/19/2015 10:44:58 AM  
 Modified: 11/19/2015 11:27:44 AM

**Electropherogram Summary Continued ...**

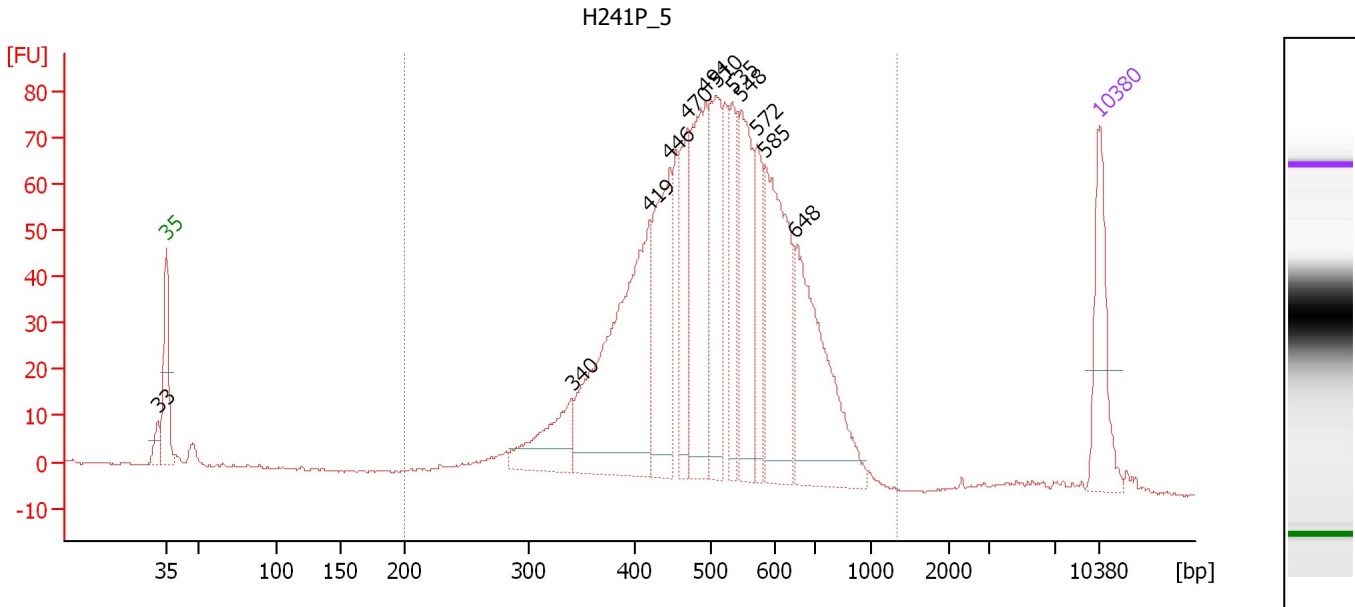

**Overall Results for sample 4 : H241P 5**

Number of peaks found: 12      Corr. Area 1: 1,441.8  
 Noise: 0.2

**Peak table for sample 4 : H241P 5**

Peak	Size [bp]	Conc. [pg/μl]	Molarity [pmol/l]	Observations	Aligned Migration Time [s]
1	33	0.00	0.0		42.35
2	35	125.00	5,411.3	Lower Marker	43.00
3	340	91.54	407.6		73.40
4	419	360.13	1,302.1		79.21
5	446	176.34	598.7		80.74
6	470	97.41	313.9		82.09
7	494	205.46	629.9		83.44
8	510	151.46	449.7		84.27
9	535	81.11	229.7		85.48
10	548	145.34	401.7		86.13
11	572	77.99	206.6		87.30
12	585	199.88	517.5		87.95
13	648	216.95	507.1		90.13
14	10,380	75.00	10.9	Upper Marker	113.00


**Region table for sample 4 : H241P 5**

From [bp]	To [bp]	Average Size [bp]	Corr. Area	Molarity [pmol/l]	Co Conc. lor [pg/μl]	% of Total	Size distribution in CV [%]
200	1,328	513	1,441.8	6,558.4	2,049.35	97	23.5

Assay Class: High Sensitivity DNA Assay  
 Data Path: C:\...ings\Bioanalyzer\2015-11-19\2015-11-19\_002\_H241P\_4\_5\_17.xad

Created: 11/19/2015 10:44:58 AM  
 Modified: 11/19/2015 11:27:44 AM

**Electropherogram Summary Continued ...**

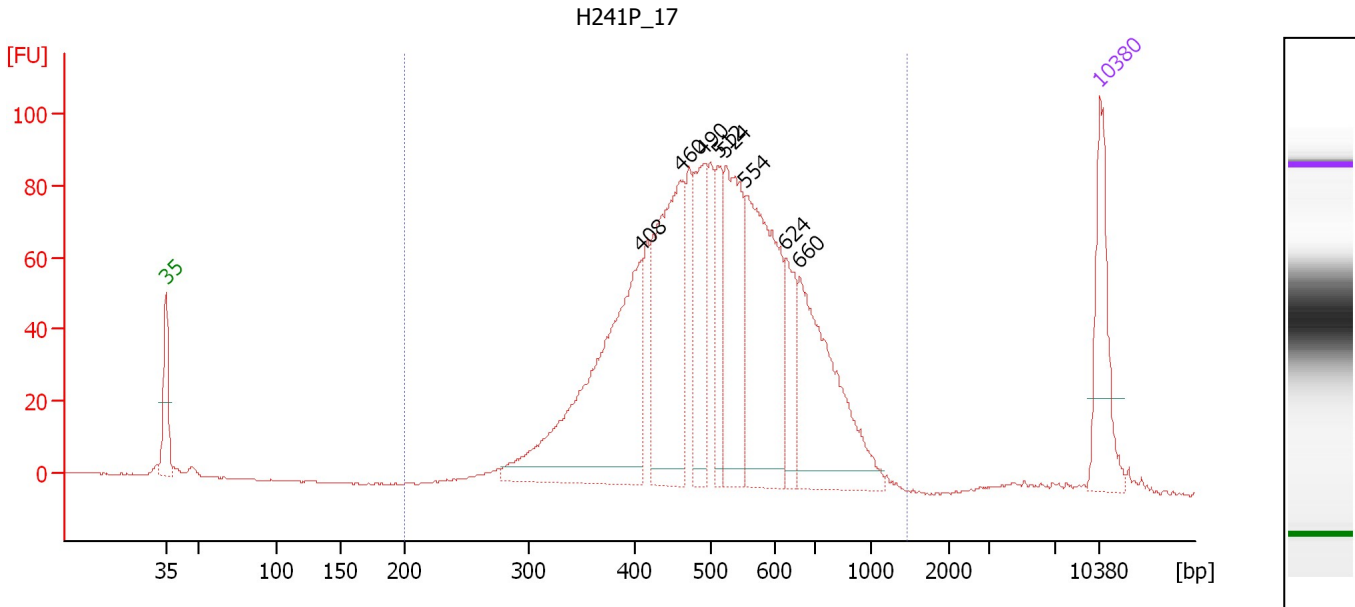

**Overall Results for sample 5 : H241P 17**

Number of peaks found: 8      Corr. Area 1: 1,701.3  
 Noise: 0.2

**Peak table for sample 5 : H241P 17**

Peak	Size [bp]	Conc. [pg/μl]	Molarity [pmol/l]	Observations	Aligned Migration Time [s]
1	35	125.00	5,411.3	Lower Marker	43.00
2	408	356.29	1,324.3		78.57
3	460	257.02	845.9		81.53
4	490	118.41	365.8		83.23
5	512	65.19	192.9		84.36
6	524	157.19	454.9		84.92
7	554	246.30	673.2		86.43
8	624	63.86	155.1		89.40
9	660	204.55	469.7		90.48
10	10,380	75.00	10.9	Upper Marker	113.00


**Region table for sample 5 : H241P 17**

From [bp]	To [bp]	Average Size [bp]	Corr. Area	Molarity [pmol/l]	Co Conc. [pg/μl]	% of Total	Size distribution in CV [%]
200	1,467	522	1,701.3	5,336.8	1,685.49	98	26.8

Assay Class: High Sensitivity DNA Assay  
Data Path: C:\...ings\Bioanalyzer\2015-11-19\2015-11-19\_002\_H241P\_4\_5\_17.xad

Created: 11/19/2015 10:44:58 AM  
Modified: 11/19/2015 11:27:44 AM

**Gel Image**

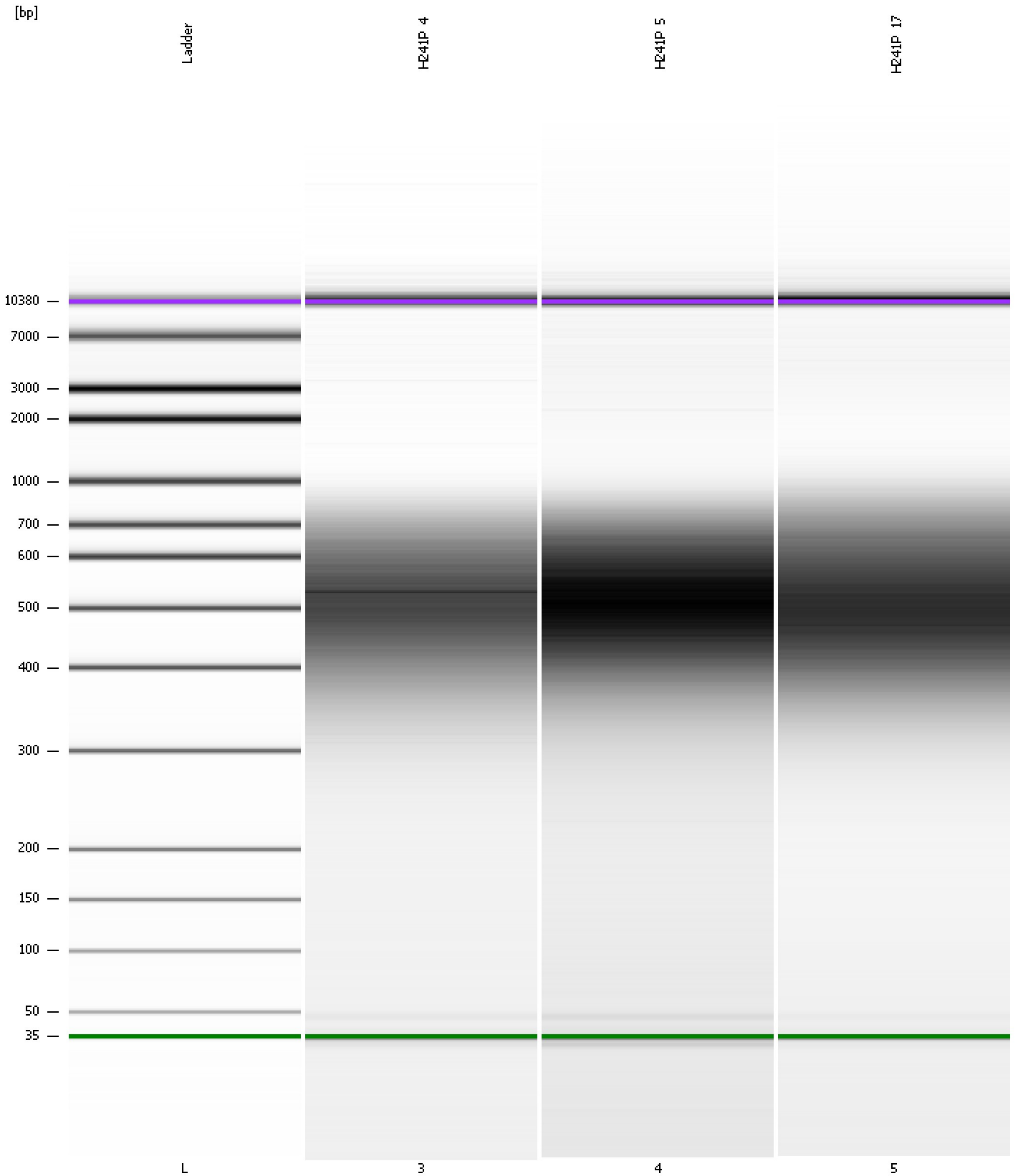




Assay Class: High Sensitivity DNA Assay  
Data Path: C:\...ings\Bioanalyzer\2015-11-19\2015-11-19\_002\_H241P\_4\_5\_17.xad

Created: 11/19/2015 10:44:58 AM  
Modified: 11/19/2015 11:27:44 AM

**Curves**

**Standard Curve**

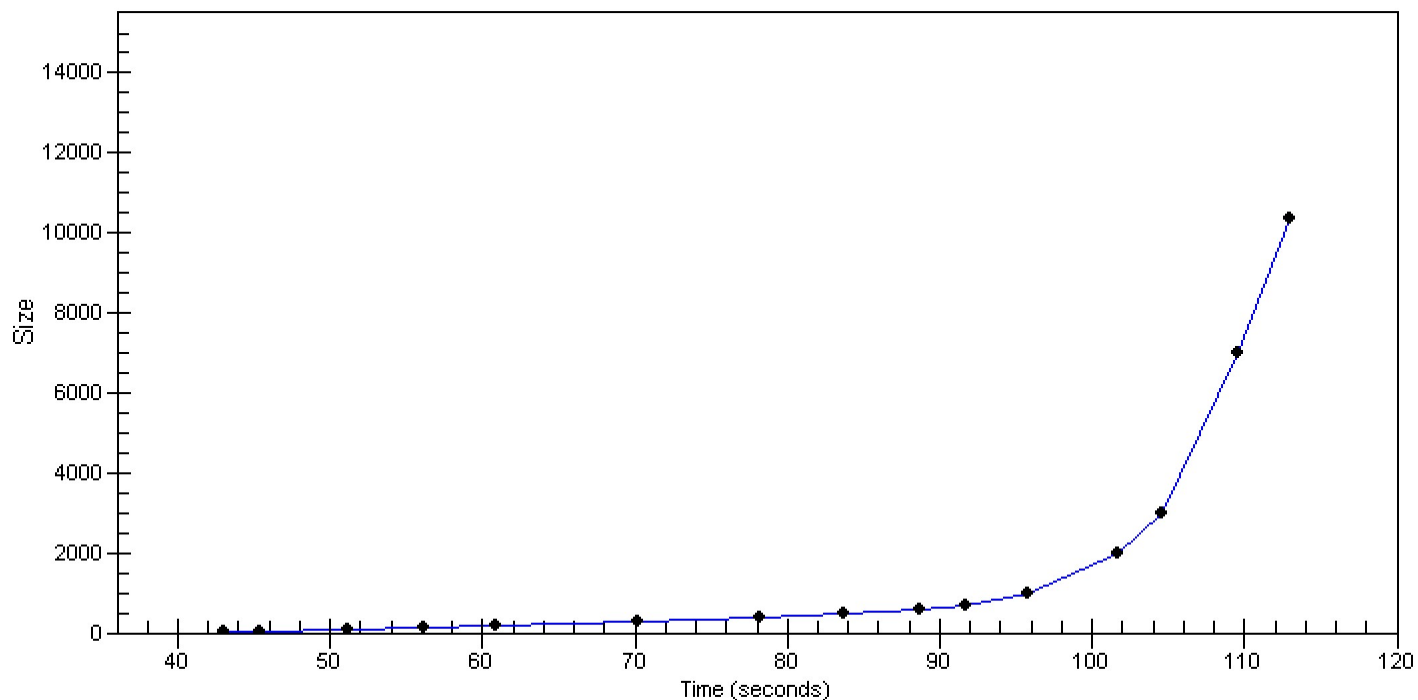

Assay Class: High Sensitivity DNA Assay  
 Data Path: C:\...ings\Bioanalyzer\2015-11-19\2015-11-19\_002\_H241P\_4\_5\_17.xad

Created: 11/19/2015 10:44:58 AM  
 Modified: 11/19/2015 11:27:44 AM

**Run Logbook**

Description	Number	Source	Category	Sub Category	Time	Time Zone	User	Host
Run ended on port 1 (Number of wells acquired: 12)		Instrument	Run		11/19/2015 11:25:22 AM	(GMT --08:00) Pacific Standard Time	UC Davis	D8XSMGH1
Run started on port 1 (File: C:\Documents and Settings\Bioanalyzer\2015-11-19\2015-11-19_002.xad)		Instrument	Run		11/19/2015 10:44:58 AM	(GMT --08:00) Pacific Standard Time	UC Davis	D8XSMGH1
Product Number : G2938B		Instrument	Run		11/19/2015 10:44:58 AM	(GMT --08:00) Pacific Standard Time	UC Davis	D8XSMGH1
Name :		Instrument	Run		11/19/2015 10:44:58 AM	(GMT --08:00) Pacific Standard Time	UC Davis	D8XSMGH1
Vendor : Agilent Technologies		Instrument	Run		11/19/2015 10:44:58 AM	(GMT --08:00) Pacific Standard Time	UC Davis	D8XSMGH1
Serial# : DE13701086		Instrument	Run		11/19/2015 10:44:58 AM	(GMT --08:00) Pacific Standard Time	UC Davis	D8XSMGH1
Firmware : C.01.069		Instrument	Run		11/19/2015 10:44:58 AM	(GMT --08:00) Pacific Standard Time	UC Davis	D8XSMGH1
Cartridge : Electrode		Instrument	Run		11/19/2015 10:44:58 AM	(GMT --08:00) Pacific Standard Time	UC Davis	D8XSMGH1