


Assay Class: High Sensitivity DNA Assay
 Data Path: C:\...2-09\2015-12-09_002_CT516S_CT5ITS3_CT5ITS2_LigationRxns.xad

Created: 12/9/2015 11:00:29 AM
 Modified: 12/9/2015 11:47:35 AM

Electrophoresis File Run Summary

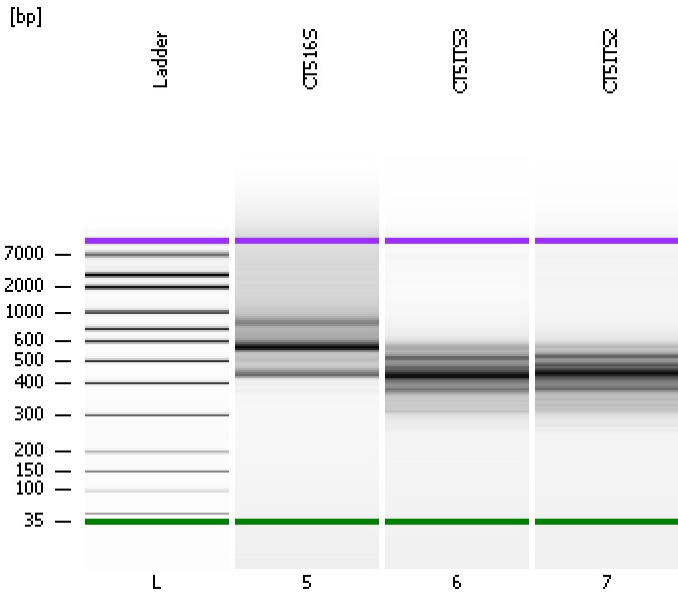

Instrument Information:

Instrument Name: DE13701086 Firmware: C.01.069
 Serial#: DE13701086 Type: G2938B

Assay Information:

Assay Origin Path: C:\Program Files\Agilent\2100 bioanalyzer\2100 expert\assays\dsDNA\High Sensitivity DNA.xsy
 Assay Class: High Sensitivity DNA Assay
 Version: 1.03
 Assay Comments: Copyright © 2003-2010 Agilent Technologies




Chip Information:

Chip Lot #:
 Reagent Kit Lot #:
 Chip Comments:

CT516S

CT5ITS3

CT5ITS2

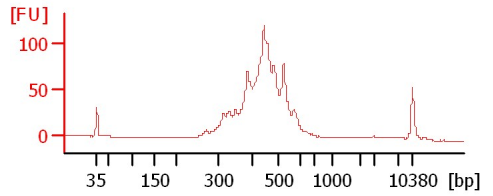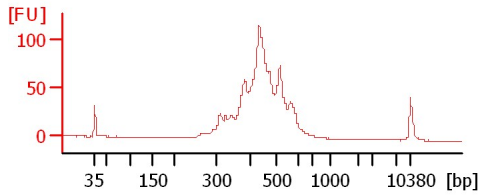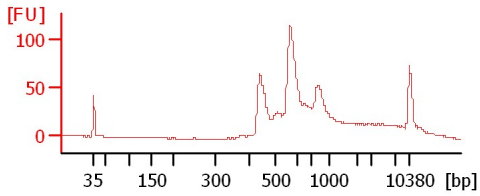

Assay Class: High Sensitivity DNA Assay
 Data Path: C:\...2-09\2015-12-09_002_CT516S_CT5ITS3_CT5ITS2_LigationRxns.xad

Created: 12/9/2015 11:00:29 AM
 Modified: 12/9/2015 11:47:35 AM

Electrophoresis File Run Summary (Chip Summary)

Sample Name	Sample Comment	Rest. Digest	Status	Observation	Result Label	Result Color
CT516S		<input type="checkbox"/>				
CT5ITS3		<input type="checkbox"/>				
CT5ITS2		<input type="checkbox"/>				
Ladder		<input type="checkbox"/>				

Chip Lot # **Reagent Kit Lot #**

Chip Comments :

Assay Class: High Sensitivity DNA Assay
 Data Path: C:\...2-09\2015-12-09_002_CT516S_CT5ITS3_CT5ITS2_LigationRxns.xad

Created: 12/9/2015 11:00:29 AM
 Modified: 12/9/2015 11:47:35 AM

Electrophoresis Assay Details

General Analysis Settings

Number of Available Sample and Ladder Wells (Max.) : 12
 Minimum Visible Range [s] : 32
 Maximum Visible Range [s] : 138
 Start Analysis Time Range [s] : 33
 End Analysis Time Range [s] : 137.5
 Ladder Concentration [pg/μl] : 1950
 Uses Standard Area for Ladder Fragments
 Lower Marker Concentration [pg/μl] : 125
 Upper Marker Concentration [pg/μl] : 75
 Used Upper Marker for Quantitation
 Standard Curve Fit is Point to Point
 Show Data Aligned to Lower and Upper Marker

Integrator Settings

Integration Start Time [s] : 33.05
 Integration End Time [s] : 137
 Slope Threshold : 0.8
 Height Threshold [FU] : 5
 Area Threshold : 0.1
 Width Threshold [s] : 0.6
 Baseline Plateau [s] : 0.5

Filter Settings

Filter Width [s] : 0.5
 Polynomial Order : 4


Ladder

Ladder Peak	Size	Area
1	35	160
2	50	210
3	100	208
4	150	221
5	200	242
6	300	270
7	400	305
8	500	306
9	600	336
10	700	321
11	1000	366
12	2000	413
13	3000	411
14	7000	400
15	10380	214

Assay Class: High Sensitivity DNA Assay
 Data Path: C:\...2-09\2015-12-09_002_CT516S_CT5ITS3_CT5ITS2_LigationRxns.xad

Created: 12/9/2015 11:00:29 AM
 Modified: 12/9/2015 11:47:35 AM

Electropherogram Summary

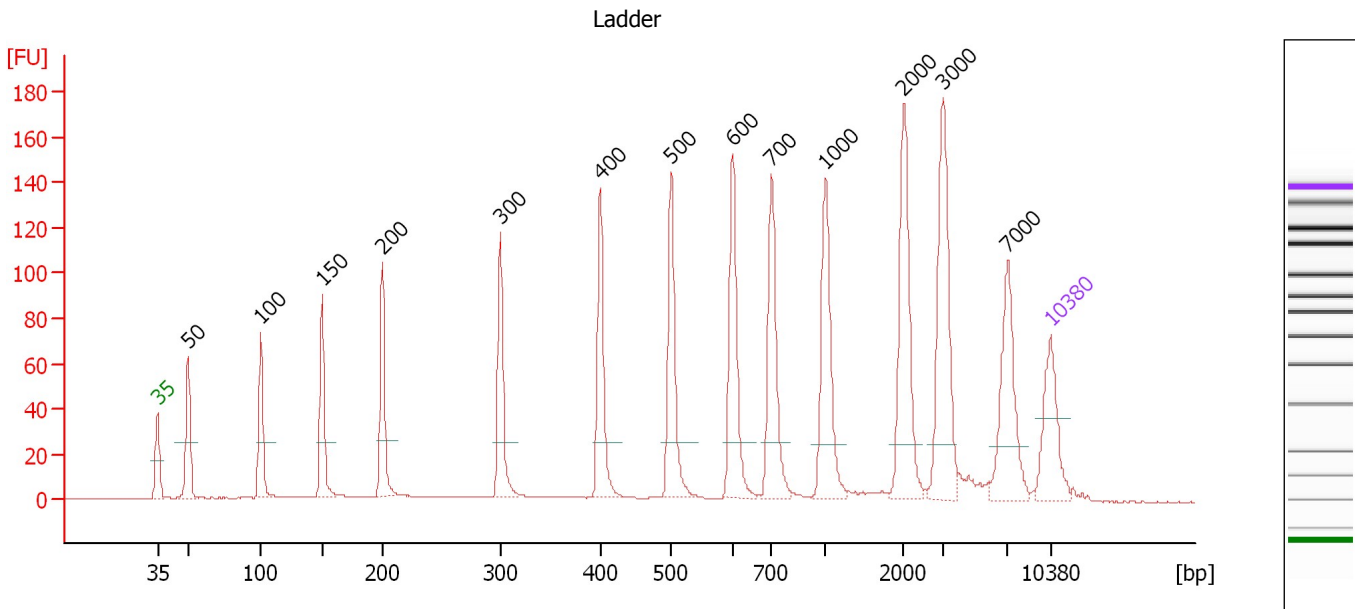

Overall Results for Ladder

Noise: 0.2


Peak table for Ladder

Peak	Size [bp]	Conc. [pg/μl]	Molarity [pmol/l]	Observations	Aligned Migration Time [s]
1	35	125.00	5,411.3	Lower Marker	43.00
2	50	150.00	4,545.5	Ladder Peak	45.38
3	100	150.00	2,272.7	Ladder Peak	51.10
4	150	150.00	1,515.2	Ladder Peak	55.91
5	200	150.00	1,136.4	Ladder Peak	60.63
6	300	150.00	757.6	Ladder Peak	69.86
7	400	150.00	568.2	Ladder Peak	77.69
8	500	150.00	454.5	Ladder Peak	83.26
9	600	150.00	378.8	Ladder Peak	88.08
10	700	150.00	324.7	Ladder Peak	91.09
11	1,000	150.00	227.3	Ladder Peak	95.37
12	2,000	150.00	113.6	Ladder Peak	101.48
13	3,000	150.00	75.8	Ladder Peak	104.54
14	7,000	150.00	32.5	Ladder Peak	109.63
15	10,380	75.00	10.9	Upper Marker	113.00

Assay Class: High Sensitivity DNA Assay
 Data Path: C:\...2-09\2015-12-09_002_CT516S_CT5ITS3_CT5ITS2_LigationRxns.xad

Created: 12/9/2015 11:00:29 AM
 Modified: 12/9/2015 11:47:35 AM

Electropherogram Summary Continued ...

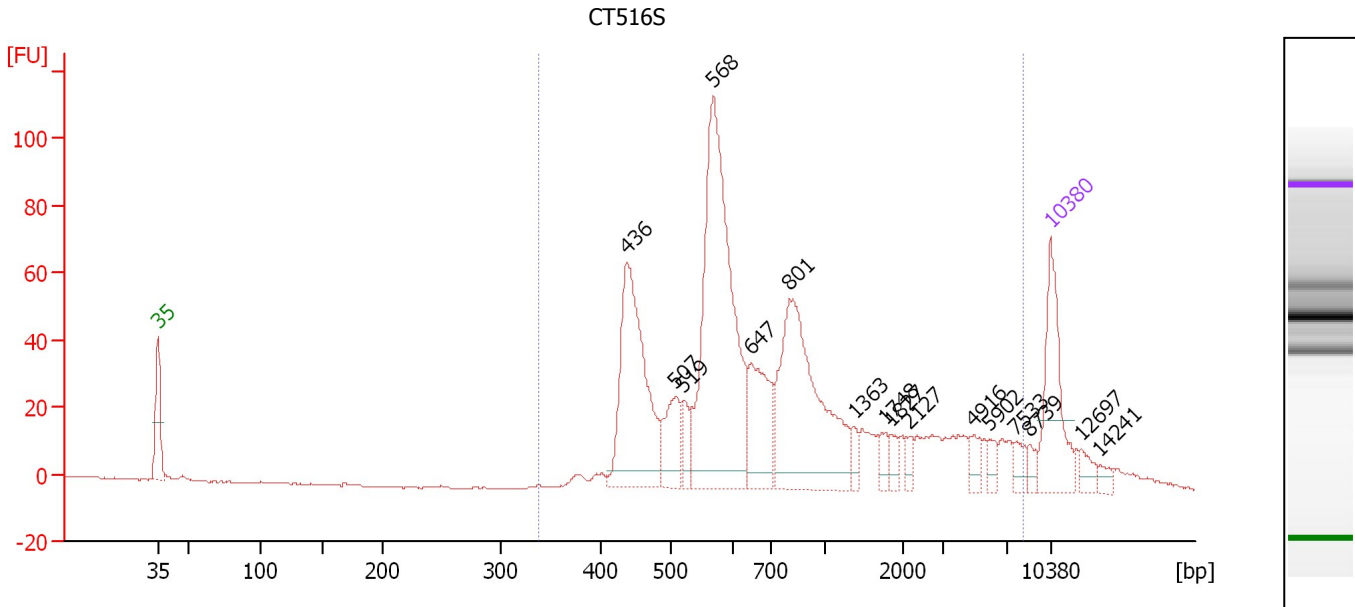

Overall Results for sample 5 : CT516S

Number of peaks found: 16 Corr. Area 1: 1,123.3
 Noise: 0.2

Peak table for sample 5 : CT516S

Peak	Size [bp]	Conc. [pg/μl]	Molarity [pmol/l]	Observations	Aligned Migration Time [s]
1	35	125.00	5,411.3	Lower Marker	43.00
2	436	219.34	761.7		79.71
3	507	55.71	166.4		83.61
4	519	21.65	63.2		84.19
5	568	375.11	1,000.1		86.55
6	647	85.77	200.9		89.49
7	801	237.49	449.4		92.53
8	1,363	11.10	12.3		97.58
9	1,748	11.67	10.1		99.94
10	1,827	11.23	9.3		100.43
11	2,127	8.69	6.2		101.87
12	4,916	12.53	3.9		106.98
13	5,902	10.39	2.7		108.23
14	7,533	11.80	2.4		110.16
15	8,739	8.12	1.4		111.36
16	10,380	75.00	10.9	Upper Marker	113.00
17	12,697	0.00	0.0		115.31
18	14,241	0.00	0.0		116.85


Region table for sample 5 : CT516S

From [bp]	To [bp]	Average Size [bp]	Corr. Area	Molarity [pmol/l]	Co Conc. lor [pg/μl]	% of Total	Size distribution in CV [%]
338	8,189	1,292	1,123.3	2,685.3	1,190.59	93	100.0

Assay Class: High Sensitivity DNA Assay
 Data Path: C:\...2-09\2015-12-09_002_CT516S_CT5ITS3_CT5ITS2_LigationRxns.xad

Created: 12/9/2015 11:00:29 AM
 Modified: 12/9/2015 11:47:35 AM

Electropherogram Summary Continued ...

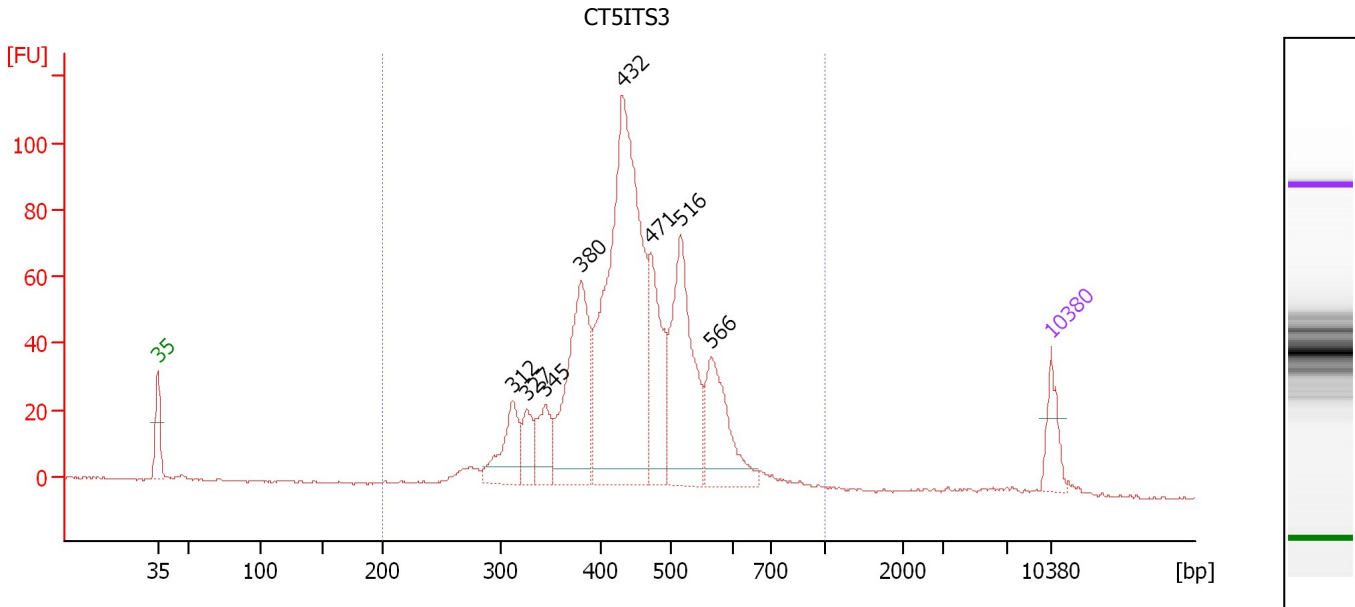

Overall Results for sample 6 : CT5ITS3

Number of peaks found: 8 Corr. Area 1: 1,152.8
 Noise: 0.3

Peak table for sample 6 : CT5ITS3

Peak	Size [bp]	Conc. [pg/μl]	Molarity [pmol/l]	Observations	Aligned Migration Time [s]
1	35	125.00	5,411.3	Lower Marker	43.00
2	312	190.82	926.4		70.80
3	327	96.47	446.9		71.98
4	345	132.03	579.4		73.40
5	380	500.15	1,992.8		76.14
6	432	1,312.14	4,602.4		79.47
7	471	318.10	1,022.3		81.67
8	516	504.61	1,482.5		84.02
9	566	272.05	728.8		86.42
10	10,380	75.00	10.9	Upper Marker	113.00


Region table for sample 6 : CT5ITS3

From [bp]	To [bp]	Average Size [bp]	Corr. Area	Molarity [pmol/l]	Co Conc. [pg/μl]	% of Total	Size distribution in CV [%]
200	1,000	449	1,152.8	12,395.9	3,495.11	98	19.4

Assay Class: High Sensitivity DNA Assay
 Data Path: C:\...2-09\2015-12-09_002_CT516S_CT5ITS3_CT5ITS2_LigationRxns.xad

Created: 12/9/2015 11:00:29 AM
 Modified: 12/9/2015 11:47:35 AM

Electropherogram Summary Continued ...

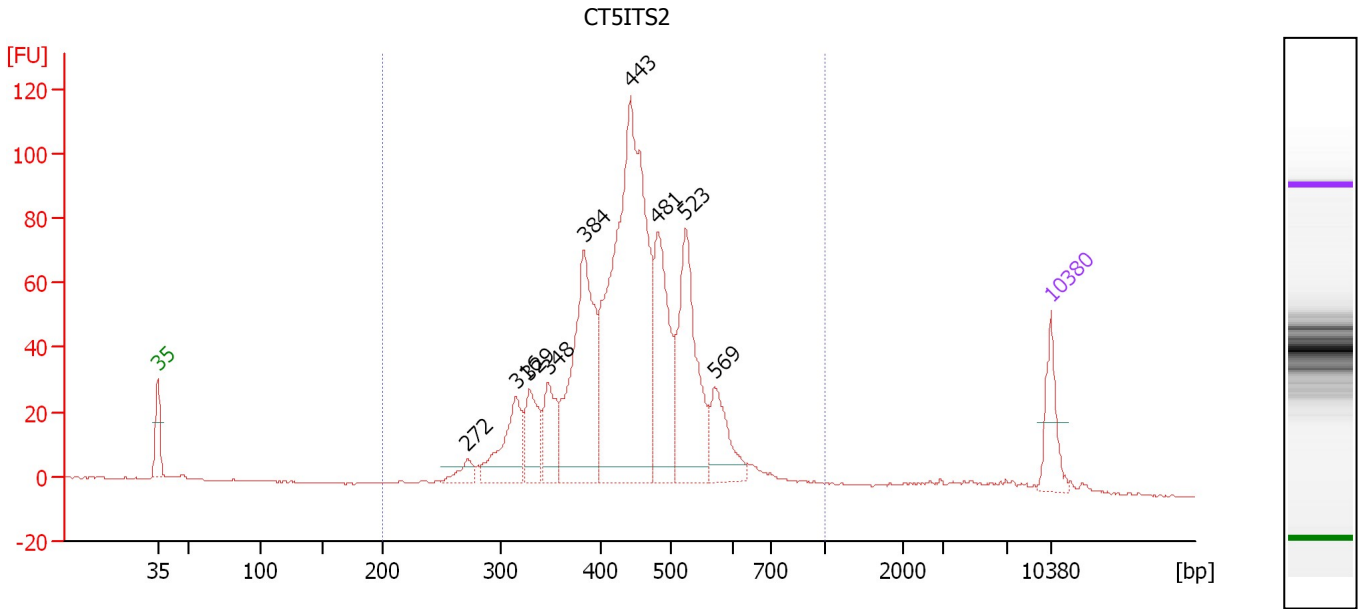

Overall Results for sample 7 : CT5ITS2

Number of peaks found: 9 Corr. Area 1: 1,236.4
 Noise: 0.2

Peak table for sample 7 : CT5ITS2

Peak	Size [bp]	Conc. [pg/μl]	Molarity [pmol/l]	Observations	Aligned Migration Time [s]
1	35	125.00	5,411.3	Lower Marker	43.00
2	272	40.99	228.5		67.26
3	316	173.99	835.3		71.08
4	329	119.08	548.5		72.12
5	348	120.66	525.4		73.61
6	384	452.23	1,784.1		76.44
7	443	1,035.87	3,546.5		80.06
8	481	290.50	915.5		82.19
9	523	345.53	1,000.8		84.37
10	569	122.39	325.6		86.61
11	10,380	75.00	10.9	Upper Marker	113.00


Region table for sample 7 : CT5ITS2

From [bp]	To [bp]	Average Size [bp]	Corr. Area	Molarity [pmol/l]	Co Conc. lor [pg/μl]	% of Total	Size distribution in CV [%]
200	1,000	448	1,236.4	10,363.7	2,901.28	95	20.4

Assay Class: High Sensitivity DNA Assay
Data Path: C:\...2-09\2015-12-09_002_CT516S_CT5ITS3_CT5ITS2_LigationRxns.xad

Created: 12/9/2015 11:00:29 AM
Modified: 12/9/2015 11:47:35 AM

Gel Image

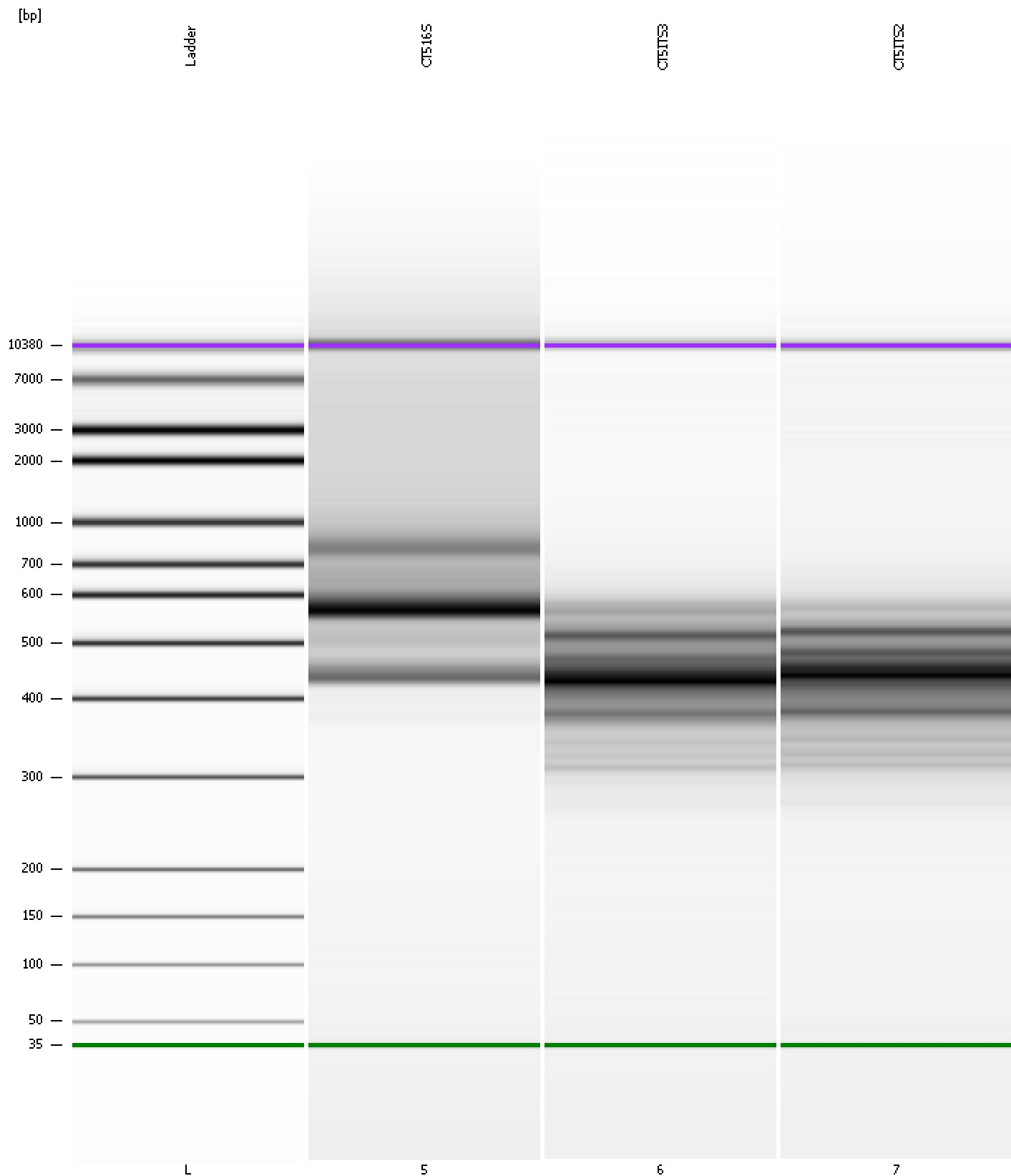

Assay Class: High Sensitivity DNA Assay
 Data Path: C:\...2-09\2015-12-09_002_CT516S_CT5ITS3_CT5ITS2_LigationRxns.xad

Created: 12/9/2015 11:00:29 AM
 Modified: 12/9/2015 11:47:35 AM

Run Logbook

Description	Number	Source	Category	Sub Category	Time	Time Zone	User	Host
Run ended on port 1 (Number of wells acquired: 12)		Instrument	Run		12/9/2015 11:40:51 AM	(GMT --08:00) Pacific Standard Time	UC Davis	D8XSMGH1
Run started on port 1 (File: C:\Documents and Settings\Bioanalyser\2015-12-09\2015-12-09_002.xad)		Instrument	Run		12/9/2015 11:00:29 AM	(GMT --08:00) Pacific Standard Time	UC Davis	D8XSMGH1
Product Number : G2938B		Instrument	Run		12/9/2015 11:00:29 AM	(GMT --08:00) Pacific Standard Time	UC Davis	D8XSMGH1
Name :		Instrument	Run		12/9/2015 11:00:29 AM	(GMT --08:00) Pacific Standard Time	UC Davis	D8XSMGH1
Vendor : Agilent Technologies		Instrument	Run		12/9/2015 11:00:29 AM	(GMT --08:00) Pacific Standard Time	UC Davis	D8XSMGH1
Serial# : DE13701086		Instrument	Run		12/9/2015 11:00:29 AM	(GMT --08:00) Pacific Standard Time	UC Davis	D8XSMGH1
Firmware : C.01.069		Instrument	Run		12/9/2015 11:00:29 AM	(GMT --08:00) Pacific Standard Time	UC Davis	D8XSMGH1
Cartridge : Electrode		Instrument	Run		12/9/2015 11:00:29 AM	(GMT --08:00) Pacific Standard Time	UC Davis	D8XSMGH1