


Assay Class: High Sensitivity DNA Assay
Data Path: C:\... and Settings\Bioanalyzer\2015-12-12\2015-12-12_001_A-G.xad

Created: 12/12/2015 2:13:29 PM
Modified: 12/12/2015 2:56:49 PM

Electrophoresis File Run Summary

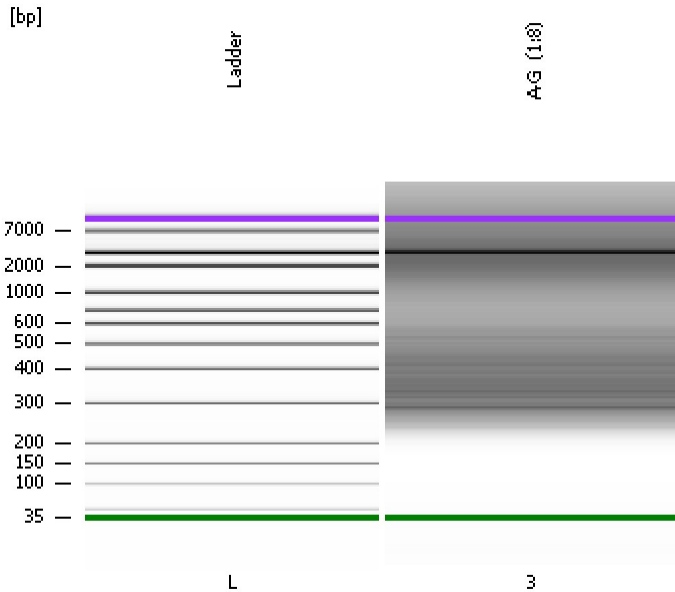

Instrument Information:

Instrument Name: DE13701086 Firmware: C.01.069
Serial#: DE13701086 Type: G2938B


Assay Information:

Assay Origin Path: C:\Program Files\Agilent\2100 bioanalyzer\2100 expert\assays\dsDNA\High Sensitivity DNA.xsy
Assay Class: High Sensitivity DNA Assay
Version: 1.03
Assay Comments: Copyright © 2003-2010 Agilent Technologies

Chip Information:

Chip Lot #:
Reagent Kit Lot #:
Chip Comments:

A-G (1:8)

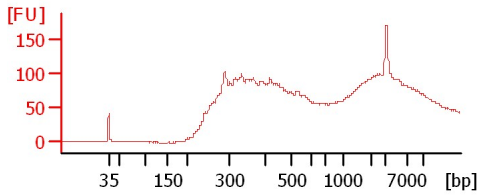

Assay Class: High Sensitivity DNA Assay
Data Path: C:\... and Settings\Bioanalyzer\2015-12-12\2015-12-12_001_A-G.xad

Created: 12/12/2015 2:13:29 PM
Modified: 12/12/2015 2:56:49 PM

Electrophoresis File Run Summary (Chip Summary)

Sample Name	Sample Comment	Rest. Digest	Status	Observation	Result Label	Result Color
A-G (1:8) Ladder		<input type="checkbox"/> <input type="checkbox"/>	 			
Chip Lot #				Reagent Kit Lot #		

Chip Comments :

Assay Class: High Sensitivity DNA Assay
Data Path: C:\... and Settings\Bioanalyzer\2015-12-12\2015-12-12_001_A-G.xad

Created: 12/12/2015 2:13:29 PM
Modified: 12/12/2015 2:56:49 PM

Electrophoresis Assay Details

General Analysis Settings

Number of Available Sample and Ladder Wells (Max.) : 12
Minimum Visible Range [s] : 32
Maximum Visible Range [s] : 138
Start Analysis Time Range [s] : 33
End Analysis Time Range [s] : 137.5
Ladder Concentration [pg/μl] : 1950
Uses Standard Area for Ladder Fragments
Lower Marker Concentration [pg/μl] : 125
Upper Marker Concentration [pg/μl] : 75
Used Upper Marker for Quantitation
Standard Curve Fit is Point to Point
Show Data Aligned to Lower and Upper Marker

Integrator Settings

Integration Start Time [s] : 33.05
Integration End Time [s] : 137
Slope Threshold : 0.8
Height Threshold [FU] : 5
Area Threshold : 0.1
Width Threshold [s] : 0.6
Baseline Plateau [s] : 0.5

Filter Settings

Filter Width [s] : 0.5
Polynomial Order : 4


Ladder

Ladder Peak	Size	Area
1	35	160
2	50	210
3	100	208
4	150	221
5	200	242
6	300	270
7	400	305
8	500	306
9	600	336
10	700	321
11	1000	366
12	2000	413
13	3000	411
14	7000	400
15	10380	214

Assay Class: High Sensitivity DNA Assay
 Data Path: C:\... and Settings\Bioanalyzer\2015-12-12\2015-12-12_001_A-G.xad

Created: 12/12/2015 2:13:29 PM
 Modified: 12/12/2015 2:56:49 PM

Electropherogram Summary

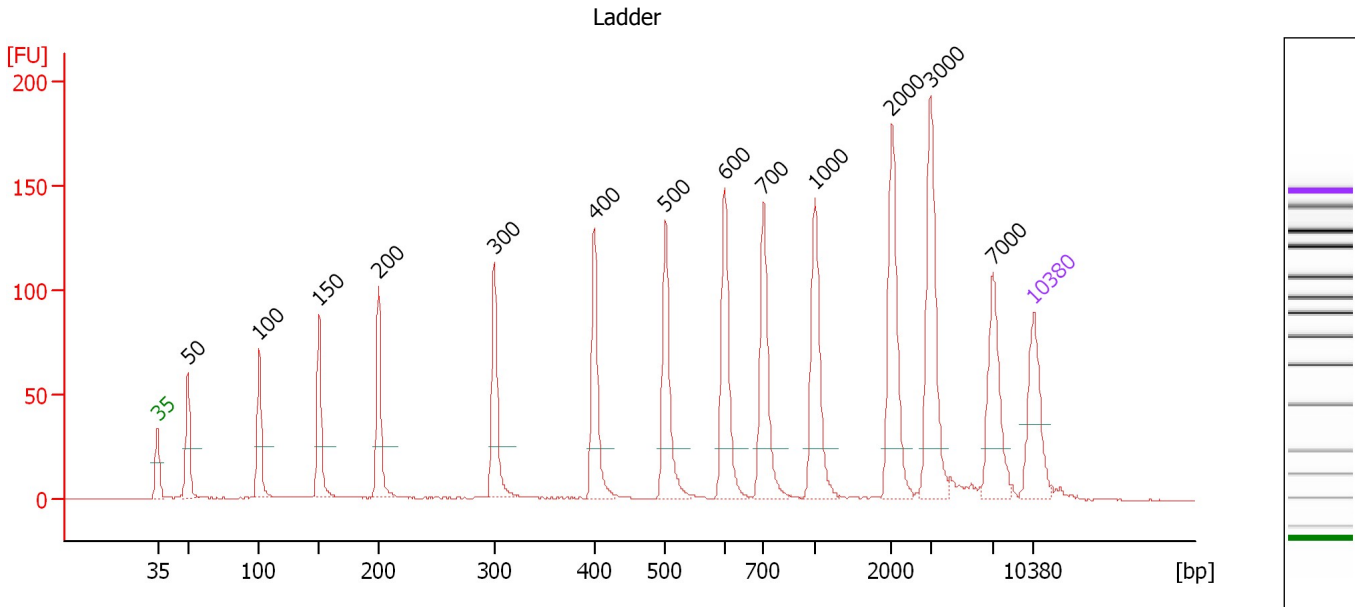

Overall Results for Ladder

Noise: 0.2


Peak table for Ladder

Peak	Size [bp]	Conc. [pg/μl]	Molarity [pmol/l]	Observations	Aligned Migration Time [s]
1	35	125.00	5,411.3	Lower Marker	43.00
2	50	150.00	4,545.5	Ladder Peak	45.42
3	100	150.00	2,272.7	Ladder Peak	51.13
4	150	150.00	1,515.2	Ladder Peak	55.91
5	200	150.00	1,136.4	Ladder Peak	60.65
6	300	150.00	757.6	Ladder Peak	69.89
7	400	150.00	568.2	Ladder Peak	77.88
8	500	150.00	454.5	Ladder Peak	83.60
9	600	150.00	378.8	Ladder Peak	88.34
10	700	150.00	324.7	Ladder Peak	91.40
11	1,000	150.00	227.3	Ladder Peak	95.49
12	2,000	150.00	113.6	Ladder Peak	101.67
13	3,000	150.00	75.8	Ladder Peak	104.78
14	7,000	150.00	32.5	Ladder Peak	109.75
15	10,380	75.00	10.9	Upper Marker	113.00

Assay Class: High Sensitivity DNA Assay
 Data Path: C:\... and Settings\Bioanalyzer\2015-12-12\2015-12-12_001_A-G.xad

Created: 12/12/2015 2:13:29 PM
 Modified: 12/12/2015 2:56:49 PM

Electropherogram Summary Continued ...

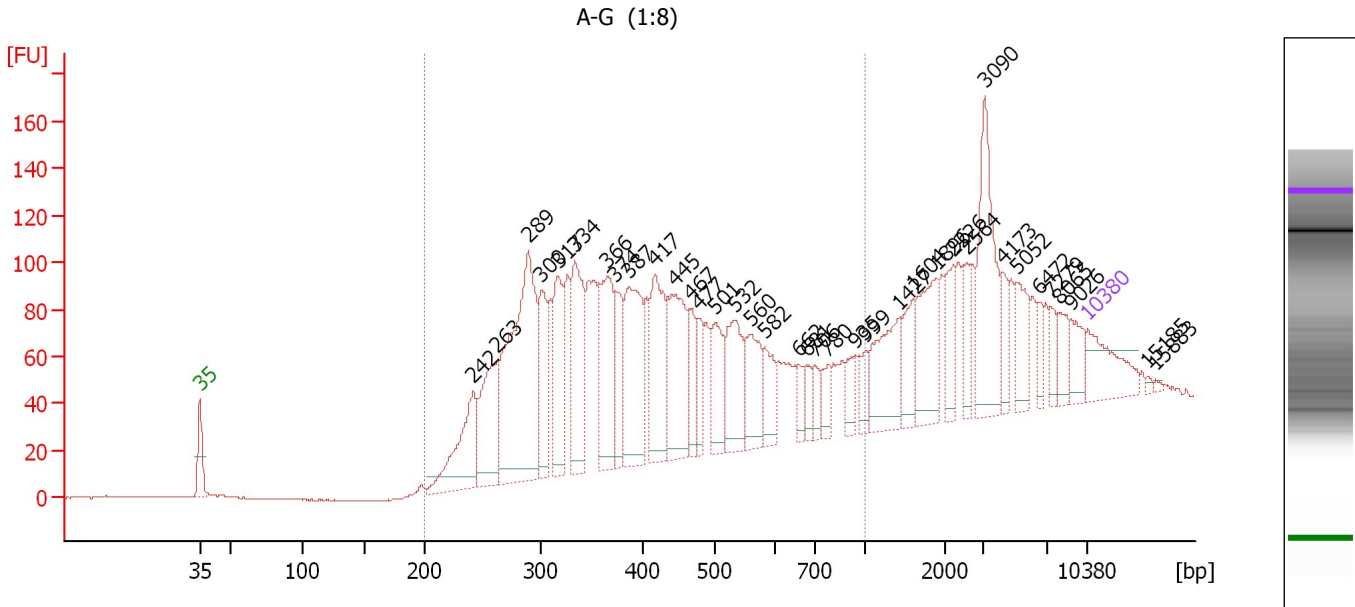

Overall Results for sample 3 : A-G (1:8)

Number of peaks found: 37 Corr. Area 1: 2,178.1
 Noise: 0.2

Peak table for sample 3 : A-G (1:8)

Peak	Size [bp]	Conc. [pg/μl]	Molarity [pmol/l]	Observations	Aligned Migration Time [s]
1	35	125.00	5,411.3	Lower Marker	43.00
2	242	198.51	1,241.7		64.55
3	263	204.06	1,177.0		66.45
4	289	536.50	2,812.3		68.88
5	300	145.47	734.8		69.89
6	317	177.30	848.0		71.24
7	334	207.18	940.9		72.58
8	366	201.28	833.9		75.14
9	374	103.85	421.2		75.77
10	387	224.71	880.4		76.82
11	417	181.56	658.9		78.88
12	445	209.21	711.7		80.48
13	467	68.37	221.9		81.70
14	477	59.41	188.7		82.29
15	501	99.59	301.4		83.63
16	532	139.07	396.3		85.10
17	560	100.19	271.0		86.45
18	582	62.14	161.7		87.50
19	662	28.36	64.9		90.23
20	681	27.81	61.9		90.82
21	706	30.48	65.4		91.49
22	780	34.21	66.4		92.50
23	935	35.35	57.3		94.60
24	999	35.59	54.0		95.48
25	1,427	118.77	126.1		98.13

Assay Class: High Sensitivity DNA Assay
 Data Path: C:\... and Settings\Bioanalyzer\2015-12-12\2015-12-12_001_A-G.xad

Created: 12/12/2015 2:13:29 PM
 Modified: 12/12/2015 2:56:49 PM

Electropherogram Summary Continued ...

... Peak table for sample 3 : A-G (1:8)

Peak	Size [bp]	Conc. [pg/μl]	Molarity [pmol/l]	Observations	Aligned Migration Time [s]
26	1,604	59.26	56.0		99.22
27	1,896	110.38	88.2		101.03
28	2,226	55.50	37.8		102.37
29	2,564	41.99	24.8		103.42
30	3,090	182.07	89.3		104.89
31	4,173	34.46	12.5		106.24
32	5,052	62.90	18.9		107.33
33	6,472	24.17	5.7		109.09
34	7,279	25.41	5.3		110.02
35	8,065	33.82	6.4		110.77
36	9,026	39.18	6.6		111.70
37	10,380	75.00	10.9	Upper Marker	113.00
38	15,185	0.00	0.0		117.62
39	15,883	0.00	0.0		118.29


Region table for sample 3 : A-G (1:8)

From [bp]	To [bp]	Average Size [bp]	Corr. Area	Molarity [pmol/l]	Co Conc. lor [pg/μl]	% of Total	Size distribution in CV [%]
200	1,000	430	2,178.1	12,533.7	3,062.84	70	36.6

Assay Class: High Sensitivity DNA Assay
Data Path: C:\... and Settings\Bioanalyzer\2015-12-12\2015-12-12_001_A-G.xad

Created: 12/12/2015 2:13:29 PM
Modified: 12/12/2015 2:56:49 PM

Gel Image

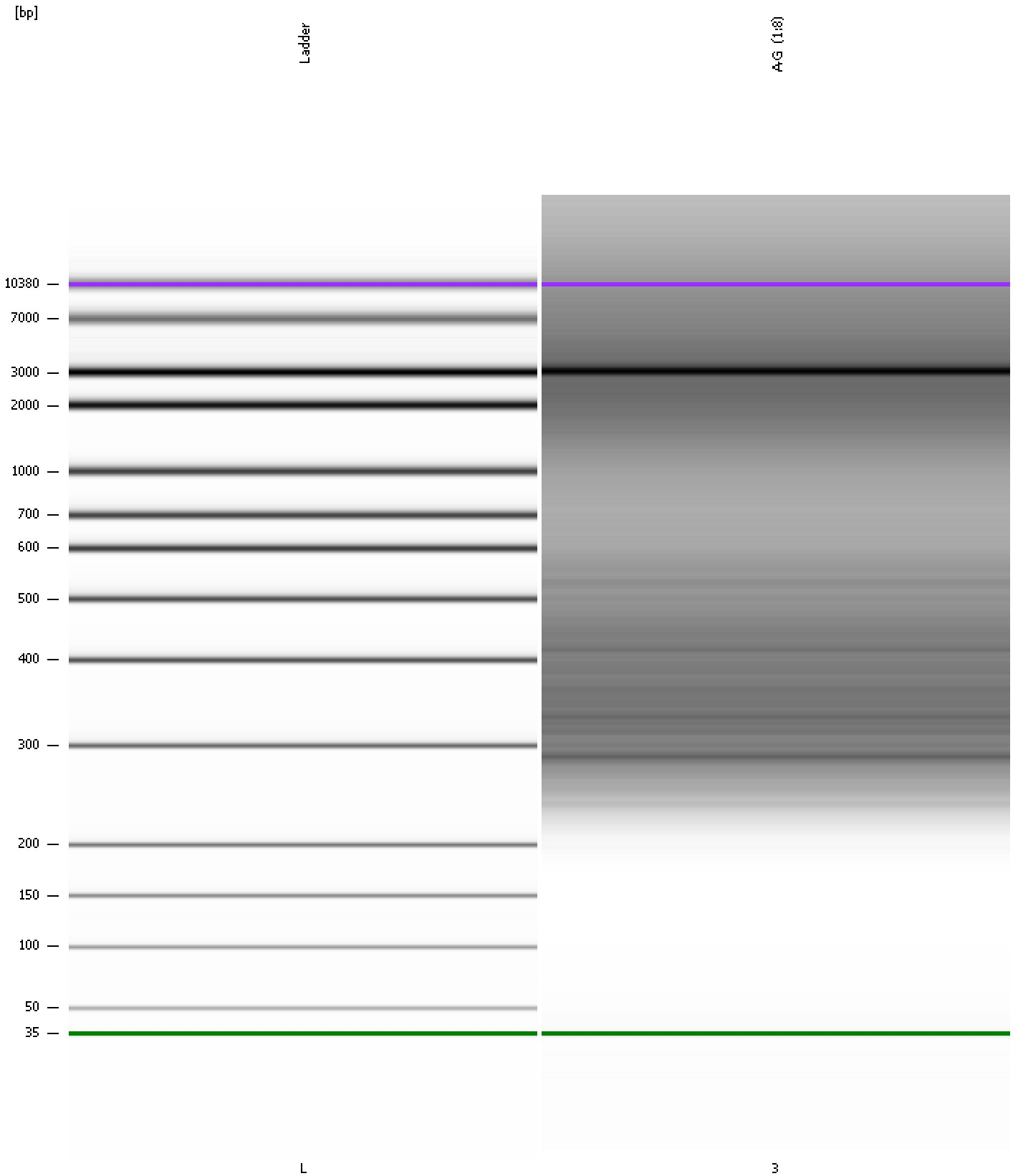

Assay Class: High Sensitivity DNA Assay
 Data Path: C:\... and Settings\Bioanalyzer\2015-12-12\2015-12-12_001_A-G.xad

Created: 12/12/2015 2:13:29 PM
 Modified: 12/12/2015 2:56:49 PM

Run Logbook

Description	Number	Source	Category	Sub Category	Time	Time Zone	User	Host
Run ended on port 1 (Number of wells acquired: 10)		Instrument	Run		12/12/2015 2:49:05 PM	(GMT --08:00) Pacific Standard Time	UC Davis	D8XSMGH1
Run started on port 1 (File: C:\Documents and Settings\Bioanalyzer\2015-12-12\2015-12-12_001.xad)		Instrument	Run		12/12/2015 2:13:34 PM	(GMT --08:00) Pacific Standard Time	UC Davis	D8XSMGH1
Product Number : G2938B		Instrument	Run		12/12/2015 2:13:34 PM	(GMT --08:00) Pacific Standard Time	UC Davis	D8XSMGH1
Name :		Instrument	Run		12/12/2015 2:13:34 PM	(GMT --08:00) Pacific Standard Time	UC Davis	D8XSMGH1
Vendor : Agilent Technologies		Instrument	Run		12/12/2015 2:13:34 PM	(GMT --08:00) Pacific Standard Time	UC Davis	D8XSMGH1
Serial# : DE13701086		Instrument	Run		12/12/2015 2:13:34 PM	(GMT --08:00) Pacific Standard Time	UC Davis	D8XSMGH1
Firmware : C.01.069		Instrument	Run		12/12/2015 2:13:34 PM	(GMT --08:00) Pacific Standard Time	UC Davis	D8XSMGH1
Cartridge : Electrode		Instrument	Run		12/12/2015 2:13:34 PM	(GMT --08:00) Pacific Standard Time	UC Davis	D8XSMGH1