


Assay Class: High Sensitivity DNA Assay
Data Path: C:\...ents and Settings\Bioanalyzer\2016-02-17\2016-02-17_001.xad

Created: 2/17/2016 10:57:37 AM
Modified: 2/17/2016 11:21:47 AM

Electrophoresis File Run Summary

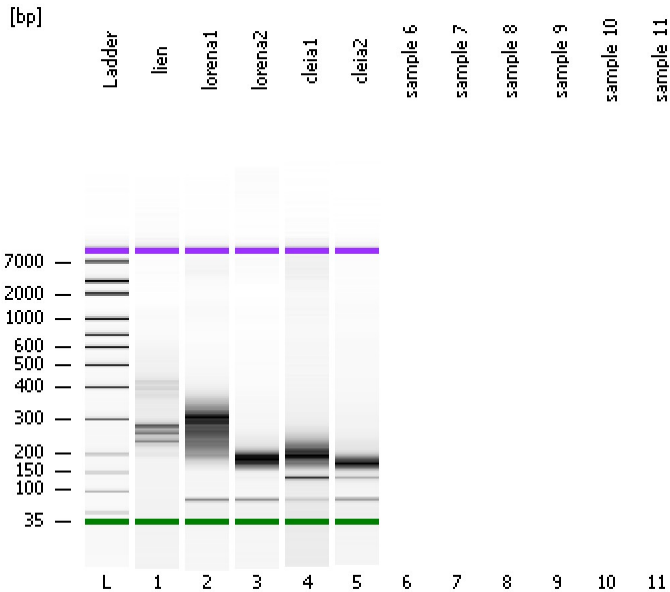

Instrument Information:

Instrument Name: DE13701086 Firmware: C.01.069
Serial#: DE13701086 Type: G2938B

Assay Information:

Assay Origin Path: C:\Program Files\Agilent\2100 bioanalyzer\2100 expert\assays\dsDNA\High Sensitivity DNA.xsy
Assay Class: High Sensitivity DNA Assay
Version: 1.03
Assay Comments: Copyright © 2003-2010 Agilent Technologies




Chip Information:

Chip Lot #: Reagent Kit Lot #: Chip Comments:

lien



lorena1

lorena2

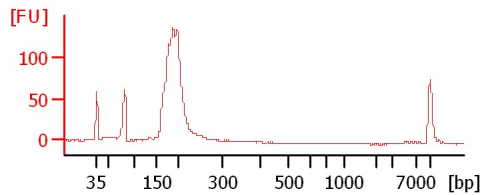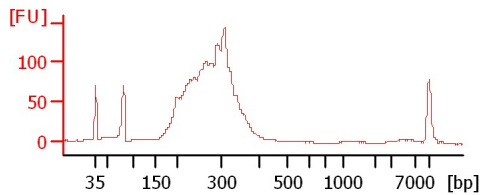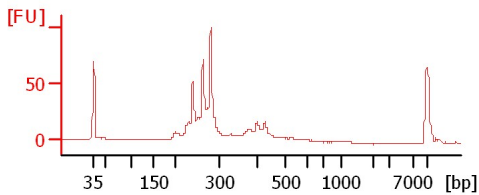

cleia1

cleia2

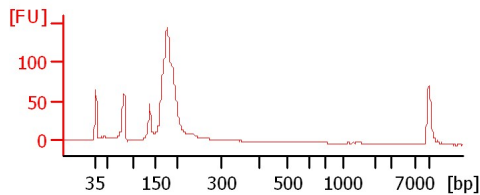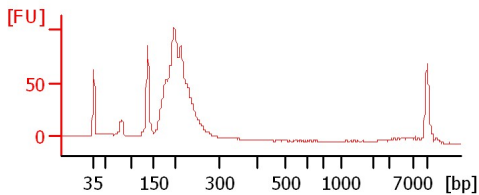

Assay Class: High Sensitivity DNA Assay
 Data Path: C:\...ents and Settings\Bioanalyzer\2016-02-17\2016-02-17_001.xad

Created: 2/17/2016 10:57:37 AM
 Modified: 2/17/2016 11:21:47 AM

Electrophoresis File Run Summary (Chip Summary)

Sample Name	Sample Comment	Rest. Digest	Status	Observation	Result Label	Result Color
lien		<input type="checkbox"/>	✓			
lorena1		<input type="checkbox"/>	✓			
lorena2		<input type="checkbox"/>	✓			
cleia1		<input type="checkbox"/>	✓			
cleia2		<input type="checkbox"/>	✓			
sample 6		<input type="checkbox"/>				
sample 7		<input type="checkbox"/>				
sample 8		<input type="checkbox"/>				
sample 9		<input type="checkbox"/>				
sample 10		<input type="checkbox"/>				
sample 11		<input type="checkbox"/>				
Ladder		<input type="checkbox"/>	✓			

Chip Lot #

Reagent Kit Lot #

Chip Comments :

Assay Class: High Sensitivity DNA Assay
Data Path: C:\...ents and Settings\Bioanalyzer\2016-02-17\2016-02-17_001.xad

Created: 2/17/2016 10:57:37 AM
Modified: 2/17/2016 11:21:47 AM

Electrophoresis Assay Details

General Analysis Settings

Number of Available Sample and Ladder Wells (Max.) : 12
Minimum Visible Range [s] : 32
Maximum Visible Range [s] : 138
Start Analysis Time Range [s] : 33
End Analysis Time Range [s] : 137.5
Ladder Concentration [pg/μl] : 1950
Uses Standard Area for Ladder Fragments
Lower Marker Concentration [pg/μl] : 125
Upper Marker Concentration [pg/μl] : 75
Used Upper Marker for Quantitation
Standard Curve Fit is Point to Point
Show Data Aligned to Lower and Upper Marker

Integrator Settings

Integration Start Time [s] : 33.05
Integration End Time [s] : 137
Slope Threshold : 0.8
Height Threshold [FU] : 5
Area Threshold : 0.1
Width Threshold [s] : 0.6
Baseline Plateau [s] : 0.5

Filter Settings

Filter Width [s] : 0.5
Polynomial Order : 4


Ladder

Ladder Peak	Size	Area
1	35	160
2	50	210
3	100	208
4	150	221
5	200	242
6	300	270
7	400	305
8	500	306
9	600	336
10	700	321
11	1000	366
12	2000	413
13	3000	411
14	7000	400
15	10380	214

Assay Class: High Sensitivity DNA Assay
 Data Path: C:\...ents and Settings\Bioanalyzer\2016-02-17\2016-02-17_001.xad

Created: 2/17/2016 10:57:37 AM
 Modified: 2/17/2016 11:21:47 AM

Electropherogram Summary

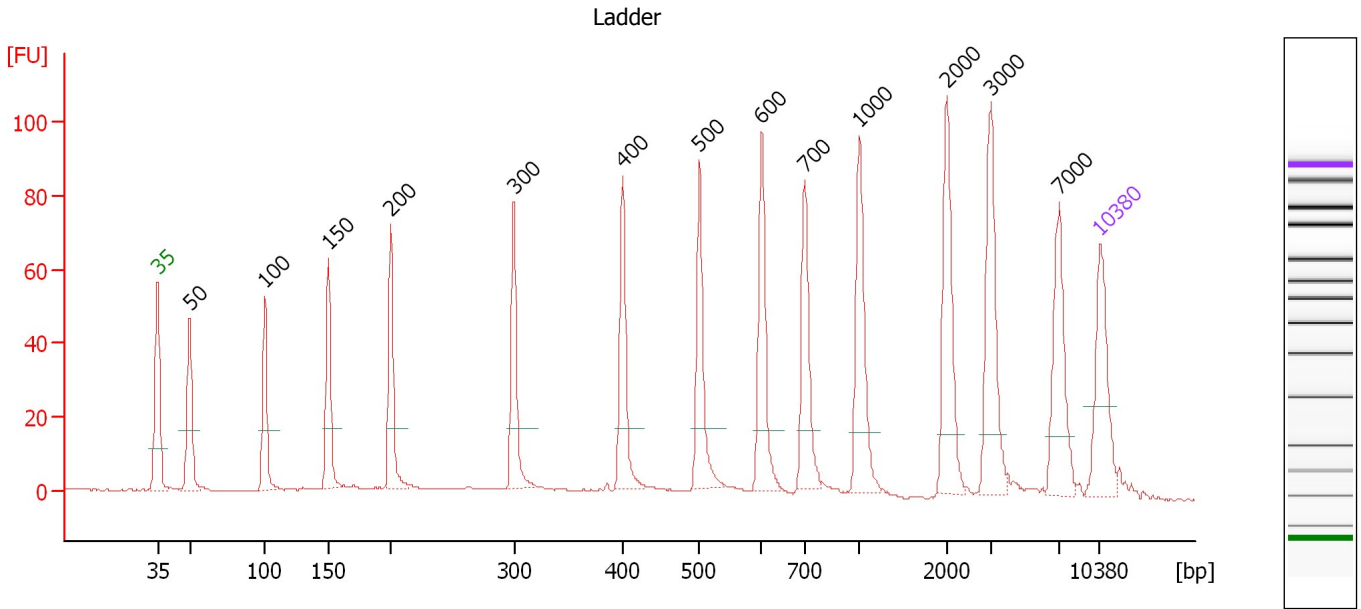

Overall Results for Ladder

Noise: 0.1


Peak table for Ladder

Peak	Size [bp]	Conc. [pg/μl]	Molarity [pmol/l]	Observations	Aligned Migration Time [s]
1	35	125.00	5,411.3	Lower Marker	43.00
2	50	150.00	4,545.5	Ladder Peak	45.39
3	100	150.00	2,272.7	Ladder Peak	50.97
4	150	150.00	1,515.2	Ladder Peak	55.69
5	200	150.00	1,136.4	Ladder Peak	60.32
6	300	150.00	757.6	Ladder Peak	69.48
7	400	150.00	568.2	Ladder Peak	77.55
8	500	150.00	454.5	Ladder Peak	83.27
9	600	150.00	378.8	Ladder Peak	87.85
10	700	150.00	324.7	Ladder Peak	91.10
11	1,000	150.00	227.3	Ladder Peak	95.15
12	2,000	150.00	113.6	Ladder Peak	101.64
13	3,000	150.00	75.8	Ladder Peak	104.89
14	7,000	150.00	32.5	Ladder Peak	109.95
15	10,380	75.00	10.9	Upper Marker	113.00

Assay Class: High Sensitivity DNA Assay
 Data Path: C:\...ents and Settings\Bioanalyzer\2016-02-17\2016-02-17_001.xad

Created: 2/17/2016 10:57:37 AM
 Modified: 2/17/2016 11:21:47 AM

Electropherogram Summary Continued ...

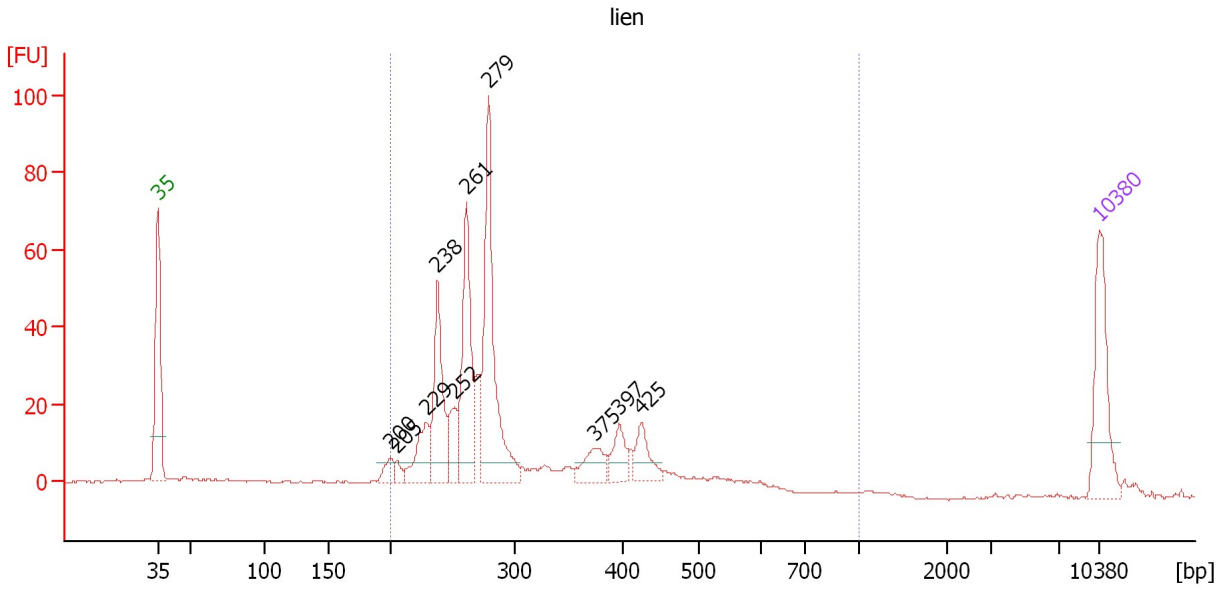

Overall Results for sample 1 : lien

Number of peaks found: 10 Corr. Area 1: 523.1
 Noise: 0.2

Peak table for sample 1 : lien

Peak	Size [bp]	Conc. [pg/μl]	Molarity [pmol/l]	Observations	Aligned Migration Time [s]
1	35	125.00	5,411.3	Lower Marker	43.00
2	200	17.81	135.0		60.31
3	205	9.89	73.1		60.79
4	229	57.25	379.5		62.94
5	238	119.12	758.5		63.80
6	252	42.20	254.2		65.04
7	261	147.09	852.4		65.95
8	279	230.84	1,252.7		67.58
9	375	35.32	142.8		75.51
10	397	31.73	121.2		77.28
11	425	36.48	130.2		78.96
12	10,380	75.00	10.9	Upper Marker	113.00


Region table for sample 1 : lien

From [bp]	To [bp]	Average Size [bp]	Corr. Area	Molarity [pmol/l]	Co Conc. lor [pg/μl]	% of Total	Size distribution in CV [%]
200	1,000	314	523.1	4,981.0	945.74	88	29.7

Assay Class: High Sensitivity DNA Assay
 Data Path: C:\...ents and Settings\Bioanalyzer\2016-02-17\2016-02-17_001.xad

Created: 2/17/2016 10:57:37 AM
 Modified: 2/17/2016 11:21:47 AM

Electropherogram Summary Continued ...

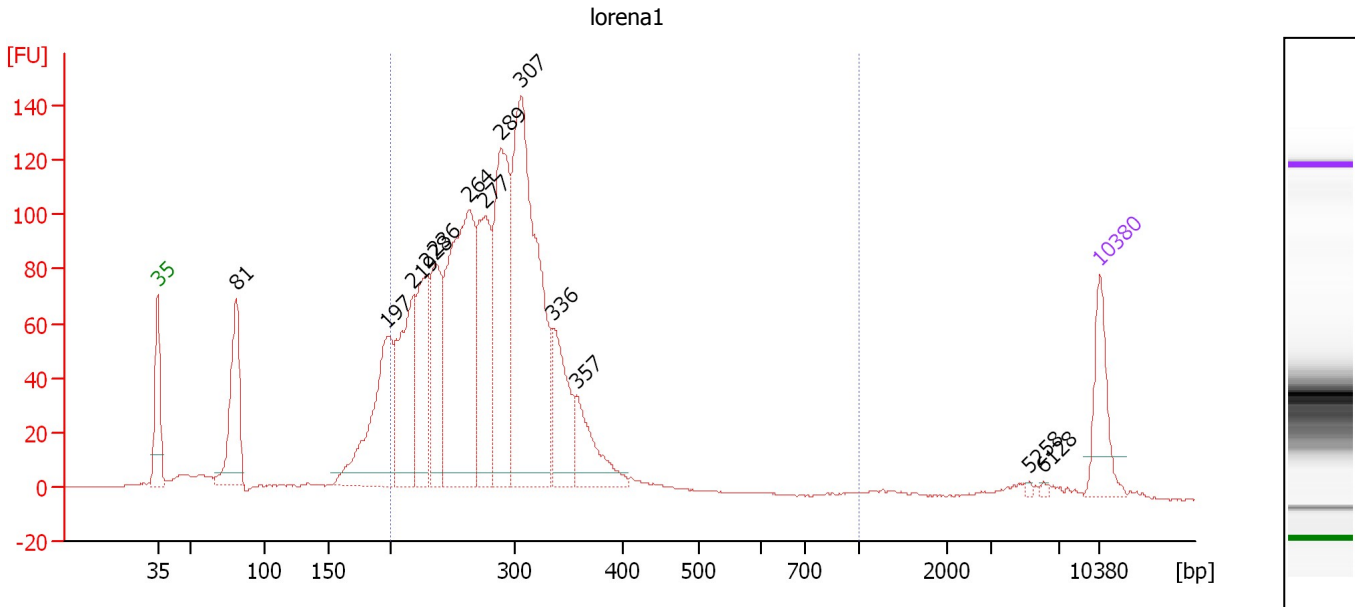

Overall Results for sample 2 : lorena1

Number of peaks found: 13 Corr. Area 1: 1,867.6
 Noise: 0.3

Peak table for sample 2 : lorena1

Peak	Size [bp]	Conc. [pg/μl]	Molarity [pmol/l]	Observations	Aligned Migration Time [s]
1	35	125.00	5,411.3	Lower Marker	43.00
2	81	225.98	4,213.2		48.88
3	197	294.65	2,264.3		60.06
4	219	253.98	1,761.1		62.02
5	228	195.57	1,297.2		62.92
6	236	203.39	1,304.5		63.64
7	264	534.58	3,069.4		66.17
8	277	258.33	1,413.4		67.37
9	289	371.87	1,950.1		68.47
10	307	654.31	3,229.6		70.04
11	336	145.65	656.8		72.39
12	357	113.84	482.7		74.11
13	5,258	2.64	0.8		107.74
14	6,128	2.73	0.7		108.84
15	10,380	75.00	10.9	Upper Marker	113.00


Region table for sample 2 : lorena1

From [bp]	To [bp]	Average Size [bp]	Corr. Area	Molarity [pmol/l]	Co Conc. lor [pg/μl]	% of Total	Size distribution in CV [%]
200	1,000	284	1,867.6	15,877.8	2,890.80	83	17.7

Assay Class: High Sensitivity DNA Assay
 Data Path: C:\...ents and Settings\Bioanalyzer\2016-02-17\2016-02-17_001.xad

Created: 2/17/2016 10:57:37 AM
 Modified: 2/17/2016 11:21:47 AM

Electropherogram Summary Continued ...

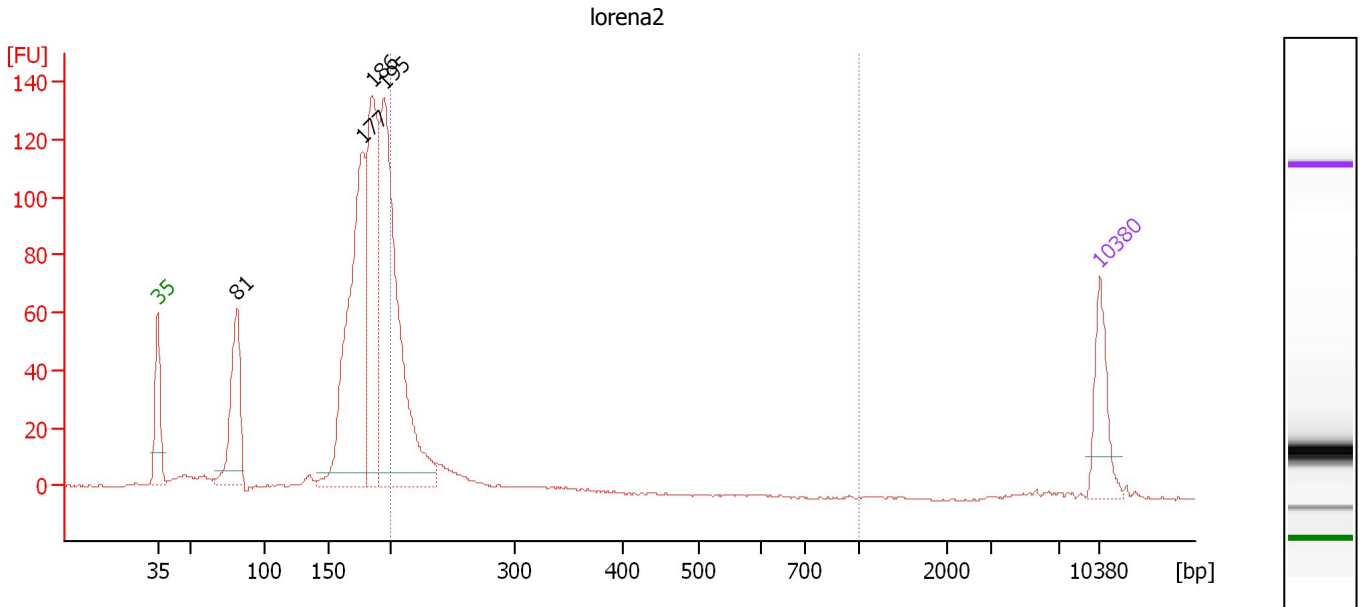

Overall Results for sample 3 : lorena2

Number of peaks found: 4 Corr. Area 1: 226.3
 Noise: 0.4

Peak table for sample 3 : lorena2

Peak	Size [bp]	Conc. [pg/μl]	Molarity [pmol/l]	Observations	Aligned Migration Time [s]
1	35	125.00	5,411.3	Lower Marker	43.00
2	81	221.58	4,121.5		48.90
3	177	564.37	4,826.1		58.21
4	186	427.75	3,491.8		58.99
5	195	677.55	5,262.2		59.87
6	10,380	75.00	10.9	Upper Marker	113.00


Region table for sample 3 : lorena2

From [bp]	To [bp]	Average Size [bp]	Corr. Area	Molarity [pmol/l]	Co Conc. lor [pg/μl]	% of Total	Size distribution in CV [%]
200	1,000	219	226.3	2,900.5	417.27	21	11.7

Assay Class: High Sensitivity DNA Assay
 Data Path: C:\...ents and Settings\Bioanalyzer\2016-02-17\2016-02-17_001.xad

Created: 2/17/2016 10:57:37 AM
 Modified: 2/17/2016 11:21:47 AM

Electropherogram Summary Continued ...

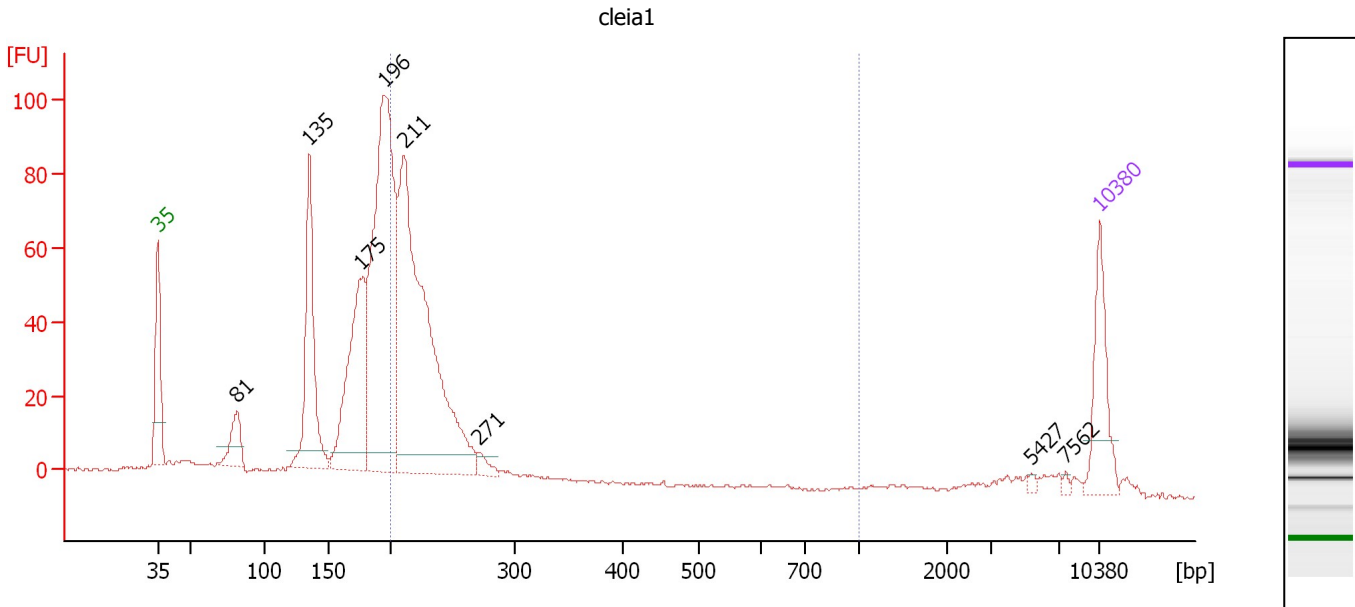

Overall Results for sample 4 : cleia1

Number of peaks found: 8 Corr. Area 1: 445.8
 Noise: 0.4

Peak table for sample 4 : cleia1

Peak	Size [bp]	Conc. [pg/μl]	Molarity [pmol/l]	Observations	Aligned Migration Time [s]
1	35	125.00	5,411.3	Lower Marker	43.00
2	81	57.41	1,068.6		48.89
3	135	238.97	2,680.3		54.28
4	175	274.17	2,370.5		58.03
5	196	612.79	4,743.6		59.93
6	211	668.41	4,801.0		61.32
7	271	16.62	92.8		66.87
8	5,427	3.43	1.0		107.96
9	7,562	3.30	0.7		110.45
10	10,380	75.00	10.9	Upper Marker	113.00


Region table for sample 4 : cleia1

From [bp]	To [bp]	Average Size [bp]	Corr. Area	Molarity [pmol/l]	Co Conc. lor [pg/μl]	% of Total	Size distribution in CV [%]
200	1,000	227	445.8	5,776.2	858.42	42	12.7

Assay Class: High Sensitivity DNA Assay
 Data Path: C:\...ents and Settings\Bioanalyzer\2016-02-17\2016-02-17_001.xad

Created: 2/17/2016 10:57:37 AM
 Modified: 2/17/2016 11:21:47 AM

Electropherogram Summary Continued ...

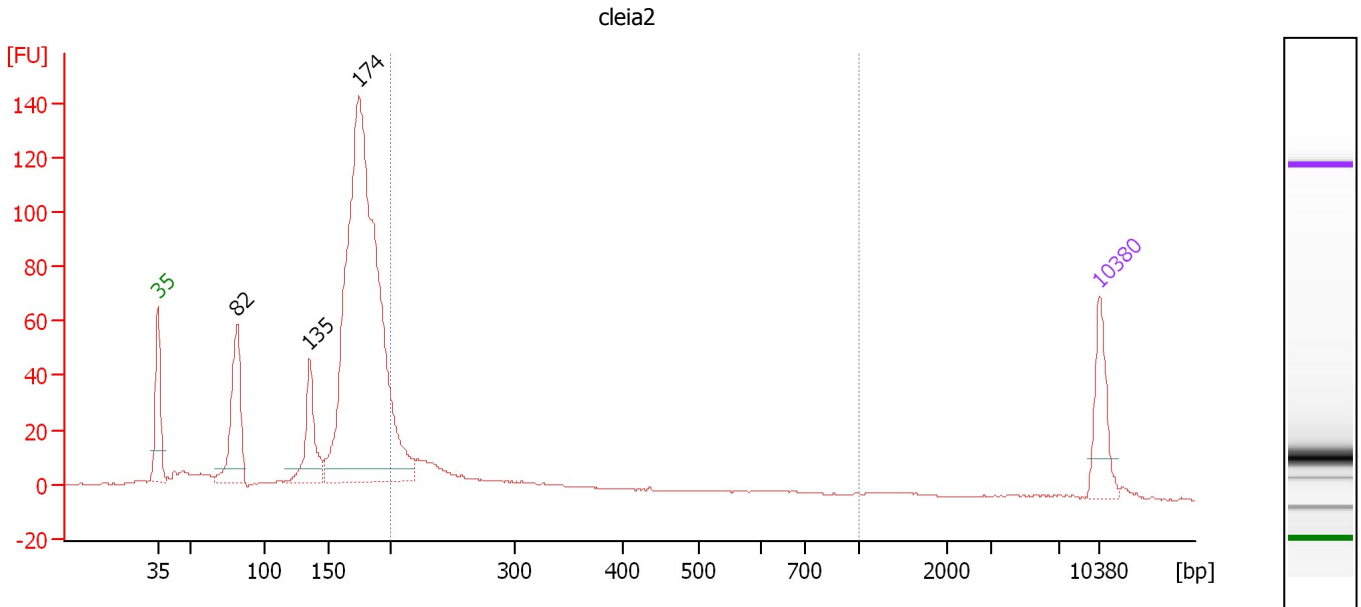

Overall Results for sample 5 : cleia2

Number of peaks found: 3 Corr. Area 1: 148.6
 Noise: 0.3

Peak table for sample 5 : cleia2

Peak	Size [bp]	Conc. [pg/μl]	Molarity [pmol/l]	Observations	Aligned Migration Time [s]
1	35	125.00	5,411.3	Lower Marker	43.00
2	82	230.96	4,284.5		48.92
3	135	142.89	1,602.6		54.28
4	174	1,411.03	12,252.2		57.96
5	10,380	75.00	10.9	Upper Marker	113.00


Region table for sample 5 : cleia2

From [bp]	To [bp]	Average Size [bp]	Corr. Area	Molarity [pmol/l]	Co Conc. lor [pg/μl]	% of Total	Size distribution in CV [%]
200	1,000	272	148.6	1,750.0	285.75	15	35.6

Assay Class: High Sensitivity DNA Assay
Data Path: C:\...ents and Settings\Bioanalyzer\2016-02-17\2016-02-17_001.xad

Created: 2/17/2016 10:57:37 AM
Modified: 2/17/2016 11:21:47 AM

Gel Image

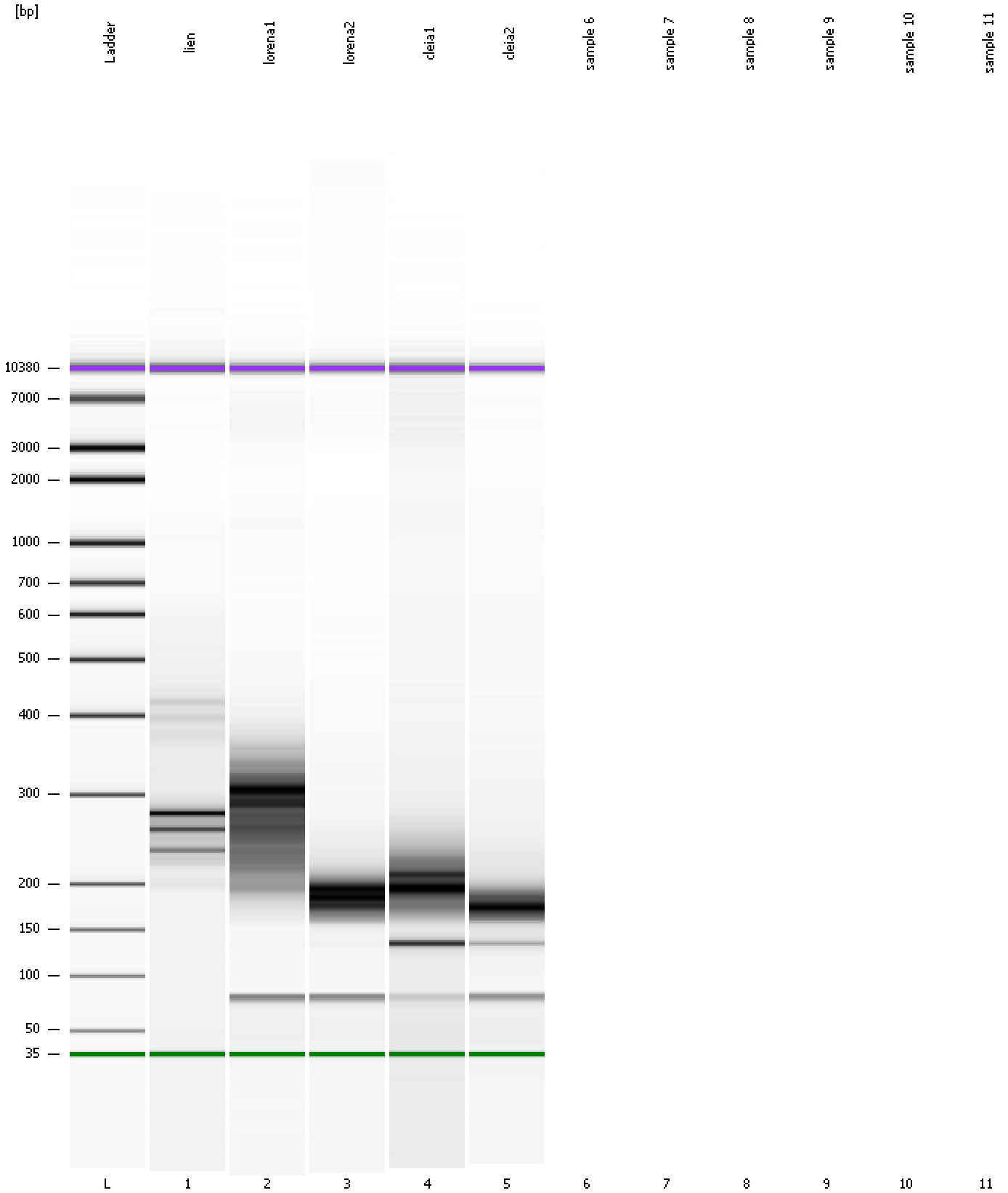

Assay Class: High Sensitivity DNA Assay
Data Path: C:\...ents and Settings\Bioanalyzer\2016-02-17\2016-02-17_001.xad

Created: 2/17/2016 10:57:37 AM
Modified: 2/17/2016 11:21:47 AM

Invalid Samples

- Sample 6 has not been run, no results available.
- Sample 7 has not been run, no results available.
- Sample 8 has not been run, no results available.
- Sample 9 has not been run, no results available.
- Sample 10 has not been run, no results available.
- Sample 11 has not been run, no results available.

Assay Class: High Sensitivity DNA Assay
 Data Path: C:\...ents and Settings\Bioanalyzer\2016-02-17\2016-02-17_001.xad

Created: 2/17/2016 10:57:37 AM
 Modified: 2/17/2016 11:21:47 AM

Run Logbook

Description	Number	Source	Category	Sub Category	Time	Time Zone	User	Host
Run ended on port 1 (Number of wells acquired: 6)		Instrument	Run		2/17/2016 11:21:46 AM	(GMT --08:00) Pacific Standard Time	UC Davis	D8XSMGH1
Run started on port 1 (File: C:\Documents and Settings\Bioanalyzer\2016-02-17\2016-02-17_001.xad)		Instrument	Run		2/17/2016 10:57:42 AM	(GMT --08:00) Pacific Standard Time	UC Davis	D8XSMGH1
Product Number : G2938B		Instrument	Run		2/17/2016 10:57:42 AM	(GMT --08:00) Pacific Standard Time	UC Davis	D8XSMGH1
Name :		Instrument	Run		2/17/2016 10:57:42 AM	(GMT --08:00) Pacific Standard Time	UC Davis	D8XSMGH1
Vendor : Agilent Technologies		Instrument	Run		2/17/2016 10:57:42 AM	(GMT --08:00) Pacific Standard Time	UC Davis	D8XSMGH1
Serial# : DE13701086		Instrument	Run		2/17/2016 10:57:42 AM	(GMT --08:00) Pacific Standard Time	UC Davis	D8XSMGH1
Firmware : C.01.069		Instrument	Run		2/17/2016 10:57:42 AM	(GMT --08:00) Pacific Standard Time	UC Davis	D8XSMGH1
Cartridge : Electrode		Instrument	Run		2/17/2016 10:57:42 AM	(GMT --08:00) Pacific Standard Time	UC Davis	D8XSMGH1