


Assay Class: High Sensitivity DNA Assay  
Data Path: C:\...ings\Bioanalyzer\2016-02-25\2016-02-25\_001\_Dag\_DY3\_DY\_6.xad

Created: 2/25/2016 2:12:58 PM  
Modified: 2/25/2016 3:32:15 PM

**Electrophoresis File Run Summary**

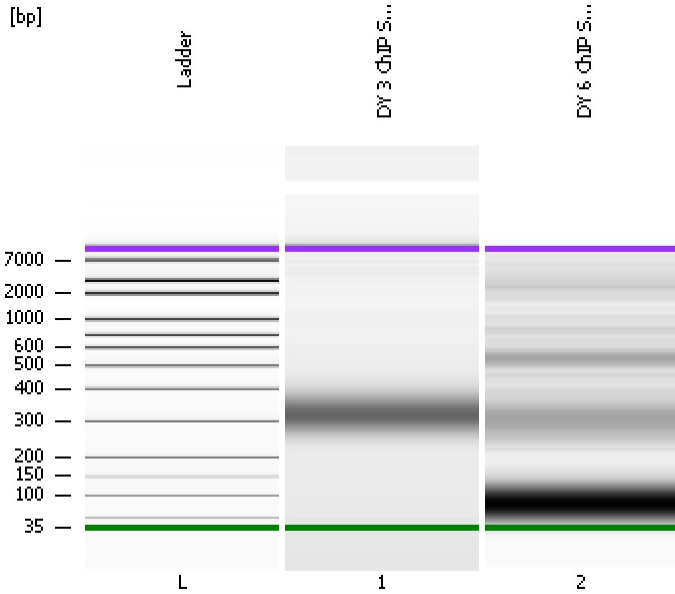

Instrument Information:

Instrument Name: DE13701086      Firmware: C.01.069  
Serial#: DE13701086      Type: G2938B

Assay Information:



Assay Origin Path: C:\Program Files\Agilent\2100 bioanalyzer\2100 expert\assays\dsDNA\High Sensitivity DNA.xsy  
Assay Class: High Sensitivity DNA Assay  
Version: 1.03  
Assay Comments: Copyright © 2003-2010 Agilent Technologies

Chip Information:

Chip Lot #:  
Reagent Kit Lot #:  
Chip Comments:

**DY 3 ChIP Seq Library**

**DY 6 ChIP Seq Library**

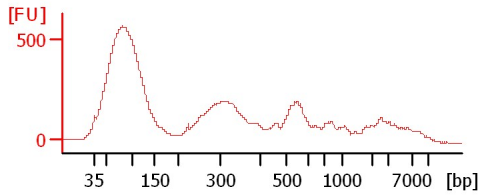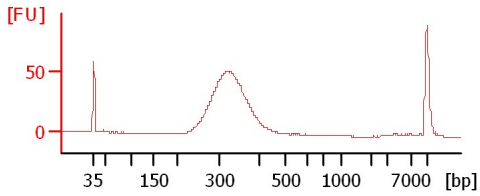

Assay Class: High Sensitivity DNA Assay  
Data Path: C:\...ings\Bioanalyzer\2016-02-25\2016-02-25\_001\_Dag\_DY3\_DY\_6.xad

Created: 2/25/2016 2:12:58 PM  
Modified: 2/25/2016 3:32:15 PM

**Electrophoresis File Run Summary (Chip Summary)**

Sample Name	Sample Comment	Rest. Digest	Status	Observation	Result Label	Result Color
DY 3 ChIP Seq Library		<input type="checkbox"/>	✓			
DY 6 ChIP Seq Library		<input type="checkbox"/>	✓			
Ladder		<input type="checkbox"/>	✓			

**Chip Lot #****Reagent Kit Lot #****Chip Comments :**

Assay Class: High Sensitivity DNA Assay  
Data Path: C:\...ings\Bioanalyzer\2016-02-25\2016-02-25\_001\_Dag\_DY3\_DY\_6.xad

Created: 2/25/2016 2:12:58 PM  
Modified: 2/25/2016 3:32:15 PM

**Electrophoresis Assay Details**

**General Analysis Settings**

Number of Available Sample and Ladder Wells (Max.) : 12  
Minimum Visible Range [s] : 32  
Maximum Visible Range [s] : 138  
Start Analysis Time Range [s] : 33  
End Analysis Time Range [s] : 137.5  
Ladder Concentration [pg/μl] : 1950  
Uses Standard Area for Ladder Fragments  
Lower Marker Concentration [pg/μl] : 125  
Upper Marker Concentration [pg/μl] : 75  
Used Upper Marker for Quantitation  
Standard Curve Fit is Point to Point  
Show Data Aligned to Lower and Upper Marker

**Integrator Settings**

Integration Start Time [s] : 33.05  
Integration End Time [s] : 137  
Slope Threshold : 0.8  
Height Threshold [FU] : 5  
Area Threshold : 0.1  
Width Threshold [s] : 0.6  
Baseline Plateau [s] : 0.5

**Filter Settings**

Filter Width [s] : 0.5  
Polynomial Order : 4


**Ladder**

Ladder Peak	Size	Area
1	35	160
2	50	210
3	100	208
4	150	221
5	200	242
6	300	270
7	400	305
8	500	306
9	600	336
10	700	321
11	1000	366
12	2000	413
13	3000	411
14	7000	400
15	10380	214

Assay Class: High Sensitivity DNA Assay  
 Data Path: C:\...ings\Bioanalyzer\2016-02-25\2016-02-25\_001\_Dag\_DY3\_DY\_6.xad

Created: 2/25/2016 2:12:58 PM  
 Modified: 2/25/2016 3:32:15 PM

**Electropherogram Summary**

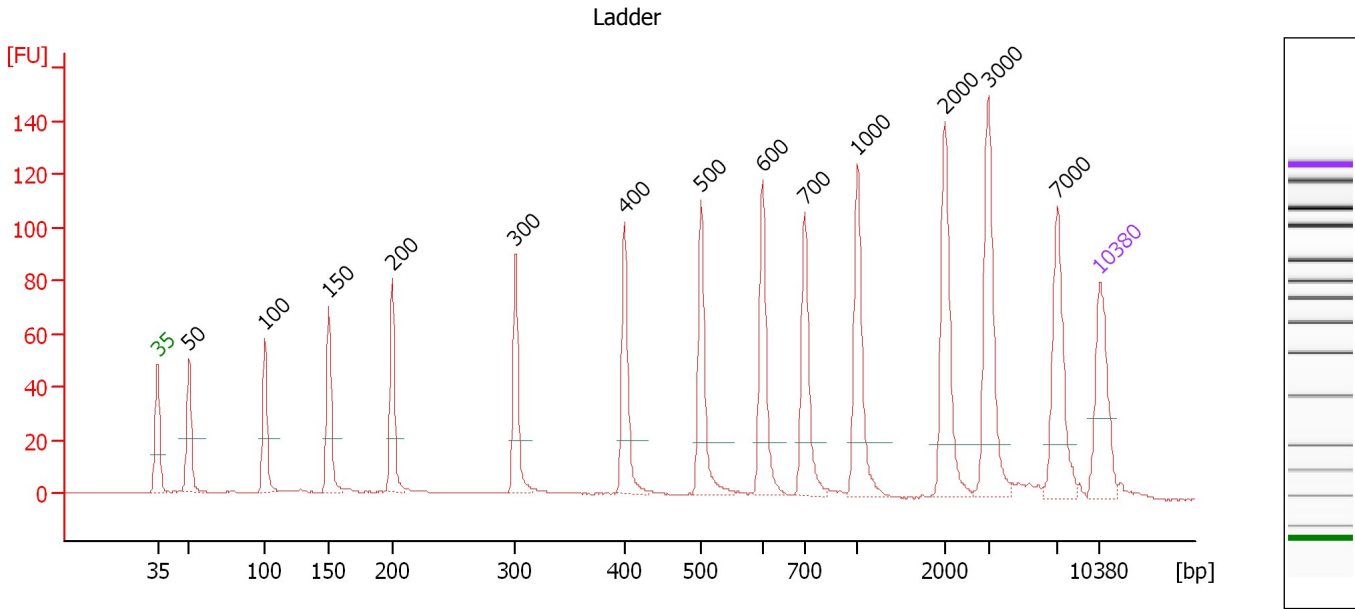

**Overall Results for Ladder**

Noise: 0.2


**Peak table for Ladder**

Peak	Size [bp]	Conc. [pg/μl]	Molarity [pmol/l]	Observations	Aligned Migration Time [s]
1	35	125.00	5,411.3	Lower Marker	43.00
2	50	150.00	4,545.5	Ladder Peak	45.34
3	100	150.00	2,272.7	Ladder Peak	50.98
4	150	150.00	1,515.2	Ladder Peak	55.72
5	200	150.00	1,136.4	Ladder Peak	60.46
6	300	150.00	757.6	Ladder Peak	69.59
7	400	150.00	568.2	Ladder Peak	77.73
8	500	150.00	454.5	Ladder Peak	83.36
9	600	150.00	378.8	Ladder Peak	87.95
10	700	150.00	324.7	Ladder Peak	91.10
11	1,000	150.00	227.3	Ladder Peak	94.99
12	2,000	150.00	113.6	Ladder Peak	101.47
13	3,000	150.00	75.8	Ladder Peak	104.72
14	7,000	150.00	32.5	Ladder Peak	109.86
15	10,380	75.00	10.9	Upper Marker	113.00

Assay Class: High Sensitivity DNA Assay  
 Data Path: C:\...ings\Bioanalyzer\2016-02-25\2016-02-25\_001\_Dag\_DY3\_DY\_6.xad

Created: 2/25/2016 2:12:58 PM  
 Modified: 2/25/2016 3:32:15 PM

**Electropherogram Summary Continued ...**

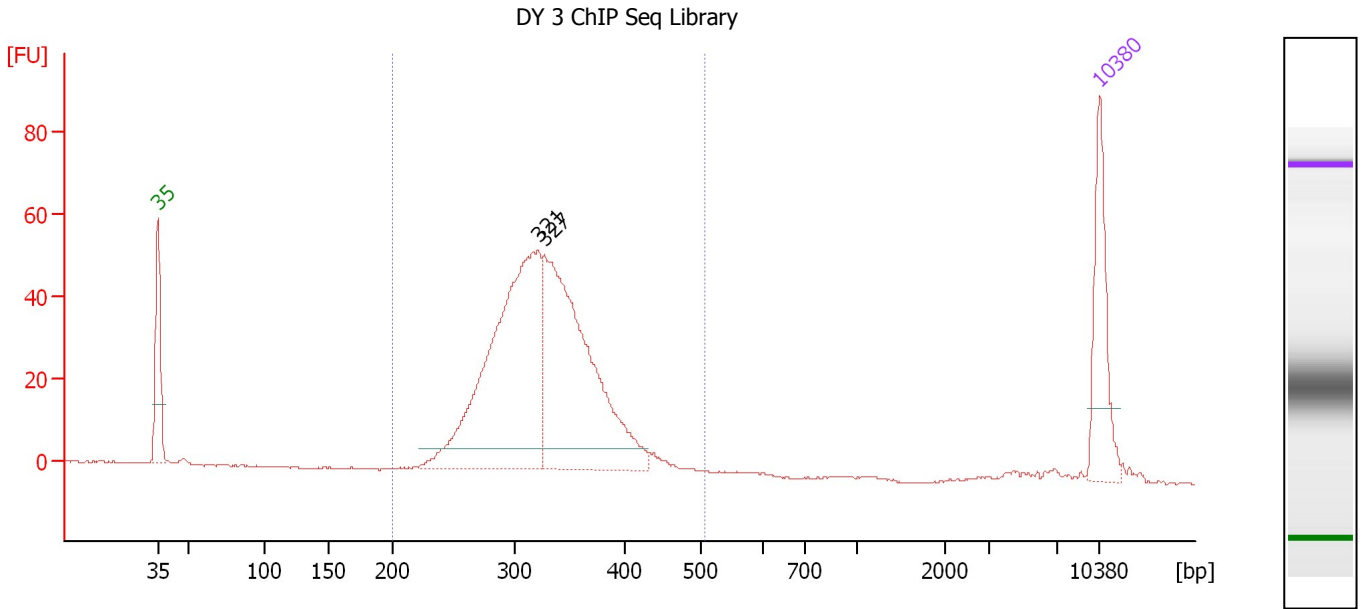

**Overall Results for sample 1 : DY 3 ChIP Seq Library**

Number of peaks found: 2                      Corr. Area 1: 612.2  
 Noise: 0.2

**Peak table for sample 1 : DY 3 ChIP Seq Library**

Peak	Size [bp]	Conc. [pg/μl]	Molarity [pmol/l]	Observations	Aligned Migration Time [s]
1	35	125.00	5,411.3	Lower Marker	43.00
2	321	439.83	2,079.3		71.26
3	327	385.55	1,788.3		71.76
4	10,380	75.00	10.9	Upper Marker	113.00


**Region table for sample 1 : DY 3 ChIP Seq Library**

From [bp]	To [bp]	Average Size [bp]	Corr. Area	Molarity [pmol/l]	Co Conc. lor [pg/μl]	% of Total	Size distribution in CV [%]
200	506	325	612.2	3,934.1	828.69	99	12.4

Assay Class: High Sensitivity DNA Assay  
 Data Path: C:\...ings\Bioanalyzer\2016-02-25\2016-02-25\_001\_Dag\_DY3\_DY\_6.xad

Created: 2/25/2016 2:12:58 PM  
 Modified: 2/25/2016 3:32:15 PM

**Electropherogram Summary Continued ...**

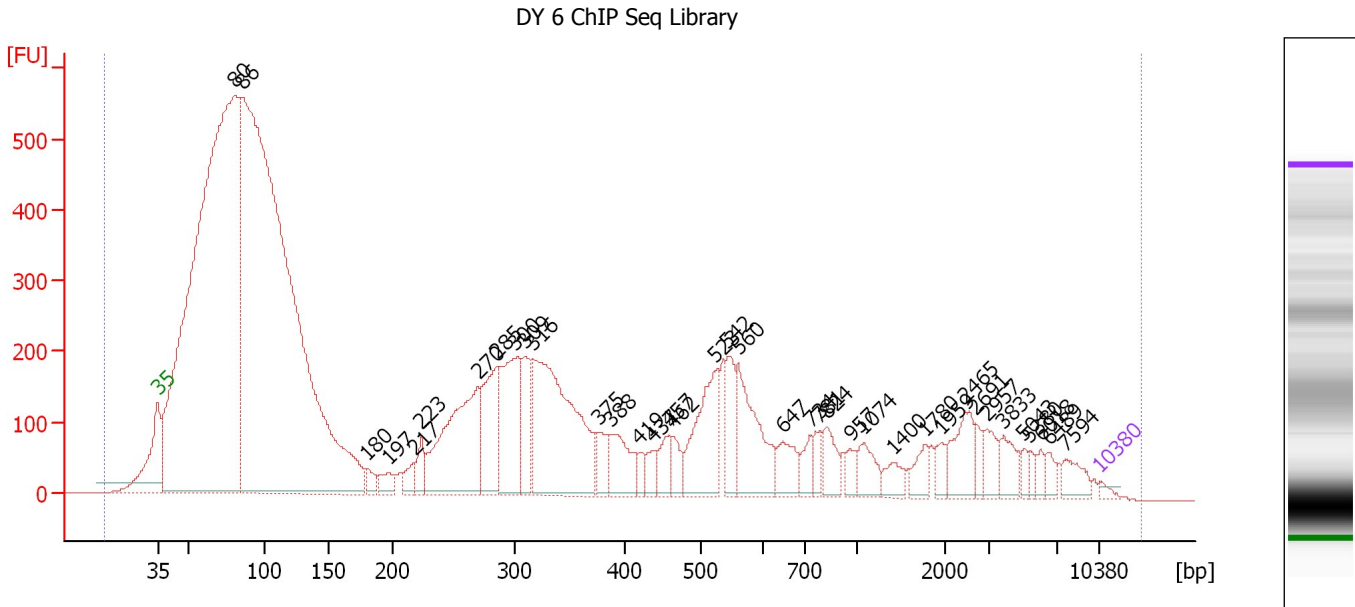

**Overall Results for sample 2 : DY 6 ChIP Seq Library**

Number of peaks found: 38                      Corr. Area 1: 17,165.1  
 Noise: 0.2

**Peak table for sample 2 : DY 6 ChIP Seq Library**

Peak	Size [bp]	Conc. [pg/μl]	Molarity [pmol/l]	Observations	Aligned Migration Time [s]
1	35	125.00	5,411.3	Lower Marker	43.00
2	80	33,139.42	626,553.4		48.74
3	86	35,877.08	635,777.6		49.35
4	180	305.61	2,575.6		58.55
5	197	361.48	2,778.9		60.19
6	217	357.31	2,491.2		62.04
7	223	475.76	3,232.6		62.56
8	270	3,754.08	21,081.0		66.84
9	285	1,855.55	9,866.4		68.22
10	300	2,474.50	12,513.6		69.56
11	309	1,111.62	5,448.5		70.33
12	316	5,213.74	24,996.6		70.90
13	375	532.53	2,148.8		75.73
14	388	1,034.39	4,042.4		76.73
15	419	227.14	821.3		78.80
16	437	347.99	1,205.4		79.84
17	457	502.45	1,664.6		80.96
18	462	369.85	1,213.1		81.22
19	523	1,868.59	5,414.8		84.41
20	542	860.21	2,406.2		85.28
21	560	1,798.20	4,861.0		86.14
22	647	632.64	1,482.2		89.42
23	734	363.15	749.6		91.54
24	781	245.41	476.3		92.14
25	824	517.26	951.3		92.70

Assay Class: High Sensitivity DNA Assay  
 Data Path: C:\...ings\Bioanalyzer\2016-02-25\2016-02-25\_001\_Dag\_DY3\_DY\_6.xad

Created: 2/25/2016 2:12:58 PM  
 Modified: 2/25/2016 3:32:15 PM

### Electropherogram Summary Continued ...

#### ... Peak table for sample 2 : DY 6 ChIP Seq Library

Peak	Size [bp]	Conc. [pg/ $\mu$ l]	Molarity [pmol/l]	Observations	Aligned Migration Time [s]
26	957	246.20	389.8		94.43
27	1,074	474.30	669.2		95.47
28	1,400	289.95	313.8		97.58
29	1,780	314.69	267.9		100.05
30	1,959	230.63	178.3		101.21
31	2,465	709.22	436.0		102.98
32	2,691	207.66	116.9		103.72
33	2,957	361.37	185.1		104.58
34	3,833	394.16	155.8		105.79
35	5,043	126.84	38.1		107.34
36	5,380	122.05	34.4		107.77
37	6,018	149.23	37.6		108.60
38	6,489	169.06	39.5		109.20
39	7,594	324.50	64.7		110.41
40	10,380	75.00	10.9	Upper Marker	113.00


#### Region table for sample 2 : DY 6 ChIP Seq Library

From [bp]	To [bp]	Average Size [bp]	Corr. Area	Molarity [pmol/l]	Co Conc. lor [pg/ $\mu$ l]	% of Total	Size distribution in CV [%]
9	13,655	671	17,165.1	1,439,926.0	99,463.78	100	100.0

Assay Class: High Sensitivity DNA Assay  
Data Path: C:\...ings\Bioanalyzer\2016-02-25\2016-02-25\_001\_Dag\_DY3\_DY\_6.xad

Created: 2/25/2016 2:12:58 PM  
Modified: 2/25/2016 3:32:15 PM

**Gel Image**

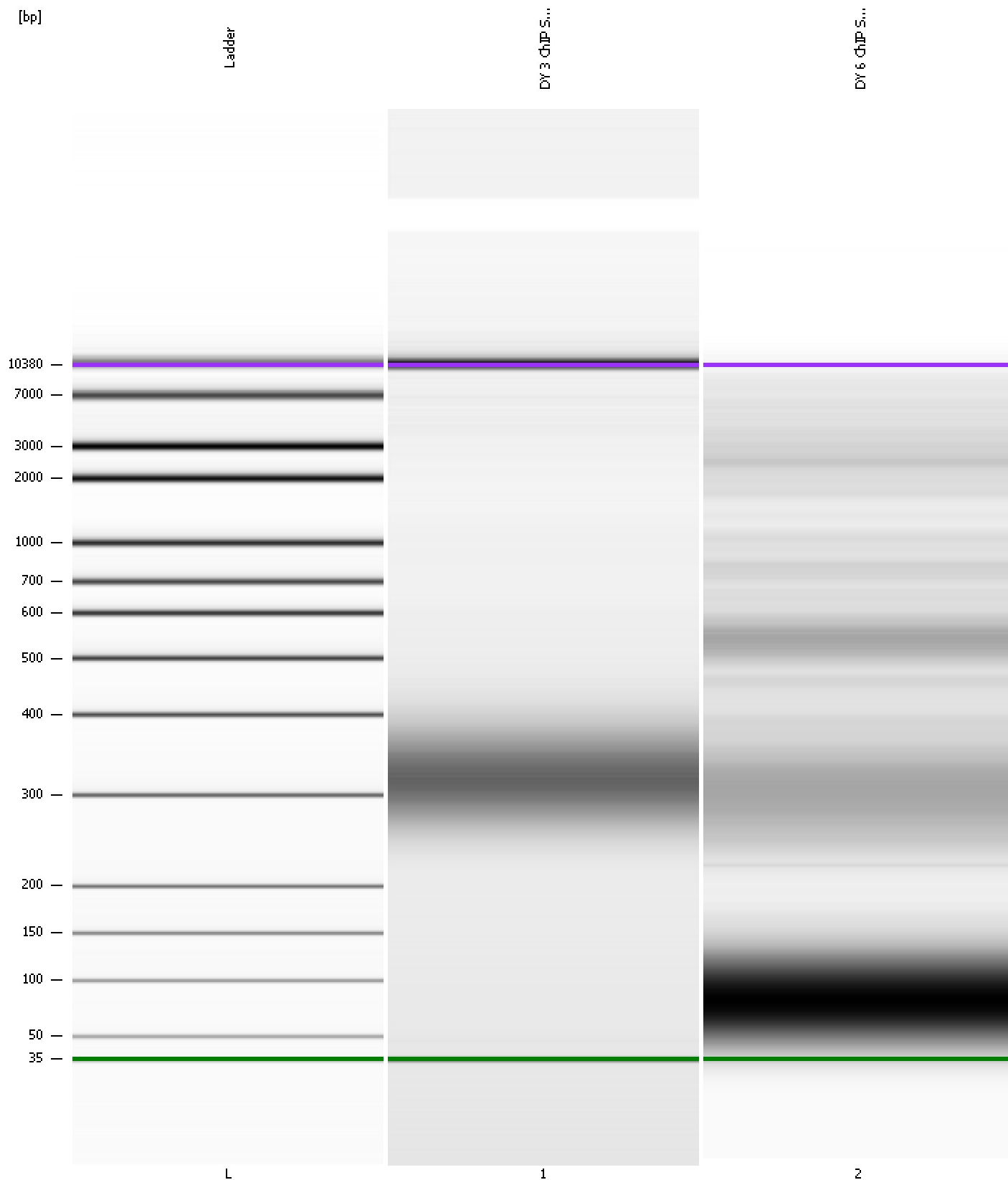



Assay Class: High Sensitivity DNA Assay  
 Data Path: C:\...ings\Bioanalyzer\2016-02-25\2016-02-25\_001\_Dag\_DY3\_DY\_6.xad

Created: 2/25/2016 2:12:58 PM  
 Modified: 2/25/2016 3:32:15 PM

**Run Logbook**

Description	Number	Source	Category	Sub Category	Time	Time Zone	User	Host
Run ended on port 1 (Number of wells acquired: 8)		Instrument	Run		2/25/2016 2:42:50 PM	(GMT --08:00) Pacific Standard Time	UC Davis	D8XSMGH1
Run started on port 1 (File: C:\Documents and Settings\Bioanalyzer\2016-02-25\2016-02-25_001.xad)		Instrument	Run		2/25/2016 2:13:03 PM	(GMT --08:00) Pacific Standard Time	UC Davis	D8XSMGH1
Product Number : G2938B		Instrument	Run		2/25/2016 2:13:03 PM	(GMT --08:00) Pacific Standard Time	UC Davis	D8XSMGH1
Name :		Instrument	Run		2/25/2016 2:13:03 PM	(GMT --08:00) Pacific Standard Time	UC Davis	D8XSMGH1
Vendor : Agilent Technologies		Instrument	Run		2/25/2016 2:13:03 PM	(GMT --08:00) Pacific Standard Time	UC Davis	D8XSMGH1
Serial# : DE13701086		Instrument	Run		2/25/2016 2:13:03 PM	(GMT --08:00) Pacific Standard Time	UC Davis	D8XSMGH1
Firmware : C.01.069		Instrument	Run		2/25/2016 2:13:03 PM	(GMT --08:00) Pacific Standard Time	UC Davis	D8XSMGH1
Cartridge : Electrode		Instrument	Run		2/25/2016 2:13:03 PM	(GMT --08:00) Pacific Standard Time	UC Davis	D8XSMGH1