


Assay Class: High Sensitivity DNA Assay
Data Path: C:\...data\2016-12-09\2016-12-09_002_HiSeq749_Libraries_72-74.xad

Created: 12/9/2016 11:02:32 AM
Modified: 12/9/2016 11:51:26 AM

Electrophoresis File Run Summary

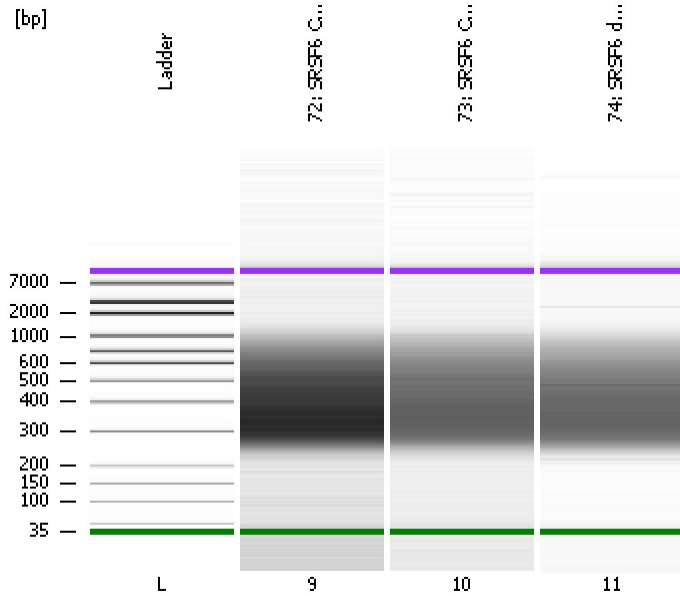

Instrument Information:




Instrument Name: DE13701086 Firmware: C.01.069
Serial#: DE13701086 Type: G2938B

Assay Information:

Assay Origin Path: C:\Program Files\Agilent\2100 bioanalyzer\2100 expert\assays\dsDNA\High Sensitivity DNA.xsy
Assay Class: High Sensitivity DNA Assay
Version: 1.03
Assay Comments: Copyright © 2003-2010 Agilent Technologies

Chip Information:

Chip Lot #:
Reagent Kit Lot #:
Chip Comments:

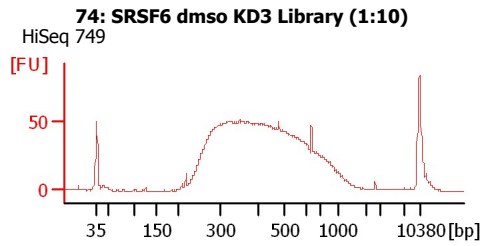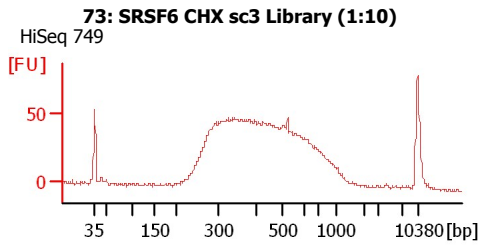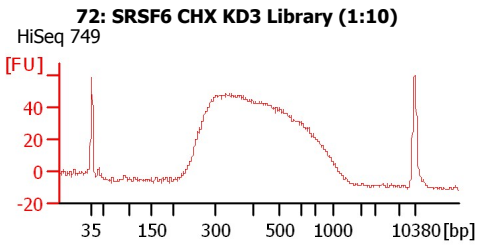

Assay Class: High Sensitivity DNA Assay
 Data Path: C:\...data\2016-12-09\2016-12-09_002_HiSeq749_Libraries_72-74.xad

Created: 12/9/2016 11:02:32 AM
 Modified: 12/9/2016 11:51:26 AM

Electrophoresis File Run Summary (Chip Summary)

Sample Name	Sample Comment	Rest. Digest	Status	Observation	Result Label	Result Color
72: SRSF6 CHX KD3 Library (1:10)	HiSeq 749	<input type="checkbox"/>	✓			
73: SRSF6 CHX sc3 Library (1:10)	HiSeq 749	<input type="checkbox"/>	✓			
74: SRSF6 dms0 KD3 Library (1:10)	HiSeq 749	<input type="checkbox"/>	✓			
Ladder		<input type="checkbox"/>	✓			

Chip Lot #

Reagent Kit Lot #

Chip Comments :

Assay Class: High Sensitivity DNA Assay
Data Path: C:\...data\2016-12-09\2016-12-09_002_HiSeq749_Libraries_72-74.xad

Created: 12/9/2016 11:02:32 AM
Modified: 12/9/2016 11:51:26 AM

Electrophoresis Assay Details

General Analysis Settings

Number of Available Sample and Ladder Wells (Max.) : 12
Minimum Visible Range [s] : 32
Maximum Visible Range [s] : 138
Start Analysis Time Range [s] : 33
End Analysis Time Range [s] : 137.5
Ladder Concentration [pg/μl] : 1950
Uses Standard Area for Ladder Fragments
Lower Marker Concentration [pg/μl] : 125
Upper Marker Concentration [pg/μl] : 75
Used Upper Marker for Quantitation
Standard Curve Fit is Point to Point
Show Data Aligned to Lower and Upper Marker

Integrator Settings

Integration Start Time [s] : 33.05
Integration End Time [s] : 137
Slope Threshold : 0.8
Height Threshold [FU] : 5
Area Threshold : 0.1
Width Threshold [s] : 0.6
Baseline Plateau [s] : 0.5

Filter Settings

Filter Width [s] : 0.5
Polynomial Order : 4


Ladder

Ladder Peak	Size	Area
1	35	160
2	50	210
3	100	208
4	150	221
5	200	242
6	300	270
7	400	305
8	500	306
9	600	336
10	700	321
11	1000	366
12	2000	413
13	3000	411
14	7000	400
15	10380	214

Assay Class: High Sensitivity DNA Assay
 Data Path: C:\...data\2016-12-09\2016-12-09_002_HiSeq749_Libraries_72-74.xad

Created: 12/9/2016 11:02:32 AM
 Modified: 12/9/2016 11:51:26 AM

Electropherogram Summary

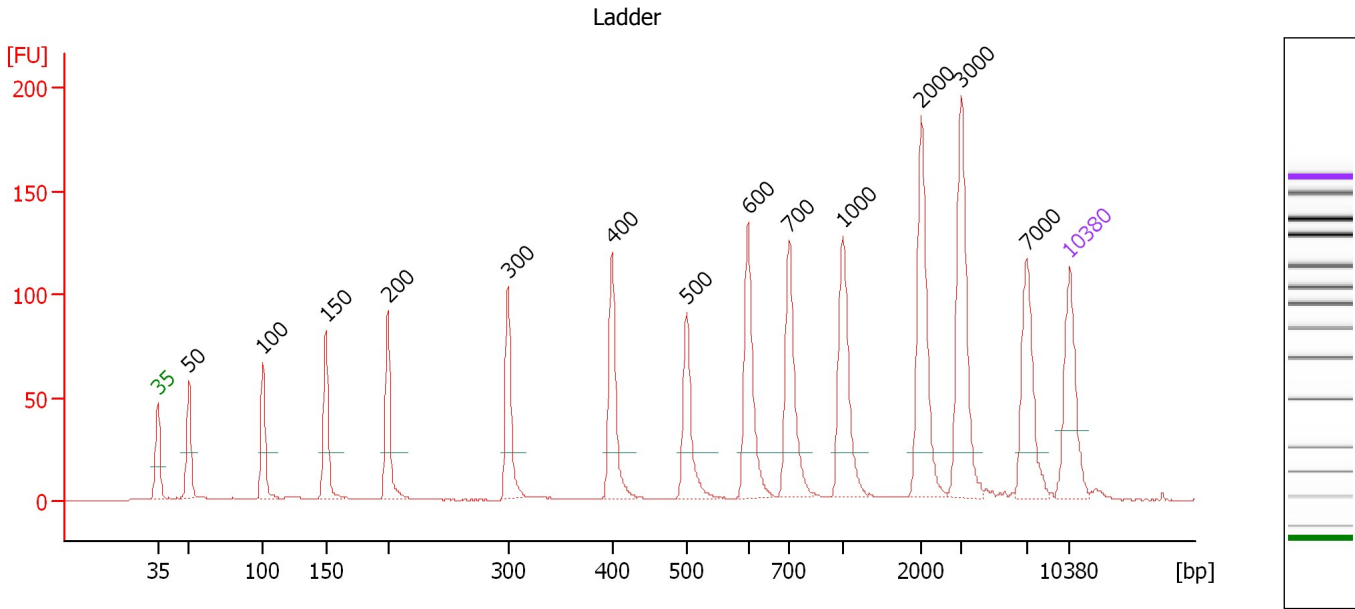

Overall Results for Ladder

Noise: 0.2

Peak table for Ladder


Peak	Size [bp]	Conc. [pg/μl]	Molarity [pmol/l]	Observations	Aligned Migration Time [s]
1	35	125.00	5,411.3	Lower Marker	43.00
2	50	150.00	4,545.5	Ladder Peak	45.40
3	100	150.00	2,272.7	Ladder Peak	51.12
4	150	150.00	1,515.2	Ladder Peak	55.91
5	200	150.00	1,136.4	Ladder Peak	60.66
6	300	150.00	757.6	Ladder Peak	69.90
7	400	150.00	568.2	Ladder Peak	77.88
8	500	150.00	454.5	Ladder Peak	83.60
9	600	150.00	378.8	Ladder Peak	88.30
10	700	150.00	324.7	Ladder Peak	91.48
11	1,000	150.00	227.3	Ladder Peak	95.59
12	2,000	150.00	113.6	Ladder Peak	101.60
13	3,000	150.00	75.8	Ladder Peak	104.64
14	7,000	150.00	32.5	Ladder Peak	109.67
15	10,380	75.00	10.9	Upper Marker	113.00

Assay Class: High Sensitivity DNA Assay
 Data Path: C:\...data\2016-12-09\2016-12-09_002_HiSeq749_Libraries_72-74.xad

Created: 12/9/2016 11:02:32 AM
 Modified: 12/9/2016 11:51:26 AM

Electropherogram Summary Continued ...

72: SRSF6 CHX KD3 Library (1:10) [HiSeq 749]

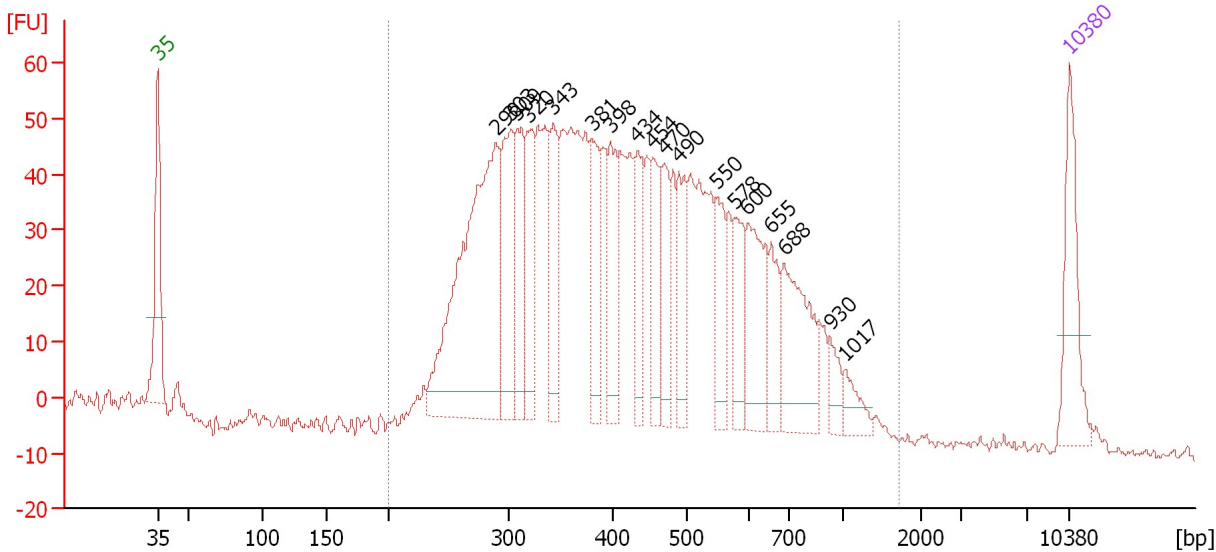

Overall Results for sample 9 : 72: SRSF6 CHX KD3 Library (1:10)

Number of peaks found: 18 Corr. Area 1: 1,592.6
 Noise: 1.2

Peak table for sample 9 : 72: SRSF6 CHX KD3 Library (1:10)

Peak	Size [bp]	Conc. [pg/μl]	Molarity [pmol/l]	Observations	Aligned Migration Time [s]
1	35	125.00	5,411.3	Lower Marker	43.00
2	290	442.33	2,314.1		68.94
3	303	138.49	693.3		70.12
4	309	104.58	513.4		70.60
5	320	88.49	419.0		71.50
6	343	91.09	401.9		73.37
7	381	80.54	320.4		76.35
8	398	86.56	329.4		77.73
9	434	67.24	234.9		79.81
10	454	64.97	216.7		80.98
11	470	62.80	202.4		81.89
12	490	68.67	212.5		83.01
13	550	57.15	157.5		85.94
14	578	57.57	150.9		87.27
15	600	95.46	241.2		88.28
16	655	49.59	114.7		90.04
17	688	122.40	269.4		91.11
18	930	21.01	34.3		94.62
19	1,017	23.59	35.2		95.69
20	10,380	75.00	10.9	Upper Marker	113.00

Region table for sample 9 : 72: SRSF6 CHX KD3 Library (1:10)


From [bp]	To [bp]	Average Size [bp]	Corr. Area	Molarity [pmol/l]	Co Conc. [pg/μl]	% of Total	Size distribution in CV [%]
200	1,714	456	1,592.6	10,501.9	2,683.91	99	39.0

Assay Class: High Sensitivity DNA Assay
 Data Path: C:\...data\2016-12-09\2016-12-09_002_HiSeq749_Libraries_72-74.xad

Created: 12/9/2016 11:02:32 AM
 Modified: 12/9/2016 11:51:26 AM

Electropherogram Summary Continued ...

73: SRSF6 CHX sc3 Library (1:10) [HiSeq 749]

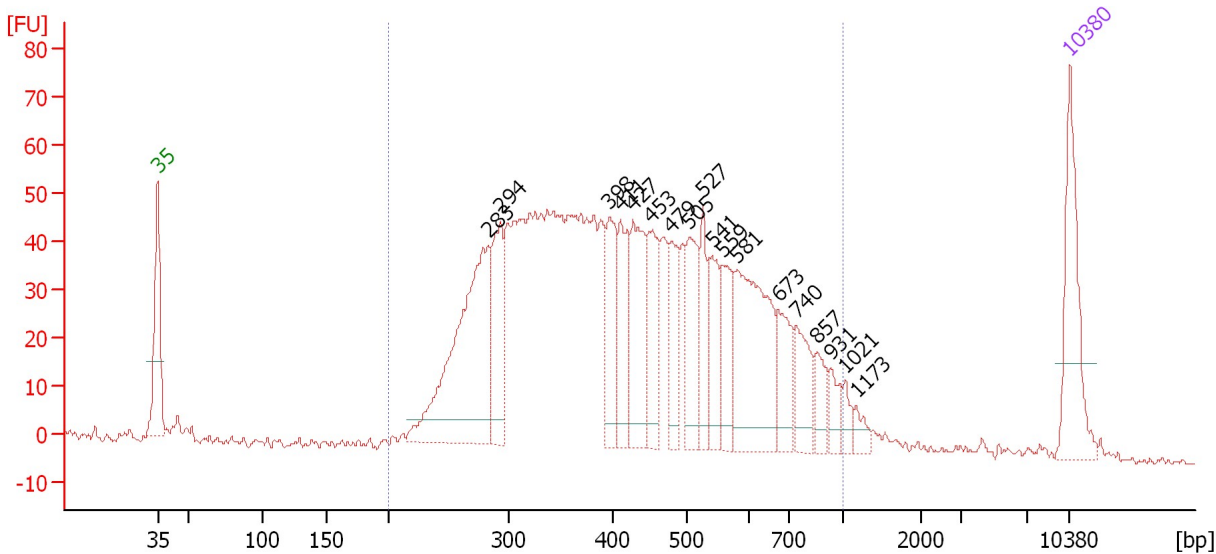

Overall Results for sample 10 : 73: SRSF6 CHX sc3 Library (1:10)

Number of peaks found: 18 Corr. Area 1: 1,457.0
 Noise: 0.9

Peak table for sample 10 : 73: SRSF6 CHX sc3 Library (1:10)

Peak	Size [bp]	Conc. [pg/μl]	Molarity [pmol/l]	Observations	Aligned Migration Time [s]
1	35	125.00	5,411.3	Lower Marker	43.00
2	283	274.60	1,468.3		68.37
3	294	95.55	492.8		69.33
4	398	67.75	258.3		77.68
5	411	69.20	254.8		78.54
6	427	100.25	355.4		79.44
7	453	61.83	207.0		80.89
8	479	50.73	160.5		82.39
9	505	68.24	204.8		83.83
10	527	49.39	142.1		84.85
11	541	52.68	147.4		85.55
12	559	52.41	142.2		86.35
13	581	154.56	402.9		87.42
14	673	44.74	100.7		90.63
15	740	47.02	96.3		92.02
16	857	24.51	43.3		93.63
17	931	18.83	30.6		94.64
18	1,021	11.92	17.7		95.71
19	1,173	10.76	13.9		96.62
20	10,380	75.00	10.9	Upper Marker	113.00

Region table for sample 10 : 73: SRSF6 CHX sc3 Library (1:10)


From [bp]	To [bp]	Average Size [bp]	Corr. Area	Molarity [pmol/l]	Co Conc. [pg/μl]	% of Total	Size distribution in CV [%]
200	1,000	449	1,457.0	7,901.1	2,016.49	97	34.6

Assay Class: High Sensitivity DNA Assay
 Data Path: C:\...data\2016-12-09\2016-12-09_002_HiSeq749_Libraries_72-74.xad

Created: 12/9/2016 11:02:32 AM
 Modified: 12/9/2016 11:51:26 AM

Electropherogram Summary Continued ...

74: SRSF6 dms0 KD3 Library (1:10) [HiSeq 749]

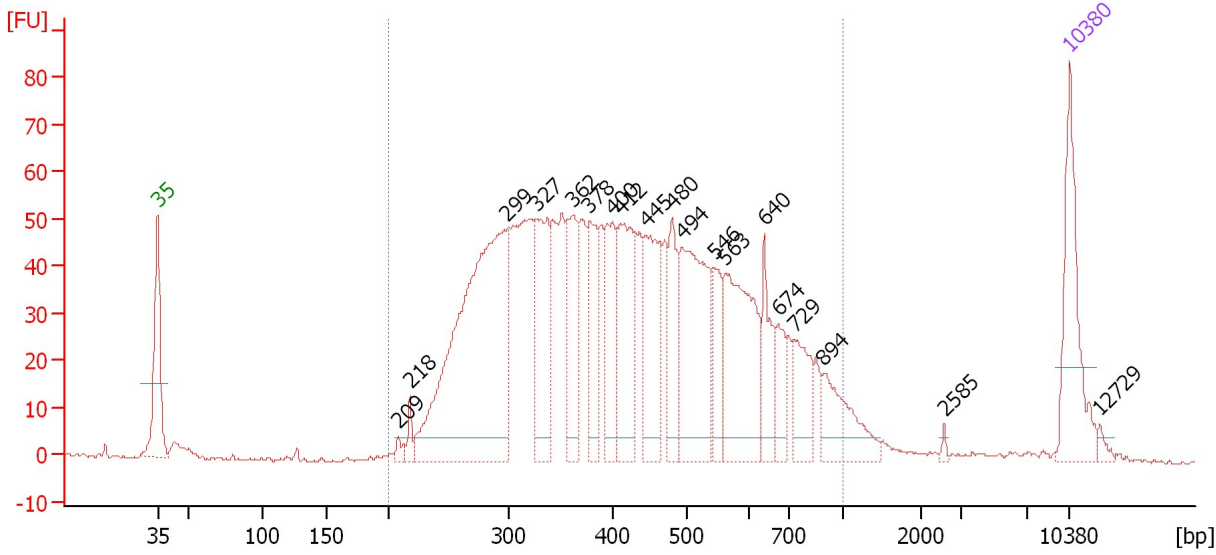

Overall Results for sample 11 : 74: SRSF6 dms0 KD3 Library (1:10)

Number of peaks found: 19 Corr. Area 1: 1,560.8
 Noise: 0.2

Peak table for sample 11 : 74: SRSF6 dms0 KD3 Library (1:10)

Peak	Size [bp]	Conc. [pg/μl]	Molarity [pmol/l]	Observations	Aligned Migration Time [s]
1	35	125.00	5,411.3	Lower Marker	43.00
2	209	7.45	54.0		61.49
3	218	12.83	89.0		62.35
4	299	428.92	2,172.2		69.83
5	327	117.92	546.9		72.03
6	362	76.26	318.8		74.88
7	378	67.20	269.4		76.12
8	400	72.45	274.3		77.89
9	412	108.84	400.7		78.54
10	445	93.32	317.5		80.47
11	480	57.15	180.4		82.46
12	494	161.29	494.5		83.27
13	546	38.78	107.7		85.74
14	563	135.57	365.0		86.55
15	640	46.05	109.1		89.56
16	674	28.81	64.8		90.63
17	729	45.34	94.3		91.87
18	894	56.00	94.9		94.13
19	2,585	2.52	1.5		103.38
20	10,380	75.00	10.9	Upper Marker	113.00
21	12,729	0.00	0.0		115.31


Region table for sample 11 : 74: SRSF6 dms0 KD3 Library (1:10)

From [bp]	To [bp]	Average Size [bp]	Corr. Area	Molarity [pmol/l]	Co Conc. [pg/μl]	% of Total	Size distribution in CV [%]
200	1,000	442	1,560.8	8,170.3	2,049.53	96	34.9

Assay Class: High Sensitivity DNA Assay
Data Path: C:\...data\2016-12-09\2016-12-09_002_HiSeq749_Libraries_72-74.xad

Created: 12/9/2016 11:02:32 AM
Modified: 12/9/2016 11:51:26 AM

Gel Image

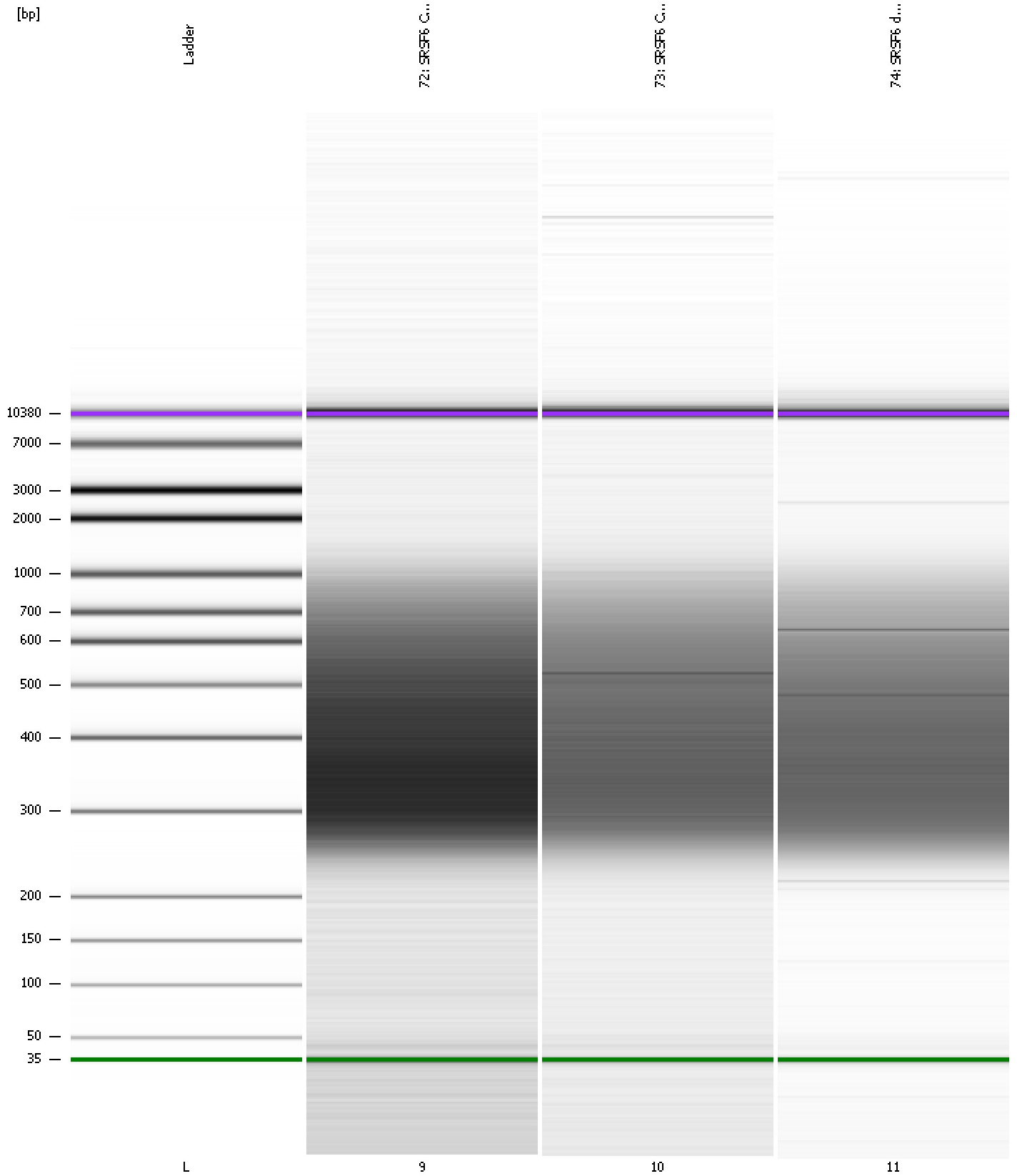

Assay Class: High Sensitivity DNA Assay
 Data Path: C:\...data\2016-12-09\2016-12-09_002_HiSeq749_Libraries_72-74.xad

Created: 12/9/2016 11:02:32 AM
 Modified: 12/9/2016 11:51:26 AM

Run Logbook

Description	Number	Source	Category	Sub Category	Time	Time Zone	User	Host
Run ended on port 1 (Number of wells acquired: 12)		Instrument	Run		12/9/2016 11:43:49 AM	(GMT --08:00) Pacific Standard Time	UC Davis	D8XSMGH1
Run started on port 1 (File: C:\Program Files\Agilent\2100 bioanalyzer\2100 expert\data\2016-12-09\2016-12-09_002.xad)		Instrument	Run		12/9/2016 11:02:38 AM	(GMT --08:00) Pacific Standard Time	UC Davis	D8XSMGH1
Product Number : G2938B		Instrument	Run		12/9/2016 11:02:38 AM	(GMT --08:00) Pacific Standard Time	UC Davis	D8XSMGH1
Name :		Instrument	Run		12/9/2016 11:02:38 AM	(GMT --08:00) Pacific Standard Time	UC Davis	D8XSMGH1
Vendor : Agilent Technologies		Instrument	Run		12/9/2016 11:02:38 AM	(GMT --08:00) Pacific Standard Time	UC Davis	D8XSMGH1
Serial# : DE13701086		Instrument	Run		12/9/2016 11:02:38 AM	(GMT --08:00) Pacific Standard Time	UC Davis	D8XSMGH1
Firmware : C.01.069		Instrument	Run		12/9/2016 11:02:38 AM	(GMT --08:00) Pacific Standard Time	UC Davis	D8XSMGH1
Cartridge : Electrode		Instrument	Run		12/9/2016 11:02:38 AM	(GMT --08:00) Pacific Standard Time	UC Davis	D8XSMGH1