


Assay Class: High Sensitivity DNA Assay  
Data Path: C:\... bioanalyzer\2100 expert\data\2018-05-03\2018-05-03\_001.xad

Created: 5/3/2018 2:55:51 PM  
Modified: 5/3/2018 3:49:23 PM

### Electrophoresis File Run Summary

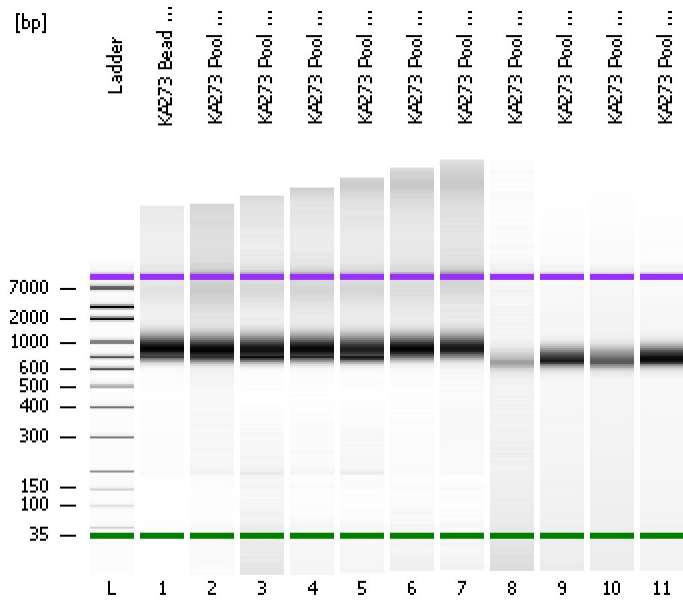

#### Instrument Information:

Instrument Name: DE34903152      Firmware: C.01.069  
Serial#: DE34903152      Type: G2938C

#### Assay Information:

Assay Origin Path: C:\Program Files\Agilent\2100 bioanalyzer\2100 expert\assays\dsDNA\High Sensitivity DNA.xsy  
Assay Class: High Sensitivity DNA Assay  
Version: 1.03  
Assay Comments: Copyright © 2003-2010 Agilent Technologies




#### Chip Information:

Chip Lot #:   
Reagent Kit Lot #:   
Chip Comments:

KA273 Bead Clean Pool

KA273 Pool 1-47




KA273 Pool 48-94

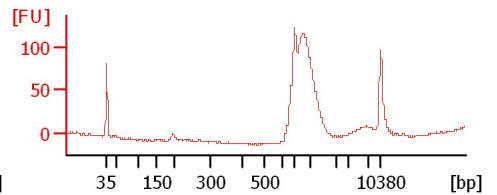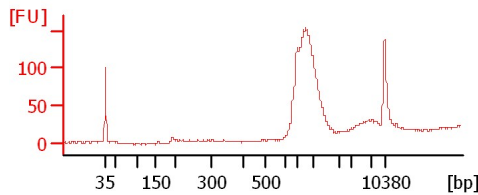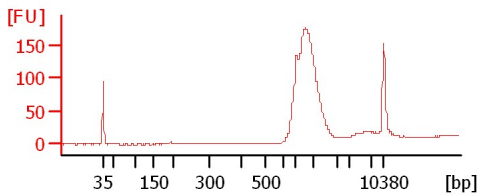

KA273 Pool 95-141

KA273 Pool 142-188




KA273 Pool 189-235

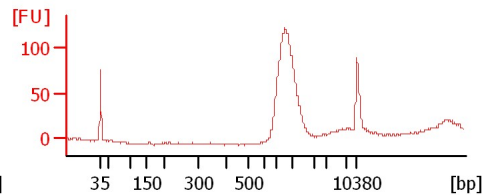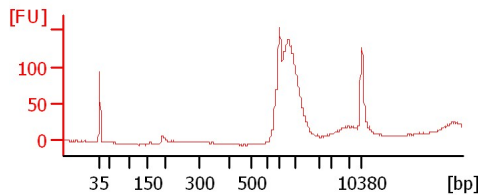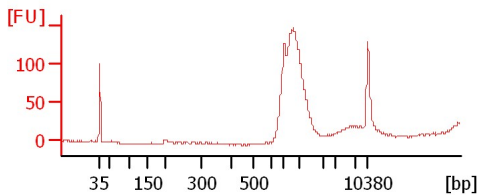

KA273 Pool 236-273



KA273 Pool Pippin 1

KA273 Pool Pippin 2

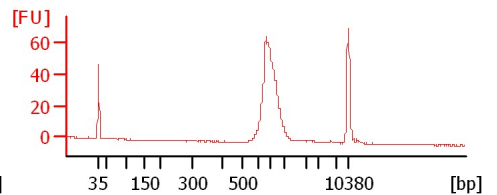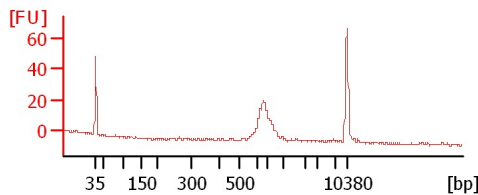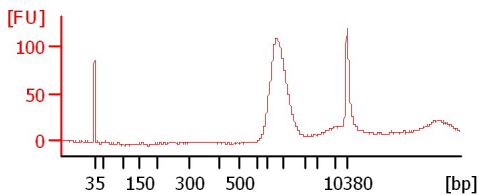

KA273 Pool Pippin 3

KA273 Pool Pippin 4

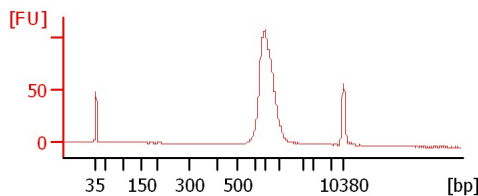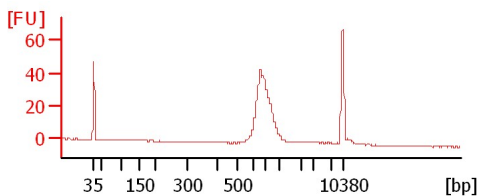

Assay Class: High Sensitivity DNA Assay  
 Data Path: C:\... bioanalyzer\2100 expert\data\2018-05-03\2018-05-03\_001.xad

Created: 5/3/2018 2:55:51 PM  
 Modified: 5/3/2018 3:49:23 PM

**Electrophoresis File Run Summary (Chip Summary)**

Sample Name	Sample Comment	Rest. Digest	Status	Observation	Result Label	Result Color
KA273 Bead Clean Pool		<input type="checkbox"/>	✓			
KA273 Pool 1-47		<input type="checkbox"/>	✓			
KA273 Pool 48-94		<input type="checkbox"/>	✓			
KA273 Pool 95-141		<input type="checkbox"/>	✓			
KA273 Pool 142-188		<input type="checkbox"/>	✓			
KA273 Pool 189-235		<input type="checkbox"/>	✓			
KA273 Pool 236-273		<input type="checkbox"/>	✓			
KA273 Pool Pippin 1		<input type="checkbox"/>	✓			
KA273 Pool Pippin 2		<input type="checkbox"/>	✓			
KA273 Pool Pippin 3		<input type="checkbox"/>	✓			
KA273 Pool Pippin 4		<input type="checkbox"/>	✓			
Ladder		<input type="checkbox"/>	✓			

**Chip Lot #**

**Reagent Kit Lot #**

**Chip Comments :**

Assay Class: High Sensitivity DNA Assay  
Data Path: C:\... bioanalyzer\2100 expert\data\2018-05-03\2018-05-03\_001.xad

Created: 5/3/2018 2:55:51 PM  
Modified: 5/3/2018 3:49:23 PM

**Electrophoresis Assay Details**

**General Analysis Settings**

Number of Available Sample and Ladder Wells (Max.) : 12  
Minimum Visible Range [s] : 32  
Maximum Visible Range [s] : 138  
Start Analysis Time Range [s] : 33  
End Analysis Time Range [s] : 137.5  
Ladder Concentration [pg/μl] : 1950  
Uses Standard Area for Ladder Fragments  
Lower Marker Concentration [pg/μl] : 125  
Upper Marker Concentration [pg/μl] : 75  
Used Upper Marker for Quantitation  
Standard Curve Fit is Point to Point  
Show Data Aligned to Lower and Upper Marker

**Integrator Settings**

Integration Start Time [s] : 33.05  
Integration End Time [s] : 137  
Slope Threshold : 0.8  
Height Threshold [FU] : 5  
Area Threshold : 0.1  
Width Threshold [s] : 0.6  
Baseline Plateau [s] : 0.5

**Filter Settings**

Filter Width [s] : 0.5  
Polynomial Order : 4


**Ladder**

Ladder Peak	Size	Area
1	35	160
2	50	210
3	100	208
4	150	221
5	200	242
6	300	270
7	400	305
8	500	306
9	600	336
10	700	321
11	1000	366
12	2000	413
13	3000	411
14	7000	400
15	10380	214

Assay Class: High Sensitivity DNA Assay  
 Data Path: C:\... bioanalyzer\2100 expert\data\2018-05-03\2018-05-03\_001.xad

Created: 5/3/2018 2:55:51 PM  
 Modified: 5/3/2018 3:49:23 PM

**Electropherogram Summary**

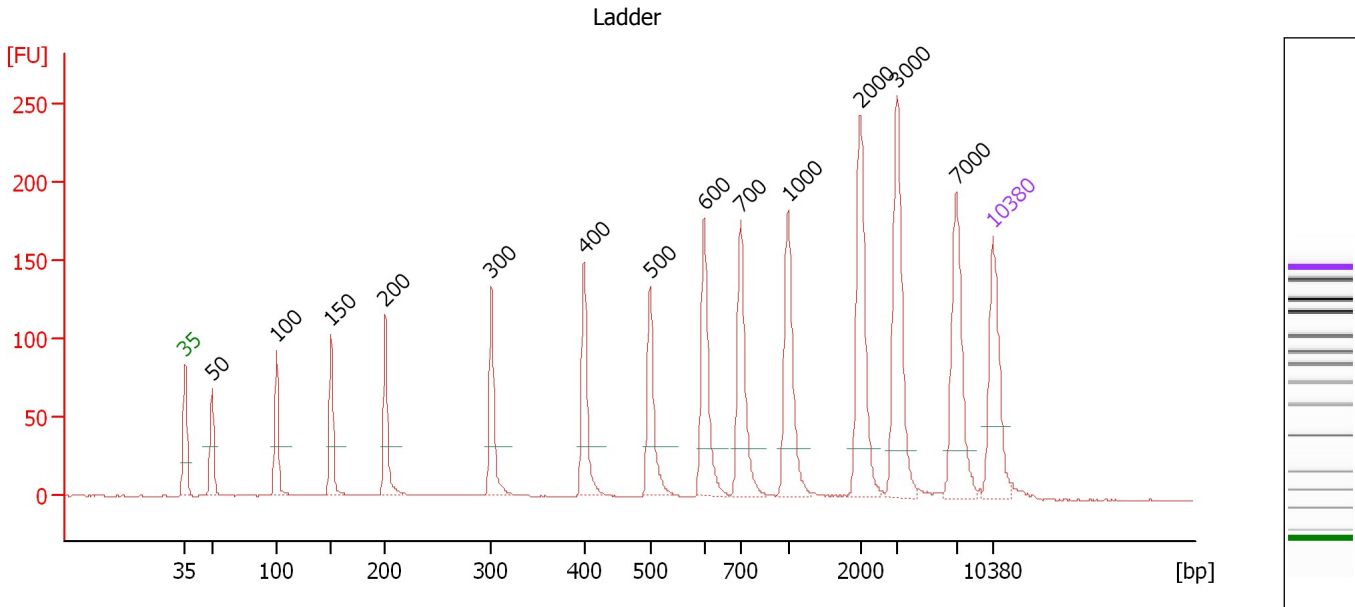

**Overall Results for Ladder**

Noise: 0.3


**Peak table for Ladder**

Peak	Size [bp]	Conc. [pg/μl]	Molarity [pmol/l]	Observations	Aligned Migration Time [s]
1	35	125.00	5,411.3	Lower Marker	43.00
2	50	150.00	4,545.5	Ladder Peak	45.36
3	100	150.00	2,272.7	Ladder Peak	50.95
4	150	150.00	1,515.2	Ladder Peak	55.66
5	200	150.00	1,136.4	Ladder Peak	60.37
6	300	150.00	757.6	Ladder Peak	69.57
7	400	150.00	568.2	Ladder Peak	77.61
8	500	150.00	454.5	Ladder Peak	83.29
9	600	150.00	378.8	Ladder Peak	88.00
10	700	150.00	324.7	Ladder Peak	91.15
11	1,000	150.00	227.3	Ladder Peak	95.26
12	2,000	150.00	113.6	Ladder Peak	101.54
13	3,000	150.00	75.8	Ladder Peak	104.68
14	7,000	150.00	32.5	Ladder Peak	109.86
15	10,380	75.00	10.9	Upper Marker	113.00

Assay Class: High Sensitivity DNA Assay  
 Data Path: C:\... bioanalyzer\2100 expert\data\2018-05-03\2018-05-03\_001.xad

Created: 5/3/2018 2:55:51 PM  
 Modified: 5/3/2018 3:49:23 PM

**Electropherogram Summary Continued ...**

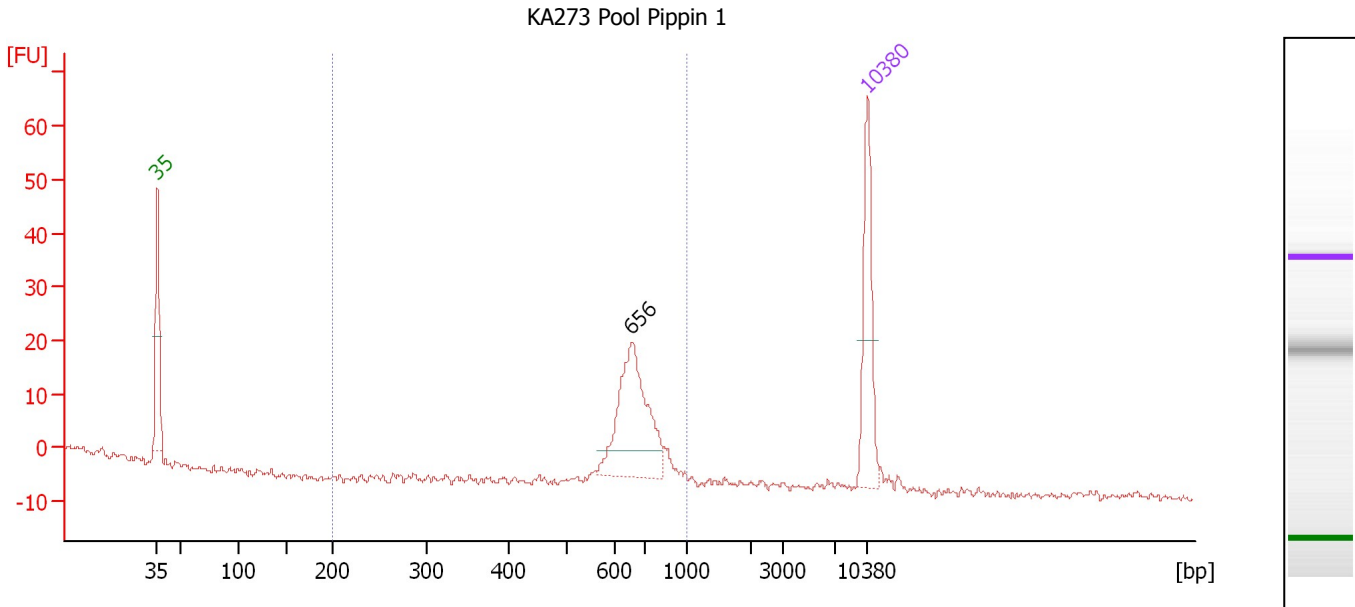

**Overall Results for sample 8 : KA273 Pool Pippin 1**

Number of peaks found: 1                      Corr. Area 1: 91.8  
 Noise: 0.6

**Peak table for sample 8 : KA273 Pool Pippin 1**

Peak	Size [bp]	Conc. [pg/μl]	Molarity [pmol/l]	Observations	Aligned Migration Time [s]
1	35	125.00	5,411.3	Lower Marker	43.00
2	656	146.44	338.1		89.77
3	10,380	75.00	10.9	Upper Marker	113.00


**Region table for sample 8 : KA273 Pool Pippin 1**

From [bp]	To [bp]	Average Size [bp]	Conc. [pg/μl]	Corr. Area	Molarity [pmol/l]	Co % of lor Total	Size distribution in CV [%]
200	1,000	680	145.16	91.8	326.1	96	10.3

Assay Class: High Sensitivity DNA Assay  
 Data Path: C:\... bioanalyzer\2100 expert\data\2018-05-03\2018-05-03\_001.xad

Created: 5/3/2018 2:55:51 PM  
 Modified: 5/3/2018 3:49:23 PM

**Electropherogram Summary Continued ...**

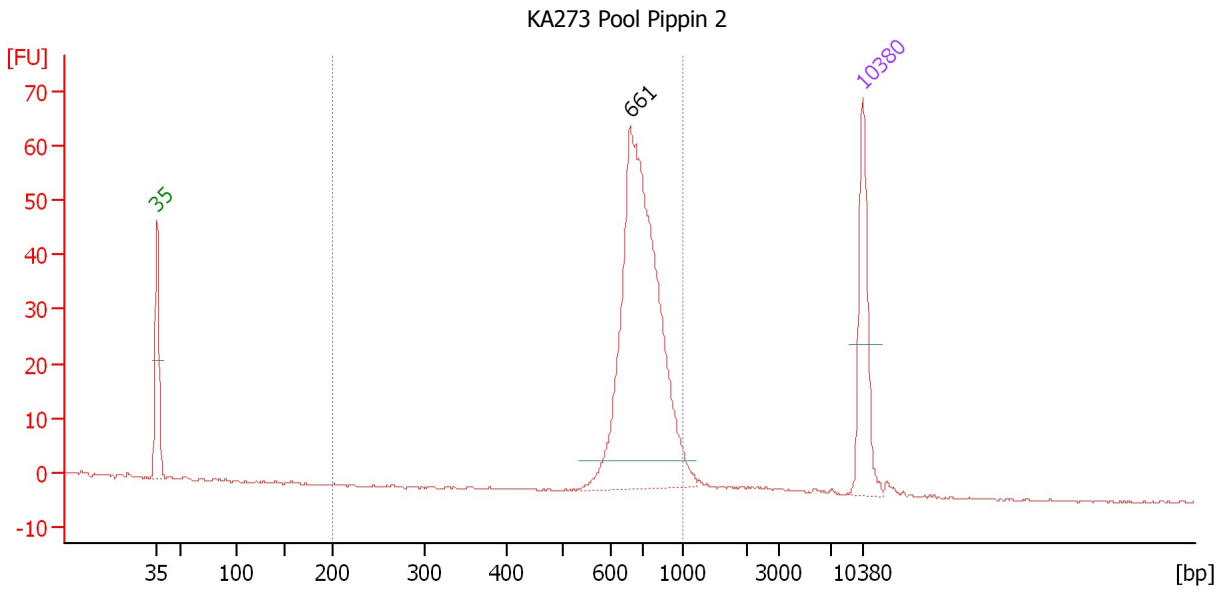

**Overall Results for sample 9 : KA273 Pool Pippin 2**

Number of peaks found: 1                      Corr. Area 1: 299.1  
 Noise: 0.2

**Peak table for sample 9 : KA273 Pool Pippin 2**

Peak	Size [bp]	Conc. [pg/μl]	Molarity [pmol/l]	Observations	Aligned Migration Time [s]
1	35	125.00	5,411.3	Lower Marker	43.00
2	661	456.00	1,044.7		89.93
3	10,380	75.00	10.9	Upper Marker	113.00


**Region table for sample 9 : KA273 Pool Pippin 2**

From [bp]	To [bp]	Average Size [bp]	Conc. [pg/μl]	Corr. Area	Molarity [pmol/l]	Co % of lor Total	Size distribution in CV [%]
200	1,000	723	434.93	299.1	924.1	95	12.7

Assay Class: High Sensitivity DNA Assay  
 Data Path: C:\... bioanalyzer\2100 expert\data\2018-05-03\2018-05-03\_001.xad

Created: 5/3/2018 2:55:51 PM  
 Modified: 5/3/2018 3:49:23 PM

**Electropherogram Summary Continued ...**

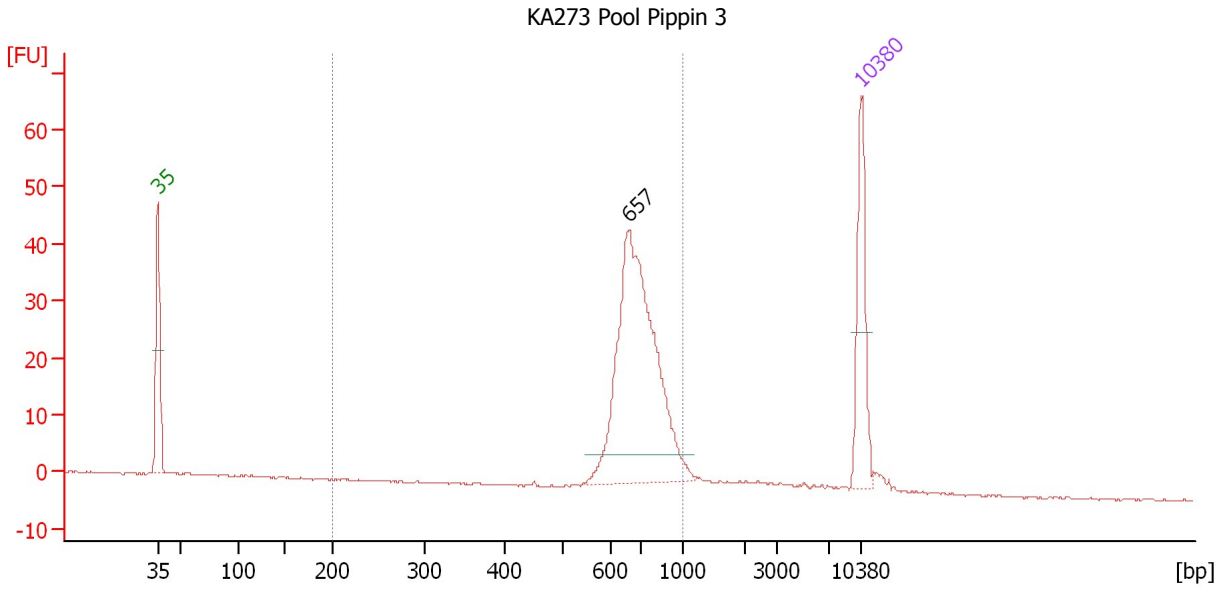

**Overall Results for sample 10 : KA273 Pool Pippin 3**

Number of peaks found: 1                      Corr. Area 1: 214.1  
 Noise: 0.1

**Peak table for sample 10 : KA273 Pool Pippin 3**

Peak	Size [bp]	Conc. [pg/μl]	Molarity [pmol/l]	Observations	Aligned Migration Time [s]
1	35	125.00	5,411.3	Lower Marker	43.00
2	657	357.14	823.1		89.81
3	10,380	75.00	10.9	Upper Marker	113.00



**Region table for sample 10 : KA273 Pool Pippin 3**

From [bp]	To [bp]	Average Size [bp]	Conc. [pg/μl]	Corr. Area	Molarity [pmol/l]	Co % of lor Total	Size distribution in CV [%]
200	1,000	714	354.00	214.1	762.7	90	13.1

Assay Class: High Sensitivity DNA Assay  
 Data Path: C:\... bioanalyzer\2100 expert\data\2018-05-03\2018-05-03\_001.xad

Created: 5/3/2018 2:55:51 PM  
 Modified: 5/3/2018 3:49:23 PM

**Electropherogram Summary Continued ...**

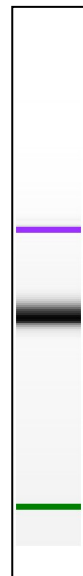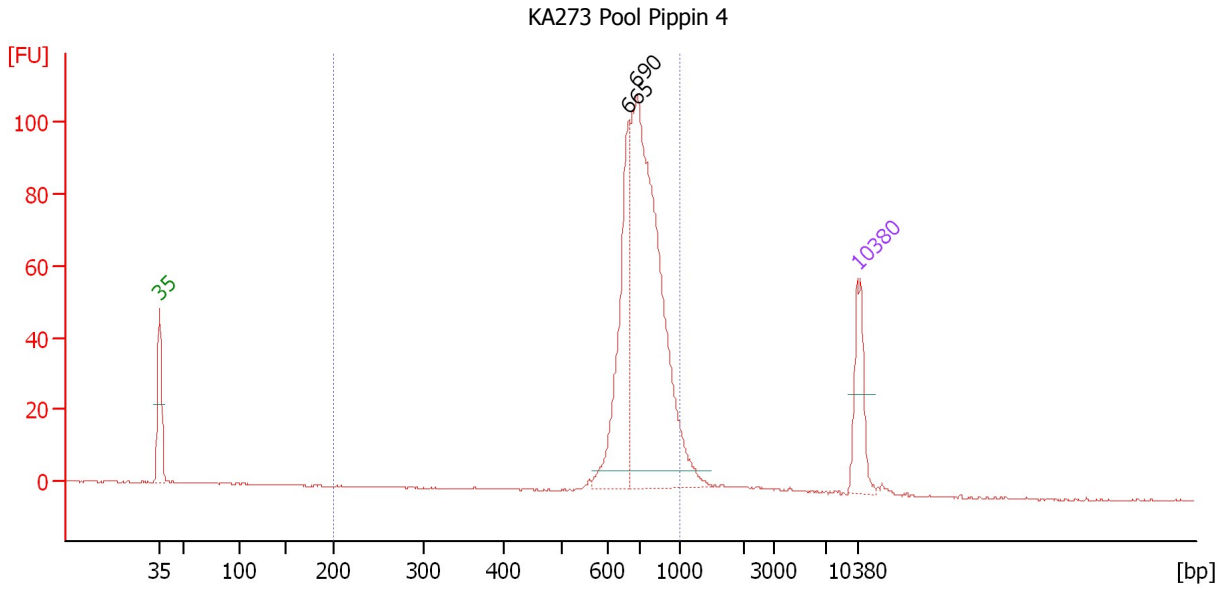

**Overall Results for sample 11 : KA273 Pool Pippin 4**

Number of peaks found: 2                      Corr. Area 1: 529.3  
 Noise: 0.2

**Peak table for sample 11 : KA273 Pool Pippin 4**

Peak	Size [bp]	Conc. [pg/μl]	Molarity [pmol/l]	Observations	Aligned Migration Time [s]
1	35	125.00	5,411.3	Lower Marker	43.00
2	665	239.99	546.9		90.04
3	690	710.27	1,558.7		90.84
4	10,380	75.00	10.9	Upper Marker	113.00

**Region table for sample 11 : KA273 Pool Pippin 4**


From [bp]	To [bp]	Average Size [bp]	Conc. [pg/μl]	Corr. Area	Molarity [pmol/l]	Co % of lor Total	Size distribution in CV [%]
200	1,000	738	898.01	529.3	1,872.1	92	13.0



Assay Class: High Sensitivity DNA Assay  
Data Path: C:\... bioanalyzer\2100 expert\data\2018-05-03\2018-05-03\_001.xad

Created: 5/3/2018 2:55:51 PM  
Modified: 5/3/2018 3:49:23 PM

**Gel Image**

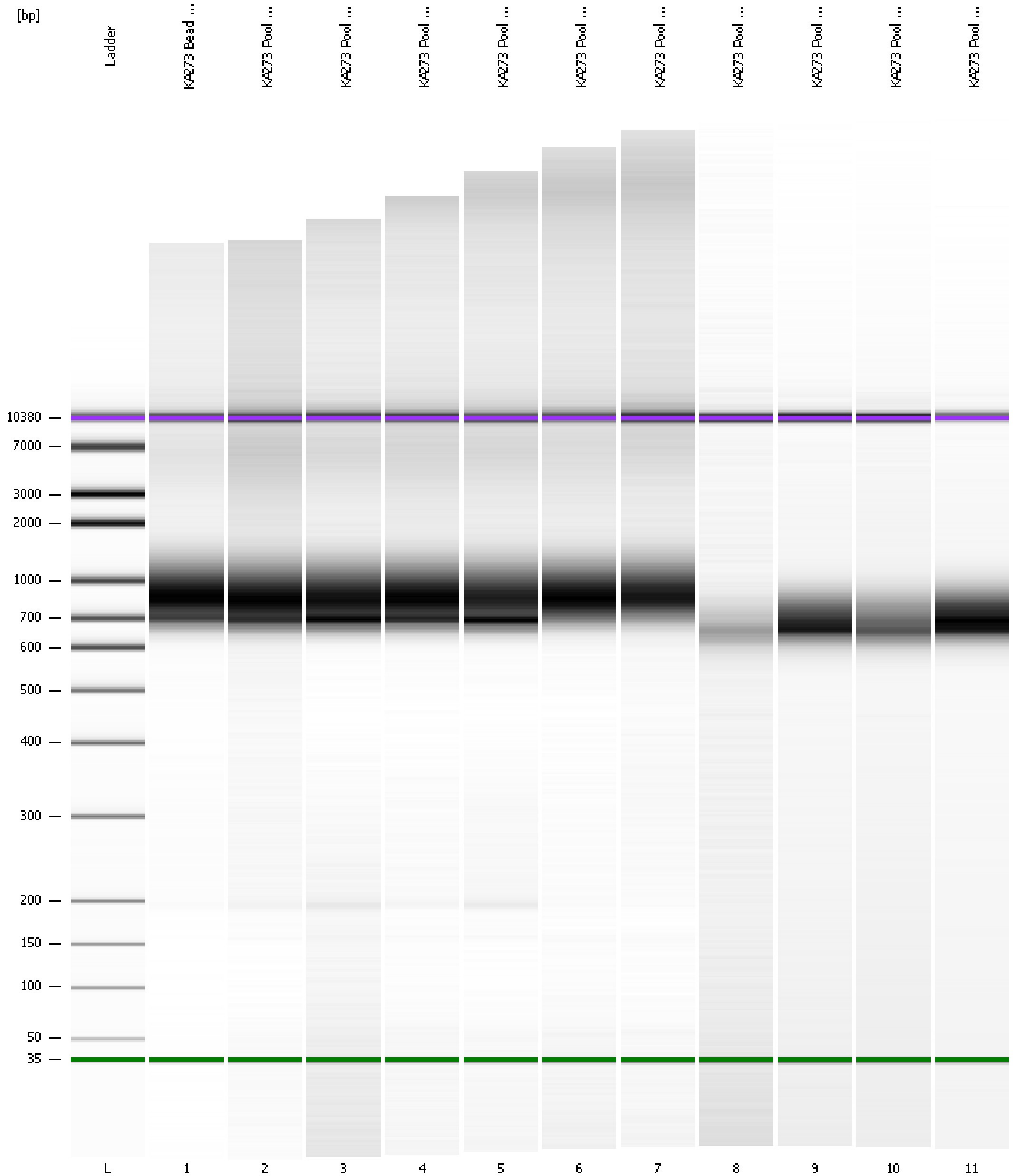

Assay Class: High Sensitivity DNA Assay  
Data Path: C:\... bioanalyzer\2100 expert\data\2018-05-03\2018-05-03\_001.xad

Created: 5/3/2018 2:55:51 PM  
Modified: 5/3/2018 3:49:23 PM

**Curves**

**Standard Curve**

