


Assay Class: High Sensitivity DNA Assay
Data Path: C:\...00 expert\data\2018-07-23\2018-07-23_002_F137_F138_LCBO.xad

Created: 7/23/2018 3:17:03 PM
Modified: 7/24/2018 2:38:27 PM

Electrophoresis File Run Summary

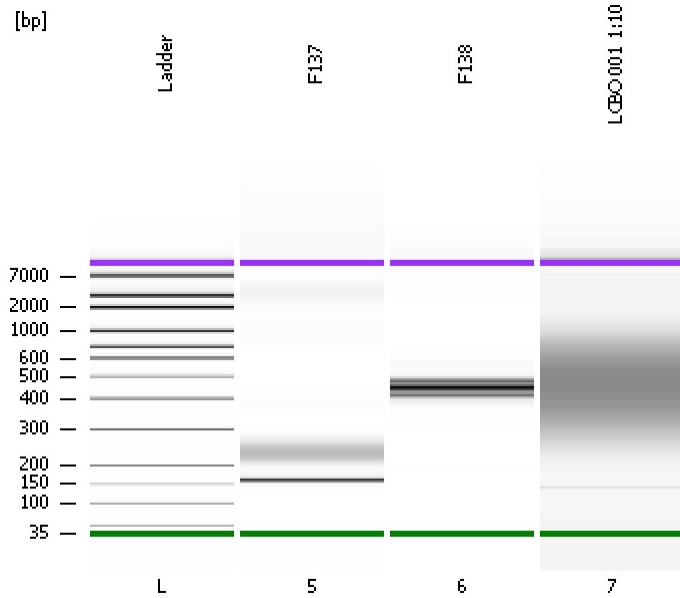

Instrument Information:




Instrument Name: DE34903152 Firmware: C.01.069
Serial#: DE34903152 Type: G2938C

Assay Information:

Assay Origin Path: C:\Program Files\Agilent\2100 bioanalyzer\2100 expert\assays\dsDNA\High Sensitivity DNA.xsy
Assay Class: High Sensitivity DNA Assay
Version: 1.03
Assay Comments: Copyright © 2003-2010 Agilent Technologies

Chip Information:

Chip Lot #:
Reagent Kit Lot #:
Chip Comments:

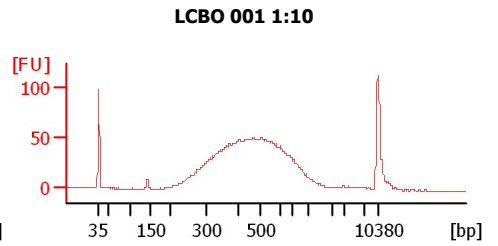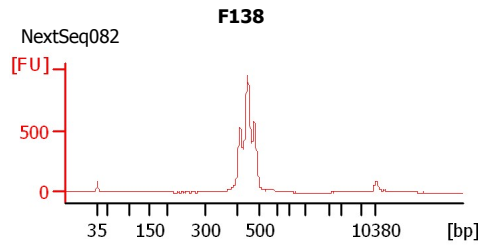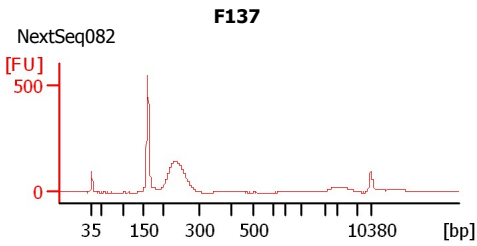

Assay Class: High Sensitivity DNA Assay
 Data Path: C:\...00 expert\data\2018-07-23\2018-07-23_002_F137_F138_LCBO.xad

Created: 7/23/2018 3:17:03 PM
 Modified: 7/24/2018 2:38:27 PM

Electrophoresis File Run Summary (Chip Summary)

Sample Name	Sample Comment	Rest. Digest	Status	Observation	Result Label	Result Color
F137	NextSeq082	<input type="checkbox"/>				
F138	NextSeq082	<input type="checkbox"/>				
LCBO 001 1:10		<input type="checkbox"/>				
Ladder		<input type="checkbox"/>				

Chip Lot # **Reagent Kit Lot #**

Chip Comments :

Assay Class: High Sensitivity DNA Assay
Data Path: C:\...00 expert\data\2018-07-23\2018-07-23_002_F137_F138_LCBO.xad

Created: 7/23/2018 3:17:03 PM
Modified: 7/24/2018 2:38:27 PM

Electrophoresis Assay Details

General Analysis Settings

Number of Available Sample and Ladder Wells (Max.) : 12
Minimum Visible Range [s] : 32
Maximum Visible Range [s] : 138
Start Analysis Time Range [s] : 33
End Analysis Time Range [s] : 137.5
Ladder Concentration [pg/μl] : 1950
Uses Standard Area for Ladder Fragments
Lower Marker Concentration [pg/μl] : 125
Upper Marker Concentration [pg/μl] : 75
Used Upper Marker for Quantitation
Standard Curve Fit is Point to Point
Show Data Aligned to Lower and Upper Marker

Integrator Settings

Integration Start Time [s] : 33.05
Integration End Time [s] : 137
Slope Threshold : 0.8
Height Threshold [FU] : 5
Area Threshold : 0.1
Width Threshold [s] : 0.6
Baseline Plateau [s] : 0.5

Filter Settings

Filter Width [s] : 0.5
Polynomial Order : 4


Ladder

Ladder Peak	Size	Area
1	35	160
2	50	210
3	100	208
4	150	221
5	200	242
6	300	270
7	400	305
8	500	306
9	600	336
10	700	321
11	1000	366
12	2000	413
13	3000	411
14	7000	400
15	10380	214

Assay Class: High Sensitivity DNA Assay
 Data Path: C:\...00 expert\data\2018-07-23\2018-07-23_002_F137_F138_LCBO.xad

Created: 7/23/2018 3:17:03 PM
 Modified: 7/24/2018 2:38:27 PM

Electropherogram Summary

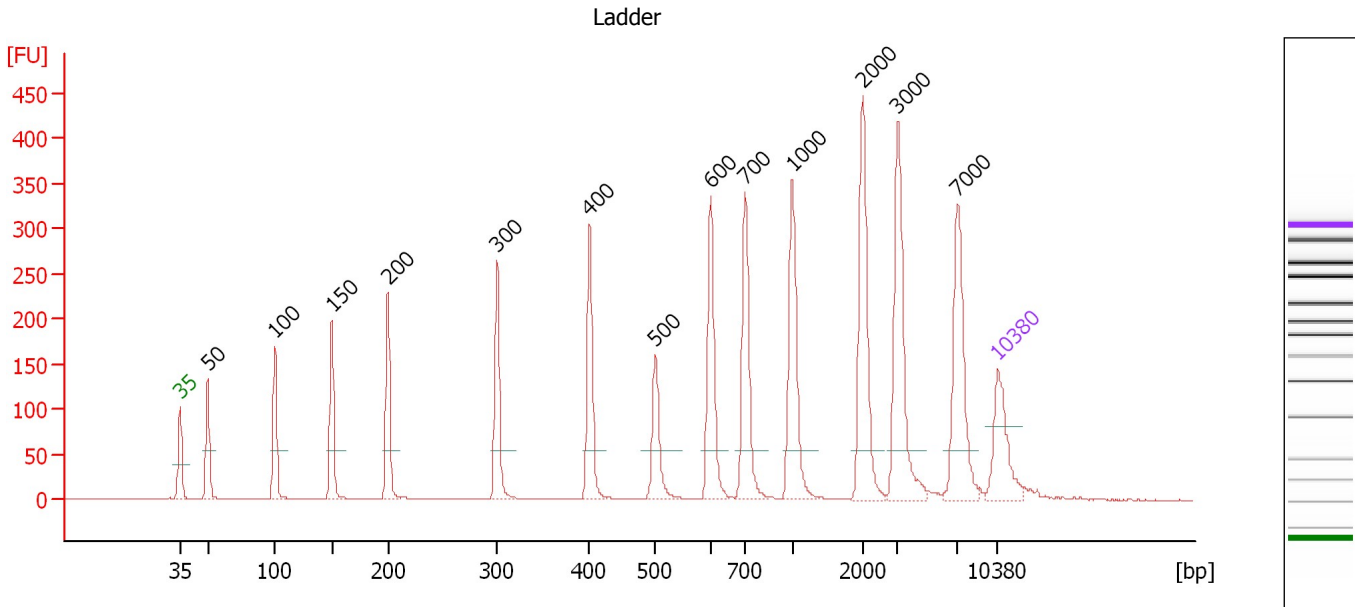

Overall Results for Ladder

Noise: 0.3


Peak table for Ladder

Peak	Size [bp]	Conc. [pg/μl]	Molarity [pmol/l]	Observations	Aligned Migration Time [s]
1	35	125.00	5,411.3	Lower Marker	43.00
2	50	150.00	4,545.5	Ladder Peak	45.37
3	100	150.00	2,272.7	Ladder Peak	51.12
4	150	150.00	1,515.2	Ladder Peak	56.01
5	200	150.00	1,136.4	Ladder Peak	60.80
6	300	150.00	757.6	Ladder Peak	70.15
7	400	150.00	568.2	Ladder Peak	78.05
8	500	150.00	454.5	Ladder Peak	83.66
9	600	150.00	378.8	Ladder Peak	88.40
10	700	150.00	324.7	Ladder Peak	91.37
11	1,000	150.00	227.3	Ladder Peak	95.39
12	2,000	150.00	113.6	Ladder Peak	101.41
13	3,000	150.00	75.8	Ladder Peak	104.42
14	7,000	150.00	32.5	Ladder Peak	109.58
15	10,380	75.00	10.9	Upper Marker	113.00

Assay Class: High Sensitivity DNA Assay
 Data Path: C:\...00 expert\data\2018-07-23\2018-07-23_002_F137_F138_LCBO.xad

Created: 7/23/2018 3:17:03 PM
 Modified: 7/24/2018 2:38:27 PM

Electropherogram Summary Continued ...

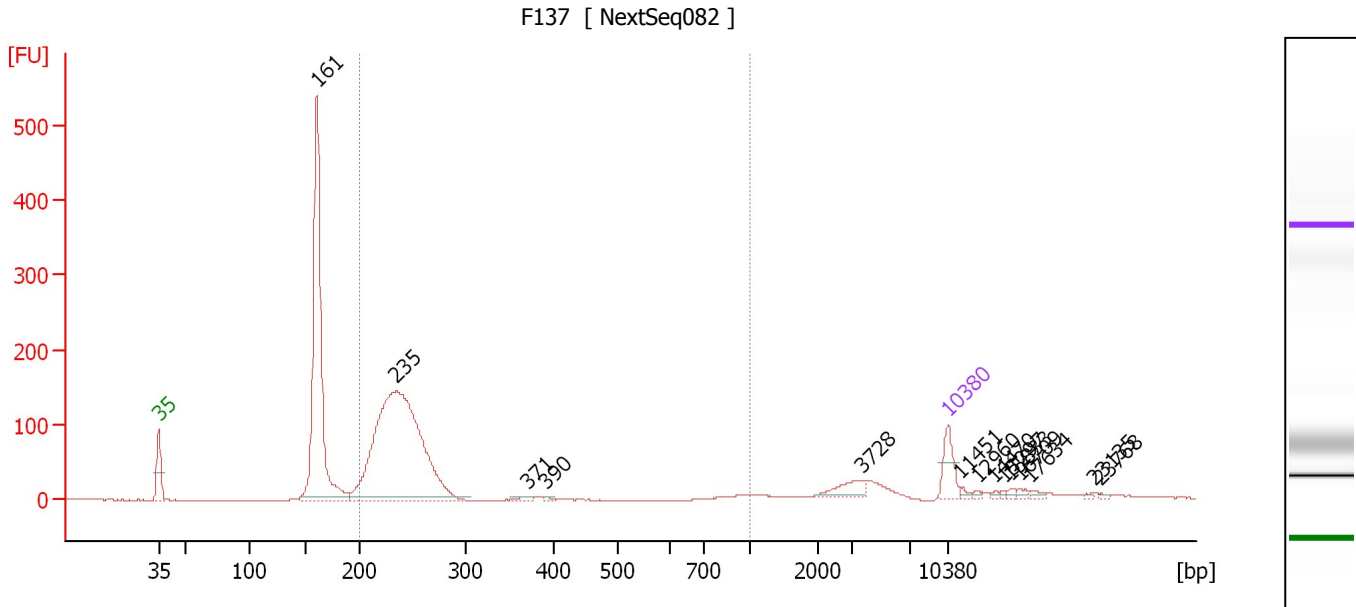

Overall Results for sample 5 : F137

Number of peaks found: 14 Corr. Area 1: 1,123.0
 Noise: 0.5

Peak table for sample 5 : F137

Peak	Size [bp]	Conc. [pg/μl]	Molarity [pmol/l]	Observations	Aligned Migration Time [s]
1	35	125.00	5,411.3	Lower Marker	43.00
2	161	1,293.70	12,182.8		57.05
3	235	1,659.64	10,686.9		64.10
4	371	15.67	63.9		75.78
5	390	7.53	29.2		77.26
6	3,728	56.71	23.1		105.36
7	10,380	75.00	10.9	Upper Marker	113.00
8	11,451	0.00	0.0		114.08
9	12,960	0.00	0.0		115.61
10	14,470	0.00	0.0		117.14
11	15,297	0.00	0.0		117.98
12	15,833	0.00	0.0		118.52
13	16,709	0.00	0.0		119.41
14	17,634	0.00	0.0		120.35
15	23,135	0.00	0.0		125.92
16	23,768	0.00	0.0		126.56


Region table for sample 5 : F137

From [bp]	To [bp]	Average Size [bp]	Conc. [pg/μl]	Corr. Area	Molarity [pmol/l]	Co % of Ior Total	Size distribution in CV [%]
200	1,000	252	1,603.58	1,123.0	10,085.8	51	35.7

Assay Class: High Sensitivity DNA Assay
 Data Path: C:\...00 expert\data\2018-07-23\2018-07-23_002_F137_F138_LCBO.xad

Created: 7/23/2018 3:17:03 PM
 Modified: 7/24/2018 2:38:27 PM

Electropherogram Summary Continued ...

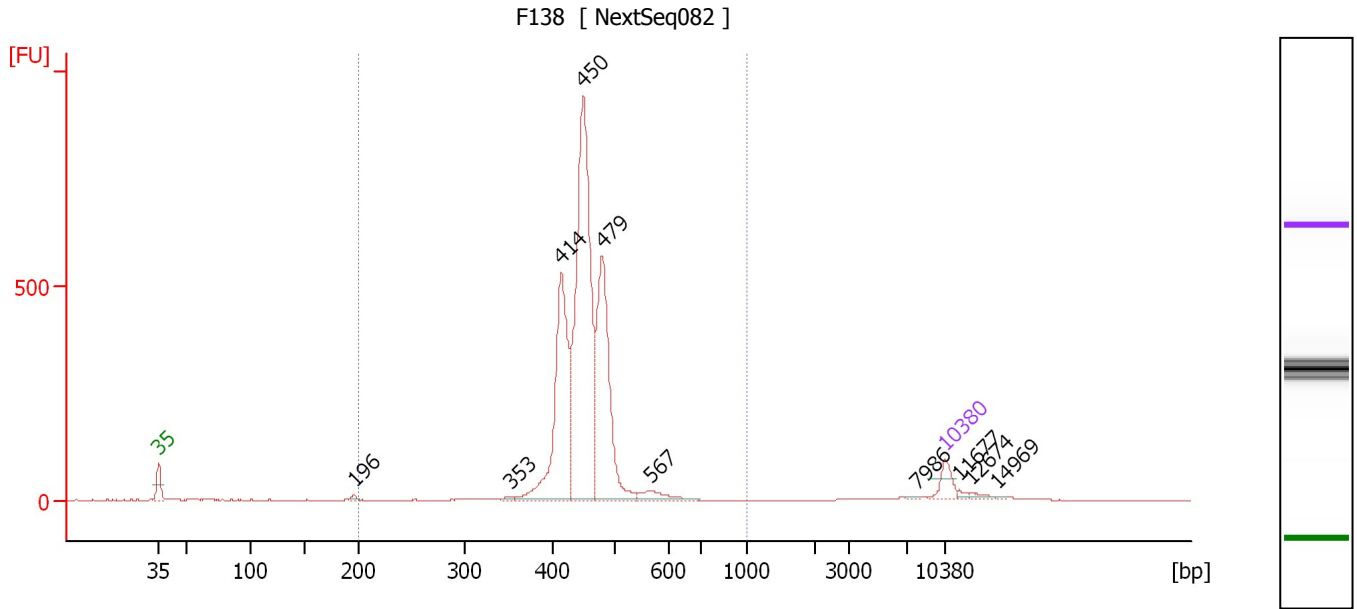

Overall Results for sample 6 : F138

Number of peaks found: 10 Corr. Area 1: 3,771.0
 Noise: 0.3

Peak table for sample 6 : F138

Peak	Size [bp]	Conc. [pg/μl]	Molarity [pmol/l]	Observations	Aligned Migration Time [s]
1	35	125.00	5,411.3	Lower Marker	43.00
2	196	14.05	108.8		60.37
3	353	11.03	47.3		74.36
4	414	1,030.37	3,766.5		78.86
5	450	1,914.54	6,452.3		80.83
6	479	989.75	3,129.0		82.49
7	567	77.16	206.2		86.84
8	7,986	4.58	0.9		110.58
9	10,380	75.00	10.9	Upper Marker	113.00
10	11,677	0.00	0.0		114.31
11	12,674	0.00	0.0		115.32
12	14,969	0.00	0.0		117.65


Region table for sample 6 : F138

From [bp]	To [bp]	Average Size [bp]	Conc. [pg/μl]	Corr. Area	Molarity [pmol/l]	Co % of lor Total	Size distribution in CV [%]
200	1,000	451	4,077.78	3,771.0	13,800.5	95	8.8

Assay Class: High Sensitivity DNA Assay
 Data Path: C:\...00 expert\data\2018-07-23\2018-07-23_002_F137_F138_LCBO.xad

Created: 7/23/2018 3:17:03 PM
 Modified: 7/24/2018 2:38:27 PM

Electropherogram Summary Continued ...

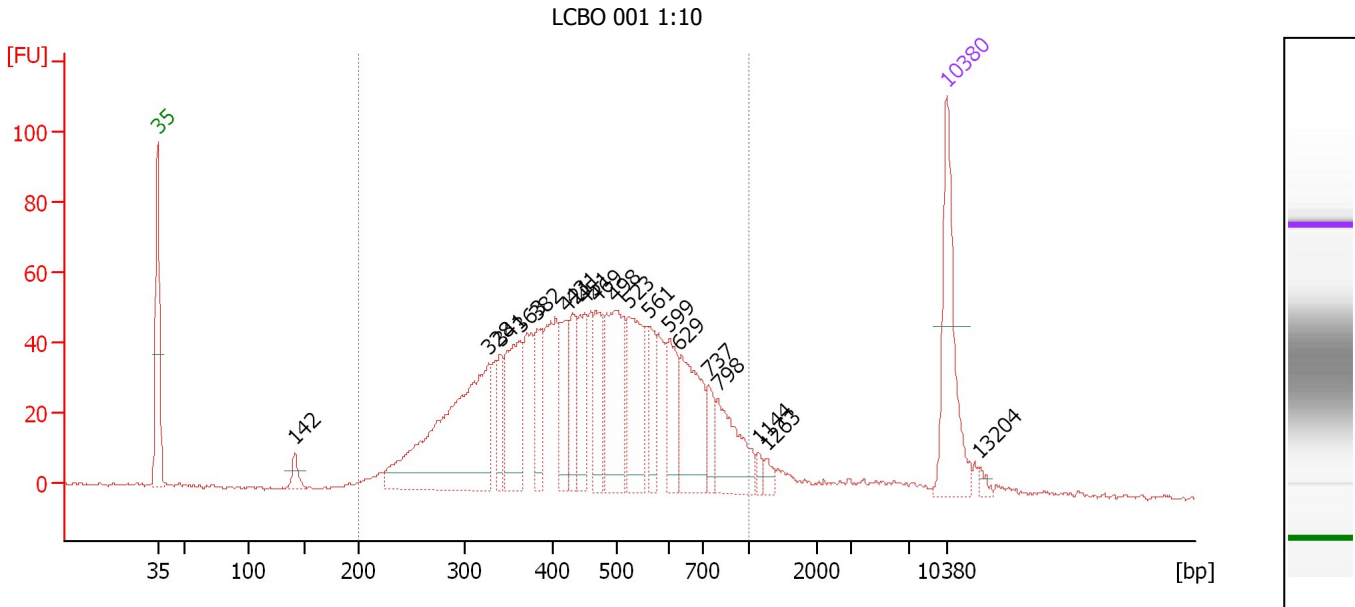

Overall Results for sample 7 : LCBO 001 1:10

Number of peaks found: 19 Corr. Area 1: 1,423.6
 Noise: 0.6

Peak table for sample 7 : LCBO 001 1:10

Peak	Size [bp]	Conc. [pg/μl]	Molarity [pmol/l]	Observations	Aligned Migration Time [s]
1	35	125.00	5,411.3	Lower Marker	43.00
2	142	14.33	153.1		55.20
3	328	247.13	1,141.7		72.36
4	341	32.50	144.4		73.38
5	363	84.72	353.8		75.12
6	382	43.96	174.2		76.65
7	421	46.84	168.8		79.20
8	431	39.33	138.4		79.76
9	451	49.94	167.9		80.88
10	469	47.40	153.2		81.91
11	498	103.25	314.2		83.54
12	523	75.45	218.5		84.77
13	561	31.74	85.7		86.55
14	599	40.27	101.9		88.34
15	629	77.53	186.8		89.26
16	737	21.16	43.5		91.86
17	798	56.73	107.7		92.68
18	1,144	5.37	7.1		96.25
19	1,263	6.51	7.8		96.97
20	10,380	75.00	10.9	Upper Marker	113.00
21	13,204	0.00	0.0		115.86

Region table for sample 7 : LCBO 001 1:10

From [bp]	To [bp]	Average Size [bp]	Conc. [pg/μl]	Corr. Area	Molarity [pmol/l]	Co % of lor Total	Size distribution in CV [%]
200	1,000	475	1,249.88	1,423.6	4,599.6	93	32.4

Assay Class: High Sensitivity DNA Assay
Data Path: C:\...00 expert\data\2018-07-23\2018-07-23_002_F137_F138_LCBO.xad

Created: 7/23/2018 3:17:03 PM
Modified: 7/24/2018 2:38:27 PM

Gel Image

