


Assay Class: High Sensitivity DNA Assay  
Data Path: C:\...xpert\data\2019-10-08\2019-10-08\_001\_HiSeq1924\_1-3\_Libs.xad

Created: 10/8/2019 10:04:50 AM  
Modified: 10/8/2019 10:58:39 AM

**Electrophoresis File Run Summary**

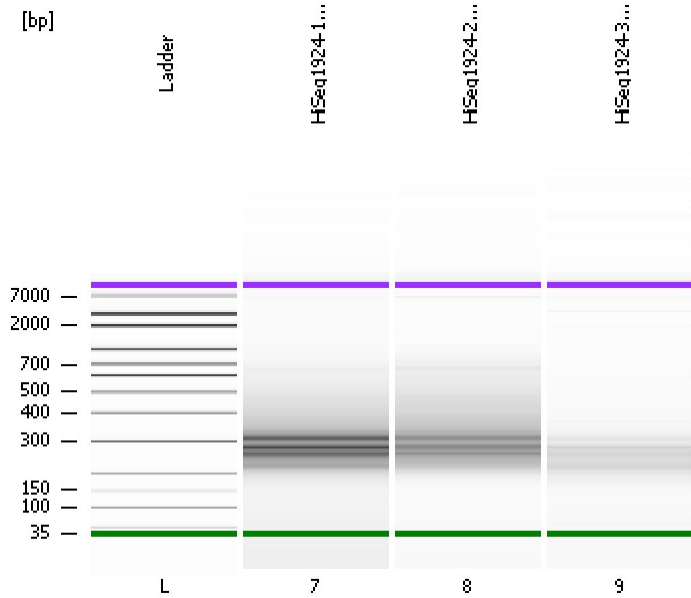

Instrument Information:

Instrument Name: DE34903152      Firmware: C.01.069  
Serial#: DE34903152      Type: G2938C

Assay Information:

Assay Origin Path: C:\Program Files\Agilent\2100 bioanalyzer\2100 expert\assays\dsDNA\High Sensitivity DNA.xsy  
Assay Class: High Sensitivity DNA Assay  
Version: 1.03  
Assay Comments: Copyright © 2003-2010 Agilent Technologies




Chip Information:

Chip Lot #:  
Reagent Kit Lot #:  
Chip Comments:

HiSeq1924-1\_Lib

HiSeq1924-2\_Lib

HiSeq1924-3\_Lib

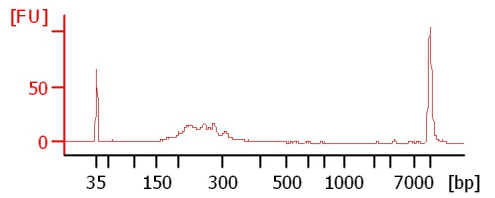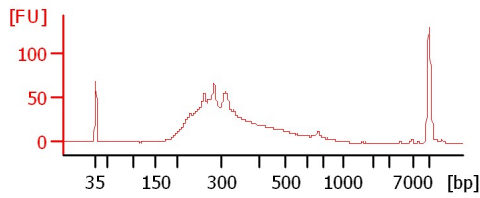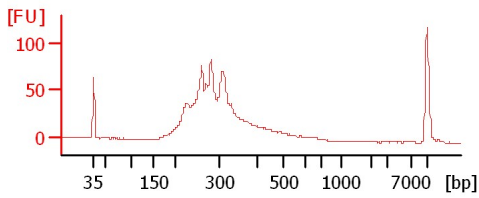

Assay Class: High Sensitivity DNA Assay  
 Data Path: C:\...xpert\data\2019-10-08\2019-10-08\_001\_HiSeq1924\_1-3\_Libs.xad

Created: 10/8/2019 10:04:50 AM  
 Modified: 10/8/2019 10:58:39 AM

**Electrophoresis File Run Summary (Chip Summary)**

Sample Name	Sample Comment	Rest. Digest	Status	Observation	Result Label	Result Color
HiSeq1924-1_Lib		<input type="checkbox"/>				
HiSeq1924-2_Lib		<input type="checkbox"/>				
HiSeq1924-3_Lib		<input type="checkbox"/>				
Ladder		<input type="checkbox"/>				

**Chip Lot #** **Reagent Kit Lot #**

**Chip Comments :**

Assay Class: High Sensitivity DNA Assay  
Data Path: C:\...xpert\data\2019-10-08\2019-10-08\_001\_HiSeq1924\_1-3\_Libs.xad

Created: 10/8/2019 10:04:50 AM  
Modified: 10/8/2019 10:58:39 AM

**Electrophoresis Assay Details**

**General Analysis Settings**

Number of Available Sample and Ladder Wells (Max.) : 12  
Minimum Visible Range [s] : 32  
Maximum Visible Range [s] : 138  
Start Analysis Time Range [s] : 33  
End Analysis Time Range [s] : 137.5  
Ladder Concentration [pg/μl] : 1950  
Uses Standard Area for Ladder Fragments  
Lower Marker Concentration [pg/μl] : 125  
Upper Marker Concentration [pg/μl] : 75  
Used Upper Marker for Quantitation  
Standard Curve Fit is Point to Point  
Show Data Aligned to Lower and Upper Marker

**Integrator Settings**

Integration Start Time [s] : 33.05  
Integration End Time [s] : 137  
Slope Threshold : 0.8  
Height Threshold [FU] : 5  
Area Threshold : 0.1  
Width Threshold [s] : 0.6  
Baseline Plateau [s] : 0.5

**Filter Settings**

Filter Width [s] : 0.5  
Polynomial Order : 4


**Ladder**

Ladder Peak	Size	Area
1	35	160
2	50	210
3	100	208
4	150	221
5	200	242
6	300	270
7	400	305
8	500	306
9	600	336
10	700	321
11	1000	366
12	2000	413
13	3000	411
14	7000	400
15	10380	214

Assay Class: High Sensitivity DNA Assay  
 Data Path: C:\...xpert\data\2019-10-08\2019-10-08\_001\_HiSeq1924\_1-3\_Libs.xad

Created: 10/8/2019 10:04:50 AM  
 Modified: 10/8/2019 10:58:39 AM

**Electropherogram Summary**

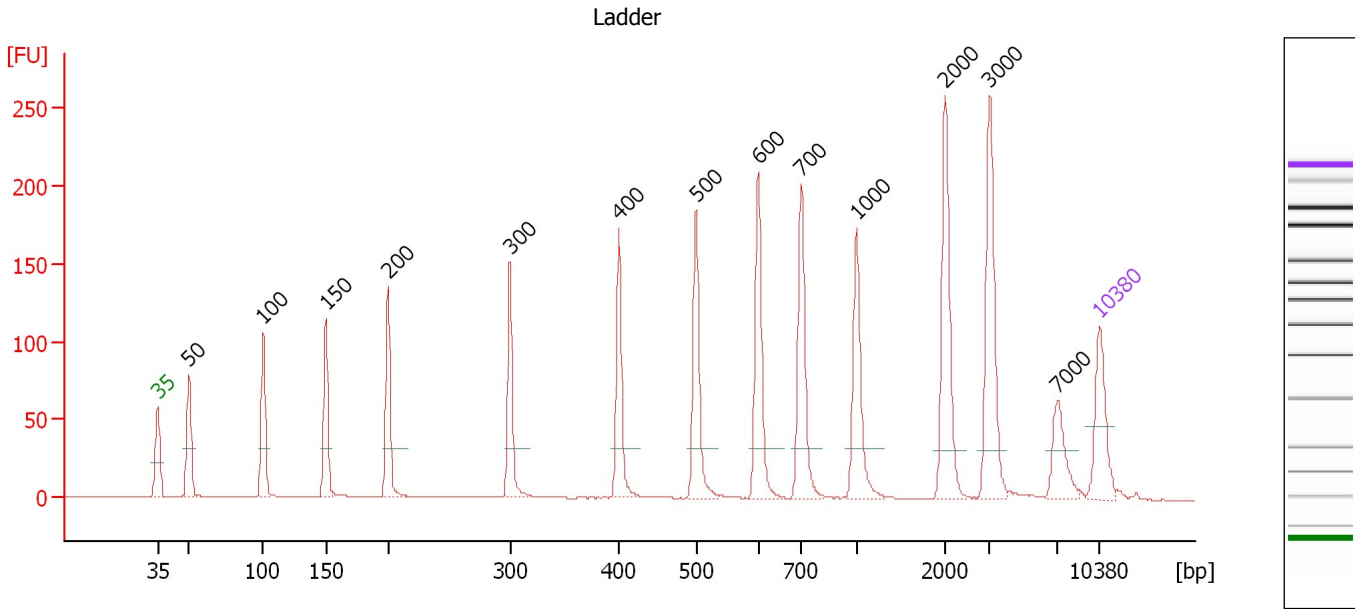

**Overall Results for Ladder**

Noise: 0.2


**Peak table for Ladder**

Peak	Size [bp]	Conc. [pg/μl]	Molarity [pmol/l]	Observations	Aligned Migration Time [s]
1	35	125.00	5,411.3	Lower Marker	43.00
2	50	150.00	4,545.5	Ladder Peak	45.36
3	100	150.00	2,272.7	Ladder Peak	50.84
4	150	150.00	1,515.2	Ladder Peak	55.50
5	200	150.00	1,136.4	Ladder Peak	60.12
6	300	150.00	757.6	Ladder Peak	69.20
7	400	150.00	568.2	Ladder Peak	77.28
8	500	150.00	454.5	Ladder Peak	83.00
9	600	150.00	378.8	Ladder Peak	87.62
10	700	150.00	324.7	Ladder Peak	90.79
11	1,000	150.00	227.3	Ladder Peak	94.92
12	2,000	150.00	113.6	Ladder Peak	101.51
13	3,000	150.00	75.8	Ladder Peak	104.83
14	7,000	150.00	32.5	Ladder Peak	109.88
15	10,380	75.00	10.9	Upper Marker	113.00

Assay Class: High Sensitivity DNA Assay  
 Data Path: C:\...xpert\data\2019-10-08\2019-10-08\_001\_HiSeq1924\_1-3\_Libs.xad

Created: 10/8/2019 10:04:50 AM  
 Modified: 10/8/2019 10:58:39 AM

**Electropherogram Summary Continued ...**

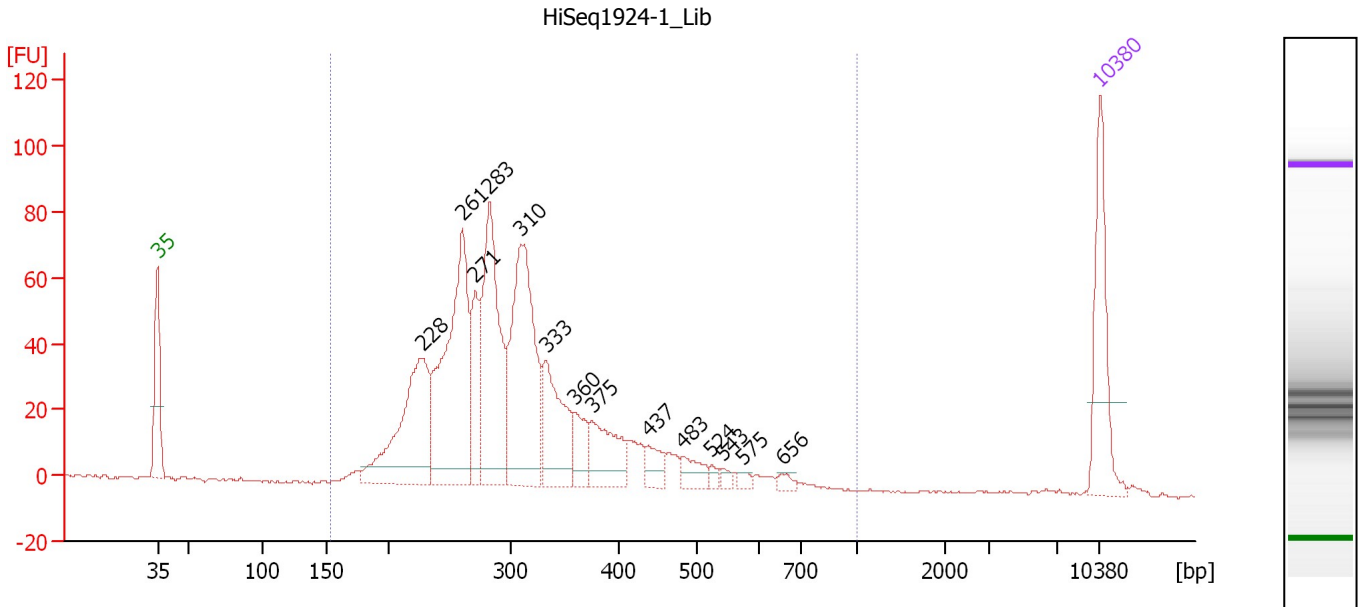

**Overall Results for sample 7 : HiSeq1924-1 Lib**

Number of peaks found: 14                      Corr. Area 1: 1,095.5  
 Noise: 0.3

**Peak table for sample 7 : HiSeq1924-1 Lib**

Peak	Size [bp]	Conc. [pg/μl]	Molarity [pmol/l]	Observations	Aligned Migration Time [s]
1	35	125.00	5,411.3	Lower Marker	43.00
2	228	197.55	1,311.4		62.68
3	261	266.49	1,545.6		65.68
4	271	77.47	432.5		66.61
5	283	206.98	1,108.8		67.64
6	310	236.73	1,158.2		69.99
7	333	101.16	460.8		71.84
8	360	33.52	141.2		74.02
9	375	66.07	266.8		75.27
10	437	21.64	75.0		79.42
11	483	21.25	66.6		82.03
12	524	5.94	17.2		84.11
13	543	6.36	17.7		84.98
14	575	6.52	17.2		86.45
15	656	6.51	15.0		89.39
16	10,380	75.00	10.9	Upper Marker	113.00


**Region table for sample 7 : HiSeq1924-1 Lib**

From [bp]	To [bp]	Average Size [bp]	Conc. [pg/μl]	Corr. Area	Molarity [pmol/l]	Co % of lor Total	Size distribution in CV [%]
155	1,000	314	1,305.12	1,095.5	6,840.7	99	28.0

Assay Class: High Sensitivity DNA Assay  
 Data Path: C:\...xpert\data\2019-10-08\2019-10-08\_001\_HiSeq1924\_1-3\_Libs.xad

Created: 10/8/2019 10:04:50 AM  
 Modified: 10/8/2019 10:58:39 AM

**Electropherogram Summary Continued ...**

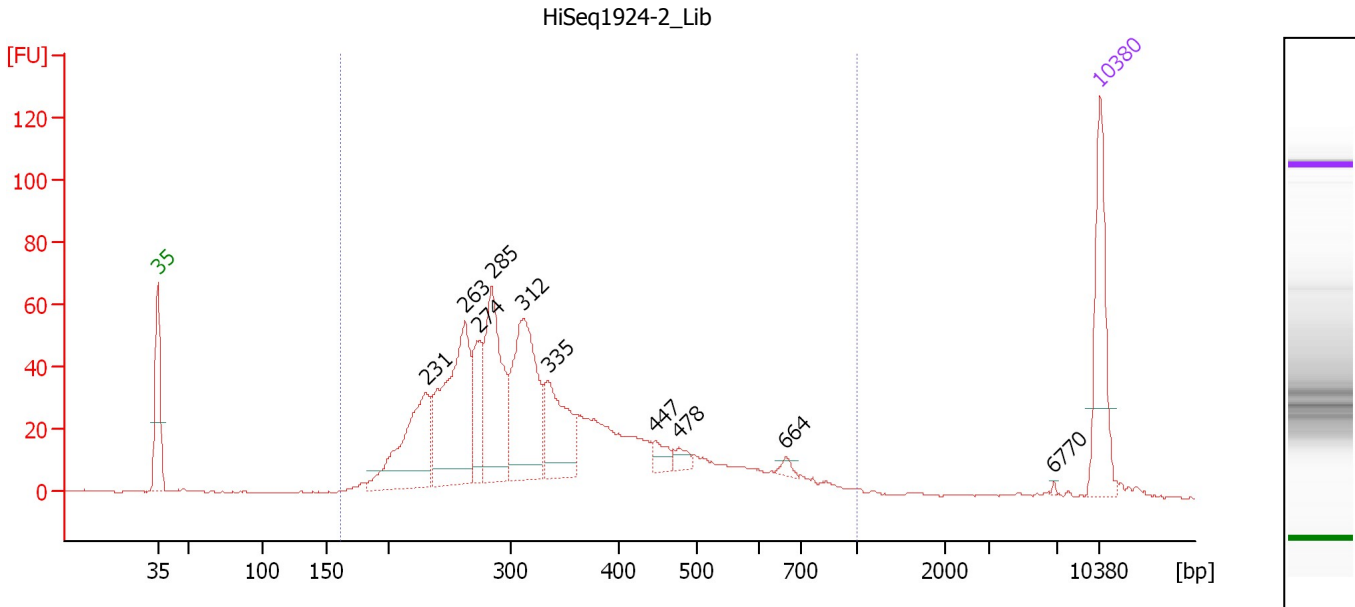

**Overall Results for sample 8 : HiSeq1924-2 Lib**

Number of peaks found: 10      Corr. Area 1: 1,088.8  
 Noise: 0.3

**Peak table for sample 8 : HiSeq1924-2 Lib**

Peak	Size [bp]	Conc. [pg/μl]	Molarity [pmol/l]	Observations	Aligned Migration Time [s]
1	35	125.00	5,411.3	Lower Marker	43.00
2	231	144.70	948.9		62.94
3	263	206.08	1,189.1		65.80
4	274	65.46	361.9		66.85
5	285	153.71	817.2		67.84
6	312	176.06	855.7		70.15
7	335	90.83	410.9		72.02
8	447	15.86	53.8		79.96
9	478	10.71	34.0		81.72
10	664	6.26	14.3		89.65
11	6,770	1.02	0.2		109.59
12	10,380	75.00	10.9	Upper Marker	113.00


**Region table for sample 8 : HiSeq1924-2 Lib**

From [bp]	To [bp]	Average Size [bp]	Conc. [pg/μl]	Corr. Area	Molarity [pmol/l]	Co % of lor Total	Size distribution in CV [%]
161	1,000	349	1,260.61	1,088.8	6,253.0	97	36.1

Assay Class: High Sensitivity DNA Assay  
 Data Path: C:\...xpert\data\2019-10-08\2019-10-08\_001\_HiSeq1924\_1-3\_Libs.xad

Created: 10/8/2019 10:04:50 AM  
 Modified: 10/8/2019 10:58:39 AM

**Electropherogram Summary Continued ...**

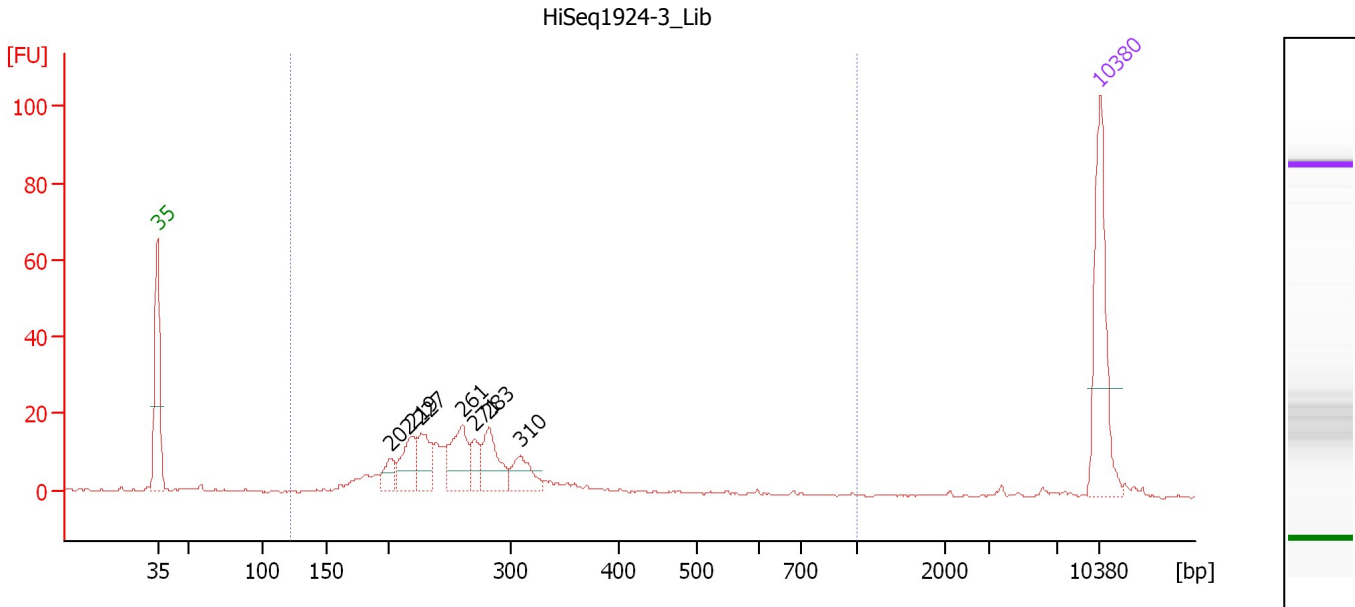

**Overall Results for sample 9 : HiSeq1924-3 Lib**

Number of peaks found: 7                      Corr. Area 1: 249.7  
 Noise: 0.3

**Peak table for sample 9 : HiSeq1924-3 Lib**

Peak	Size [bp]	Conc. [pg/μl]	Molarity [pmol/l]	Observations	Aligned Migration Time [s]
1	35	125.00	5,411.3	Lower Marker	43.00
2	202	19.22	144.3		60.28
3	219	39.29	271.3		61.88
4	227	34.75	231.6		62.60
5	261	49.85	289.5		65.65
6	271	21.18	118.5		66.54
7	283	42.75	229.0		67.64
8	310	27.90	136.6		69.97
9	10,380	75.00	10.9	Upper Marker	113.00


**Region table for sample 9 : HiSeq1924-3 Lib**

From [bp]	To [bp]	Average Size [bp]	Conc. [pg/μl]	Corr. Area	Molarity [pmol/l]	Co % of Ior Total	Size distribution in CV [%]
122	1,000	274	350.08	249.7	2,107.5	88	32.8

Assay Class: High Sensitivity DNA Assay  
Data Path: C:\...xpert\data\2019-10-08\2019-10-08\_001\_HiSeq1924\_1-3\_Libs.xad

Created: 10/8/2019 10:04:50 AM  
Modified: 10/8/2019 10:58:39 AM

**Gel Image**

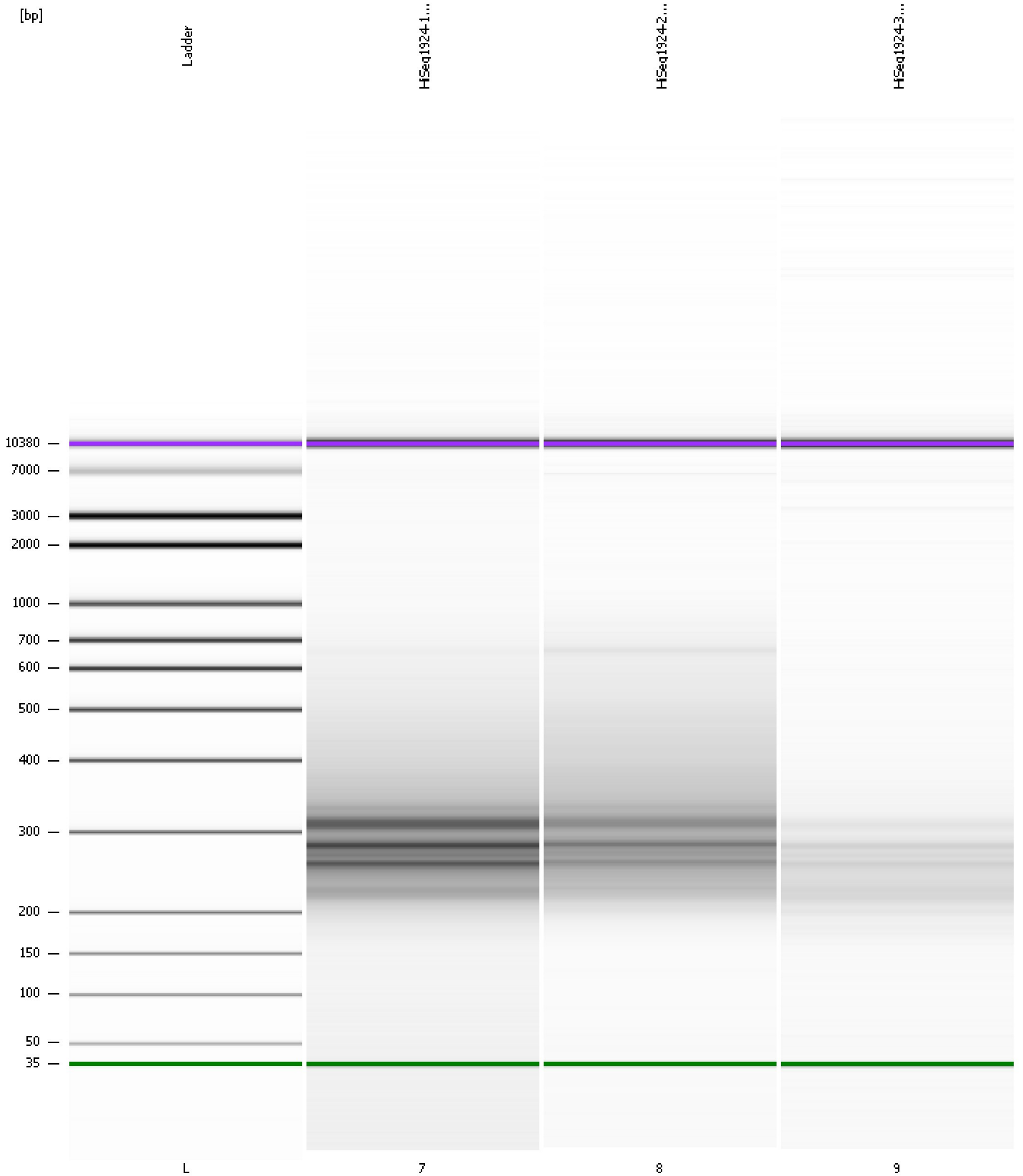




Assay Class: High Sensitivity DNA Assay  
Data Path: C:\...xpert\data\2019-10-08\2019-10-08\_001\_HiSeq1924\_1-3\_Libs.xad

Created: 10/8/2019 10:04:50 AM  
Modified: 10/8/2019 10:58:39 AM

**Curves**

**Standard Curve**

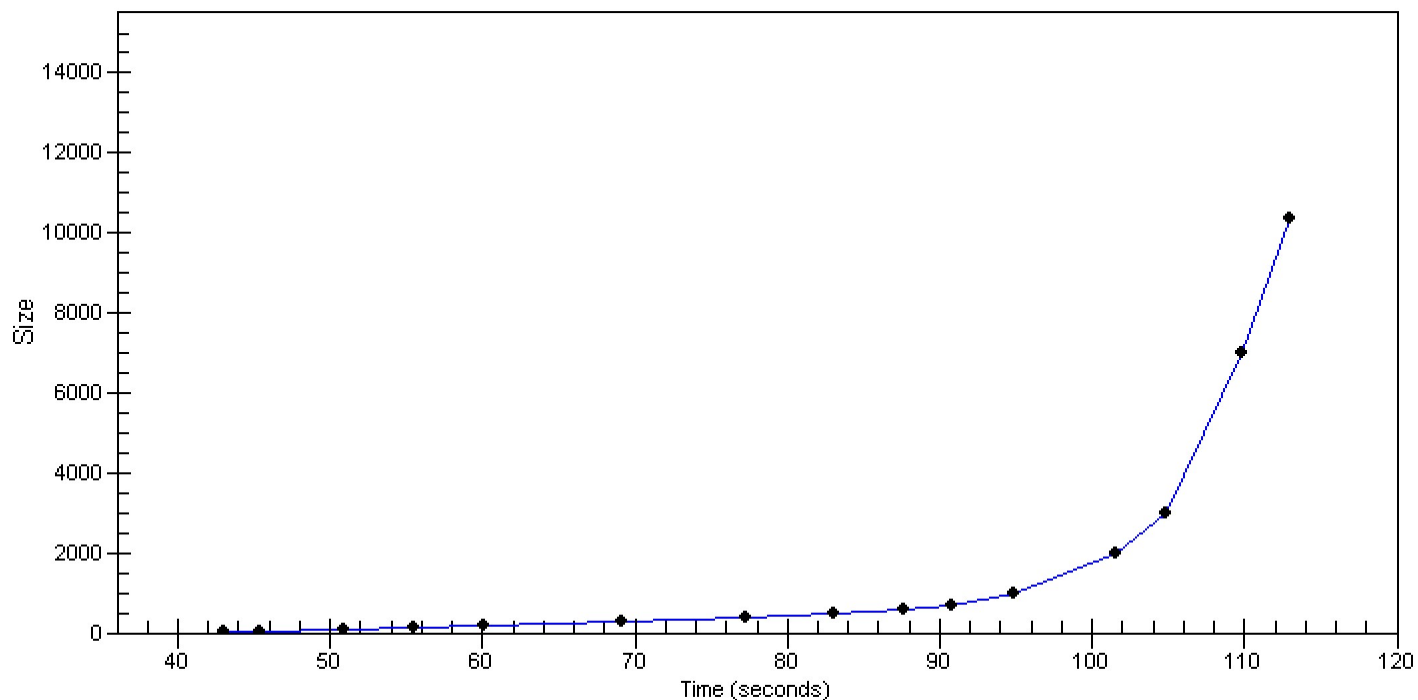

Assay Class: High Sensitivity DNA Assay  
 Data Path: C:\...xpert\data\2019-10-08\2019-10-08\_001\_HiSeq1924\_1-3\_Libs.xad

Created: 10/8/2019 10:04:50 AM  
 Modified: 10/8/2019 10:58:39 AM

**Run Logbook**

Description	Number	Source	Category	Sub Category	Time	Time Zone	User	Host
Run ended on port 1 (Number of wells acquired: 12)		Instrument	Run		10/8/2019 10:46:08 AM	(GMT --07:00) Pacific Standard Time	UC Davis	D8XSMGH1
Run started on port 1 (File: C:\Program Files\Agilent\2100 bioanalyzer\2100 expert\data\2019-10-08\2019-10-08_001.xad)		Instrument	Run		10/8/2019 10:04:55 AM	(GMT --07:00) Pacific Standard Time	UC Davis	D8XSMGH1
Product Number : G2938C		Instrument	Run		10/8/2019 10:04:55 AM	(GMT --07:00) Pacific Standard Time	UC Davis	D8XSMGH1
Name :		Instrument	Run		10/8/2019 10:04:55 AM	(GMT --07:00) Pacific Standard Time	UC Davis	D8XSMGH1
Vendor : Agilent Technologies		Instrument	Run		10/8/2019 10:04:55 AM	(GMT --07:00) Pacific Standard Time	UC Davis	D8XSMGH1
Serial# : DE34903152		Instrument	Run		10/8/2019 10:04:55 AM	(GMT --07:00) Pacific Standard Time	UC Davis	D8XSMGH1
Firmware : C.01.069		Instrument	Run		10/8/2019 10:04:55 AM	(GMT --07:00) Pacific Standard Time	UC Davis	D8XSMGH1
Cartridge : Electrode		Instrument	Run		10/8/2019 10:04:55 AM	(GMT --07:00) Pacific Standard Time	UC Davis	D8XSMGH1