


Assay Class: DNA 12000
Data Path: C:\...ents and Settings\Bioanalyzer\2013-11-20\2013-11-20_004.xad

Created: 11/20/2013 3:03:07 PM
Modified: 11/20/2013 3:33:40 PM

Electrophoresis File Run Summary

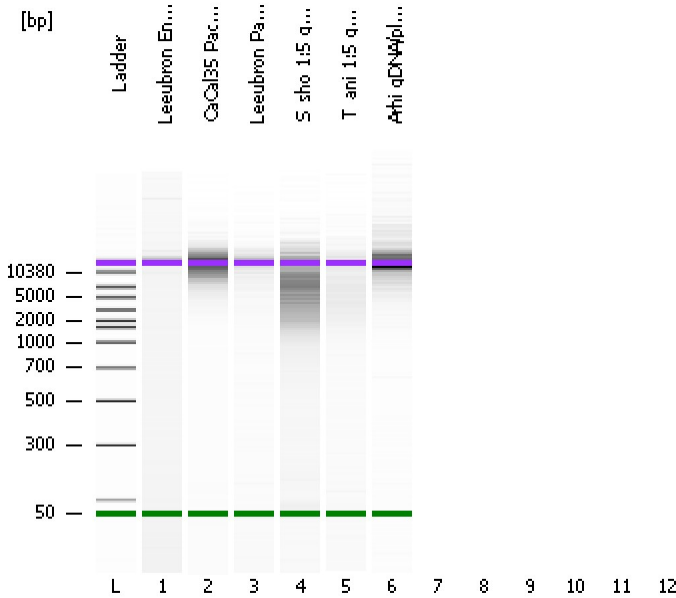

Instrument Information:

Instrument Name: DE13701086 Firmware: C.01.069
Serial#: DE13701086 Type: G2938B

Assay Information:

Assay Origin Path: C:\Program Files\Agilent\2100 bioanalyzer\2100 expert\assays\dsDNA\DNA 12000 Series II.xsy


Assay Class: DNA 12000
Version: 2.4
Assay Comments: DNA Analysis 100 -12000 bp

© Copyright 2003-2009 Agilent Technologies, Inc.


Chip Information:

Chip Lot #:
Reagent Kit Lot #:
Chip Comments:


Leeubron End Repair (1:2)

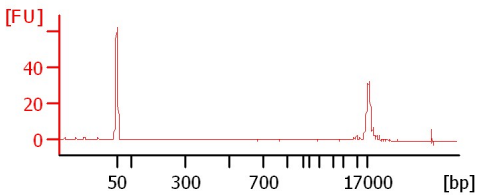


CaCal35 PacBio Lib 1:3

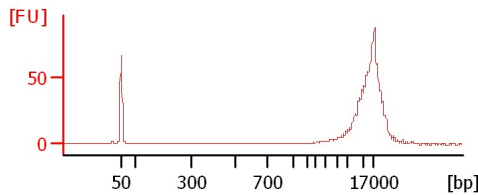


Leeubron PacBio Lib

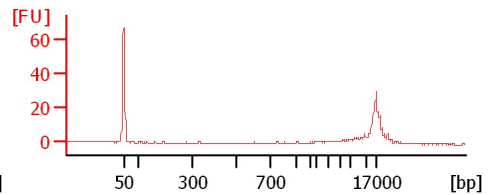


S_sho 1:5 gDNA

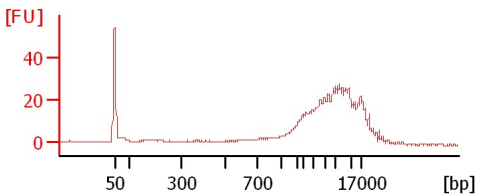

T_ani 1:5 gDNA

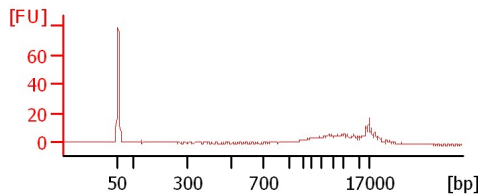

Arhi gDNA/plasmids 1:5

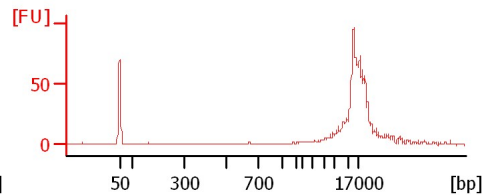

Assay Class: DNA 12000
Data Path: C:\...ents and Settings\Bioanalyzer\2013-11-20\2013-11-20_004.xad

Created: 11/20/2013 3:03:07 PM
Modified: 11/20/2013 3:33:40 PM

Electrophoresis File Run Summary (Chip Summary)

Sample Name	Sample Comment	Rest. Digest	Status	Observation	Result Label	Result Color
Leeubron End Repair (1:2)		<input type="checkbox"/>	✓			
CaCal35 PacBio Lib 1:3		<input type="checkbox"/>	✓			
Leeubron PacBio Lib		<input type="checkbox"/>	✓			
S_sho 1:5 gDNA		<input type="checkbox"/>	✓			
T_ani 1:5 gDNA		<input type="checkbox"/>	✓			
Arhi gDNA/plasmids 1:5		<input type="checkbox"/>	✓			
		<input type="checkbox"/>				
		<input type="checkbox"/>				
		<input type="checkbox"/>				
		<input type="checkbox"/>				
		<input type="checkbox"/>				
Ladder		<input type="checkbox"/>	✓			
Chip Lot #				Reagent Kit Lot #		

Chip Comments :

Assay Class: DNA 12000
Data Path: C:\...ents and Settings\Bioanalyzer\2013-11-20\2013-11-20_004.xad

Created: 11/20/2013 3:03:07 PM
Modified: 11/20/2013 3:33:40 PM

Electrophoresis Assay Details

General Analysis Settings

Number of Available Sample and Ladder Wells (Max.) : 13
Minimum Visible Range [s] : 20
Maximum Visible Range [s] : 99
Start Analysis Time Range [s] : 20
End Analysis Time Range [s] : 98.95
Ladder Concentration [ng/μl] : 44
Uses Standard Area for Ladder Fragments
Lower Marker Concentration [ng/μl] : 8.3
Upper Marker Concentration [ng/μl] : 4.2
Used Upper Marker for Quantitation
This is a Qualitative Assay Only
Standard Curve Fit is Point to Point
Show Data Aligned to Lower and Upper Marker

Integrator Settings

Integration Start Time [s] : 20
Integration End Time [s] : 98.95
Slope Threshold : 0.8
Height Threshold [FU] : 20
Area Threshold : 0.1
Width Threshold [s] : 0.5
Baseline Plateau [s] : 0.5

Filter Settings

Filter Width [s] : 0.5
Polynomial Order : 4


Ladder

Ladder Peak	Size	Area
1	50	120
2	100	48
3	300	61
4	500	77
5	700	81
6	1000	86
7	1500	92
8	2000	92
9	3000	96
10	5000	100
11	7000	98
12	10380	102
13	17000	110

Assay Class: DNA 12000
 Data Path: C:\...ents and Settings\Bioanalyzer\2013-11-20\2013-11-20_004.xad

Created: 11/20/2013 3:03:07 PM
 Modified: 11/20/2013 3:33:40 PM

Electropherogram Summary

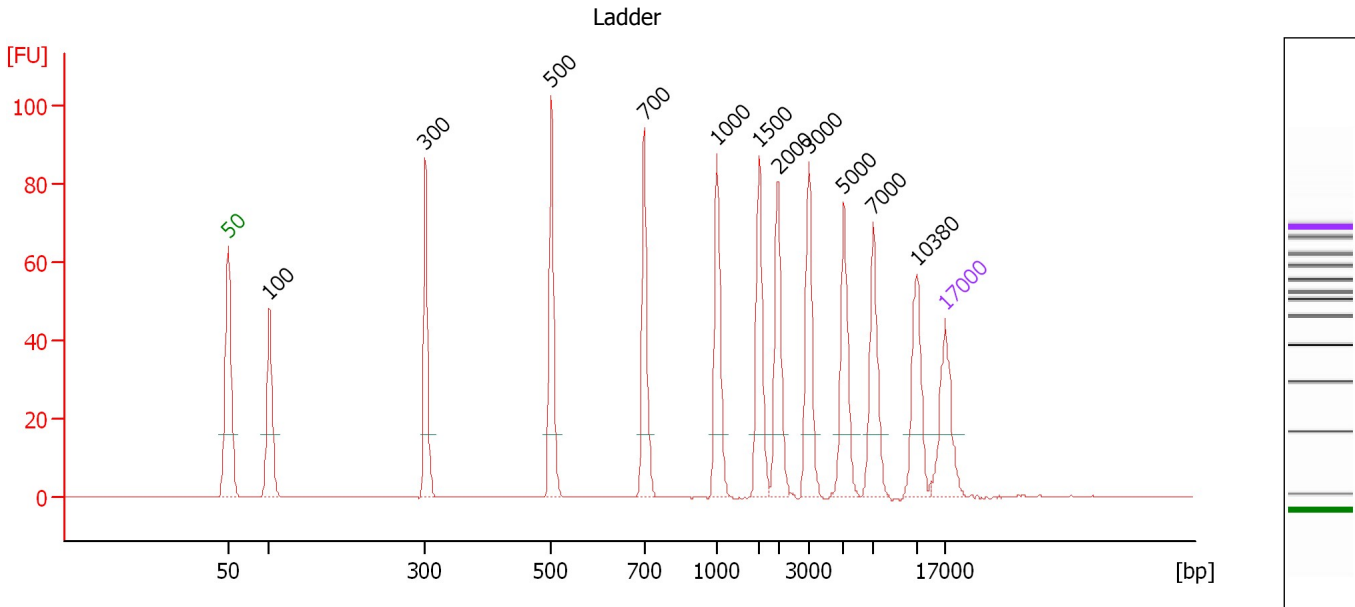

Peak table for Ladder


Peak	Size [bp]	Conc. [ng/μl]	Molarity [nmol/l]	Observations
1	50	8.30	251.5	Lower Marker
2	100	4.00	60.6	Ladder Peak
3	300	4.00	20.2	Ladder Peak
4	500	4.00	12.1	Ladder Peak
5	700	4.00	8.7	Ladder Peak
6	1,000	4.00	6.1	Ladder Peak
7	1,500	4.00	4.0	Ladder Peak
8	2,000	4.00	3.0	Ladder Peak
9	3,000	4.00	2.0	Ladder Peak
10	5,000	4.00	1.2	Ladder Peak
11	7,000	4.00	0.9	Ladder Peak
12	10,380	4.00	0.6	Ladder Peak
13	17,000	4.20	0.4	Upper Marker

Assay Class: DNA 12000
Data Path: C:\...ents and Settings\Bioanalyzer\2013-11-20\2013-11-20_004.xad

Created: 11/20/2013 3:03:07 PM
Modified: 11/20/2013 3:33:40 PM

Electropherogram Summary Continued ...

Leeubron End Repair (1:2)

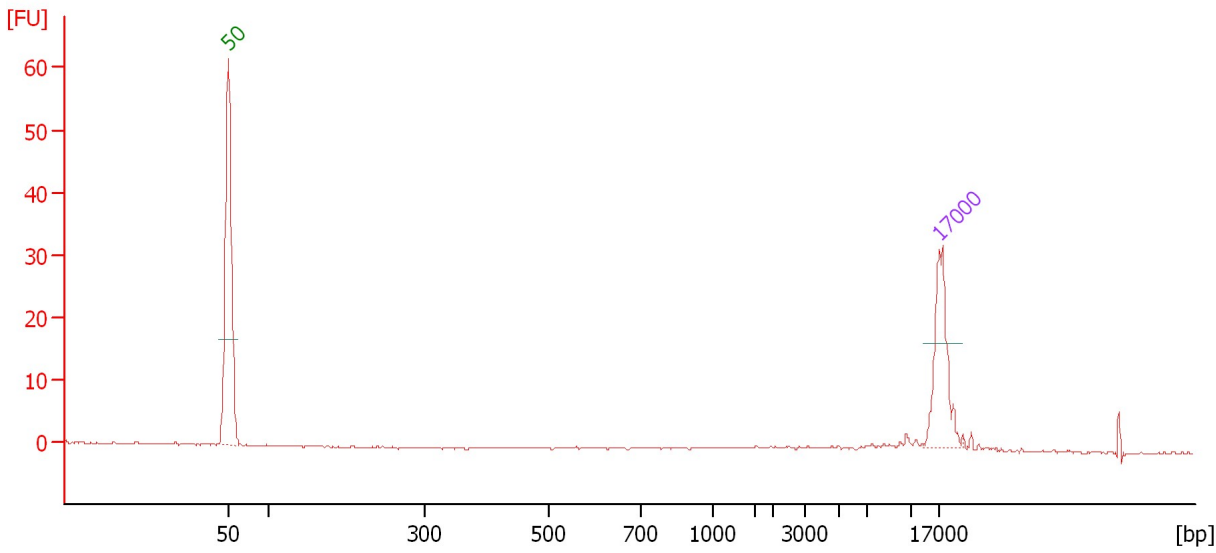

Overall Results for sample 1 : Leeubron End Repair (1:2)

Number of peaks found: 0

Peak table for sample 1 : Leeubron End Repair (1:2)


Peak	Size [bp]	Conc. [ng/μl]	Molarity [nmol/l]	Observations
1	50	8.30	251.5	Lower Marker
2	17,000	4.20	0.4	Upper Marker

Assay Class: DNA 12000
Data Path: C:\...ents and Settings\Bioanalyzer\2013-11-20\2013-11-20_004.xad

Created: 11/20/2013 3:03:07 PM
Modified: 11/20/2013 3:33:40 PM

Electropherogram Summary Continued ...

CaCal35 PacBio Lib 1:3

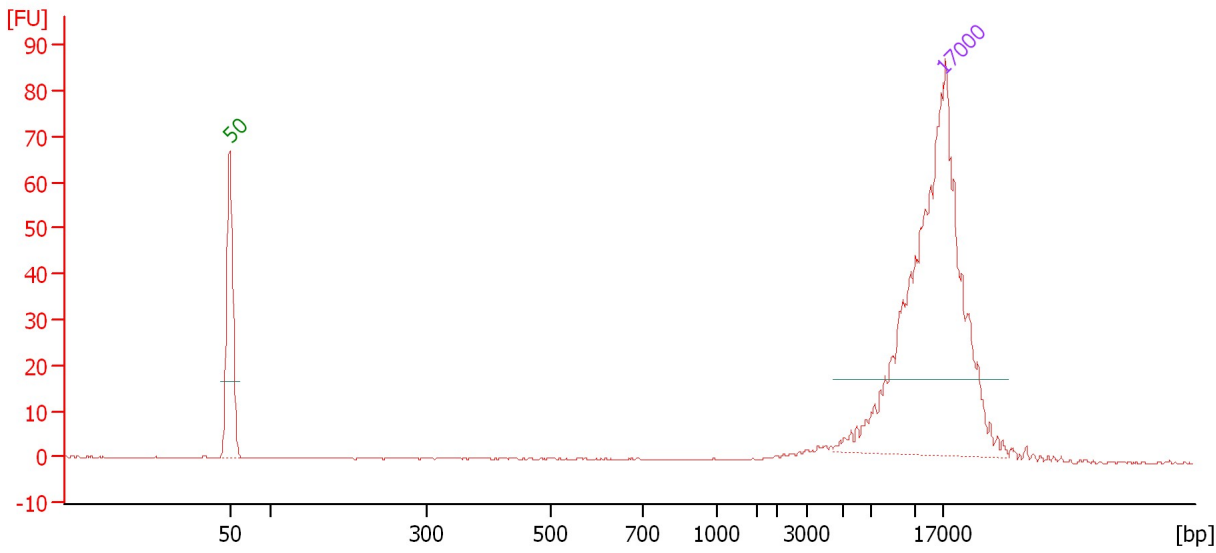

Overall Results for sample 2 : CaCal35 PacBio Lib 1:3

Number of peaks found: 0

Peak table for sample 2 : CaCal35 PacBio Lib 1:3


Peak	Size [bp]	Conc. [ng/μl]	Molarity [nmol/l]	Observations
1	50	8.30	251.5	Lower Marker
2	17,000	4.20	0.4	Upper Marker

Assay Class: DNA 12000
Data Path: C:\...ents and Settings\Bioanalyzer\2013-11-20\2013-11-20_004.xad

Created: 11/20/2013 3:03:07 PM
Modified: 11/20/2013 3:33:40 PM

Electropherogram Summary Continued ...

Leeubron PacBio Lib

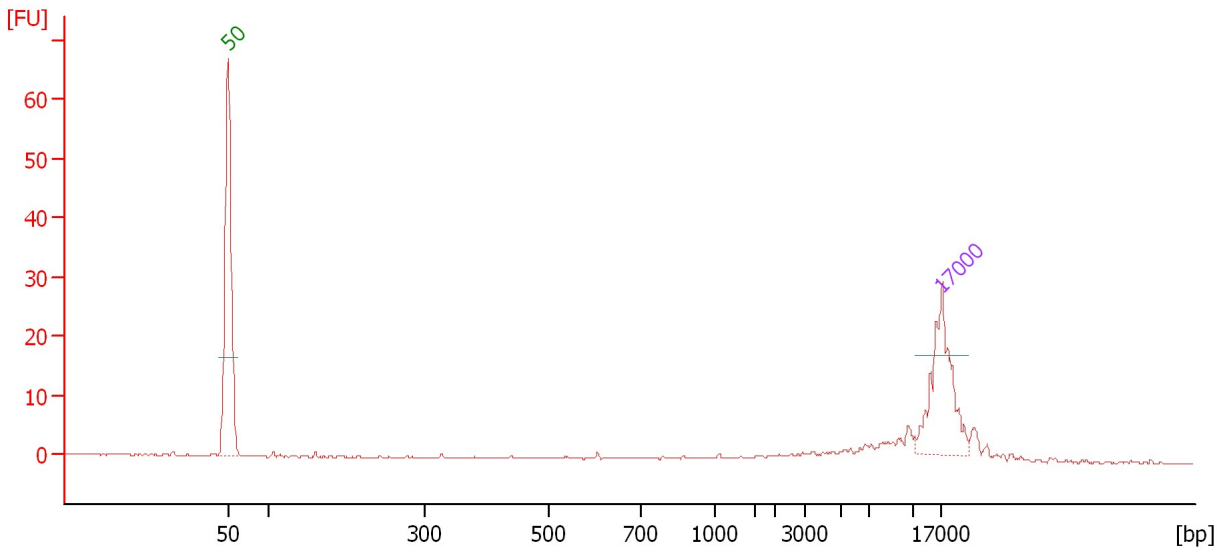

Overall Results for sample 3 : Leeubron PacBio Lib

Number of peaks found: 0

Peak table for sample 3 : Leeubron PacBio Lib


Peak	Size [bp]	Conc. [ng/μl]	Molarity [nmol/l]	Observations
1	50	8.30	251.5	Lower Marker
2	17,000	4.20	0.4	Upper Marker

Assay Class: DNA 12000
 Data Path: C:\...ents and Settings\Bioanalyzer\2013-11-20\2013-11-20_004.xad

Created: 11/20/2013 3:03:07 PM
 Modified: 11/20/2013 3:33:40 PM

Electropherogram Summary Continued ...

S_sho 1:5 gDNA

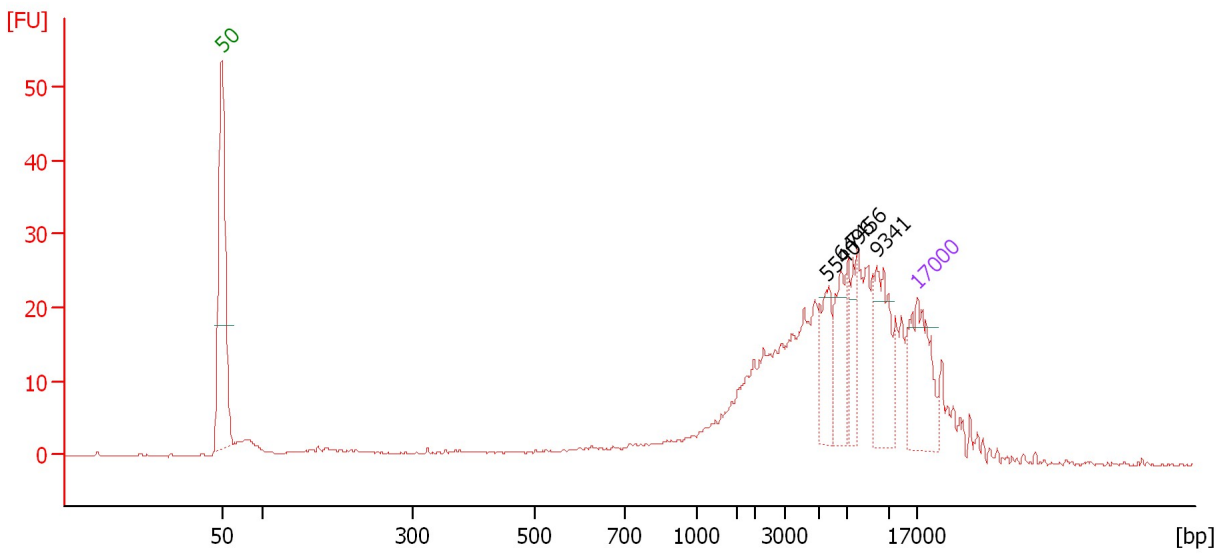

Overall Results for sample 4 : S_sho 1:5 gDNA

Number of peaks found: 4

Peak table for sample 4 : S_sho 1:5 gDNA


Peak	Size [bp]	Conc. [ng/μl]	Molarity [nmol/l]	Observations
1	50	8.30	251.5	Lower Marker
2	5,540	2.38	0.7	
3	6,496	2.54	0.6	
4	7,456	1.52	0.3	
5	9,341	3.84	0.6	
6	17,000	4.20	0.4	Upper Marker

Assay Class: DNA 12000
Data Path: C:\...ents and Settings\Bioanalyzer\2013-11-20\2013-11-20_004.xad

Created: 11/20/2013 3:03:07 PM
Modified: 11/20/2013 3:33:40 PM

Electropherogram Summary Continued ...

T_ani 1:5 gDNA

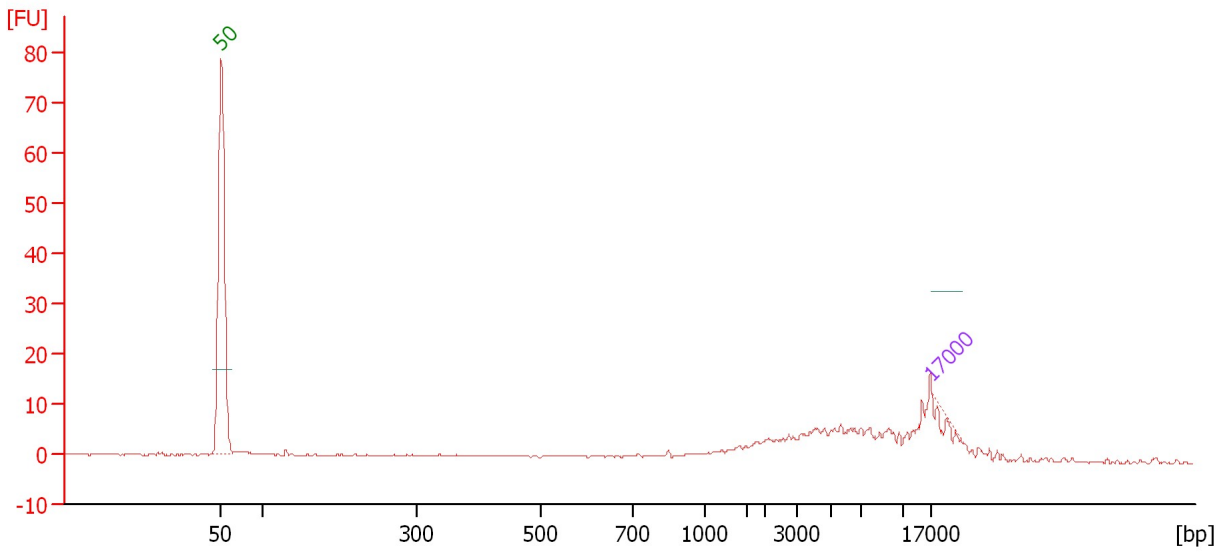

Overall Results for sample 5 : T_ani 1:5 gDNA

Number of peaks found: 0

Peak table for sample 5 : T_ani 1:5 gDNA


Peak	Size [bp]	Conc. [ng/μl]	Molarity [nmol/l]	Observations
1	50	8.30	251.5	Lower Marker
2	17,000	4.20	0.4	Upper Marker

Assay Class: DNA 12000
Data Path: C:\...ents and Settings\Bioanalyzer\2013-11-20\2013-11-20_004.xad

Created: 11/20/2013 3:03:07 PM
Modified: 11/20/2013 3:33:40 PM

Electropherogram Summary Continued ...

Arhi gDNA/plasmids 1:5

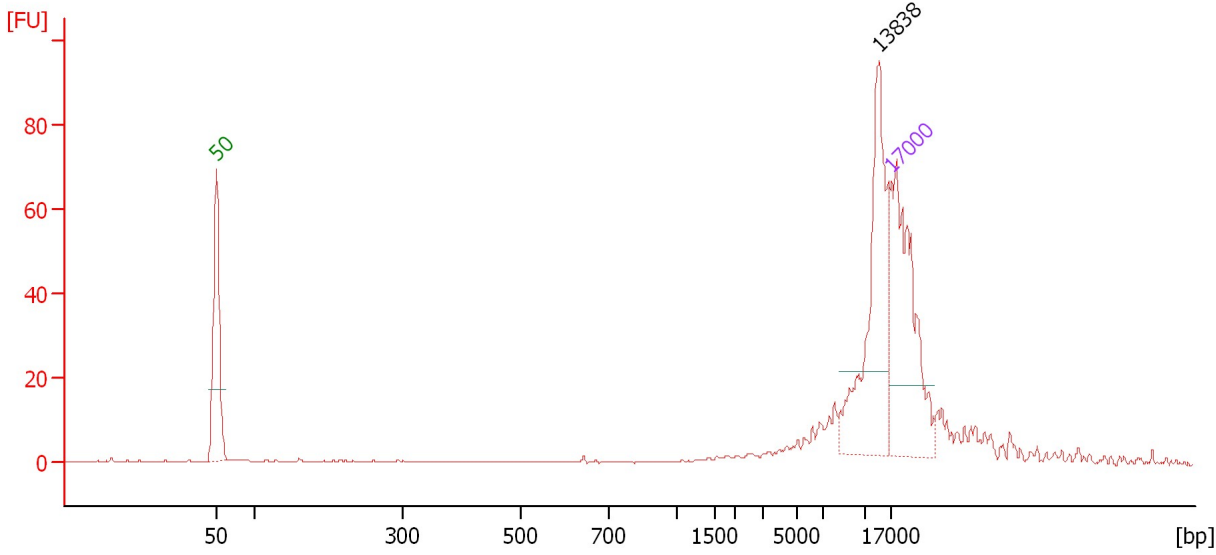

Overall Results for sample 6 : Arhi gDNA/plasmids 1:5

Number of peaks found: 1


Peak table for sample 6 : Arhi gDNA/plasmids 1:5

Peak	Size [bp]	Conc. [ng/μl]	Molarity [nmol/l]	Observations
1	50	8.30	251.5	Lower Marker
2	13,838	4.17	0.5	
3	17,000	4.20	0.4	Upper Marker

Assay Class: DNA 12000
Data Path: C:\...ents and Settings\Bioanalyzer\2013-11-20\2013-11-20_004.xad

Created: 11/20/2013 3:03:07 PM
Modified: 11/20/2013 3:33:40 PM

Gel Image

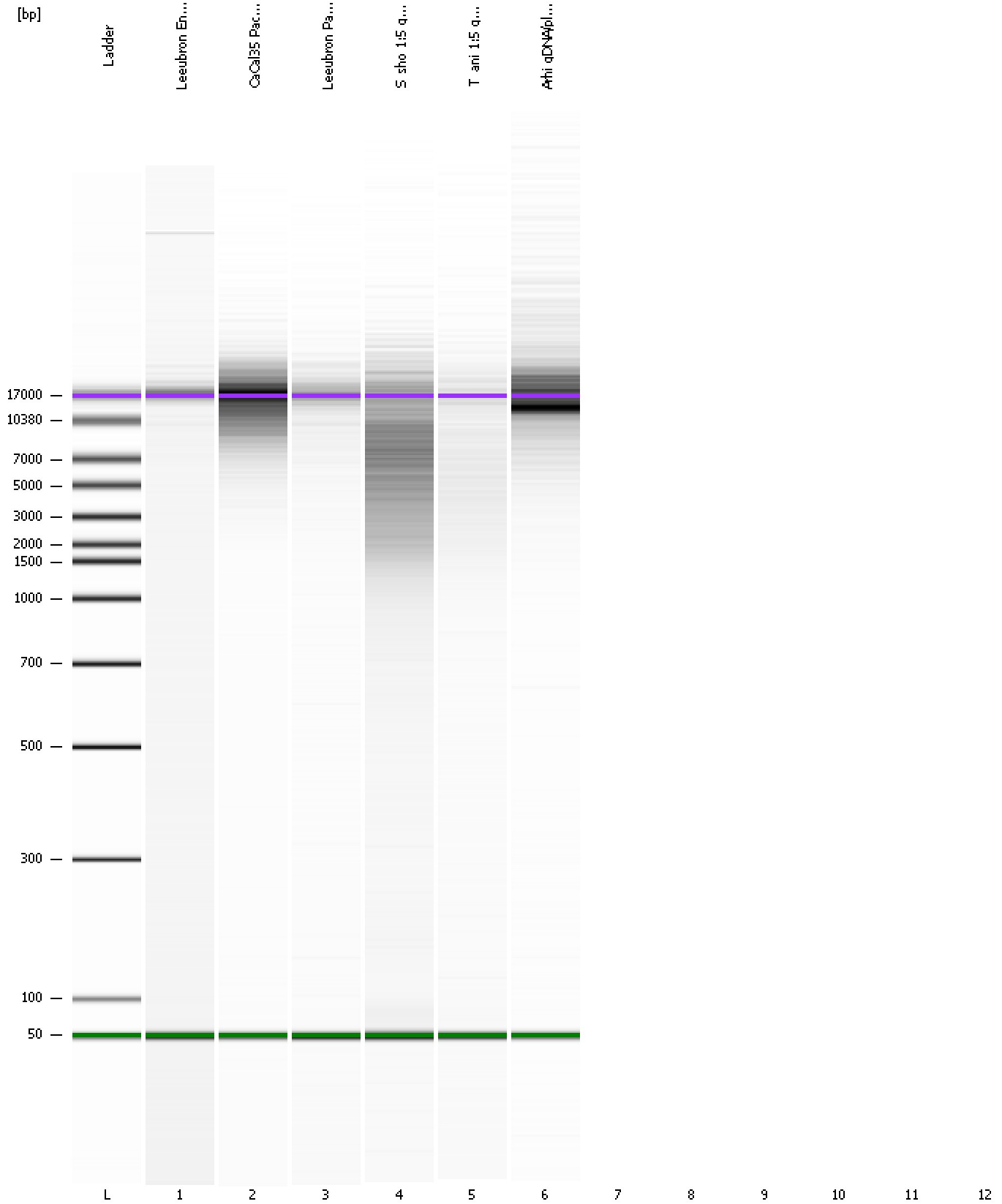

Assay Class: DNA 12000
Data Path: C:\...ents and Settings\Bioanalyzer\2013-11-20\2013-11-20_004.xad

Created: 11/20/2013 3:03:07 PM
Modified: 11/20/2013 3:33:40 PM

Invalid Samples

Sample 7 has not been run, no results available.

Sample 8 has not been run, no results available.

Sample 9 has not been run, no results available.

Sample 10 has not been run, no results available.

Sample 11 has not been run, no results available.

Sample 12 has not been run, no results available.