


Assay Class: High Sensitivity DNA Assay  
Data Path: C:\...ings\Bioanalyzer\2014-01-22\2014-01-22\_007\_D10-D13\_Pool.xad

Created: 1/22/2014 5:29:09 PM  
Modified: 1/22/2014 6:14:08 PM

**Electrophoresis File Run Summary**

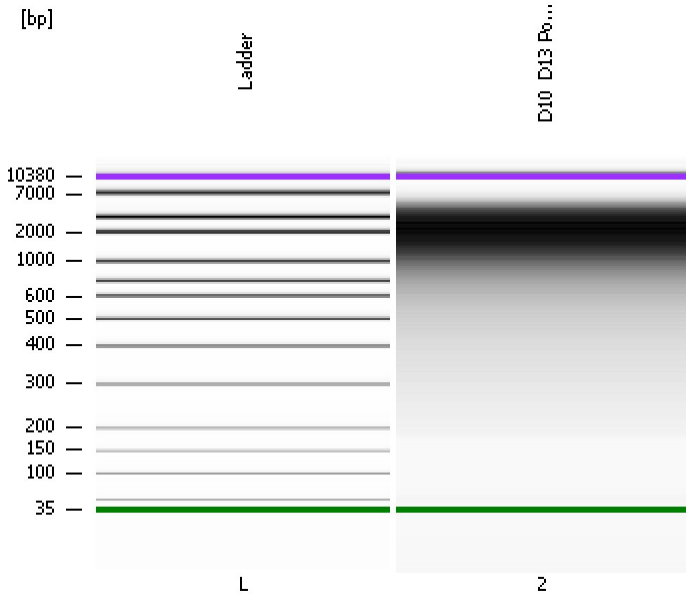

Instrument Information:

Instrument Name: DE13701086      Firmware: C.01.069  
Serial#: DE13701086      Type: G2938B


Assay Information:

Assay Origin Path: C:\Program Files\Agilent\2100 bioanalyzer\2100 expert\assays\dsDNA\High Sensitivity DNA.xsy  
Assay Class: High Sensitivity DNA Assay  
Version: 1.03  
Assay Comments: Copyright © 2003-2010 Agilent Technologies

Chip Information:

Chip Lot #:  
Reagent Kit Lot #:  
Chip Comments:

**D10\_D13 Pool 1:2**

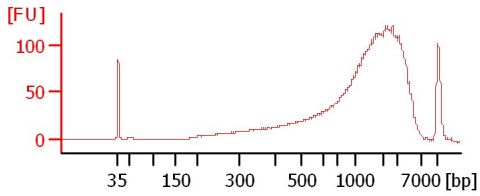

Assay Class: High Sensitivity DNA Assay  
Data Path: C:\...ings\Bioanalyzer\2014-01-22\2014-01-22\_007\_D10-D13\_Pool.xad

Created: 1/22/2014 5:29:09 PM  
Modified: 1/22/2014 6:14:08 PM

**Electrophoresis File Run Summary (Chip Summary)**

Sample Name	Sample Comment	Rest. Digest	Status	Observation	Result Label	Result Color
D10_D13 Pool 1:2		<input type="checkbox"/>	✓			
Ladder		<input type="checkbox"/>	✓			
<b>Chip Lot #</b>				<b>Reagent Kit Lot #</b>		

**Chip Comments :**

Assay Class: High Sensitivity DNA Assay  
Data Path: C:\...ings\Bioanalyzer\2014-01-22\2014-01-22\_007\_D10-D13\_Pool.xad

Created: 1/22/2014 5:29:09 PM  
Modified: 1/22/2014 6:14:08 PM

**Electrophoresis Assay Details**

**General Analysis Settings**

Number of Available Sample and Ladder Wells (Max.) : 12  
Minimum Visible Range [s] : 32  
Maximum Visible Range [s] : 138  
Start Analysis Time Range [s] : 33  
End Analysis Time Range [s] : 137.5  
Ladder Concentration [pg/μl] : 1950  
Uses Standard Area for Ladder Fragments  
Lower Marker Concentration [pg/μl] : 125  
Upper Marker Concentration [pg/μl] : 75  
Used Upper Marker for Quantitation  
Standard Curve Fit is Point to Point  
Show Data Aligned to Lower and Upper Marker

**Integrator Settings**

Integration Start Time [s] : 33.05  
Integration End Time [s] : 137  
Slope Threshold : 0.8  
Height Threshold [FU] : 5  
Area Threshold : 0.1  
Width Threshold [s] : 0.6  
Baseline Plateau [s] : 0.5

**Filter Settings**

Filter Width [s] : 0.5  
Polynomial Order : 4


**Ladder**

Ladder Peak	Size	Area
1	35	160
2	50	210
3	100	208
4	150	221
5	200	242
6	300	270
7	400	305
8	500	306
9	600	336
10	700	321
11	1000	366
12	2000	413
13	3000	411
14	7000	400
15	10380	214

Assay Class: High Sensitivity DNA Assay  
 Data Path: C:\...ings\Bioanalyzer\2014-01-22\2014-01-22\_007\_D10-D13\_Pool.xad

Created: 1/22/2014 5:29:09 PM  
 Modified: 1/22/2014 6:14:08 PM

**Electropherogram Summary**

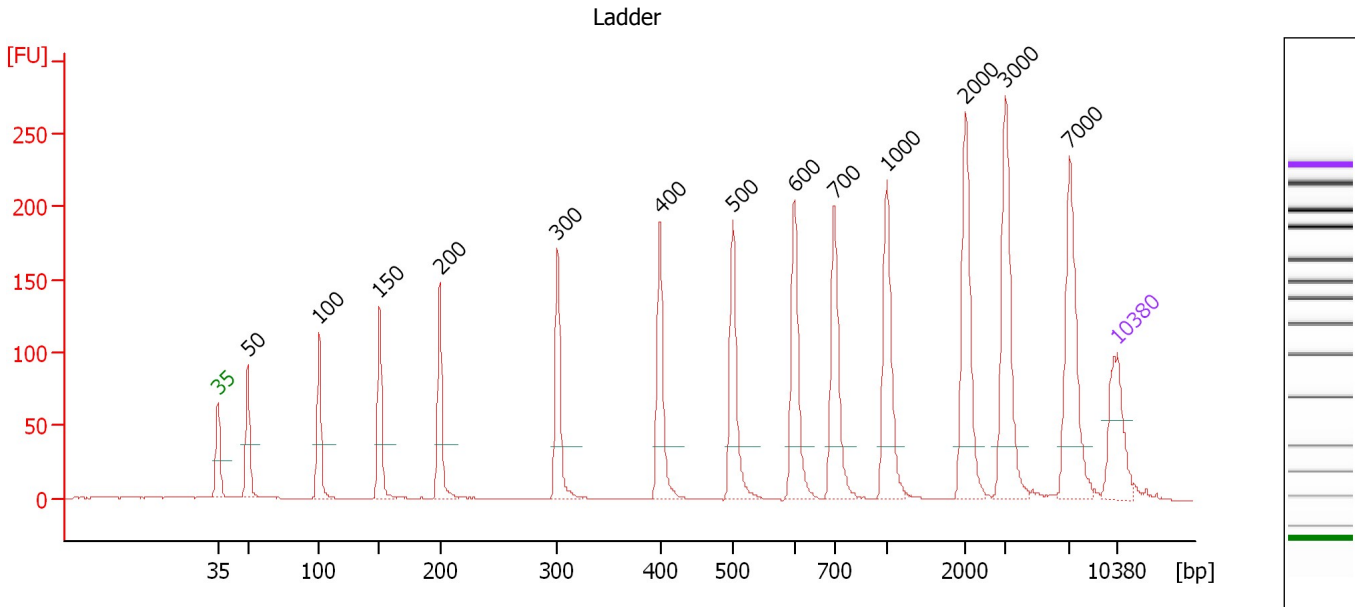

**Overall Results for Ladder**

Noise: 0.2


**Peak table for Ladder**

Peak	Aligned Migration Time [s]	Size [bp]	Conc. [pg/μl]	Molarity [pmol/l]	Observations
1	43.00	35	125.00	5,411.3	Lower Marker
2	45.33	50	150.00	4,545.5	Ladder Peak
3	50.89	100	150.00	2,272.7	Ladder Peak
4	55.58	150	150.00	1,515.2	Ladder Peak
5	60.27	200	150.00	1,136.4	Ladder Peak
6	69.45	300	150.00	757.6	Ladder Peak
7	77.42	400	150.00	568.2	Ladder Peak
8	83.11	500	150.00	454.5	Ladder Peak
9	87.88	600	150.00	378.8	Ladder Peak
10	91.00	700	150.00	324.7	Ladder Peak
11	95.11	1,000	150.00	227.3	Ladder Peak
12	101.25	2,000	150.00	113.6	Ladder Peak
13	104.32	3,000	150.00	75.8	Ladder Peak
14	109.30	7,000	150.00	32.5	Ladder Peak
15	113.00	10,380	75.00	10.9	Upper Marker

Assay Class: High Sensitivity DNA Assay  
 Data Path: C:\...ings\Bioanalyzer\2014-01-22\2014-01-22\_007\_D10-D13\_Pool.xad

Created: 1/22/2014 5:29:09 PM  
 Modified: 1/22/2014 6:14:08 PM

**Electropherogram Summary Continued ...**

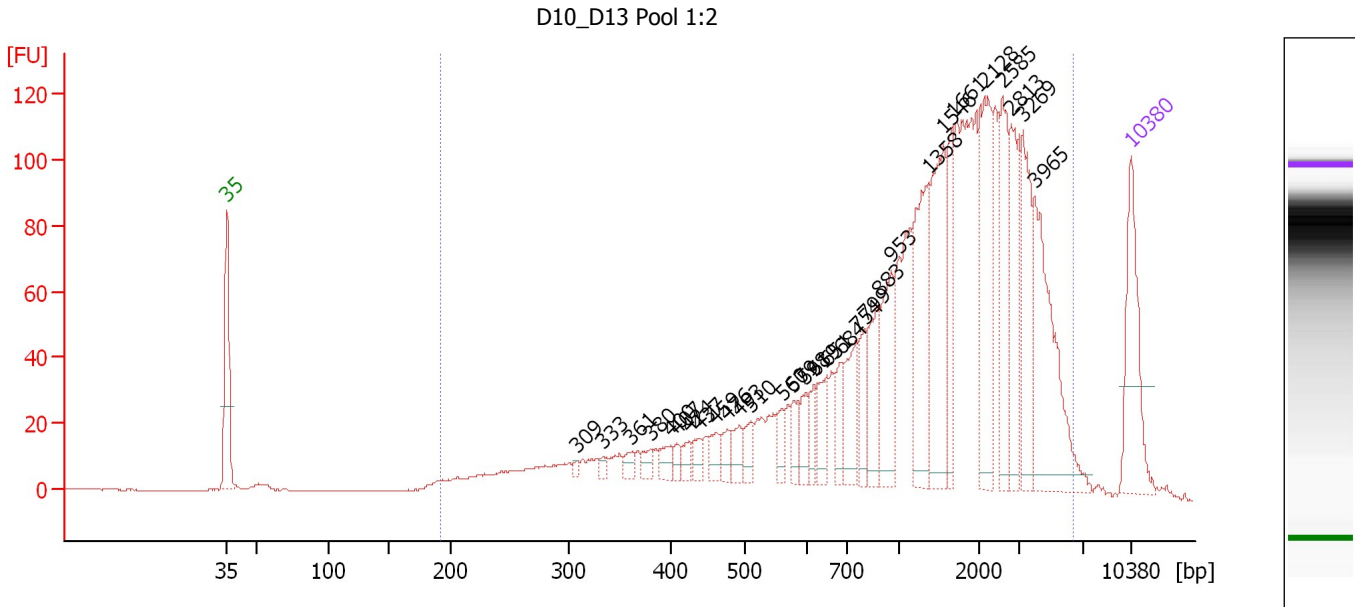

**Overall Results for sample 2 : D10 D13 Pool 1:2**

Number of peaks found: 30                      Corr. Area 1: 2,203.0  
 Noise: 0.2

**Peak table for sample 2 : D10 D13 Pool 1:2**


Peak	Aligned Migration Time [s]	Size [bp]	Conc. [pg/μl]	Molarity [pmol/l]	Observations
1	43.00	35	125.00	5,411.3	Lower Marker
2	70.13	309	4.69	23.0	
3	72.12	333	6.07	27.6	
4	74.35	361	10.58	44.4	
5	75.83	380	11.81	47.1	
6	77.40	400	15.87	60.2	
7	77.81	407	8.26	30.8	
8	78.81	424	13.11	46.8	
9	79.51	437	12.90	44.8	
10	80.79	459	17.69	58.4	
11	81.74	476	16.04	51.1	
12	82.69	493	20.56	63.2	
13	83.60	510	19.04	56.5	
14	85.99	560	16.51	44.6	
15	86.90	579	18.10	47.3	
16	87.77	598	21.43	54.3	
17	88.47	619	17.75	43.5	
18	89.46	651	31.32	72.9	
19	90.49	684	24.23	53.7	
20	91.73	754	50.18	100.9	
21	92.35	799	33.44	63.4	
22	93.51	883	48.74	83.6	
23	94.46	953	66.85	106.3	
24	97.31	1,358	104.50	116.6	
25	98.46	1,546	110.07	107.8	

Assay Class: High Sensitivity DNA Assay  
 Data Path: C:\...ings\Bioanalyzer\2014-01-22\2014-01-22\_007\_D10-D13\_Pool.xad


Created: 1/22/2014 5:29:09 PM  
 Modified: 1/22/2014 6:14:08 PM

### Electropherogram Summary Continued ...

#### ... Peak table for sample 2 : D10 D13 Pool 1:2

Peak	Aligned Migration Time [s]	Size [bp]	Conc. [pg/ $\mu$ l]	Molarity [pmol/l]	Observations
26	99.17	1,661	48.94	44.7	
27	101.64	2,128	103.94	74.0	
28	103.05	2,585	74.68	43.8	
29	103.75	2,813	59.88	32.3	
30	104.66	3,269	73.92	34.3	
31	105.53	3,965	126.26	48.2	
32	 113.00	10,380	75.00	10.9	Upper Marker

#### Region table for sample 2 : D10 D13 Pool 1:2

From [bp]	To [bp]	Average Size [bp]	Molarity [pmol/l]	Conc. [pg/ $\mu$ l]	Corr. Area	% of Total	Size distribution in CV [%]	Color
192	6,377	1,745	3,851.0	1,989.46	2,203.0	98	73.2	

Assay Class: High Sensitivity DNA Assay  
Data Path: C:\...ings\Bioanalyzer\2014-01-22\2014-01-22\_007\_D10-D13\_Pool.xad

Created: 1/22/2014 5:29:09 PM  
Modified: 1/22/2014 6:14:08 PM

**Gel Image**

