


Assay Class: DNA 12000  
Data Path: C:\...ents and Settings\Bioanalyzer\2014-10-17\2014-10-17\_002.xad

Created: 10/17/2014 12:21:20 PM  
Modified: 10/17/2014 12:42:30 PM

**Electrophoresis File Run Summary**

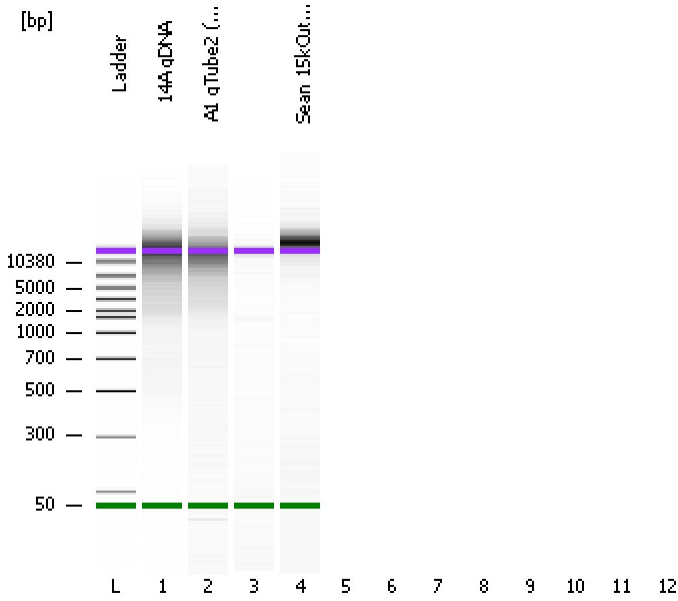

Instrument Information:

Instrument Name: DE13701086      Firmware: C.01.069  
Serial#: DE13701086      Type: G2938B

Assay Information:

Assay Origin Path: C:\Program Files\Agilent\2100 bioanalyzer\2100 expert\assays\dsDNA\DNA 12000 Series II.xsy

Assay Class: DNA 12000

Version: 2.4

Assay Comments: DNA Analysis 100 -12000 bp

© Copyright 2003-2009 Agilent Technologies, Inc.


Chip Information:

Chip Lot #:



Reagent Kit Lot #:

Chip Comments:


**14A gDNA**

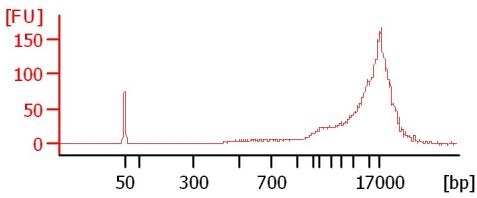

**A1 gTube2 (4000 rpm)**

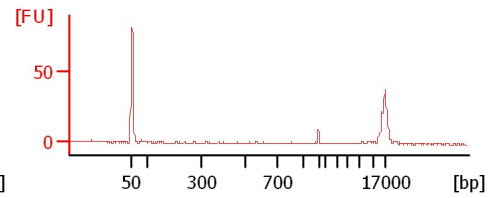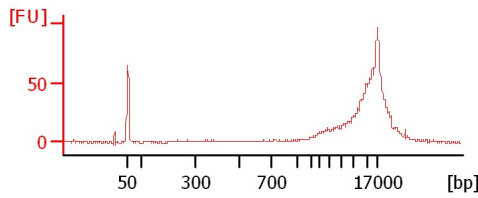

**Sean 15kCut lib**

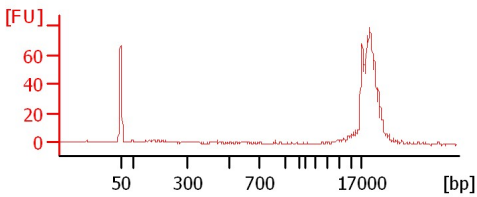

Assay Class: DNA 12000  
 Data Path: C:\...ents and Settings\Bioanalyzer\2014-10-17\2014-10-17\_002.xad

Created: 10/17/2014 12:21:20 PM  
 Modified: 10/17/2014 12:42:30 PM

**Electrophoresis File Run Summary (Chip Summary)**

Sample Name	Sample Comment	Rest. Digest	Status	Observation	Result Label	Result Color
14A gDNA		<input type="checkbox"/>	✓			
A1 gTube2 (4000 rpm)		<input type="checkbox"/>	✓			
		<input type="checkbox"/>	✓			
Sean 15kCut lib		<input type="checkbox"/>	✓			
		<input type="checkbox"/>				
		<input type="checkbox"/>				
		<input type="checkbox"/>				
		<input type="checkbox"/>				
		<input type="checkbox"/>				
		<input type="checkbox"/>				
Ladder		<input type="checkbox"/>	✓			
<b>Chip Lot #</b>				<b>Reagent Kit Lot #</b>		

Chip Comments :

Assay Class: DNA 12000  
Data Path: C:\...ents and Settings\Bioanalyzer\2014-10-17\2014-10-17\_002.xad

Created: 10/17/2014 12:21:20 PM  
Modified: 10/17/2014 12:42:30 PM

**Electrophoresis Assay Details**

**General Analysis Settings**

Number of Available Sample and Ladder Wells (Max.) : 13  
Minimum Visible Range [s] : 20  
Maximum Visible Range [s] : 99  
Start Analysis Time Range [s] : 20  
End Analysis Time Range [s] : 98.95  
Ladder Concentration [ng/μl] : 44  
Uses Standard Area for Ladder Fragments  
Lower Marker Concentration [ng/μl] : 8.3  
Upper Marker Concentration [ng/μl] : 4.2  
Used Upper Marker for Quantitation  
This is a Qualitative Assay Only  
Standard Curve Fit is Point to Point  
Show Data Aligned to Lower and Upper Marker

**Integrator Settings**

Integration Start Time [s] : 20  
Integration End Time [s] : 98.95  
Slope Threshold : 0.8  
Height Threshold [FU] : 20  
Area Threshold : 0.1  
Width Threshold [s] : 0.5  
Baseline Plateau [s] : 0.5

**Filter Settings**

Filter Width [s] : 0.5  
Polynomial Order : 4


**Ladder**

Ladder Peak	Size	Area
1	50	120
2	100	48
3	300	61
4	500	77
5	700	81
6	1000	86
7	1500	92
8	2000	92
9	3000	96
10	5000	100
11	7000	98
12	10380	102
13	17000	110

Assay Class: DNA 12000  
 Data Path: C:\...ents and Settings\Bioanalyzer\2014-10-17\2014-10-17\_002.xad

Created: 10/17/2014 12:21:20 PM  
 Modified: 10/17/2014 12:42:30 PM

**Electropherogram Summary**

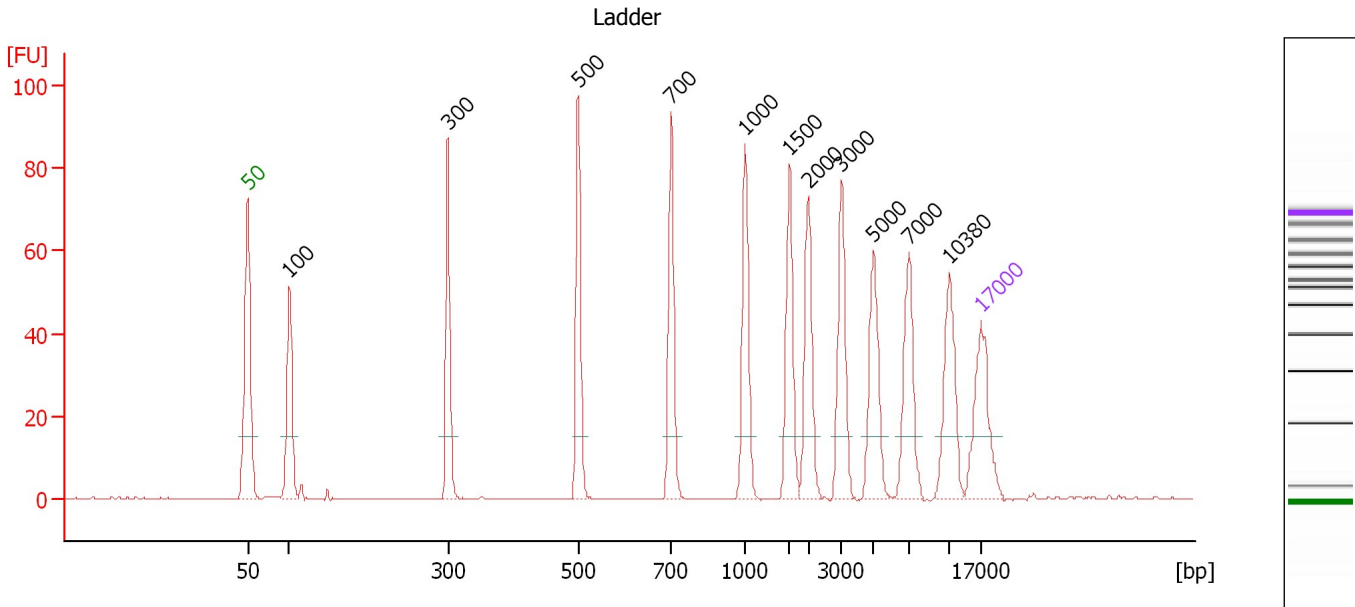


**Peak table for Ladder**

Peak	Size [bp]	Conc. [ng/μl]	Molarity [nmol/l]	Observations
1	50	8.30	251.5	Lower Marker
2	100	4.00	60.6	Ladder Peak
3	300	4.00	20.2	Ladder Peak
4	500	4.00	12.1	Ladder Peak
5	700	4.00	8.7	Ladder Peak
6	1,000	4.00	6.1	Ladder Peak
7	1,500	4.00	4.0	Ladder Peak
8	2,000	4.00	3.0	Ladder Peak
9	3,000	4.00	2.0	Ladder Peak
10	5,000	4.00	1.2	Ladder Peak
11	7,000	4.00	0.9	Ladder Peak
12	10,380	4.00	0.6	Ladder Peak
13	17,000	4.20	0.4	Upper Marker

Assay Class: DNA 12000  
 Data Path: C:\...ents and Settings\Bioanalyzer\2014-10-17\2014-10-17\_002.xad

Created: 10/17/2014 12:21:20 PM  
 Modified: 10/17/2014 12:42:30 PM

**Electropherogram Summary Continued ...**

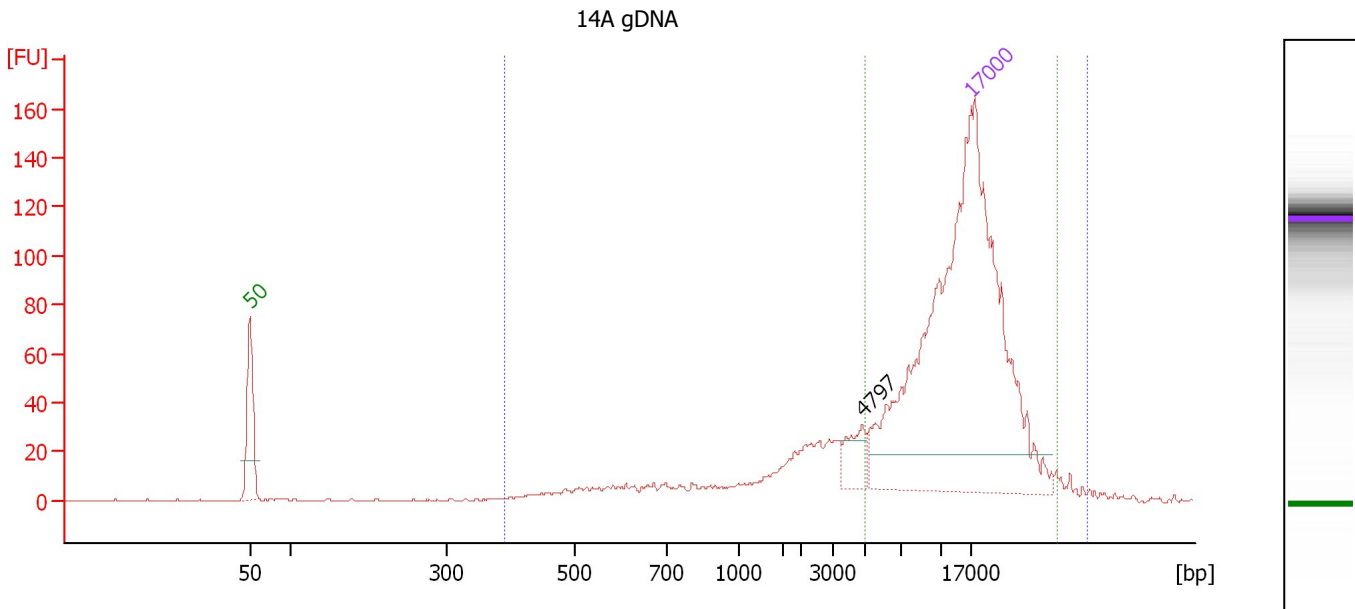

**Overall Results for sample 1 : 14A gDNA**

Number of peaks found: 1      Area 2: 10.8  
 Area 1: 243.9

**Peak table for sample 1 : 14A gDNA**

Peak	Size [bp]	Conc. [ng/μl]	Molarity [nmol/l]	Observations
1	50	8.30	251.5	Lower Marker
2	4,797	0.22	0.1	
3	17,000	4.20	0.4	Upper Marker


**Region table for sample 1 : 14A gDNA**

From [bp]	To [bp]	Average Size [bp]	Conc. [ng/μl]	Area	% of Total	Size distribution in CV [%]	Color
391	41,697	13,250	1.42	243.9	97	64.5	Blue
4,903	35,224	15,711	0.06	10.8	4	43.4	Green

Assay Class: DNA 12000  
 Data Path: C:\...ents and Settings\Bioanalyzer\2014-10-17\2014-10-17\_002.xad

Created: 10/17/2014 12:21:20 PM  
 Modified: 10/17/2014 12:42:30 PM

**Electropherogram Summary Continued ...**

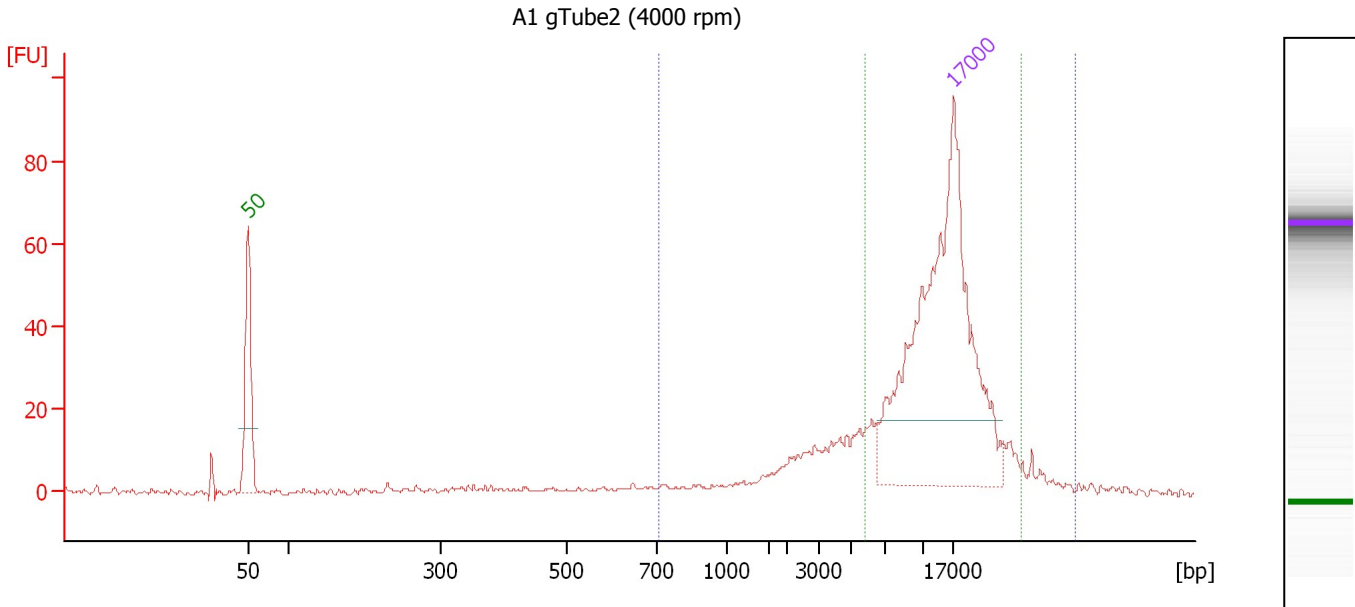

**Overall Results for sample 2 : A1 gTube2 (4000 rpm)**

Number of peaks found: 0 Area 2: 27.4  
 Area 1: 123.9

**Peak table for sample 2 : A1 gTube2 (4000 rpm)**

Peak	Size [bp]	Conc. [ng/μl]	Molarity [nmol/l]	Observations
1	50	8.30	251.5	Lower Marker
2	17,000	4.20	0.4	Upper Marker


**Region table for sample 2 : A1 gTube2 (4000 rpm)**

From [bp]	To [bp]	Average Size [bp]	Conc. [ng/μl]	Area	% of Total	Size distribution in CV [%]	Color
709	43,341	13,621	1.64	123.9	81	58.7	Blue
5,813	31,774	15,049	0.35	27.4	18	39.5	Green

Assay Class: DNA 12000  
Data Path: C:\...ents and Settings\Bioanalyzer\2014-10-17\2014-10-17\_002.xad

Created: 10/17/2014 12:21:20 PM  
Modified: 10/17/2014 12:42:30 PM

**Electropherogram Summary Continued ...**

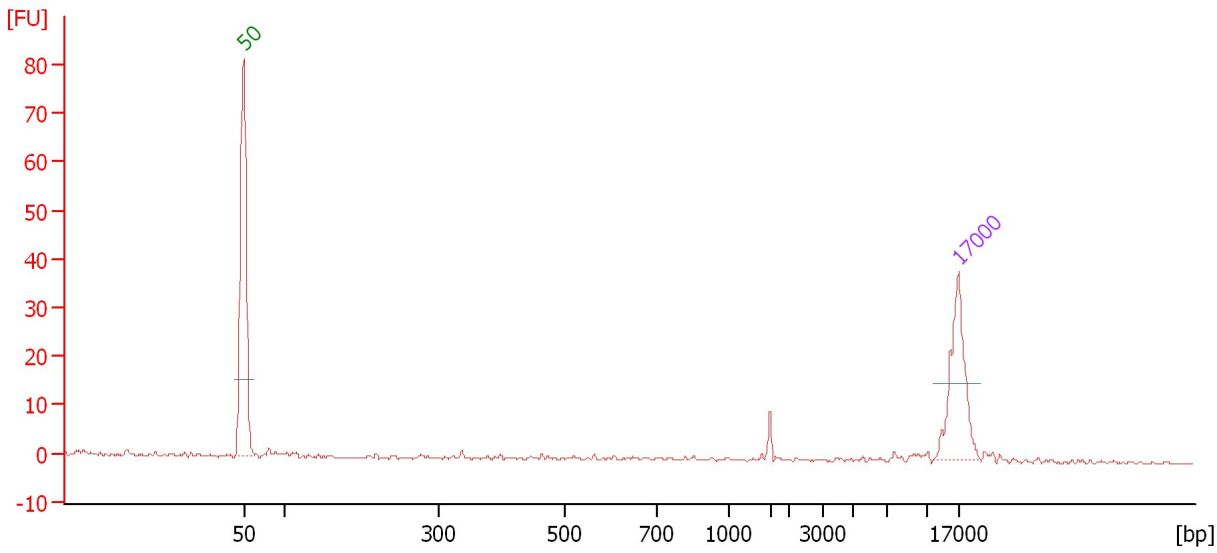

**Overall Results for sample 3 :**

Number of peaks found: 0


**Peak table for sample 3 :**

Peak	Size [bp]	Conc. [ng/μl]	Molarity [nmol/l]	Observations
1	50	8.30	251.5	Lower Marker
2	17,000	4.20	0.4	Upper Marker

Assay Class: DNA 12000  
 Data Path: C:\...ents and Settings\Bioanalyzer\2014-10-17\2014-10-17\_002.xad

Created: 10/17/2014 12:21:20 PM  
 Modified: 10/17/2014 12:42:30 PM

**Electropherogram Summary Continued ...**

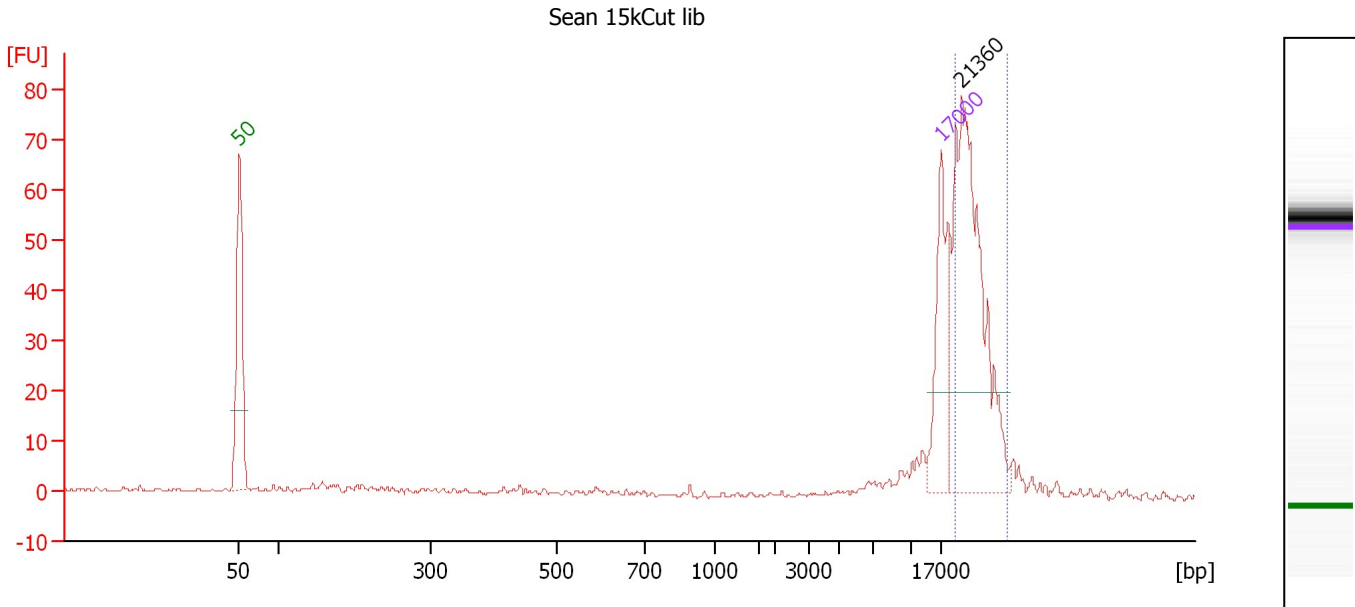

**Overall Results for sample 4 : Sean 15kCut lib**

Number of peaks found: 0 Area 1: 165.1

**Peak table for sample 4 : Sean 15kCut lib**

Peak	Size [bp]	Conc. [ng/μl]	Molarity [nmol/l]	Observations
1	50	8.30	251.5	Lower Marker
2	17,000	4.20	0.4	Upper Marker
3	21,360	0.00	0.0	

**Region table for sample 4 : Sean 15kCut lib**


From [bp]	To [bp]	Average Size [bp]	Conc. [ng/μl]	Area	% of Total	Size distribution in CV [%]	Color
19,800	31,221	23,724	10.95	165.1	71	12.0	Blue



Assay Class: DNA 12000  
Data Path: C:\...ents and Settings\Bioanalyzer\2014-10-17\2014-10-17\_002.xad

Created: 10/17/2014 12:21:20 PM  
Modified: 10/17/2014 12:42:30 PM

**Gel Image**

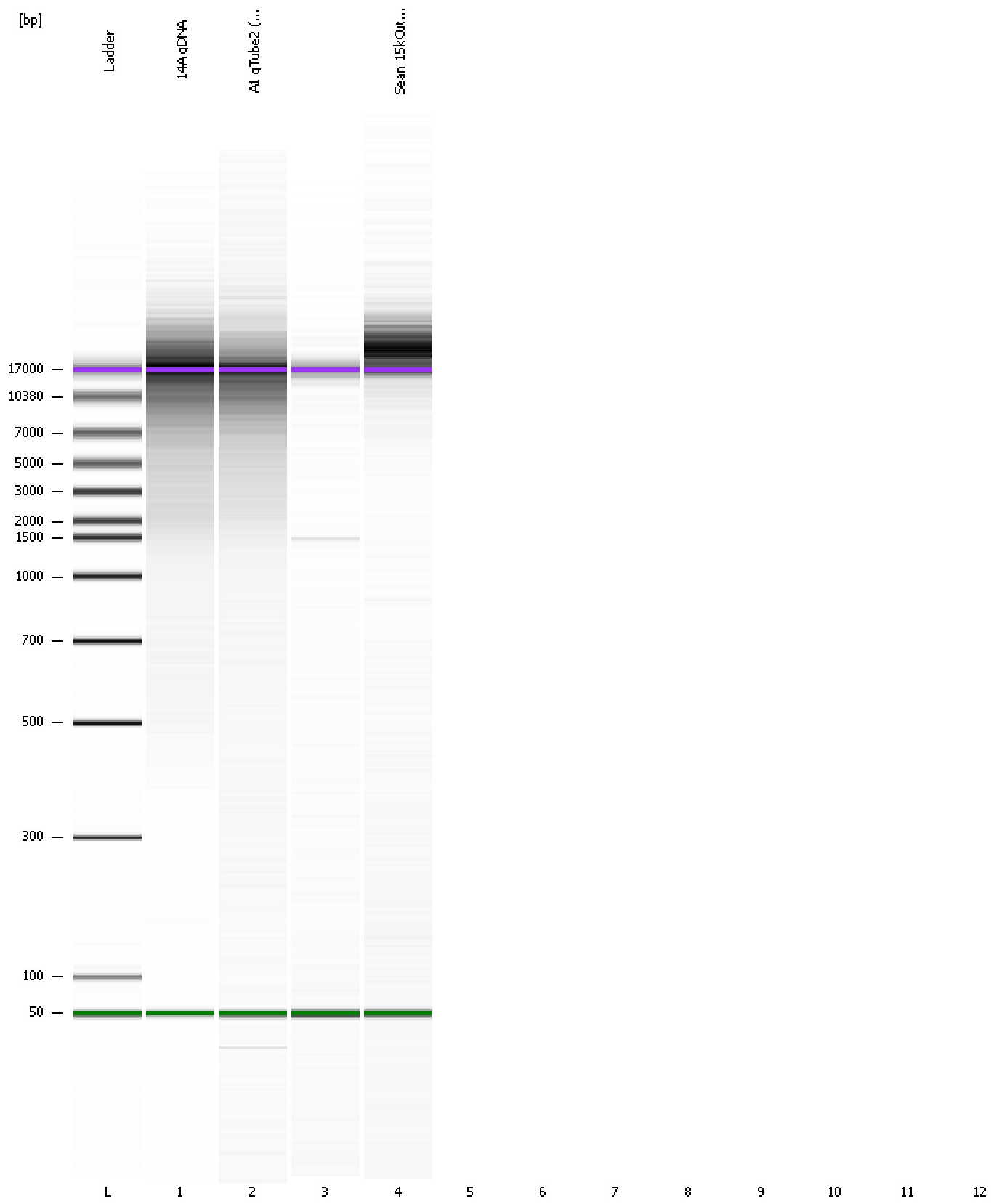


Assay Class: DNA 12000  
Data Path: C:\...ents and Settings\Bioanalyzer\2014-10-17\2014-10-17\_002.xad

Created: 10/17/2014 12:21:20 PM  
Modified: 10/17/2014 12:42:30 PM

### Curves

#### Standard Curve

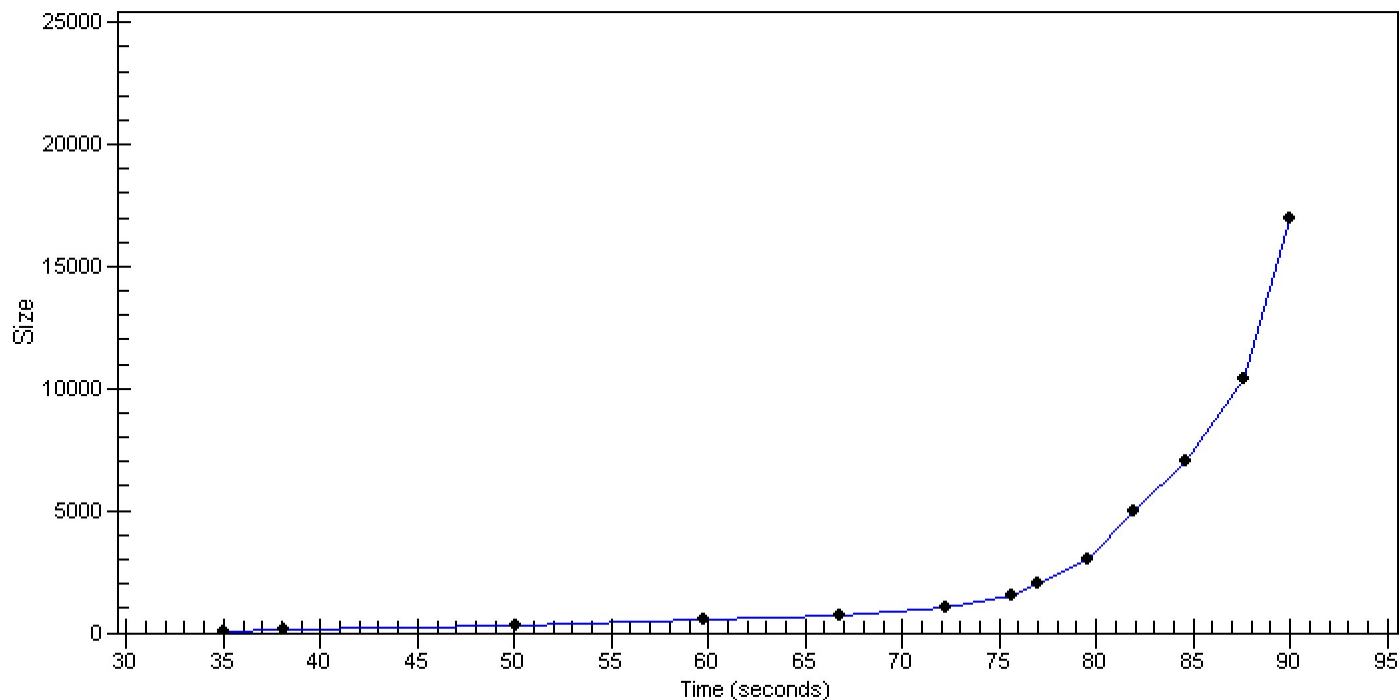

Assay Class: DNA 12000  
Data Path: C:\...ents and Settings\Bioanalyzer\2014-10-17\2014-10-17\_002.xad

Created: 10/17/2014 12:21:20 PM  
Modified: 10/17/2014 12:42:30 PM

**Invalid Samples**

Sample 5 has not been run, no results available.

Sample 6 has not been run, no results available.

Sample 7 has not been run, no results available.

Sample 8 has not been run, no results available.

Sample 9 has not been run, no results available.

Sample 10 has not been run, no results available.

Sample 11 has not been run, no results available.

Sample 12 has not been run, no results available.

Assay Class: DNA 12000  
 Data Path: C:\...ents and Settings\Bioanalyzer\2014-10-17\2014-10-17\_002.xad

Created: 10/17/2014 12:21:20 PM  
 Modified: 10/17/2014 12:42:30 PM

**Run Logbook**

Description	Number	Source	Category	Sub Category	Time	Time Zone	User	Host
Run ended on port 1 (Number of wells acquired: 5)		Instrument	Run		10/17/2014 12:39:21 PM	(GMT --07:00) Pacific Standard Time	UC Davis	D8XSMGH1
Run started on port 1 (File: C:\Documents and Settings\Bioanalyzer\2014-10-17\2014-10-17_002.xad)		Instrument	Run		10/17/2014 12:21:25 PM	(GMT --07:00) Pacific Standard Time	UC Davis	D8XSMGH1
Product Number : G2938B		Instrument	Run		10/17/2014 12:21:25 PM	(GMT --07:00) Pacific Standard Time	UC Davis	D8XSMGH1
Name :		Instrument	Run		10/17/2014 12:21:25 PM	(GMT --07:00) Pacific Standard Time	UC Davis	D8XSMGH1
Vendor : Agilent Technologies		Instrument	Run		10/17/2014 12:21:25 PM	(GMT --07:00) Pacific Standard Time	UC Davis	D8XSMGH1
Serial# : DE13701086		Instrument	Run		10/17/2014 12:21:25 PM	(GMT --07:00) Pacific Standard Time	UC Davis	D8XSMGH1
Firmware : C.01.069		Instrument	Run		10/17/2014 12:21:25 PM	(GMT --07:00) Pacific Standard Time	UC Davis	D8XSMGH1
Cartridge : Electrode		Instrument	Run		10/17/2014 12:21:25 PM	(GMT --07:00) Pacific Standard Time	UC Davis	D8XSMGH1